


# Runthrough® NS

PTCA Guidewire



Runthrough® NS is a PTCA guidewire intended to be used in interventional cardiology procedures.

Runthrough® NS is an innovative guidewire designed for workhorse lesions, as well as challenging lesions.

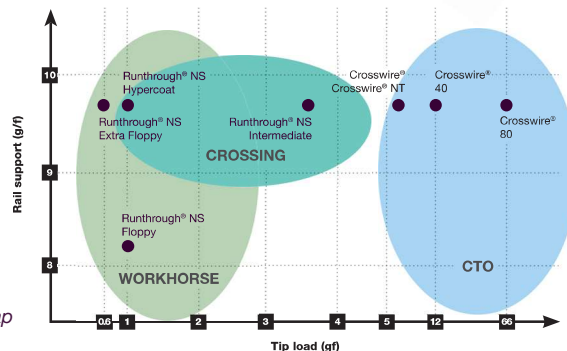

Guidewire segmentation map

## Product Characteristics

- Nitinol alloy distal core offers kink resistance and durability
- High strength stainless steel shaft designed to offer enhanced steerability and trackability
- Innovative DuoCore combines a great torque response with a durable tip
- Durable Terumo M Coat™ hydrophilic coating enables smooth manipulation
- The distal tip presents 2 mm silicone coating for an improved tactile feed-back

## General Specifications

Diameter	0.014" (0.36 mm)
Wire length	180 cm
Radiopaque length	3 cm
Distal shape	Straight
Distal coating	Hydrophilic
Proximal coating	PTFE

## Item Specifications

Item	Clinical application	Tip load (gf)	Tip shaping part	Rail support	Coating	Distal tip coating	Item Reference Code
Runthrough® NS Extra Floppy	Workhorse	0.6 gf	10 mm	Stiffer	Hydrophilic	2 mm silicone	TW-AS418XA
Runthrough® NS Floppy	Workhorse	1 gf	10 mm	Moderate	Hydrophilic	2 mm silicone	TW-AS418FA
Runthrough® NS Hypercoat	Crossing	1 gf	14 mm	Stiffer	Enhanced hydrophilic	-	TW-DS418FH
Runthrough® NS Intermediate	Crossing	3.6 gf	10 mm	Stiffer	Hydrophilic	2 mm silicone	TW-DS418IA
Runthrough® NS Extension Wire	Extension	—	—	—	—	—	TW-JJ415

Please quote above item reference code when placing an order