



Care for Life



www.yx-yiliao.com
www.yaxinmed.com

YaXin-Medical


SuZhou YaXin Medical Products Co.,Ltd



SuZhou YaXin Medical Products Co.,Ltd

Add : NO.12 , Zhongta Road , Mudu Town , Suzhou , China
Mail: yaxin@yx-yiliao.com
Tel: +0086-512-66613547
Fax: +0086-512-66264505
web : www.yaxinmed.com



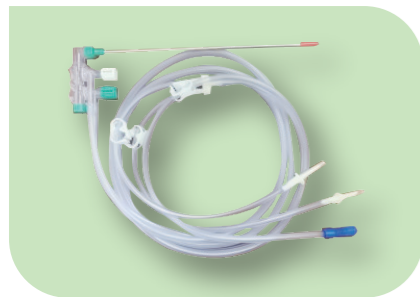

Disposable Suction Irrigation System

- With dual functions of suction and irrigation, with one-handed performance
- Single irrigation channel and double irrigation channel are available
- Pistol and Billard type are available.
- Multiple lengths(280/330/420mm) and diameters(3/5/10) are available

Item No.	Type
031001	Pistol
031002	Billiard



Ventricular External Drainage Device

Components: Silicone Ventricle Drainage Catheter, Guide Pin, Speed Regulator, Press Regulation Ball, Drainage Bag.

Features

Silicone Ventricle Drainage Catheter

- Series sizes, there are sizes for Fr8, Fr10, Fr12, Fr14, Fr16, Fr18
- No artifacts under X-ray, operate with CT guide
- There is scale from 3–15 cm
- Large side hole design to avoid obstruct
- Press regulation Ball
- Regulate the patient' s intracranial pressure by change the height of the ball.
- Top with air filter to prevent siphon action
- To avoid infection in Cranial Cavity by backflow.

Make accurate drainage with scale

Drainage Bag

- With scale to quantify the volume.
- With air filter to avoid accumulation.
- Trocar are available

Item No.	Style
031101	without adjustable plate
031102	with adjustable plate

Others



Irrigation Syringe

- Medical grade PVC material
- Capacity : 60ml
- For wound flushing

Item No.	040101
----------	--------



Irrigation Syringe

- Medical grade PVC material
- Capacity : 80ml
- Used together with gastric tube for patients to feed


Item No.	040102
----------	--------



Feeding Bag

- Medical grade PVC material
- Size: 1000ml, 1200ml
- Type: Gravity, Pump
- Rigid neck for easy filling and handing
- Bottom exit port allows complete drainage
- With flow regulator

Item No.	Size	Style
04020101	1000 ml	Gravity
04020102	1200 ml	Gravity
04020103	1000 ml	Pump
04020104	1200 ml	Pump

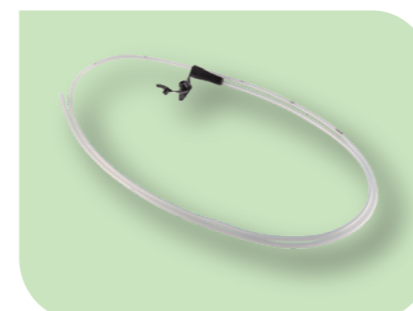

PVC Stomach Tube

- Medical grade PVC material
- X-ray opaque line throughout the length
- Size: Fr6 – Fr18
- Color-coded for visualization of size
- Length: 45cm / 100cm

Item No.	Size	Length
05040106	Fr 06	45 cm
05040108	Fr 08	45 cm
05040110	Fr 10	45 cm
05040112	Fr 12	100 cm
05040114	Fr 14	100 cm
05040116	Fr 16	100 cm
05040118	Fr 18	100 cm

Outstanding Clinical Performances Faster and Safe Postoperation Healing



Disposable Feeding Bag



YAXIN MEDICAL

苏州市亚新医疗用品有限公司 制造
Suzhou Yaxin Medical Supplies Co., Ltd.

Introduction:

Feeding bag should be considered for malnourished patients, or for those who are at the risk of malnutrition and have functional gastrointestinal tract but are unable to maintain an adequate or safe oral intake.

Usually used together with feeding tube and feeding syringe.

Features:

There are TWO main styles of feeding bag, **Pump type** and **Gravity type**.

Pump style is designed to be used with feeding pumps.


Gravity style does not require an enteral feeding pump to be used.

We have different volume bag, 500ml, 1000ml, 1200ml, 1500ml and 2000ml. The most popular size is 1000ml and 1200ml.

Details:

Disposable Feeding Bag consists of PVC bag, PVC tube, flow control valves, flow observation window and connector.

Robert clamp is available but not necessary.

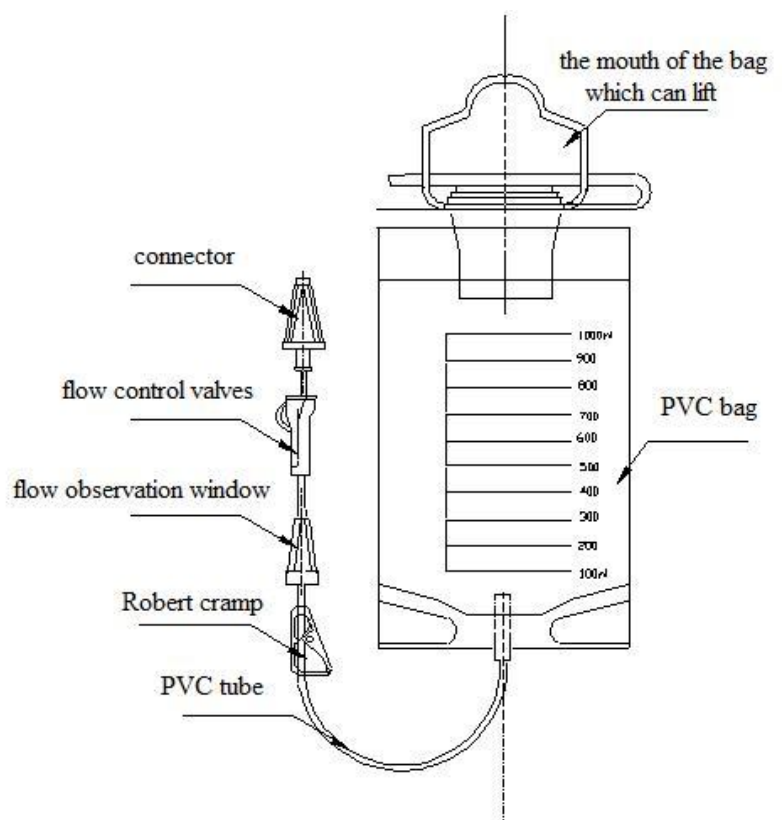

Directions:

- 1.Close control clamp completely
- 2.Fill bag with desired amount of formula
- 3.Close top cap
- 4.Hang bag or refrigerate for later use
- 5.Open control clamp. Fill entire line with fluid and adjust drops speed using the roller clamp
- 6.After a few minutes re-adjust control drops speed (If necessary)
- 7.If desired,you may fill cooling additives to ice pouch

PRECAUTIONS:

- 1.Use within 24 hours after filling.
- 2.Single patient use only
- 3.Refrigerate filled bags during storage
- 4.Not for I.V. use
- 5.Do not reuse
- 6.Use of this device is limited to those qualified in its use.

Feeding Bag Model



Size	Ref	Style
1000ml	04020101	Gravity
1200ml	04020102	Gravity
1000ml	04020103	Pump
1200ml	04020104	Pump

Caution:

1. Read instruction for use before using the device.
2. This product is sterilized for single use only, discard after single use. Please do not repeat sterilization, do not clean, do not use it once again.
3. Should be used the product immediately after opening the individual package, please destroy when it was used.
4. Do not use it if the inner package has been previously opened or damaged, or package be affected with damp, or the product went moldy.
5. Do not use after expiry date. If over the expiry date, please discard it.
6. Sterile by ETO. Single use only.

Storage:

Store in a cool, dry and dark place. Avoid direct sunlight.


Actual Photos



04020102

1200ml, Gravity

Packing Photos

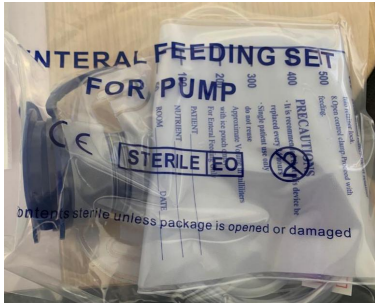

1 pc / Poly Bag, 100 bags / Carton(52*40*35cm, 9kg)



Disposable Feeding Bag



YAXIN MEDICAL

Suzhou Yaxin Medical Products Co., Ltd.

Ad:No.12, Zhongta Road, Mudu Town, Suzhou 215000, Jiangsu province, China

Mob:0086- 15250178777

E-mail:jimwu@yx-yiliao.com

Tel: 0086-512-66613547

Fax:0086-512-66264505

Web: yxmedical.en.alibaba.com, www.yx-yiliao.com