



Mosaic P 8C

Smart Demo Mosaic P 8C

Made for
iPhone | iPad | iPod

Digitally signed by Chirtoacă Ludmila
Date: 2023.12.08 01:24:24 EET
Reason: MoldSign Signature
Location: Moldova



Data Sheet

www.rexton.com

REXTON


pag. 1

Mosaic P 8C · Technical Data

Type
Output sound pressure level
OSPL 90 at 1.6 kHz
OSPL 90 (Peak)
HFA-OSPL 90
Gain
FOG at 1.6 kHz
FOG (peak)
HFA-FOG
Reference test gain
Frequency, noise and directivity
Frequency range
Equivalent input noise
Total harmonic distortion at 500 / 800 / 1600 / 3200 Hz
Tinnitus Function broadband
AI-DI
Latency
Inductive coil sensitivity
MASL (1 mA/m) at 1.6 kHz
HFA MASL (1 mA/m)
HFA SPLITS (left/right)
RSETS (left/right)
HFA SPLIV
Battery
Battery voltage
Battery current drain
Battery life (cell zinc air)
Battery life (rechargeable)
IRIL IEC 60118-13:2016 Ed. 4.0
700-960 MHz (rating)
1400-2000 MHz (rating)
2000-2700 MHz (rating)
ANSI C63.19-2011
800-950 MHz (rating)
1600-2500 MHz (rating)

Earhook		ThinTube	
			
2 ccm coupler	Ear simulator	2 ccm coupler	Ear simulator
–	138 dB SPL	–	128 dB SPL
135 dB SPL	141 dB SPL	130 dB SPL	133 dB SPL
129 dB SPL	–	118 dB SPL	–
–	76 dB	–	61 dB
77 dB	81 dB	65 dB	68 dB
70 dB	–	55 dB	–
52 dB	62 dB	41 dB	52 dB
100 - 6000 Hz	120 - 5900 Hz	100 - 5800 Hz	100 - 5900 Hz
18 dB SPL	17 dB SPL	21 dB SPL	19 dB SPL
4 / 3 / 1 / 1 %	5 / 4 / 1 / – %	2 / 2 / 1 / 1 %	2 / 2 / 2 / – %
80 dB SPL	–	80 dB SPL	–
4.0 dB		4.0 dB	
< 15 ms		< 15 ms	
–	105 dB SPL	–	93 dB SPL
99 dB SPL	–	85 dB SPL	–
109 / 109 dB SPL	–	99 / 99 dB SPL	–
–3 / –3 dB	–	–2 / –2 dB	–
111 dB SPL	–	101 dB SPL	–
1.3 V		1.3 V	
2.0 mA	2.0 mA	1.7 mA	1.7 mA
~ 120 h		~ 130 h	
–		–	
user		user	
user		user	
user		user	
M4 / T4		M4 / T4	
M4 / T4		M4 / T4	


Mosaic P 8C · Fitting Range

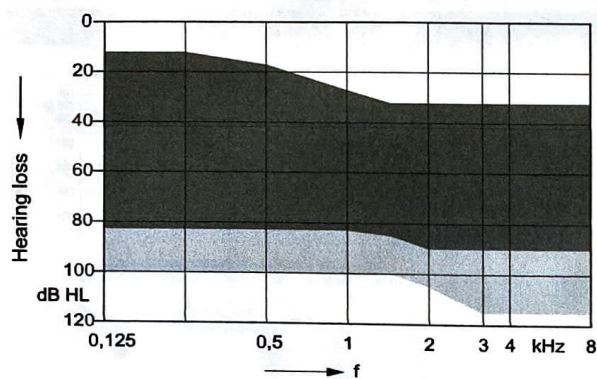

ThinTube double tip
Earhook


Earhook · Basic Data

2 ccm coupler

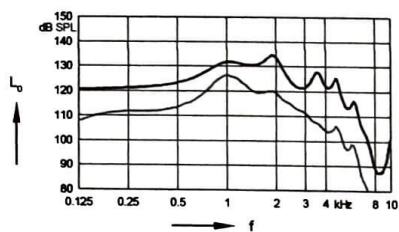


Max. Output sound pressure level
($L_i = 90$ dB)

Full on gain
($L_i = 50$ dB)

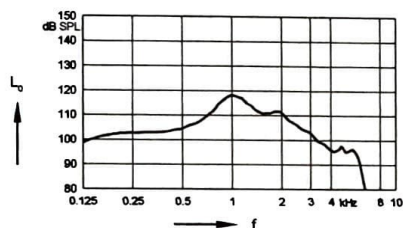


Frequency response
($L_i = 60$ dB)

Ear simulator

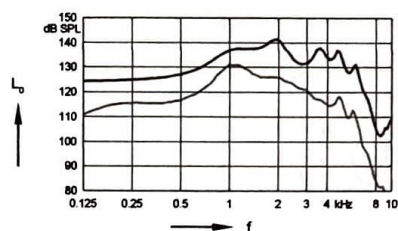


Max. Output sound pressure level
($L_i = 90$ dB)

Full on gain
($L_i = 50$ dB)

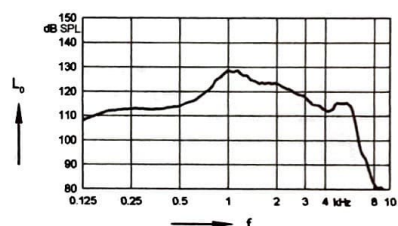


Basic acoustic response
($L_i = 60$ dB)

Inductive response

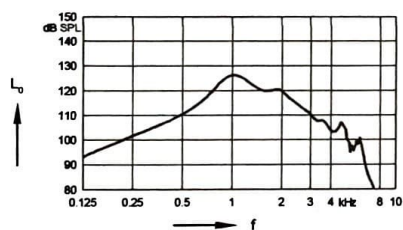


Inductive response
($H = 10$ mA/m)

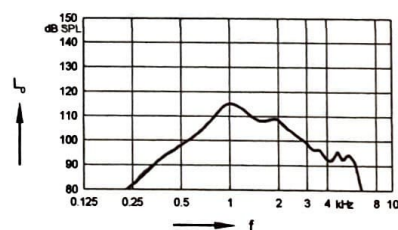

SPLITS curve left
($H = 31.6$ mA/m)



SPLITS curve right
($H = 31.6$ mA/m)

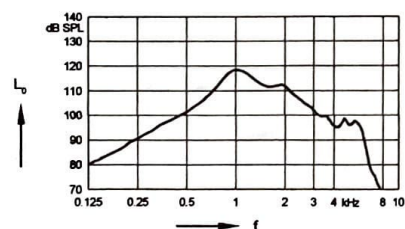



SPLIV curve
($H = 31.6$ mA/m)

ThinTube · Basic Data



2 ccm coupler

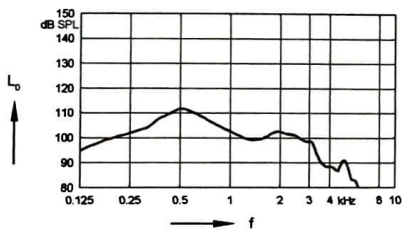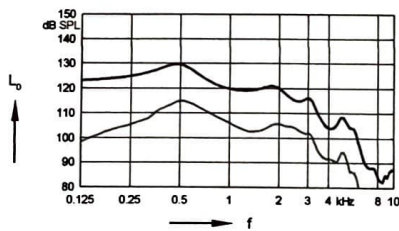


Ear simulator

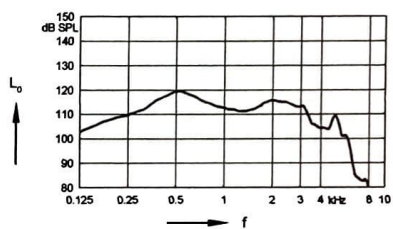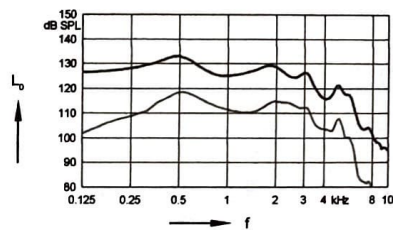

Inductive response

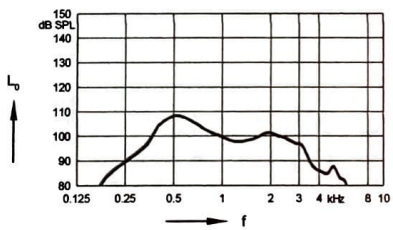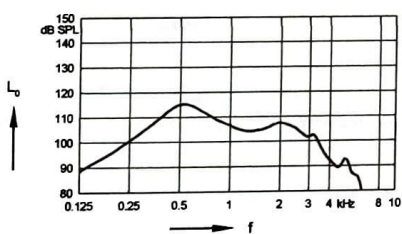

SPLITS curve right
($H = 31.6$ mA/m)

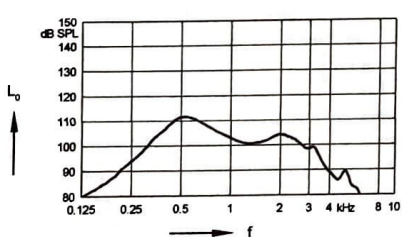

Mosaic P 8C · Features and Accessories

MyCore Platform	80 8C	60 8C	40 8C	30 8C	20 8C
Signal processing (channels) / Gain/MPO (handles)	48 / 20	32 / 16	24 / 12	16 / 8	16 / 8
Hearing programs	6	6	6	4	4
Direct Audio Streaming / Made for iPhone	●	●	●	●	●
My Voice ¹⁾	●	●	●	—	—
Wireless Sync ¹⁾	●	●	●	●	●
Volume and control coupling ¹⁾	●	●	●	●	●
MyCore Speech					
HD Bandwidth (up to 10 kHz)	●	—	—	—	—
iFocus 360 ¹⁾	automatic	manual	—	—	—
Focus 360	●	●	—	—	—
HD Directionality	●	●	—	—	—
Stereo iLock ¹⁾	●	—	—	—	—
Directional iLock ¹⁾	Premium	Premium	Standard	—	—
Voice Ranger	●	●	●	●	●
XPhone ¹⁾	●	●	●	—	—
Multichannel Adaptive Directional Microphone	●	●	●	●	●
Automatic Directional Microphone	●	●	●	●	●
Fixed Directional Microphone	●	●	●	●	●
Bandwidth Compression	●	●	●	●	●
Intelligent Feedback Preventer	●	●	●	●	●
MyCore Sound Quality and Comfort					
Dynamic Extender	●	●	●	●	●
Auto Volume ²⁾	●	●	●	●	●
Microphone-pattern adjustment ^{1) 3)}	Premium	High	—	—	—
Reverb Reducer	●	—	—	—	—
Music Enhancer	Premium	High	—	—	—
iOmni	●	●	●	●	●
Sound Smoothing (settings)	3	3	on/off	on/off	—
Intelligent Wind Noise Cancellation ¹⁾	Premium	High	—	—	—
Wind Noise Cancellation	●	●	●	●	—
Noise Management	●	●	●	●	●
Tinnitus Sound Function ⁴⁾ (presets)	6	6	6	1	—
Tinnitus Notch Function ⁴⁾	●	●	●	—	—
MyCore Automatic Optimization					
Smart Automatic Equalizer	Premium	Premium	High	—	—
Smart Automatic Acclimatization	Premium	Premium	High	High	Standard
Automatic Classifier	●	●	●	●	●
Data Logging	●	●	●	●	●

¹⁾ Bilateral fitting required

²⁾ Streaming only

³⁾ requires Connexx Smart Direct App

⁴⁾ Availability is country-dependent

● available — not available

Performance levels: Premium High Standard

pag. 6

Mosaic P 8C · Features and Accessories

Style specific features	80 8C / 60 8C / 40 8C	30 8C / 20 8C
SecureTec protection	IP68	IP68
Charging contacts	—	—
Battery Size	13	13
Battery door on/off function	●	●
Nanocoated housing	●	●
Wireless programming	●	●
Instrument configurations		
Flat cover	—	—
Rotary volume control	—	—
Push button	—	—
Rocker switch	●	●
Color conversion kit	○	○
Battery door - integrated telecoil	○	○
Battery door - child lock	—	—
Small earhook	○	○
Programming accessories		
ConnexxAir / ConnexxLink	— / —	— / —
Noahlink Wireless	●	●
Programming adapter / cable	size 13	size 13
Accessories		
Connexx Smart Key	○	○
CROS inoX 8C	—	—
CROS RIC 8C	○	—
CROS Li RIC 8C	—	—
Connexx Smart Transmitter 2,4	○	○
Connexx Smart Mic	○	○
Apps		
Connexx Smart Direct App	○	○

● available ○ optional — not available

Mosaic P 8C

Abbreviations and Standards

Abbreviations

The following abbreviations are used in this datasheet:

OSPL	Output Sound Pressure Level
HFA	High Frequency Average
FOG	Full-On Gain
MASL	Magneto Acoustical Sensitivity Level
SPLITS	Coupler SPL for an Inductive Telephone Simulator
RSETS	Relative Equivalent Telephone Sensitivity
SPLIV	SPL In a Vertical magnetic field
AI-DI	Articulation Index - Directivity Index
IRIL	Input Related Interference Level
RTF	Reference Test Frequency

Standards

- ▶ All measurements with the 2 ccm coupler were performed according to ANSI S3.22-2014 and IEC 60118-0:2015 if applicable.
- ▶ All measurements with an ear simulator were performed according to IEC 118-0/A1:1994 and to DIN 45605 (frequency range) if applicable.
- ▶ Curves and figures representing FOG are measured with 20 dB reduction and 70 dB SPL input level.
- ▶ Figures representing Equivalent Input Noise incorporate a moderate expansion.
- ▶ Inductive coil sensitivity values, inductive response curves and T ratings apply for instruments with telecoil battery door only.
- ▶ Tinnitus noiser measurement conditions: all tinnitus single frequency sliders in max position, master volume slider in default position (0 dB) and local volume control in default position.
- ▶ The current consumption is measured in reference test setting (RTS) according to the applicable standards. Due to the settling behaviour of hearing instruments supporting RF (radio frequency), the battery current is measured 3 minutes after turning on (note: no pairing).
- ▶ The battery life is based on first fit settings using 60% of the fitting range and an ISTS (International Speech Test Signal) input signal at 65 dB SPL (note: pairing established). The actual battery life is determined by battery quality, hearing loss, sound environment, usage and activated feature set.
- ▶ The following acoustic connections / ear pieces were used:
 - Earhook
 - ThinTube
- ▶ HD Bandwidth up to 10 kHz for 80 8C devices only.

Made for

iPhone | iPad | iPod

"Made for iPod", "Made for iPhone", and "Made for iPad" mean that an electronic accessory has been designed to connect specifically to iPod, iPhone, or iPad, respectively, and has been certified by the developer to meet Apple performance standards. Apple is not responsible for the operation of this device or its compliance with safety and regulatory standards. Please note that the use of this accessory with iPod, iPhone, or iPad may affect wireless performance.

The information in this document contains general descriptions of the technical options available, which do not always have to be present in individual cases and are subject to change without prior notice. The required features should therefore be specified in each individual case at the time of conclusion of the respective contract.

pro-P