Content

Dental implantation system	01-04
Endo motors	05-14
Endo Ultrasonic Activator	15-16
Endodontic obturation system	17-20
Root apex locators	21-23
Pulp tester, Caries detector devices	24
Niti file system	25-30
Gutta cuts	31
Interproximal enamel kit	32
Operating microscope	33-34
Dental electrical motors	35-42
LED curing lights	43-49
Teeth whitening accelerator	50
Dental alginate mixers, amalgamators	51-52
Air abrasion system with spray	53-54
Accessory of electric equipment	55-57
4 way spray handpiece	F0 C4
4 way spray nanupiece	58-61
Led handpiece with generator	62-63
Led handpiece with generator	62-63
Led handpiece with generator Optical handpiece	62-63 64-65
Led handpiece with generator Optical handpiece 1way spray handpiece	62-63 64-65 66-67
Led handpiece with generator Optical handpiece 1way spray handpiece 3way spray handpiece	62-63 64-65 66-67 68-69
Led handpiece with generator Optical handpiece 1way spray handpiece 3way spray handpiece Color handpiece	62-63 64-65 66-67 68-69 70-71
Led handpiece with generator Optical handpiece 1way spray handpiece 3way spray handpiece Color handpiece Nature handpiece	62-63 64-65 66-67 68-69 70-71 72-73
Led handpiece with generator Optical handpiece 1way spray handpiece 3way spray handpiece Color handpiece Nature handpiece Nature handpiece	62-63 64-65 66-67 68-69 70-71 72-73 76-77
Led handpiece with generator Optical handpiece 1way spray handpiece 3way spray handpiece Color handpiece Nature handpiece Nature handpiece Speed increased contra angle	62-63 64-65 66-67 68-69 70-71 72-73 76-77 78-83
Led handpiece with generator Optical handpiece 1way spray handpiece 3way spray handpiece Color handpiece Nature handpiece Nature handpiece Speed increased contra angle Reduction contra angle	62-63 64-65 66-67 68-69 70-71 72-73 76-77 78-83 84-85
Led handpiece with generator Optical handpiece 1way spray handpiece 3way spray handpiece Color handpiece Nature handpiece Nature handpiece Speed increased contra angle Reduction contra angle Low speed handpiece set Contra angle handpiece, air motors,	62-63 64-65 66-67 68-69 70-71 72-73 76-77 78-83 84-85 86-89
Led handpiece with generator Optical handpiece 1way spray handpiece 3way spray handpiece Color handpiece Nature handpiece Nature handpiece Speed increased contra angle Reduction contra angle Low speed handpiece set Contra angle handpiece, air motors, stright handpiece	62-63 64-65 66-67 68-69 70-71 72-73 76-77 78-83 84-85 86-89 90-93
Led handpiece with generator Optical handpiece 1way spray handpiece 3way spray handpiece Color handpiece Nature handpiece Nature handpiece Speed increased contra angle Reduction contra angle Low speed handpiece set Contra angle handpiece, air motors, stright handpiece Economical complete kit	62-63 64-65 66-67 68-69 70-71 72-73 76-77 78-83 84-85 86-89 90-93 94





C-SAILOR PRO Dental implantation system Intuitive Powerful Safe & Reliable





Patent design, delicate contracted

This product has been applied layout-design patent, Counterfeit must be investigated.



The large size screen ensures that the treatment parameters are visible at any time



Fiber optic contra angle



Built-in water supply system 4 shifts water supply up to 150ml/mir



-Full touch 7 inch screen-Large size screen ensures that the treatment parameters are visible at any time



Full function pedal controller Aseptic operation can be achieved giving the maximum freedom of the handle





C-SMART-I PILOT Endo motor

Large screen, but also wireless Simple, but also rich functional Split Endo motor & Apex locator





Order name: C-SMART-I PILOT

Technical:

Adaptor working voltage: Input: AC100~240V, 50Hz Output: DC10V, 1.5A

Rotation speed range: 150~1000rpm

Torque: 0.6~3.9N.cm

File system: 6 mainstream brands 30 endo file systems


200+ file parameters

Volume: 24x24x7.5cm

Weight: 1.5kg

Packing material: carton

QTY per master carton: 10 pc/carton

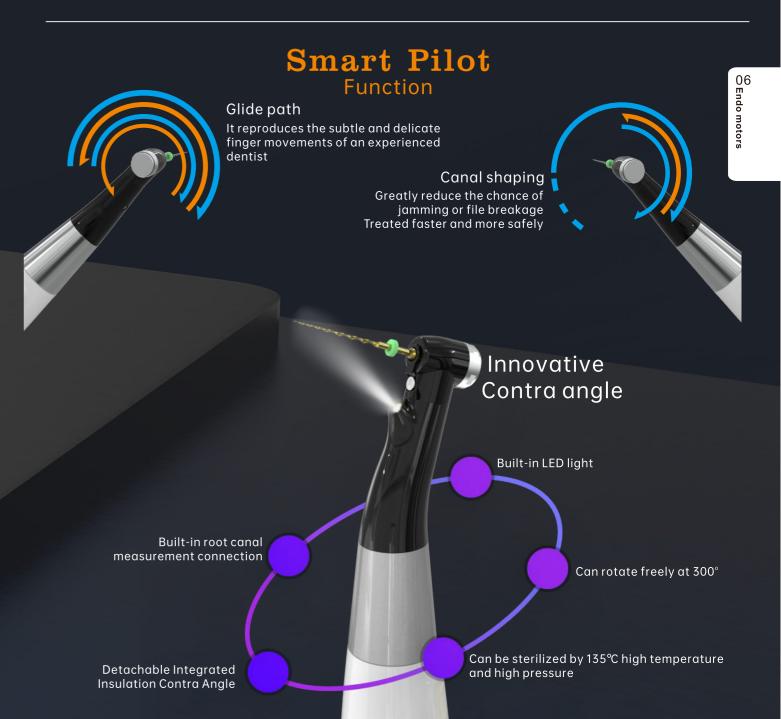

With Apex locator function

Model: C-SMART-I PRO Endo motors




Order name: C-SMART-I PRO Technical: AC100-240V,50/60HZ; Adaptor working voltage Adjustable torque,range 0.6-5.2N.cm(6-52mN.m) 5 setting for adjustable 150rpm-650rpm rotational speed ranging: 23.5X20X14.5 Volume(cm) Weight(kg)/PCS Packing material carton QTY per master carton 6 pc/ctn

Built-in mainstream system



Endo motor & Apex locator Dual mode

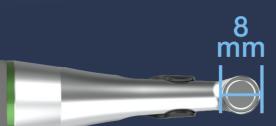

16:1 mini head contra angle











Wireless endo motor & apex locator Powerful functions Precise measurements Efficient and safe





Order name: C-SMART-MINI AP Technical: AC100-240V 50/60HZ; Adaptor working voltage Adjustable torque,range 0.6-3.9N.cm(6-39mN.m) 5 setting for adjustable 150-600rpm rotational speed ranging: 19X13X6 Volume(cm) Weight(kg)/PCS Packing material carton QTY per master carton 12 pc/ctn





With 6 mainstream brands, 31 muilt file system Almost 200 Endo files parameters



16:1 Mini Head Design 360° Rotation, Meet the needs of various treatment angles



3 working mode Auto-switch work mode Colorful and sound prompt







Non-contact recharge Faster and safer












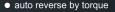


• auto reverse by torque • auto reverse by apical position


Order name: C-SMART-MINI

Technical:

 Adaptor working voltage	100-240Vac,50/60HZ;
 Adjustable torque,range	0.6-3.9N.cm(6-39mN.m)
 5 setting for adjustable rotational speed ranging:	125rpm-625rpm
Volume(cm)	19X13X6
Weight(kg)/PCS	0.7
Packing material	carton
QTY per master carton	12 pc/ctn



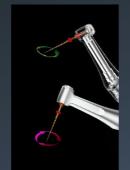

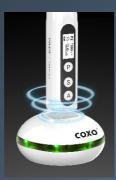







Technical:

Adaptor working voltage	AC100-240V,50/60HZ;
Adjustable torque,range	0.6-3.9N.cm(6-39mN.m)
5 setting for adjustable rotational speed ranging:	125rpm-625rpm
Volume(cm)	19X13X6
Weight(kg)/PCS	0.7
Packing material	carton
QTY per master carton	12 pc/ctn

Order name: C-SMART-MINI (LED)















C-SMART-I Endo motor Comprehensive upgrade Both apex locate and endo at the same time



Endo motors Model: C-SMART-I



Order name: C-SMART-I


Feature:


- 1- Detachable 1:1 mini contra angle, multi-angle
- 2-Color wide LCD screen
- 3-Original switzerland electrical-motor
- 4-Rotation speed and torque can be programmed
- 5-Dual-frequency apex locator
- 6-4 working modes
- 7-Auto apical reverse
- 8- Auto apical slow down
- 9-Reciprocating cutting mode
- 10-High capacity Li-ion chargeable battery
- 11-Auto power off and memory function

Technical:

AC100-240V,50/60HZ;
0.3-5.2N.cm(3-52mN.m)
125rpm-625rpm
24X21X9.5
1.8
carton
8 pc/ctn


This product has been applied layout-design patent, Counterfeit must be investigated; NO.ZL201230082365.5

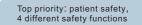

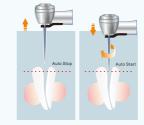


Auto Start/Stop

The file starts rotating auto-matically when inserted into the canal, and stops automati-

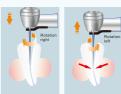

Auto Torque Reverse

The file stops rotating automatically if the torque exceeds the specified limit of one of 11 torque levels.

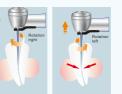

Auto Apical Reverse



When the file reaches at a specified measurement. The file rotation stops and reverses automatically.



Auto Apical Slow Down

When the file is approaching to apical foramen, the motor auto-matically slows down for safer trearment.

4 working modes















ULTRA SMART Endo Ultrasonic Activator Safer Cleaner Professional



Model: Ultra Smart

Endo Ultrasonic Activator

Order name: ULTRA-SMART

Feature:

- 1. High frequency ultrasound up to 38kHz.
- 2. Variable frequency control, self adjusting frequency in treatment.
- 3. Activation sequences of 20 sec.
- 4. With LED lighting system.

Technical:

Volume: 17x17x6 cm Weight: 0.65 kg

Packing material: carton


QTY per master carton: 8 pcs/carton









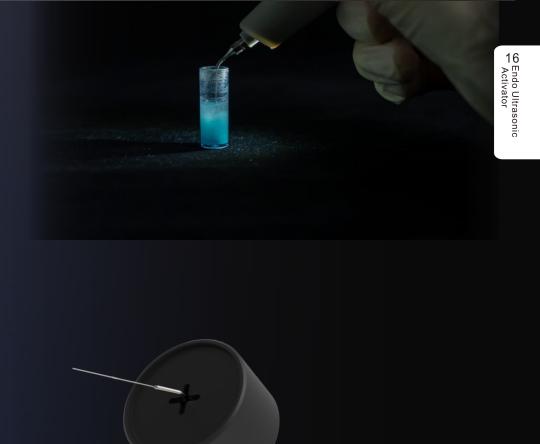

Force-limiting wrench
Safer and more convienient

Working temperature: 150 ° C, 180 ° C, 200 ° C, 230 ° C

Adapter: input AC 100~240V, output DC5V 1.5A

Use time 1.5H

Needle specifications: 23G, 25G



New wrench design Help you to push more accurately



Accurate temperature control
Automatically cool down to 150°C on the base
Save electricity and prevent glue overflow

Order name: C-FILL Mini P type

Battery: DC3.7V 2000mAH

Heating time: 4S(With heat protection function, the single heating time is up to 4 S)

Working temperature: 150 ° C, 180 ° C, 200 ° C, 230 ° C Adapter: input AC 100~240V, output DC5V 1.5A

Use time 1.5H

Plunger specifications: F 45/02, FM 50/04, M 55/06, ML 55/08

Different angle Easy to reach any treatment places





Intelligent reminder Sound and flashing prompt With protection function, single heating time up to 4s



Precise temperature control Good performance give you the most smooth cutting of gutta-percha



Needle		
Gauge	Diameter	Length
23ga	0.65mm	23mm
23ga	0.65mm	23mm
25ga	0.55mm	23mm
25ga	0.55mm	23mm

Heat plunger
Size | Color code | Taper/Size

#XF Green 0.02#35

#F White 0.02#45

0.04#50

0.06#55

0.08#55

#FM

#M #ML





C-FILL MINI Endodontic Obturation system

Precise temperature control

Perfect filling

Order name: C-FILL Mini



The tip of plunger with minimum diameter 0.45mm, 02 taper, can press the effective filling position that is 3 mm from the Apex



The gun head can be rotated by ±100° for e asy use.



Double-sided button operation, convenient and fast



360 ° dynamic charging display, more humane.



Model: C-FILL



Order name: C-FILL apack



Order name: C-FILL β back



Size	Color code	Taper/Size
#F	White	0.02 #45
#FM	Yellow	0.04 #50
#M	Red	0.06#55
#ML	Red	0.08 #55

Heat plunger



Needle		
Gauge	Diameter	Length
23ga	0.65mm	23mm
23ga	0.65mm	23mm
25ga	0.55mm	23mm
25ga	0.55mm	23mm

Endodontic obturation system

Feature:

- 1-Wireless design, reduce wrist fatigue
- 2-Holding comfortable, easy to use
- 3-Heat plunger has easy –lock structure. and can rotate 360 $^{\circ}$
- 4-Can be choose heating method according type of heat plunger.
- 5-Fast heating, quick operation
- 6-Large capacity battery, double battery charging and spare.
- 1-Wireless handpiece, easy to use
- 2-Holding and operating very comfortable
- 3-Exactly temperature-control and thermal protector cap can be avoiding scald
- 4-Gun needle can be 360 rotation, more easy to get the fill position
- 5-Gun needle screw design, effectively prevent the glue;
- 6-It can be at 200 degree in 30s, can be melting all kinds of autta percha
- 7-Large capacity for battery ,which can use long time

Technical:

Li-ion chargeable Battery:		DC3.7V 2200mAh	DC3.7V 2200mAh
Heating time:		5s	30s to200°C
Working temperature:		150°C,180°C, 200°C,230°C	160°C,180°C, 200°C,230°C
Adaptor	Input:	AC100-240V	AC100-240V
, aupioi	Output:	DV5V,1.5A	DV5V,1.5A
Heat plunger		F,FM,M,ML	23G,25G
Operating time		continuous	for 1.5A
Volume(cm)		20.5X18	X10
Weight(kg)/PCS		2.0	
Packing material		carton	
QTY per master carton		10 pc/cti	1



#80 Stainless steel #100 Stainless steel
Light and practical,
They are used in conjunction with
the obturation system

#50 NiTi


Order name: CX265-59




Order name: C-ROOT I

Technical:

Adaptor working voltage:	AC100-240V,50/60HZ;
Battery:	DC 7.4V 1800mAh
	changeable Li-ion battery
Display:	color TFT-LCD touch screen
Display specifications:	size:73.44x48.96mm(3.5inch)
Display specifications.	resolution:480x320
Charge time:	3hours
Volume(cm):	16X14X6.5
Weight(kg)/PCS	0.7
Packing material	carton
QTY per master carton	12 pc/ctn

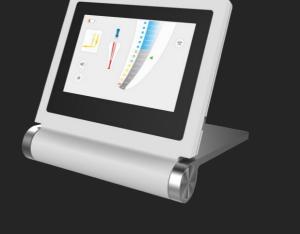

Colorful touch screen
Easy to adjust various parameters

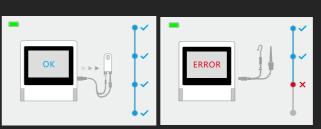

Automatic check function
Ensure you the accuracy of each measurement



Simple and light
The thickness is only 2cm, but it
brings unimaginable accuracy

Pulp tester



Feature:

- 1-Root cannal apex finder and pulp vital measurement;
- 2-Color LCD screen display and can adjust the angle
- 3-Multi-frequency design, exact accurate;
- 4-Accurate current-stabilized generate technology
- 5-Self-test function in apex locator;
- 6-The narrow zone of the apex can be adjust if the doctor need;
- 7-Sound prompt for the apex in the narrow zone of the pulp vital:
- 8-Rechargeable Li-ion battery and charger

Technical


- 1.Adaptor voltage: AC100-240V, 50/60HZ
- 2.Battery:7.4V/2000mAh changeable Li-ion battery
- 3.Charge time: 3 hours
- 4. Display: color LCD Screen, buzzer

Volume/Weight


- 1.Volume(cm):24X19X9 2.Weight(KG)/PCS:2.0KG
- 3. Packing material: carton
- 4.QTY per master carton: 6 pc/ctn

Two working modes

Root apex locator with pulp tester

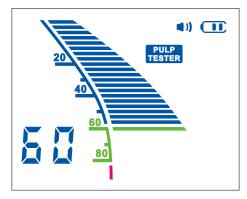

Pulp tester mode

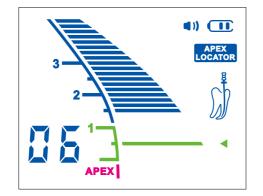

Apex locator mode

Accessories



Order name: C-PULSE

Feature:

- 1-Preset speed mode(high-mid-low speed)
- 2-Agentle, pulsed stimulus begins to increase at a rate specifically designed for patient comfort.
- 3-If the patient indicates perception, simly release the button, the stimulus stops immediately, but the numbers remain frozen on the screen for about 3 minutes.
- 4-Specifically designed for patient comfort.
- 5-Convenient to operate.
- 6-Maximum stimulus -- 80. between 0-40, the patient feels ache and anesthesia, alive teeth nerve.
- 7-Between 40-80, with above-mentioned reaction, part of teeth nerve dead; 80, no above-mentioned reaction, dead
- 8-The device turns itself off in three minutes after operation (No wasted batteries).

Caries detection device

This product has been applied layout-design patent, Counterfeit must be investigated; NO.ZL200730333925.9

Model: C-HUNTER

Model: C-PULSE

Feature: 1-Effective function for decayed tooth- detection, quickly, easily, painless 2-Wireless operation, Lightness using 3-Support minimally invasive surgery, safely, lighten patient's fear

Technical:

treatment operation

Adaptor working voltage:	AC100-240V,50/60HZ;
Lightwave Data:	LED intensities:40-280mw
	LED wavelength:380-420nm
Volume(cm):	16X14X6.5
Weight(kg)/PCS	0.5
Packing material	carton
QTY per master carton	6 pc/ctn

4-Mini slender mainframe designing, provide the best

Innovative caries detection Competitive with the products of

Simple, Effctive,

Order name: C-HUNTER

Accessories





Adapter

Function:

Bacteria in dentine infected by caries leave behing metabolic products(porphyrins).when an exposed cavity is illuminated with violet light, these porphyrins show a red fluorescence, thereby indicating the significant areas for caries excavation.

C-hunter aids detection and monitoring of carious tooth substance, in particular in dentin, during the excavation of previously opened cavities.

According to the particular fluorescence properties of the degradation product which is produced by the cariogenic bacteria, substances with teeth for removing infected, to remove the infected material on the teeth.