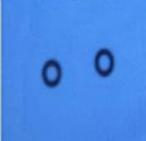


RENO-S130D PM (Preventive Maintenance) KIT

PM (Preventive Maintenance) should be performed replacing consumable parts and regular inspection according to PM Schedule for maintaining RENO Plasma sterilizer in good working order. PM parts has important role to extend RENO Plasma sterilizer's life and prevent Error occurring for customer satisfaction. RENOSEM Co., Ltd. supply PM Kit 1,2,3 set for performing proper replacement by the detail below.

S130D PM KIT-1 (Every 600 cycles)						
No.	Part	Picture	Q'ty	Spec.	Cycle	Year
	PM06-0010-00 (Option 1)					
	PM06-0020-00 (Option 2)					
1-1	Rotary Pump Oil (Laybold)		1	1 Liter	600	0.5
1-2	Rotary Pump Oil (WooSung)		3	1 Liter	600	0.5
2	Activated Carbon Filter		1		600	0.5
3-1	O-ring (Inlet&Outlet) (Laybold)		1		600	0.5
3-2	O-ring (Inlet&Outlet) (WooSung)		1		600	0.5
4	Oil Paper & Glove		1		600	0.5
5	Oil container		1	3 Liter	-	-
6	Funnel		1		-	-
7	Garbage Vinyl		1		-	-

※ **Remark**
Replace PM KIT-1 every 6 month or 600 cycles, whichever comes first.

※ **Cleaning**
Injection port hole cleaning should be performed every 600 cycles.

S130D PM KIT-2 (Every 1,200 cycles)							※ Remark Replace PM KIT-2 every 1 year or 1,200 cycles, whichever comes first. ※ Cleaning Injection ball inside of SFS Nozzle cleaning should be performed every 1,200 cycles.
No.	Part	Picture	Q'ty	Spec.	Cycle	Year	
	PM06-0030-00 (Option 1)						
	PM06-0040-00 (Option 2)						
	PM06-0060-00 (Option 3)						
1-1	Option 1 Oil Mist Trap (Laybold)		1		1,200	1.0	
1-2	Option 2 Oil Mist Trap (WooSung)		1	SUS	1,200	1.0	
1-3	Option 3 Oil Mist Trap-Inside Filter Set (WooSung)		1	Set	1,200	1.0	
2	Heat-registing O-ring		1		1,200	1.0	
3	NW25 Clamp		1		-	-	
S130D PM KIT-3 (Every 2,400 cycles)							※ Remark Replace PM KIT-3 every 2 years or 2,400 cycles, whichever comes first.
No.	Part	Picture	Q'ty	Spec.	Cycle	Year	
	PM06-0050-00						
1	Catalyst Filter		1		2,400	2.0	
2	Injection Valve		2	Shintec	2,400	2.0	
3	HEPA Filter		1		2,400	2.0	
4	Door Gasket		2	Hardness 40	2,400	2.0	

※ Notice

PM schedule can be different consumable parts and spare parts depending on the Operating condition of RENO Plasma Sterilizer.

S130D

No.	Model	New Part No.	Description	Qty	Unit	IMAGE
44	S130D	S011-0020-00	Regulator(0.5~7 bar)	1	ea	
45	S130D	S011-0050-00	Regulator (0.2~2bar)	1	ea	
46	S130D	S004-0020-00	CYLINDER	1	ea	
47	S130D	S004-0010-00	Cylinder(H2O2)	1	ea	
48	S130D	E030-0050-00	Niddle Cylinder Detection	1	ea	
49	S130D	E030-0060-00	Pressure Gauge for compressor	1	ea	
50	S130D	S060-0011-00	NEEDLE ASS'Y	1	ea	
51	S130D	E001-0010-00	Circuit Breaker (MCCB)	1	ea	
52	S130D	E009-0040-00	Data Logger	1	ea	
53	S130D	S040-0010-00	CASTER(LOCK TYPE)	1	ea	
54	S130D	S015-0010-00	VACUUM CENTERING	1	ea	
55	S130D	S015-0040-00	Vacuum Centering (NW40)	1	ea	
56	S130D	S015-0031-00	Centering (Metal, NW25)	1	ea	
57	S130D	S014-0010-00	VACUUM CLAMP (NW25)	1	ea	
58	S130D	S014-0020-00	VACUUM CLAMP (NW40)	1	ea	

Lot 120