


# Chemical Indication Sticker

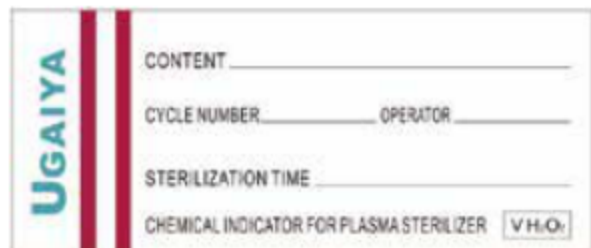

The image shows a template for a chemical indication sticker. On the left, there is a vertical red bar with the word "UGAIYA" written vertically in green. To the right of the bar, there are four horizontal lines for text entry, each preceded by a label: "CONTENT", "CYCLE NUMBER", "OPERATOR", and "STERILIZATION TIME". Below these lines, the text "CHEMICAL INDICATOR FOR PLASMA STERILIZER" is printed, followed by a small box containing the chemical formula  $VH_2O_2$ .

Chemical indication sticker is PET with 2cm width and 50m length. It is with chemical indication ink and would turn from blue to yellow when exposing to hydrogen peroxide low temperature plasma. It could monitor if plasma reach required area and it is used for non-woven fabrics packaging.

Code	Model	Description	Package
211.003	UGCI0201	Single layer	50 Sheet/Bag
211.004	UGCI0201B	Double layer	50 Sheet/Bag