



# MULTI-STAGE DILATION BALLOONS

## THREE DIAMETERS, ONE BALLOON

The multi-stage dilation balloon from MICRO-TECH facilitates simple and reliable dilation of stenoses. Each balloon can be expanded to three different, fixed diameters for progressive dilation of stenosis in three pre-defined stages. The pressure is distributed evenly over the entire working length of the balloon at every diameter. The

rounded tips enable a good endoscopic view of the stenosis, even during dilation, ensuring controlled, safe intervention and good therapeutic success. Thanks to the preloaded stiffed wire with its flexible, hydrophilic, radiopaque tip, the balloons are readily introduced and are then rapidly deflated after successful dilation.

### SPECIFIC CHARACTERISTICS


- Three different diameters with three separate pressure levels
- 2 X-ray markers on opposite ends to indicate the exact working length of the balloon
- Rounded Balloon shoulders for clear endoscopic view of the balloon's distal end
- Conical, atraumatic catheter tip
- Preloaded stiffed wire with hydrophilic, radiopaque tip



Separate openings for stiffed wire and insufflation

## SINGLE-USE INFLATION-INSTRUMENT

Sturdy construction and an ergonomic handle for easy control. With quick-lock mechanism for one-handed deflation, even under full pressure. The clear cylinder and the easy to read scale allow for a close look at the content and the amount. The exact gauge guarantees reliable pressure measurement.

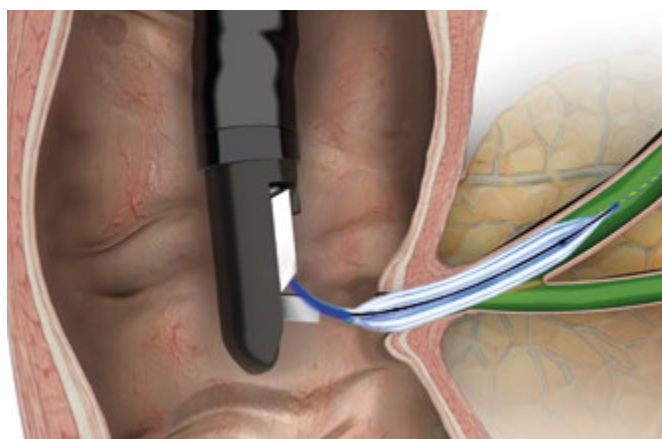

## SPECIFICATIONS

REF	Content ml air/liquid	Max. pressure ATM	Feature
<b>SINGLE-USE INFLATION-INSTRUMENT</b>			
BID-15-60-01	60	15	3-way-stopcock

REF	Balloon Ø mm	Balloon 3-stagedilation pressure (atm)	Balloon length mm	Catheter Ø Fr	Catheter length mm
<b>MULTI-STAGE DILATION BALLOONS</b>					
MBD-AS-BB-2	6 / 7 / 8	3 / 6 / 10	30	7	2300
MBD-AM-BA-2	6 / 7 / 8	3 / 6 / 10	55	7	1800
MBD-AM-BB-2	6 / 7 / 8	3 / 6 / 10	55	7	2300
MBD-AL-BA-2	6 / 7 / 8	3 / 6 / 10	80	7	1800
MBD-BS-BB-2	8 / 9 / 10	3 / 5.5 / 9	30	7	2300
MBD-BM-BA-2	8 / 9 / 10	3 / 5.5 / 9	55	7	1800
MBD-BM-BB-2	8 / 9 / 10	3 / 5.5 / 9	55	7	2300
MBD-BL-BA-2	8 / 9 / 10	3 / 5.5 / 9	80	7	1800
MBD-CS-BB-2	10 / 11 / 12	3 / 5 / 8	30	7	2300
MBD-CM-BA-2	10 / 11 / 12	3 / 5 / 8	55	7	1800
MBD-CM-BB-2	10 / 11 / 12	3 / 5 / 8	55	7	2300
MBD-CL-BA-2	10 / 11 / 12	3 / 5 / 8	80	7	1800
MBD-DS-BB-2	12 / 13.5 / 15	3 / 4.5 / 8	30	7	2300
MBD-DM-BA-2	12 / 13.5 / 15	3 / 4.5 / 8	55	7	1800
MBD-DM-BB-2	12 / 13.5 / 15	3 / 4.5 / 8	55	7	2300
MBD-DL-BA-2	12 / 13.5 / 15	3 / 4.5 / 8	80	7	1800
MBD-EM-BA-2	15 / 16.5 / 18	3 / 4.5 / 7	55	7	1800
MBD-EM-BB-2	15 / 16.5 / 18	3 / 4.5 / 7	55	7	2300
MBD-EL-BA-2	15 / 16.5 / 18	3 / 4.5 / 7	80	7	1800
MBD-FM-BA-2	18 / 19 / 20	3 / 4.5 / 6	55	7	1800
MBD-FM-BB-2	18 / 19 / 20	3 / 4.5 / 6	55	7	2300
MBD-FL-BA-2	18 / 19 / 20	3 / 4.5 / 6	80	7	1800

Packaging unit: 1 piece / suitable for guide wires with 0.035"



# DILATION BALLOONS

## SIMPLY EXPAND RELIABLY AND EFFECTIVELY

The dilation balloons from MICRO-TECH stand for maximum quality. MICRO-TECH produces all the dilation balloons itself – using state-of-the-art European and American production technology. Two x-ray markings indicate the exact working length of the balloon, thereby ensuring its

precise positioning. Under application of pressure, the balloon is rapidly filled to acquire its predefined expansion width. Furthermore, the high-pressure design ensures the balloon retains its shape over the entire working length, thereby allowing optimum dilation success.

### SPECIFIC CHARACTERISTICS

- Dimensionally stable high pressure design
- 2 X-ray markers on opposite ends to indicate the exact working length of the balloon
- Separate guide wire and insufflation lumens
- Conical shaped catheter tip
- Rapid expansion
- Constant dilation over the entire length
- Guide-wire compatible up to 0.035 inches




X-ray marking



Conical tip

## SUCCESSFUL IN PRACTICAL USE

Stenoses are effectively dilated using the MICRO-TECH dilation balloons. The balloon dilates the stenosis safely and gently and prepares the area to be treated in an optimum manner. The uniform distribution of pressure over the entire working length of the balloon and the outstanding dimensional stability during dilation guarantee excellent therapeutic success over the entire length of the stenosis.

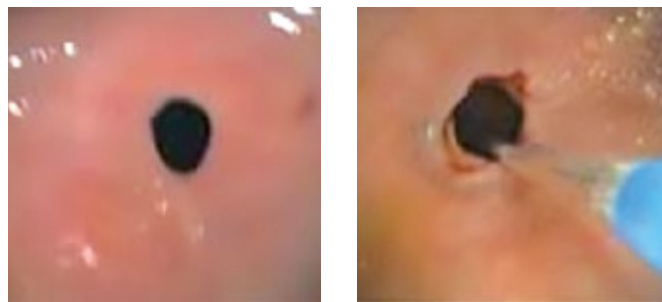

Endoscopic control during dilation therapy

## SPECIFICATIONS

REF	Balloon Ø mm	Balloon length mm	Catheter Ø Fr	Catheter length mm	Channel Ø mm	Application
<b>ACHALASIE-BALLOONS</b>						
BDC-26/80-7/10	26	80	7	1000	–	Achalasie
BDC-30/80-7/10	30	80	7	1000	–	Achalasie
BDC-34/80-7/10	34	80	7	1000	–	Achalasie
<b>FIXED WIRE BALLOONS</b>						
FWB-06/80-6/18-M	6-7-8	80	6	1800	2.8	Oesophagus
FWB-08/80-6/18-M	8-9-10	80	6	1800	2.8	Oesophagus
FWB-10/80-6/18-M	10-11-12	80	6	1800	2.8	Oesophagus
FWB-12/80-6/18-M	12-13.5-15	80	6	1800	2.8	Oesophagus
FWB-15/80-6/18-M	15-16.5-18	80	6	1800	2.8	Oesophagus
FWB-18/80-6/18-M	18-19-20	80	6	1800	2.8	Oesophagus

Packaging unit: 1 piece / suitable for guide wires with 0.035"