

Institute of Virology

Prof. Dr. Dr. Thomas W. Vahlenkamp
OIE Reference Laboratory for EBL
Email: thomas.vahlenkamp@uni-leipzig.de
An den Tierkliniken 29
04103 Leipzig (Germany)

+49(0) 0341-97 38200
Fax: +49(0) 341-97 38219
Leipzig 08.03.2017

Report: INgezim BLV COMPAC 2.0 ELISA
Delivery address: INGENASA
Hermanos García Noblejas, 39
28037 MADRID SPANIEN
Order number: -
Material: INgezim BLV COMPAC 2.0 ELISA (Prod Ref: 12.BLV.K3)
Lot: 110416; Expiry: 10/2017; Registration: 0808 RD
Request: BLV-ELISA with OIE reference serum for EBL 'E05'
Results:


1) INGEZIM BLV COMPAC 2.0 ELISA (Prod Ref: 12.BLV.K3)

Sample	Mean OD	Result
Positive control	0,12	positive
Negative control	1,96	negative
E05 dilution 1:2	0,07	positive
E05 dilution 1:10	0,10	positive
E05 dilution 1:50	0,15	positive
E05 dilution 1:100	0,17	positive
E05 dilution 1:1000	0,35	positive
E05 dilution 1:10.000	0,69	positive
E05 dilution 1:20.000	1,13	doubtful
E05 dilution 1:40.000	1,84	negative
E05 dilution 1:80.000	2,16	negative

Cut off negative: 1,2246; Cut off positive: 1,0408

The test is valid. Using a dilution of 1:10.000 of OIE reference serum for EBL 'E05', the INGEZIM BLV COMPAC 2.0 ELISA (Prod Ref: 12.BLV.K3) ELISA displays a positive result.

Kind regards

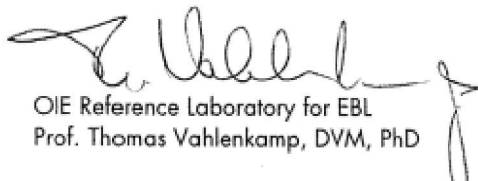

OIE Reference Laboratory for EBL
Prof. Thomas Vahlenkamp, DVM, PhD