The materials and information provided via this brochure are distributed internationally, but not all the products or services found on this brochure are necessarily available in your country or area. Please contact local office or agents to make sure of the availability of the relevant products or services.

Mindray Medical Netherlands B.V.

Drs.W.van Rovenstraat 8, 3871 AN Hoevelaken. The Netherlands Voice: (31-33) 254-4911

Mindray Medical Italy S.r.l.

Via Leonardo DA Vinci, 158, 20090 Trezzano Sul Naviglio (MI), Italy Voice: (39-02) 5737-401 Fax: (39-02) 5523-2018

Avenida Manoteras 38, Edificio B, Bajo B008. 28050-Madrid Spain Voice: (34-91) 392-3754 Fax: (34-91) 392-3756

Mindray Medical Germany GmbH

Goebelstrasse 21, 3. OG, 64293 Darmstadt Voice: (49 6151) 3910 0 Fax: (49 6151) 3910 300

Mindray House, Kingfisher Way, Hinchingbrooke Business Park, Huntingdon, Cambs, PE29 6FN Voice: (44) 1480-416-840 Fax: (44) 1480-436-588

Mindray Medical France SARL Europarc Creteil 1 allée des cerisiers 94000 CRETEIL, France Voice: (33-1) 4513-9150 Fax: (33-1) 4513-9151

Mindray Medical Sweden AB

Rissneleden 136, 174 57 Sundbyberg, Sweden Voice: (46) 8 555 54 121

Mindray Medical Rus Co.Ltd

2 Zvenigorodskaya Street, 13, bld 41, 123022 Moscow, Russia Voice: (7-499) 553-6036

Mindray Medical Colombia S.A.S

Room No.702, Torre UNIKA, Carrera 9 No.77-67, Bogota, Colombia Voice: (57-1) 3130892 / 3210916

Mindray do Brasil - Comércio e Distribuição de Equipamentos Médicos Ltda.

Rua Tavares Bastos. 329. Perdizes. São Paulo - SP. CEP 05012-020 Voice: (55-11) 3124-8026 Fax: (55-11) 3078-8035

Mindray Medical Mexico S. de R.L. de C.V. Félix Parra # 175 Colonia. San José Insurgentes Sur Delegación Benito Juárez 03900 México, D. F. Voice: (52-55) 5661-9450 / (52-55) 5662-6620

PT.Mindray Medical Indonesia

The East 11th Floor, Unit 08 Jl. Dr. Ide Anak Agung Gde Agung Kav. E 3.2 No.1 Jakarta 12950 Voice: (6221) 2902-7280 Fax: (6221) 2902-7283

Mindray Medical Egypt

Ground Floor -11 El Imam Metwally El Sharawi Street, Sheraton Heliopolis - El Nozha, Cairo, Egypt Voice: (2) 02 22691047 Fax: (2) 02 22691046

Mindray Medical Australia

ADD: 19 Harker St, Burwood VIC 3125. Australia Voice: (61-3) 9830 6610 Fax: (61-2) 9550 2841

Mindray Medical Nigeria

Mindray Medical India Pvt. Ltd. (Delhi Office)

16th floor, Building No. 9-B, DLF Cyber City, DLF Phase- 3, Gurgaon- 122002, Harvana Voice: (91-124) 4632-488 Fax: (91-124) 4632499

Mindray Medical India Pvt. Ltd. (Mumbai office) B-404, City Point, Andheri Kurla Road, Andheri East, Mumbai, 400 059, India Voice: (91-22) 4020-0000

Mindray Medikal Teknoloji Istanbul Ltd. Sti.

KÖYALTI MEVKİİ 29 EKİM CAD.N.1 KUYUMCUKENT 2.NO.PLZ.3.KT. D.1- 9 34520 BAHÇELİEVLER/YENİBOSNA, İSTANBUL, TÜRKİYE Voice: (90-212) 482 08 77 Fax: (90-212) 482 08 78


Mindray Medikal Teknoloji Istanbul Ltd. Sti. Branch Office Ankara

TURAN GÜNEŞ BULVARI ALEKSANDER DÜBÇEK CAD. NO:18/A CANKAYA/ANKARA, TÜRKİYE Fax: (90-312) 441 96 93

Mindray Medical Vietnam Co.,Ltd.
Unit 2, 4th Floor, "B" Building, Broadway Office Park,
102 Nguyen Luong Bang Blvd., Tan Phu Ward, Dist. 7, HCMC, Vietnam Voice: (84) 8 5413-6970 Fax: (84) 8 5413-6971

Mindray Medical Thailand Limited 1768 Thai summit Tower Room 407, 4th floor, New petchaburee Road, Kwaeng Bangkapi, Khet Huai Khwang, Bangkok 10320, Thailand Voice: (66) 2652-7288 Fax: (66) 2652-7266

Mindray Medical Korea KNI Tower 9F 823-19 Yeoksam dong, Seoul-ganum, Korea Fax: 008225688042

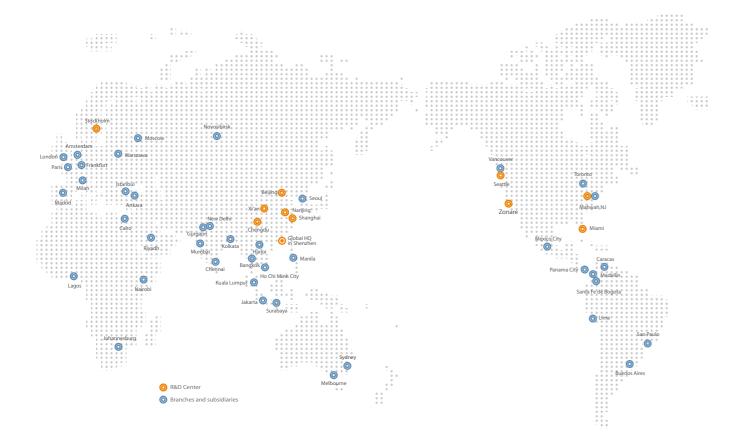

Mindray Building, Keji 12th Road South, High-tech Industrial Park, Nanshan, Shenzhen 518057, P.R. China Tel: +86 755 8188 8998 Fax: +86 755 26582680 E-mail: intl-market@mindray.com www.mindray.com

Mindray is listed on the NYSE under the symbol "MR"

© 2014 Shenzhen Mindray Bio-Medical Electronics Co., Ltd. All rights reserved. Specifications subject to changes without prior notice. Pls refer to user manual for manufacturer information P/N:ENG- new WATO EX-35-210X285X8PX20151120





new WATO EX-35 Anesthesia Machine

better performance than performance



Since 2006 Mindray has successfully installed over 20,000 anesthetic machines to satisfied customers all over the world.

Mindray takes immense pride knowing that every few seconds a clinician, somewhere in the world, will be switching on a Mindray anesthesia machine and using it with complete confidence and satisfaction. Over the last decade Mindray has continued to work closely with clinicians, across the world, to recognize and understand the clinical challenges encountered everyday and overcome them with new innovative and intuitive solutions. With this in mind, Mindray is now pleased to offer the new improved WATO EX-35 with advanced features, functionality and fresh new look for 2106.



new WATO EX-35

better performance than performance



Touch from Now...

Featuring a new 10.4" touchscreen with quick access parameter controls and flat menu design. The new WATO EX-35 UI enables the user to make adjustments to mode and monitor settings easier than ever witnessed before at this product level.

The touchable Smart Alarm system provides real-time detailed information that quickly and clearly alerts the clinician to adverse situations supporting early resolution.

Better Performance...

Employing the latest technologies the new WATO EX-35 supports an extensive range of features and functionality ensuring all anesthetic procedures are delivered safely, efficiently and effectively. The enhanced features including but not limited to:

Integrated Plug-and-play AG Module (auto ID of 5 agent and N_2 O), integrated multi-mode ventilator supports controlled and support modes of ventilation. The unique Mindray Ventilator Driven Gas Auto Switch function ensures the patient is ventilated without interruption at all times...

Better Usability...

Designed with essential usability in mind, the state-of-the-art design helps your to reduce your complex daily workload, gives you a fresh user experience.

The ergonomic handles, the foot pedal, make you to move the machine easily and use the machine for a long time without extra physical fatigue.



Easy Setting by Touch

The new WATO EX-35 supports a 10.4" high resolution touch-screen which supports a maximum of 3 ventilation waveforms (EtCO₂, Pressure, Flow & Volume), Spirometry Loop and full gas monitoring on one screen whilst allowing direct access to all parameters and settings. The new tabular, flat level, menu structure with quick touch control is a global first for devices at this market level.

Touchable Smart Alarm Management

The Smart Alarm system provides real-time detailed information that quickly and clearly alerts the clinician to adverse situations supporting early resolution.

The comprehensive alarm log includes physical alarm listings with displayed parameter snapshot and a technical alarm diagnosis software that offers possible causes for alarms which assists immediate trouble-shooting in the immediate situation.







With conveniently located handles and smooth running castors the new WATO EX-35 can easily be positioned effortlessly within the workplace.

The large easy-clean work surface provides plenty of space for note writing, placement of drugs and intubation equipment and an ergonomically designed footrest offers the user total comfort throughout the longest of cases.





Innovative P-n-P Gas Monitoring

Our improved Plug and Play single slot CO_2 module or double slot anesthetic agent module provides comprehensive, visual, breath-by-breath analysis for FiO_2 , EtO_2 , CO_2 , N_2O and auto detection of 5 anesthetic agents. What's more the modules can be moved between functionally similar Mindray monitors and anesthetic devices further enhancing device configuration and flexibility to suit your needs.

Driven Gas Auto-switch


Safety is at the heart of all Mindray's devices and the new WATO EX-35 is no exception. The all new ventilator gas Auto-switch function ensures the patient is mechanically ventilated at all times even if the primary ventilator drive gas supply fails.





High Performance Ventilator

he new WATO EX-35 supports a full range of controlled and support modes of ventilation and with dynamic tidal volume compensation, from as low as Tve 20ml, the new WATO EX-35 is suitable for all patient categories from neonates to adult.

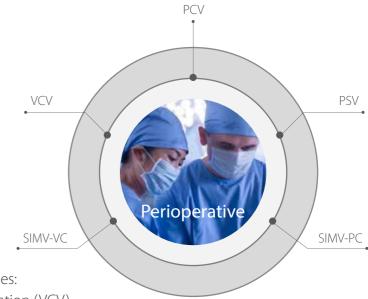

Available ventilation modes:

- Volume controlled Ventilation (VCV)
- Pressure Controlled Ventilation (PCV)
- Pressure Support Ventilation + apnea ventilation (PSV + apnea ventilation)
- Synchronized Intermittent Mandatory Ventilation VC (SIMV -VC)
- Synchronized Intermittent Mandatory Ventilation PC (SIMV -PC)