

Central Venous Catheter (CVC) Kit



TUORen

CE₂₇₉₇

Regulation (EU)
2017/745 (MDR)





CVC Kit

The Central Venous Catheter (CVC) Kit is used to insert into the central venous system for the infusions of fluids or medications, blood sampling or pressure measurement.

Central Venous Catheter

- Adopting top quality TPU material with good biocompatibility, reducing the infection rate.
- It automatically softens to adapt to the human body temperature after placing the catheter.
- The inner and outer walls are smooth, which reduces thrombosis and catheter placement resistance.
- Visualized under X-ray, accurate positioning.

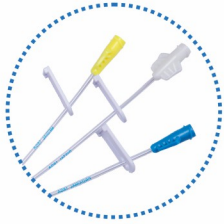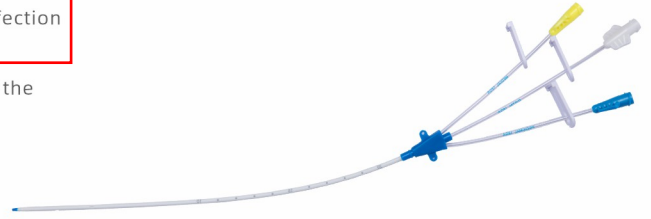

Luer connector

- Unique Luer connector design, anti-split and safer.
- Different length extension tubes with different color connectors, more convenient for identification.
- Comply with the ISO 80369-7 standard.



Extension line

- Transparent tube body makes it easy to observe, soft and elastic, anti-kinking, safer for clinical application.



Distal tip

- Soft and reduces the risk of damage to the blood vessel wall.





Side holes of catheter tip

- The side holes of multi-lumen catheter are on different directions, which suitable for a variety of drugs, and avoid contraindications to drug compatibility.

Component


Guide wire

- Anti-kinking;
- J-shaped end, soft tip, Nickel-titanium alloy, reduces damage.

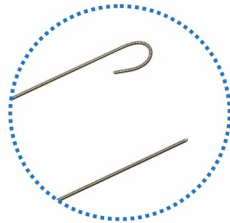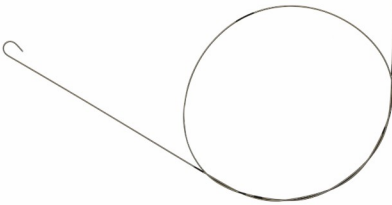

Dilator

- Smooth and round, low expansion resistance, moderate hardness and toughness, and more convenient to operate.

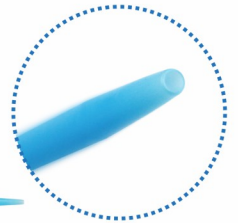

Introducer Needle

- The straight introducer needle and Y model introducer needle are available, and both are equipped with a check valve, reduce the possibility of bleeding and air entering the blood vessels.
- The thin needle reduces puncture resistance and damage.




Fastener

- The unique inner buckle structure design makes the clip and catheter fit more closely, reduces the the occurrence of fastener loosening and catheter replacement, relieve patients' pain and medical care troubles.



Scalpel

- Humanized design with sliding protective cap, anti-damage to the medical staff.

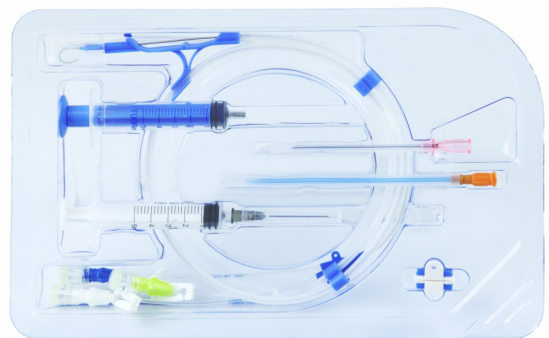

Specification

single lumen

Catheter OD	Model Specifications	Catheter length(cm)	Dialator	Guidewire	Cat. No
22Ga	1-22GA-08CM-Y	8	4FR*5CM	0.018*60CM	CVMY12208
	1-22GA-13CM-Y	13	4FR*5CM	0.018*60CM	CVMY12213
	1-22GA-16CM-Y	16	4FR*5CM	0.018*60CM	CVMY12216
	1-22GA-20CM-Y	20	4FR*5CM	0.018*60CM	CVMY12220
	1-22GA-30CM-Y	30	4FR*5CM	0.018*60CM	CVMY12230
20Ga	1-20GA-13CM-Y	13	4FR*5CM	0.021*60CM	CVMY12013
	1-20GA-16CM-Y	16	4FR*5CM	0.021*60CM	CVMY12016
	1-20GA-20CM-Y	20	4FR*5CM	0.021*60CM	CVMY12020
18Ga	1-18GA-13CM-Y	13	5.5FR*6CM	0.021*60CM	CVMY11813
	1-18GA-16CM-Y	16	5.5FR*6CM	0.021*60CM	CVMY11816
	1-18GA-20CM-Y	20	5.5FR*6CM	0.021*60CM	CVMY11820
16Ga	1-16GA-16CM-Y	16	6.5FR*10CM	0.035*60CM	CVMY11616
	1-16GA-20CM-Y	20	6.5FR*10CM	0.035*60CM	CVMY11620
	1-16GA-30CM-Y	30	6.5FR*10CM	0.035*60CM	CVMY11630
14Ga	1-14GA-16CM-Y	16	8FR*10CM	0.035*60CM	CVMY11416
	1-14GA-20CM-Y	20	8FR*10CM	0.035*60CM	CVMY11420
	1-14GA-30CM-Y	30	8FR*10CM	0.035*60CM	CVMY11430
12Ga	1-12GA-16CM-Y	16	8.5FR*10CM	0.035*60CM	CVMY11216
	1-12GA-20CM-Y	20	8.5FR*10CM	0.035*60CM	CVMY11220
	1-12GA-30CM-Y	30	8.5FR*10CM	0.035*60CM	CVMY11230

double lumen

Catheter OD	Model Specifications	Catheter length(cm)	Dialator	Guidewire	Cat. No
4Fr	2-4FR-05CM-Y	5	5.5FR*6CM	0.018*60CM	CVMY24005
	2-4FR-08CM-Y	8	5.5FR*6CM	0.018*60CM	CVMY24008
	2-4FR-13CM-Y	13	5.5FR*6CM	0.018*60CM	CVMY24013
	2-4FR-16CM-Y	16	5.5FR*6CM	0.018*60CM	CVMY24016
	2-4FR-20CM-Y	20	5.5FR*6CM	0.018*60CM	CVMY24020
5Fr	2-5FR-13CM-Y	13	5.5FR*6CM	0.021*60CM	CVMY25013
	2-5FR-16CM-Y	16	5.5FR*6CM	0.021*60CM	CVMY25016
	2-5FR-20CM-Y	20	5.5FR*6CM	0.021*60CM	CVMY25020
6Fr	2-6FR-13CM-Y	13	6.5FR*8CM	0.021*60CM	CVMY26013
	2-6FR-16CM-Y	16	6.5FR*8CM	0.021*60CM	CVMY26016
	2-6FR-20CM-Y	20	6.5FR*8CM	0.021*60CM	CVMY26020
7Fr	2-7FR-16CM-Y	16	8.5FR*10CM	0.035*60CM	CVMY27016
	2-7FR-20CM-Y	20	8.5FR*10CM	0.035*60CM	CVMY27020
	2-7FR-30CM-Y	30	8.5FR*10CM	0.035*60CM	CVMY27030
8Fr	2-8FR-16CM-Y	16	8.5FR*10CM	0.035*60CM	CVMY28016
	2-8FR-20CM-Y	20	8.5FR*10CM	0.035*60CM	CVMY28020
	2-8FR-30CM-Y	30	8.5FR*10CM	0.035*60CM	CVMY28030

triple lumen

Catheter OD	Model Specifications	Catheter length(cm)	Dialator	Guidewire	Cat. No
4Fr	3-4FR-08CM-Y	8	5.5FR*6CM	0.018*60CM	CVMY34008
	3-4FR-13CM-Y	13	5.5FR*6CM	0.018*60CM	CVMY34013
	3-4FR-16CM-Y	16	5.5FR*6CM	0.018*60CM	CVMY34016
	3-4FR-20CM-Y	20	5.5FR*6CM	0.018*60CM	CVMY34020
5.5Fr	3-5.5FR-08CM-Y	8	6.5FR*8CM	0.021*60CM	CVMY35508
	3-5.5FR-13CM-Y	13	6.5FR*8CM	0.021*60CM	CVMY35513
	3-5.5FR-16CM-Y	16	6.5FR*8CM	0.021*60CM	CVMY35516
	3-5.5FR-20CM-Y	20	6.5FR*8CM	0.021*60CM	CVMY35520
7Fr	3-7FR-16CM-Y	16	8.5FR*10CM	0.035*60CM	CVMY37016
	3-7FR-20CM-Y	20	8.5FR*10CM	0.035*60CM	CVMY37020
	3-7FR-30CM-Y	30	8.5FR*10CM	0.035*60CM	CVMY37030
8.5Fr	3-8.5FR-16CM-Y	16	9FR*10CM	0.035*60CM	CVMY38516
	3-8.5FR-20CM-Y	20	9FR*10CM	0.035*60CM	CVMY38520
	3-8.5FR-30CM-Y	30	9FR*10CM	0.035*60CM	CVMY38530

quadruple lumen

Catheter OD	Model Specifications	Catheter length(cm)	Dialator	Guidewire	Cat. No
8.5Fr	4-8.5FR-20CM-Y	20	9FR*10CM	0.035*60CM	CVMY48520
	4-8.5FR-30CM-Y	30	9FR*10CM	0.035*60CM	CVMY48530



Y model Introducer needle

Specification

single lumen

Catheter OD	Model Specifications	Catheter length(cm)	Dialator	Guidewire	Cat. No
22Ga	1-22GA-08CM-I	8	4FR*5CM	0.018*60CM	ICVM12208
	1-22GA-13CM-I	13	4FR*5CM	0.018*60CM	ICVM12213
	1-22GA-16CM-I	16	4FR*5CM	0.018*60CM	ICVM12216
	1-22GA-20CM-I	20	4FR*5CM	0.018*60CM	ICVM12220
	1-22GA-30CM-I	30	4FR*5CM	0.018*60CM	ICVM12230
20Ga	1-20GA-13CM-I	13	4FR*5CM	0.021*60CM	ICVM12013
	1-20GA-16CM-I	16	4FR*5CM	0.021*60CM	ICVM12016
	1-20GA-20CM-I	20	4FR*5CM	0.021*60CM	ICVM12020
18Ga	1-18GA-13CM-I	13	5.5FR*6CM	0.021*60CM	ICVM11813
	1-18GA-16CM-I	16	5.5FR*6CM	0.021*60CM	ICVM11816
	1-18GA-20CM-I	20	5.5FR*6CM	0.021*60CM	ICVM11820
16Ga	1-16GA-16CM-I	16	6.5FR*10CM	0.035*60CM	ICVM11616
	1-16GA-20CM-I	20	6.5FR*10CM	0.035*60CM	ICVM11620
	1-16GA-30CM-I	30	6.5FR*10CM	0.035*60CM	ICVM11630
14Ga	1-14GA-16CM-I	16	8FR*10CM	0.035*60CM	ICVM11416
	1-14GA-20CM-I	20	8FR*10CM	0.035*60CM	ICVM11420
	1-14GA-30CM-I	30	8FR*10CM	0.035*60CM	ICVM11430
12Ga	1-12GA-16CM-I	16	8.5FR*10CM	0.035*60CM	ICVM11216
	1-12GA-20CM-I	20	8.5FR*10CM	0.035*60CM	ICVM11220
	1-12GA-30CM-I	30	8.5FR*10CM	0.035*60CM	ICVM11230

double lumen

Catheter OD	Model Specifications	Catheter length(cm)	Dialator	Guidewire	Cat. No
4Fr	2-4FR-05CM-I	5	5.5FR*6CM	0.018*60CM	ICVM24005
	2-4FR-08CM-I	8	5.5FR*6CM	0.018*60CM	ICVM24008
	2-4FR-13CM-I	13	5.5FR*6CM	0.018*60CM	ICVM24013
	2-4FR-16CM-I	16	5.5FR*6CM	0.018*60CM	ICVM24016
	2-4FR-20CM-I	20	5.5FR*6CM	0.018*60CM	ICVM20420
5Fr	2-5FR-13CM-I	13	5.5FR*6CM	0.021*60CM	ICVM25013
	2-5FR-16CM-I	16	5.5FR*6CM	0.021*60CM	ICVM25016
	2-5FR-20CM-I	20	5.5FR*6CM	0.021*60CM	ICVM25020
6Fr	2-6FR-13CM-I	13	6.5FR*8CM	0.021*60CM	ICVM26013
	2-6FR-16CM-I	16	6.5FR*8CM	0.021*60CM	ICVM26016
	2-6FR-20CM-I	20	6.5FR*8CM	0.021*60CM	ICVM26020
7Fr	2-7FR-16CM-I	16	8.5FR*10CM	0.035*60CM	ICVM27016
	2-7FR-20CM-I	20	8.5FR*10CM	0.035*60CM	ICVM27020
	2-7FR-30CM-I	30	8.5FR*10CM	0.035*60CM	ICVM27030
8Fr	2-8FR-16CM-I	16	8.5FR*10CM	0.035*60CM	ICVM28016
	2-8FR-20CM-I	20	8.5FR*10CM	0.035*60CM	ICVM28020
	2-8FR-30CM-I	30	8.5FR*10CM	0.035*60CM	ICVM28030

triple lumen

Catheter OD	Model Specifications	Catheter length(cm)	Dialator	Guidewire	Cat. No
4Fr	3-4FR-08CM-I	8	5.5FR*6CM	0.018*60CM	ICVM34008
	3-4FR-13CM-I	13	5.5FR*6CM	0.018*60CM	ICVM34013
	3-4FR-16CM-I	16	5.5FR*6CM	0.018*60CM	ICVM34016
	3-4FR-20CM-I	20	5.5FR*6CM	0.018*60CM	ICVM34020
5.5Fr	3-5.5FR-08CM-I	8	6.5FR*8CM	0.021*60CM	ICVM35508
	3-5.5FR-13CM-I	13	6.5FR*8CM	0.021*60CM	ICVM35513
	3-5.5FR-16CM-I	16	6.5FR*8CM	0.021*60CM	ICVM35516
	3-5.5FR-20CM-I	20	6.5FR*8CM	0.021*60CM	ICVM35520
7Fr	3-7FR-16CM-I	16	8.5FR*10CM	0.035*60CM	ICVM37016
	3-7FR-20CM-I	20	8.5FR*10CM	0.035*60CM	ICVM37020
	3-7FR-30CM-I	30	8.5FR*10CM	0.035*60CM	ICVM37030
8.5Fr	3-8.5FR-16CM-I	16	9FR*10CM	0.035*60CM	ICVM38516
	3-8.5FR-20CM-I	20	9FR*10CM	0.035*60CM	ICVM38520
	3-8.5FR-30CM-I	30	9FR*10CM	0.035*60CM	ICVM38530

quadruple lumen

Catheter OD	Model Specifications	Catheter length(cm)	Dialator	Guidewire	Cat. No
8.5Fr	4-8.5FR-20CM-I	20	9FR*10CM	0.035*60CM	ICVM48520
	4-8.5FR-30CM-I	30	9FR*10CM	0.035*60CM	ICVM48530



Straight introducer needle



TUOREN MEDICAL

+86-373-8605444
info@bjtuoren.com; info@etuoren.com
www.tuoren.com



Annex 1: List of Model Specifications of Central Venous Catheter (CVC) Kit														
Product Code	Catalogue Number	Specification												
		Central Venous Catheter	Guide wire	Dilator	Introducer Needle		Fastener	Hypodermic Syringe		Blue introducer syringe	Hypodermic Needle	Scalpel	Injection site cap	Extension line clamp
					Y Model	I Model		I Model	II Model					
CVMY12208	1-22Ga-08cm-Y	1-22Ga-08cm	0.018"*60cm	4Fr*5cm	20Ga	-	22Ga	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM12208	1-22Ga-08cm-I	1-22Ga-08cm	0.018"*60cm	4Fr*5cm	-	20Ga	22Ga	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS12208	1-22Ga-08cm-S	1-22Ga-08cm	0.018"*60cm	4Fr*5cm	-	20Ga	22Ga	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY12213	1-22Ga-13cm-Y	1-22Ga-13cm	0.018"*60cm	4Fr*5cm	20Ga	-	22Ga	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM12213	1-22Ga-13cm-I	1-22Ga-13cm	0.018"*60cm	4Fr*5cm	-	20Ga	22Ga	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS12213	1-22Ga-13cm-S	1-22Ga-13cm	0.018"*60cm	4Fr*5cm	-	20Ga	22Ga	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY12216	1-22Ga-16cm-Y	1-22Ga-16cm	0.018"*60cm	4Fr*5cm	20Ga	-	22Ga	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM12216	1-22Ga-16cm-I	1-22Ga-16cm	0.018"*60cm	4Fr*5cm	-	20Ga	22Ga	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS12216	1-22Ga-16cm-S	1-22Ga-16cm	0.018"*60cm	4Fr*5cm	-	20Ga	22Ga	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY12220	1-22Ga-20cm-Y	1-22Ga-20cm	0.018"*60cm	4Fr*5cm	20Ga	-	22Ga	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM12220	1-22Ga-20cm-I	1-22Ga-20cm	0.018"*60cm	4Fr*5cm	-	20Ga	22Ga	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS12220	1-22Ga-20cm-S	1-22Ga-20cm	0.018"*60cm	4Fr*5cm	-	20Ga	22Ga	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY12230	1-22Ga-30cm-Y	1-22Ga-30cm	0.018"*60cm	4Fr*5cm	20Ga	-	22Ga	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM12230	1-22Ga-30cm-I	1-22Ga-30cm	0.018"*60cm	4Fr*5cm	-	20Ga	22Ga	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS12230	1-22Ga-30cm-S	1-22Ga-30cm	0.018"*60cm	4Fr*5cm	-	20Ga	22Ga	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY12013	1-20Ga-13cm-Y	1-20Ga-13cm	0.021"*60cm	4Fr*5cm	20Ga	-	20Ga	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM12013	1-20Ga-13cm-I	1-20Ga-13cm	0.021"*60cm	4Fr*5cm	-	20Ga	20Ga	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS12013	1-20Ga-13cm-S	1-20Ga-13cm	0.021"*60cm	4Fr*5cm	-	20Ga	20Ga	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY12016	1-20Ga-16cm-Y	1-20Ga-16cm	0.021"*60cm	4Fr*5cm	20Ga	-	20Ga	5mL	5mL	-	0.7mm	11#	Regular	1.7#

Annex 1: List of Model Specifications of Central Venous Catheter (CVC) Kit														
Product Code	Catalogue Number	Specification												
		Central Venous Catheter	Guide wire	Dilator	Introducer Needle		Fastener	Hypodermic Syringe		Blue introducer syringe	Hypodermic Needle	Scalpel	Injection site cap	Extension line clamp
					Y Model	I Model		I Model	II Model					
ICVM12016	1-20Ga-16cm-I	1-20Ga-16cm	0.021"*60cm	4Fr*5cm	-	20Ga	20Ga	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS12016	1-20Ga-16cm-S	1-20Ga-16cm	0.021"*60cm	4Fr*5cm	-	20Ga	20Ga	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY12020	1-20Ga-20cm-Y	1-20Ga-20cm	0.021"*60cm	4Fr*5cm	20Ga	-	20Ga	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM12020	1-20Ga-20cm-I	1-20Ga-20cm	0.021"*60cm	4Fr*5cm	-	20Ga	20Ga	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS12020	1-20Ga-20cm-S	1-20Ga-20cm	0.021"*60cm	4Fr*5cm	-	20Ga	20Ga	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY11813	1-18Ga-13cm-Y	1-18Ga-13cm	0.021"*60cm	5.5Fr*6cm	20Ga	-	18Ga	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM11813	1-18Ga-13cm-I	1-18Ga-13cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS11813	1-18Ga-13cm-S	1-18Ga-13cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY11816	1-18Ga-16cm-Y	1-18Ga-16cm	0.021"*60cm	5.5Fr*6cm	20Ga	-	18Ga	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM11816	1-18Ga-16cm-I	1-18Ga-16cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS11816	1-18Ga-16cm-S	1-18Ga-16cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY11820	1-18Ga-20cm-Y	1-18Ga-20cm	0.021"*60cm	5.5Fr*6cm	20Ga	-	18Ga	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM11820	1-18Ga-20cm-I	1-18Ga-20cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS11820	1-18Ga-20cm-S	1-18Ga-20cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY11616	1-16Ga-16cm-Y	1-16Ga-16cm	0.035"*60cm	6.5Fr*10cm	18Ga	-	16Ga	5mL	5mL	-	0.7mm	11#	Regular	2.2#
ICVM11616	1-16Ga-16cm-I	1-16Ga-16cm	0.035"*60cm	6.5Fr*10cm	-	18Ga	16Ga	5mL	-	5mL	0.7mm	11#	Regular	2.2#
CVMS11616	1-16Ga-16cm-S	1-16Ga-16cm	0.035"*60cm	6.5Fr*10cm	-	18Ga	16Ga	5mL	-	-	0.7mm	11#	Regular	2.2#
CVMY11620	1-16Ga-20cm-Y	1-16Ga-20cm	0.035"*60cm	6.5Fr*10cm	18Ga	-	16Ga	5mL	5mL	-	0.7mm	11#	Regular	2.2#
ICVM11620	1-16Ga-20cm-I	1-16Ga-20cm	0.035"*60cm	6.5Fr*10cm	-	18Ga	16Ga	5mL	-	5mL	0.7mm	11#	Regular	2.2#
CVMS11620	1-16Ga-20cm-S	1-16Ga-20cm	0.035"*60cm	6.5Fr*10cm	-	18Ga	16Ga	5mL	-	-	0.7mm	11#	Regular	2.2#
CVMY11630	1-16Ga-30cm-Y	1-16Ga-30cm	0.035"*60cm	6.5Fr*10cm	18Ga	-	16Ga	5mL	5mL	-	0.7mm	11#	Regular	2.2#
ICVM11630	1-16Ga-30cm-I	1-16Ga-30cm	0.035"*60cm	6.5Fr*10cm	-	18Ga	16Ga	5mL	-	5mL	0.7mm	11#	Regular	2.2#

Annex 1: List of Model Specifications of Central Venous Catheter (CVC) Kit														
Product Code	Catalogue Number	Specification												
		Central Venous Catheter	Guide wire	Dilator	Introducer Needle		Fastener	Hypodermic Syringe		Blue introducer syringe	Hypodermic Needle	Scalpel	Injection site cap	Extension line clamp
					Y Model	I Model		I Model	II Model					
CVMS11630	1-16Ga-30cm-S	1-16Ga-30cm	0.035"*60cm	6.5Fr*10cm	-	18Ga	16Ga	5mL	-	-	0.7mm	11#	Regular	2.2#
CVMY11416	1-14Ga-16cm-Y	1-14Ga-16cm	0.035"*60cm	8Fr*10cm	18Ga	-	14Ga	5mL	5mL	-	0.7mm	11#	Regular	2.5#
ICVM11416	1-14Ga-16cm-I	1-14Ga-16cm	0.035"*60cm	8Fr*10cm	-	18Ga	14Ga	5mL	-	5mL	0.7mm	11#	Regular	2.5#
CVMS11416	1-14Ga-16cm-S	1-14Ga-16cm	0.035"*60cm	8Fr*10cm	-	18Ga	14Ga	5mL	-	-	0.7mm	11#	Regular	2.5#
CVMY11420	1-14Ga-20cm-Y	1-14Ga-20cm	0.035"*60cm	8Fr*10cm	18Ga	-	14Ga	5mL	5mL	-	0.7mm	11#	Regular	2.5#
ICVM11420	1-14Ga-20cm-I	1-14Ga-20cm	0.035"*60cm	8Fr*10cm	-	18Ga	14Ga	5mL	-	5mL	0.7mm	11#	Regular	2.5#
CVMS11420	1-14Ga-20cm-S	1-14Ga-20cm	0.035"*60cm	8Fr*10cm	-	18Ga	14Ga	5mL	-	-	0.7mm	11#	Regular	2.5#
CVMY11430	1-14Ga-30cm-Y	1-14Ga-30cm	0.035"*60cm	8Fr*10cm	18Ga	-	14Ga	5mL	5mL	-	0.7mm	11#	Regular	2.5#
ICVM11430	1-14Ga-30cm-I	1-14Ga-30cm	0.035"*60cm	8Fr*10cm	-	18Ga	14Ga	5mL	-	5mL	0.7mm	11#	Regular	2.5#
CVMS11430	1-14Ga-30cm-S	1-14Ga-30cm	0.035"*60cm	8Fr*10cm	-	18Ga	14Ga	5mL	-	-	0.7mm	11#	Regular	2.5#
CVMY11216	1-12Ga-16cm-Y	1-12Ga-16cm	0.035"*60cm	9Fr*10cm	18Ga	-	8.5Fr	5mL	5mL	-	0.7mm	11#	Regular	3.0#
ICVM11216	1-12Ga-16cm-I	1-12Ga-16cm	0.035"*60cm	9Fr*10cm	-	18Ga	8.5Fr	5mL	-	5mL	0.7mm	11#	Regular	3.0#
CVMS11216	1-12Ga-16cm-S	1-12Ga-16cm	0.035"*60cm	9Fr*10cm	-	18Ga	8.5Fr	5mL	-	-	0.7mm	11#	Regular	3.0#
CVMY11220	1-12Ga-20cm-Y	1-12Ga-20cm	0.035"*60cm	9Fr*10cm	18Ga	-	8.5Fr	5mL	5mL	-	0.7mm	11#	Regular	3.0#
ICVM11220	1-12Ga-20cm-I	1-12Ga-20cm	0.035"*60cm	9Fr*10cm	-	18Ga	8.5Fr	5mL	-	5mL	0.7mm	11#	Regular	3.0#
CVMS11220	1-12Ga-20cm-S	1-12Ga-20cm	0.035"*60cm	9Fr*10cm	-	18Ga	8.5Fr	5mL	-	-	0.7mm	11#	Regular	3.0#
CVMY11230	1-12Ga-30cm-Y	1-12Ga-30cm	0.035"*60cm	9Fr*10cm	18Ga	-	8.5Fr	5mL	5mL	-	0.7mm	11#	Regular	3.0#
ICVM11230	1-12Ga-30cm-I	1-12Ga-30cm	0.035"*60cm	9Fr*10cm	-	18Ga	8.5Fr	5mL	-	5mL	0.7mm	11#	Regular	3.0#
CVMS11230	1-12Ga-30cm-S	1-12Ga-30cm	0.035"*60cm	9Fr*10cm	-	18Ga	8.5Fr	5mL	-	-	0.7mm	11#	Regular	3.0#
CVMY24005	2-4Fr-05cm-Y	2-4Fr-05cm	0.021"*60cm	5.5Fr*6cm	20Ga	-	18Ga	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM24005	2-4Fr-05cm-I	2-4Fr-05cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS24005	2-4Fr-05cm-S	2-4Fr-05cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	-	0.7mm	11#	Regular	1.7#

Annex 1: List of Model Specifications of Central Venous Catheter (CVC) Kit														
Product Code	Catalogue Number	Specification												
		Central Venous Catheter	Guide wire	Dilator	Introducer Needle		Fastener	Hypodermic Syringe		Blue introducer syringe	Hypodermic Needle	Scalpel	Injection site cap	Extension line clamp
					Y Model	I Model		I Model	II Model					
CVMY24008	2-4Fr-08cm-Y	2-4Fr-08cm	0.021"*60cm	5.5Fr*6cm	20Ga	-	18Ga	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM24008	2-4Fr-08cm-I	2-4Fr-08cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS24008	2-4Fr-08cm-S	2-4Fr-08cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY24013	2-4Fr-13cm-Y	2-4Fr-13cm	0.021"*60cm	5.5Fr*6cm	20Ga	-	18Ga	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM24013	2-4Fr-13cm-I	2-4Fr-13cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS24013	2-4Fr-13cm-S	2-4Fr-13cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY24016	2-4Fr-16cm-Y	2-4Fr-16cm	0.021"*60cm	5.5Fr*6cm	20Ga	-	18Ga	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM24016	2-4Fr-16cm-I	2-4Fr-16cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS24016	2-4Fr-16cm-S	2-4Fr-16cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY24020	2-4Fr-20cm-Y	2-4Fr-20cm	0.021"*60cm	5.5Fr*6cm	20Ga	-	18Ga	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM24020	2-4Fr-20cm-I	2-4Fr-20cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS24020	2-4Fr-20cm-S	2-4Fr-20cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY25013	2-5Fr-13cm-Y	2-5Fr-13cm	0.021"*60cm	5.5Fr*6cm	20Ga	-	16Ga	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM25013	2-5Fr-13cm-I	2-5Fr-13cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	16Ga	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS25013	2-5Fr-13cm-S	2-5Fr-13cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	16Ga	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY25016	2-5Fr-16cm-Y	2-5Fr-16cm	0.021"*60cm	5.5Fr*6cm	20Ga	-	16Ga	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM25016	2-5Fr-16cm-I	2-5Fr-16cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	16Ga	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS25016	2-5Fr-16cm-S	2-5Fr-16cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	16Ga	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY25020	2-5Fr-20cm-Y	2-5Fr-20cm	0.021"*60cm	5.5Fr*6cm	20Ga	-	16Ga	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM25020	2-5Fr-20cm-I	2-5Fr-20cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	16Ga	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS25020	2-5Fr-20cm-S	2-5Fr-20cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	16Ga	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY26013	2-6Fr-13cm-Y	2-6Fr-13cm	0.021"*60cm	6.5Fr*8cm	20Ga	-	14Ga	5mL	5mL	-	0.7mm	11#	Regular	2.2#

Annex 1: List of Model Specifications of Central Venous Catheter (CVC) Kit														
Product Code	Catalogue Number	Specification												
		Central Venous Catheter	Guide wire	Dilator	Introducer Needle		Fastener	Hypodermic Syringe		Blue introducer syringe	Hypodermic Needle	Scalpel	Injection site cap	Extension line clamp
					Y Model	I Model		I Model	II Model					
ICVM26013	2-6Fr-13cm-I	2-6Fr-13cm	0.021"*60cm	6.5Fr*8cm	-	20Ga	14Ga	5mL	-	5mL	0.7mm	11#	Regular	2.2#
CVMS26013	2-6Fr-13cm-S	2-6Fr-13cm	0.021"*60cm	6.5Fr*8cm	-	20Ga	14Ga	5mL	-	-	0.7mm	11#	Regular	2.2#
CVMY26016	2-6Fr-16cm-Y	2-6Fr-16cm	0.021"*60cm	6.5Fr*8cm	20Ga	-	14Ga	5mL	5mL	-	0.7mm	11#	Regular	2.2#
ICVM26016	2-6Fr-16cm-I	2-6Fr-16cm	0.021"*60cm	6.5Fr*8cm	-	20Ga	14Ga	5mL	-	5mL	0.7mm	11#	Regular	2.2#
CVMS26016	2-6Fr-16cm-S	2-6Fr-16cm	0.021"*60cm	6.5Fr*8cm	-	20Ga	14Ga	5mL	-	-	0.7mm	11#	Regular	2.2#
CVMY26020	2-6Fr-20cm-Y	2-6Fr-20cm	0.021"*60cm	6.5Fr*8cm	20Ga	-	14Ga	5mL	5mL	-	0.7mm	11#	Regular	2.2#
ICVM26020	2-6Fr-20cm-I	2-6Fr-20cm	0.021"*60cm	6.5Fr*8cm	-	20Ga	14Ga	5mL	-	5mL	0.7mm	11#	Regular	2.2#
CVMS26020	2-6Fr-20cm-S	2-6Fr-20cm	0.021"*60cm	6.5Fr*8cm	-	20Ga	14Ga	5mL	-	-	0.7mm	11#	Regular	2.2#
CVMY27016	2-7Fr-16cm-Y	2-7Fr-16cm	0.035"*60cm	8.5Fr*10cm	18Ga	-	7Fr	5mL	5mL	-	0.7mm	11#	Regular	2.2#
ICVM27016	2-7Fr-16cm-I	2-7Fr-16cm	0.035"*60cm	8.5Fr*10cm	-	18Ga	7Fr	5mL	-	5mL	0.7mm	11#	Regular	2.2#
CVMS27016	2-7Fr-16cm-S	2-7Fr-16cm	0.035"*60cm	8.5Fr*10cm	-	18Ga	7Fr	5mL	-	-	0.7mm	11#	Regular	2.2#
CVMY27020	2-7Fr-20cm-Y	2-7Fr-20cm	0.035"*60cm	8.5Fr*10cm	18Ga	-	7Fr	5mL	5mL	-	0.7mm	11#	Regular	2.2#
ICVM27020	2-7Fr-20cm-I	2-7Fr-20cm	0.035"*60cm	8.5Fr*10cm	-	18Ga	7Fr	5mL	-	5mL	0.7mm	11#	Regular	2.2#
CVMS27020	2-7Fr-20cm-S	2-7Fr-20cm	0.035"*60cm	8.5Fr*10cm	-	18Ga	7Fr	5mL	-	-	0.7mm	11#	Regular	2.2#
CVMY27030	2-7Fr-30cm-Y	2-7Fr-30cm	0.035"*60cm	8.5Fr*10cm	18Ga	-	7Fr	5mL	5mL	-	0.7mm	11#	Regular	2.2#
ICVM27030	2-7Fr-30cm-I	2-7Fr-30cm	0.035"*60cm	8.5Fr*10cm	-	18Ga	7Fr	5mL	-	5mL	0.7mm	11#	Regular	2.2#
CVMS27030	2-7Fr-30cm-S	2-7Fr-30cm	0.035"*60cm	8.5Fr*10cm	-	18Ga	7Fr	5mL	-	-	0.7mm	11#	Regular	2.2#
CVMY28016	2-8Fr-16cm-Y	2-8Fr-16cm	0.035"*60cm	8.5Fr*10cm	18Ga	-	7Fr	5mL	5mL	-	0.7mm	11#	Regular	2.5#
ICVM28016	2-8Fr-16cm-I	2-8Fr-16cm	0.035"*60cm	8.5Fr*10cm	-	18Ga	7Fr	5mL	-	5mL	0.7mm	11#	Regular	2.5#
CVMS28016	2-8Fr-16cm-S	2-8Fr-16cm	0.035"*60cm	8.5Fr*10cm	-	18Ga	7Fr	5mL	-	-	0.7mm	11#	Regular	2.5#
CVMY28020	2-8Fr-20cm-Y	2-8Fr-20cm	0.035"*60cm	8.5Fr*10cm	18Ga	-	7Fr	5mL	5mL	-	0.7mm	11#	Regular	2.5#
ICVM28020	2-8Fr-20cm-I	2-8Fr-20cm	0.035"*60cm	8.5Fr*10cm	-	18Ga	7Fr	5mL	-	5mL	0.7mm	11#	Regular	2.5#

Annex 1: List of Model Specifications of Central Venous Catheter (CVC) Kit														
Product Code	Catalogue Number	Specification												
		Central Venous Catheter	Guide wire	Dilator	Introducer Needle		Fastener	Hypodermic Syringe		Blue introducer syringe	Hypodermic Needle	Scalpel	Injection site cap	Extension line clamp
					Y Model	I Model		I Model	II Model					
CVMS28020	2-8Fr-20cm-S	2-8Fr-20cm	0.035"*60cm	8.5Fr*10cm	-	18Ga	7Fr	5mL	-	-	0.7mm	11#	Regular	2.5#
CVMY28030	2-8Fr-30cm-Y	2-8Fr-30cm	0.035"*60cm	8.5Fr*10cm	18Ga	-	7Fr	5mL	5mL	-	0.7mm	11#	Regular	2.5#
ICVM28030	2-8Fr-30cm-I	2-8Fr-30cm	0.035"*60cm	8.5Fr*10cm	-	18Ga	7Fr	5mL	-	5mL	0.7mm	11#	Regular	2.5#
CVMS28030	2-8Fr-30cm-S	2-8Fr-30cm	0.035"*60cm	8.5Fr*10cm	-	18Ga	7Fr	5mL	-	-	0.7mm	11#	Regular	2.5#
CVMY34008	3-4Fr-08cm-Y	3-4Fr-08cm	0.021"*60cm	5.5Fr*6cm	20Ga	-	18Ga	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM34008	3-4Fr-08cm-I	3-4Fr-08cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS34008	3-4Fr-08cm-S	3-4Fr-08cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY34013	3-4Fr-13cm-Y	3-4Fr-13cm	0.021"*60cm	5.5Fr*6cm	20Ga	-	18Ga	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM34013	3-4Fr-13cm-I	3-4Fr-13cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS34013	3-4Fr-13cm-S	3-4Fr-13cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY34016	3-4Fr-16cm-Y	3-4Fr-16cm	0.021"*60cm	5.5Fr*6cm	20Ga	-	18Ga	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM34016	3-4Fr-16cm-I	3-4Fr-16cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS34016	3-4Fr-16cm-S	3-4Fr-16cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY34020	3-4Fr-20cm-Y	3-4Fr-20cm	0.021"*60cm	5.5Fr*6cm	20Ga	-	18Ga	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM34020	3-4Fr-20cm-I	3-4Fr-20cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS34020	3-4Fr-20cm-S	3-4Fr-20cm	0.021"*60cm	5.5Fr*6cm	-	20Ga	18Ga	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY35508	3-5.5Fr-08cm-Y	3-5.5Fr-08cm	0.021"*60cm	6.5Fr*8cm	20Ga	-	5.5Fr	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM35508	3-5.5Fr-08cm-I	3-5.5Fr-08cm	0.021"*60cm	6.5Fr*8cm	-	20Ga	5.5Fr	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS35508	3-5.5Fr-08cm-S	3-5.5Fr-08cm	0.021"*60cm	6.5Fr*8cm	-	20Ga	5.5Fr	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY35513	3-5.5Fr-13cm-Y	3-5.5Fr-13cm	0.021"*60cm	6.5Fr*8cm	20Ga	-	5.5Fr	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM35513	3-5.5Fr-13cm-I	3-5.5Fr-13cm	0.021"*60cm	6.5Fr*8cm	-	20Ga	5.5Fr	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS35513	3-5.5Fr-13cm-S	3-5.5Fr-13cm	0.021"*60cm	6.5Fr*8cm	-	20Ga	5.5Fr	5mL	-	-	0.7mm	11#	Regular	1.7#

Annex 1: List of Model Specifications of Central Venous Catheter (CVC) Kit														
Product Code	Catalogue Number	Specification												
		Central Venous Catheter	Guide wire	Dilator	Introducer Needle		Fastener	Hypodermic Syringe		Blue introducer syringe	Hypodermic Needle	Scalpel	Injection site cap	Extension line clamp
					Y Model	I Model		I Model	II Model					
CVMY35516	3-5.5Fr-16cm-Y	3-5.5Fr-16cm	0.021"*60cm	6.5Fr*8cm	20Ga	-	5.5Fr	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM35516	3-5.5Fr-16cm-I	3-5.5Fr-16cm	0.021"*60cm	6.5Fr*8cm	-	20Ga	5.5Fr	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS35516	3-5.5Fr-16cm-S	3-5.5Fr-16cm	0.021"*60cm	6.5Fr*8cm	-	20Ga	5.5Fr	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY35520	3-5.5Fr-20cm-Y	3-5.5Fr-20cm	0.021"*60cm	6.5Fr*8cm	20Ga	-	5.5Fr	5mL	5mL	-	0.7mm	11#	Regular	1.7#
ICVM35520	3-5.5Fr-20cm-I	3-5.5Fr-20cm	0.021"*60cm	6.5Fr*8cm	-	20Ga	5.5Fr	5mL	-	5mL	0.7mm	11#	Regular	1.7#
CVMS35520	3-5.5Fr-20cm-S	3-5.5Fr-20cm	0.021"*60cm	6.5Fr*8cm	-	20Ga	5.5Fr	5mL	-	-	0.7mm	11#	Regular	1.7#
CVMY37016	3-7Fr-16cm-Y	3-7Fr-16cm	0.035"*60cm	8.5Fr*10cm	18Ga	-	7Fr	5mL	5mL	-	0.7mm	11#	Regular	2.2#
ICVM37016	3-7Fr-16cm-I	3-7Fr-16cm	0.035"*60cm	8.5Fr*10cm	-	18Ga	7Fr	5mL	-	5mL	0.7mm	11#	Regular	2.2#
CVMS37016	3-7Fr-16cm-S	3-7Fr-16cm	0.035"*60cm	8.5Fr*10cm	-	18Ga	7Fr	5mL	-	-	0.7mm	11#	Regular	2.2#
CVMY37020	3-7Fr-20cm-Y	3-7Fr-20cm	0.035"*60cm	8.5Fr*10cm	18Ga	-	7Fr	5mL	5mL	-	0.7mm	11#	Regular	2.2#
ICVM37020	3-7Fr-20cm-I	3-7Fr-20cm	0.035"*60cm	8.5Fr*10cm	-	18Ga	7Fr	5mL	-	5mL	0.7mm	11#	Regular	2.2#
CVMS37020	3-7Fr-20cm-S	3-7Fr-20cm	0.035"*60cm	8.5Fr*10cm	-	18Ga	7Fr	5mL	-	-	0.7mm	11#	Regular	2.2#
CVMY37030	3-7Fr-30cm-Y	3-7Fr-30cm	0.035"*60cm	8.5Fr*10cm	18Ga	-	7Fr	5mL	5mL	-	0.7mm	11#	Regular	2.2#
ICVM37030	3-7Fr-30cm-I	3-7Fr-30cm	0.035"*60cm	8.5Fr*10cm	-	18Ga	7Fr	5mL	-	5mL	0.7mm	11#	Regular	2.2#
CVMS37030	3-7Fr-30cm-S	3-7Fr-30cm	0.035"*60cm	8.5Fr*10cm	-	18Ga	7Fr	5mL	-	-	0.7mm	11#	Regular	2.2#
CVMY38516	3-8.5Fr-16cm-Y	3-8.5Fr-16cm	0.035"*60cm	9Fr*10cm	18Ga	-	8.5Fr	5mL	5mL	-	0.7mm	11#	Regular	2.5#
ICVM38516	3-8.5Fr-16cm-I	3-8.5Fr-16cm	0.035"*60cm	9Fr*10cm	-	18Ga	8.5Fr	5mL	-	5mL	0.7mm	11#	Regular	2.5#
CVMS38516	3-8.5Fr-16cm-S	3-8.5Fr-16cm	0.035"*60cm	9Fr*10cm	-	18Ga	8.5Fr	5mL	-	-	0.7mm	11#	Regular	2.5#
CVMY38520	3-8.5Fr-20cm-Y	3-8.5Fr-20cm	0.035"*60cm	9Fr*10cm	18Ga	-	8.5Fr	5mL	5mL	-	0.7mm	11#	Regular	2.5#
ICVM38520	3-8.5Fr-20cm-I	3-8.5Fr-20cm	0.035"*60cm	9Fr*10cm	-	18Ga	8.5Fr	5mL	-	5mL	0.7mm	11#	Regular	2.5#
CVMS38520	3-8.5Fr-20cm-S	3-8.5Fr-20cm	0.035"*60cm	9Fr*10cm	-	18Ga	8.5Fr	5mL	-	-	0.7mm	11#	Regular	2.5#
CVMY38530	3-8.5Fr-30cm-Y	3-8.5Fr-30cm	0.035"*60cm	9Fr*10cm	18Ga	-	8.5Fr	5mL	5mL	-	0.7mm	11#	Regular	2.5#

Annex 1: List of Model Specifications of Central Venous Catheter (CVC) Kit														
Product Code	Catalogue Number	Specification												
		Central Venous Catheter	Guide wire	Dilator	Introducer Needle		Fastener	Hypodermic Syringe		Blue introducer syringe	Hypodermic Needle	Scalpel	Injection site cap	Extension line clamp
					Y Model	I Model		I Model	II Model					
ICVM38530	3-8.5Fr-30cm-I	3-8.5Fr-30cm	0.035"*60cm	9Fr*10cm	-	18Ga	8.5Fr	5mL	-	5mL	0.7mm	11#	Regular	2.5#
CVMS38530	3-8.5Fr-30cm-S	3-8.5Fr-30cm	0.035"*60cm	9Fr*10cm	-	18Ga	8.5Fr	5mL	-	-	0.7mm	11#	Regular	2.5#
CVMY48520	4-8.5Fr-20cm-Y	4-8.5Fr-20cm	0.035"*60cm	9Fr*10cm	18Ga	-	8.5Fr	5mL	5mL	-	0.7mm	11#	Regular	2.2#
ICVM48520	4-8.5Fr-20cm-I	4-8.5Fr-20cm	0.035"*60cm	9Fr*10cm	-	18Ga	8.5Fr	5mL	-	5mL	0.7mm	11#	Regular	2.2#
CVMS48520	4-8.5Fr-20cm-S	4-8.5Fr-20cm	0.035"*60cm	9Fr*10cm	-	18Ga	8.5Fr	5mL	-	-	0.7mm	11#	Regular	2.2#
CVMY48530	4-8.5Fr-30cm-Y	4-8.5Fr-30cm	0.035"*60cm	9Fr*10cm	18Ga	-	8.5Fr	5mL	5mL	-	0.7mm	11#	Regular	2.2#
ICVM48530	4-8.5Fr-30cm-I	4-8.5Fr-30cm	0.035"*60cm	9Fr*10cm	-	18Ga	8.5Fr	5mL	-	5mL	0.7mm	11#	Regular	2.2#
CVMS48530	4-8.5Fr-30cm-S	4-8.5Fr-30cm	0.035"*60cm	9Fr*10cm	-	18Ga	8.5Fr	5mL	-	-	0.7mm	11#	Regular	2.2#
Remarks	The number of Injection site caps and Extension line clamps should be consistent with the number of lumens of the catheter.													

Annex 2: List of flow rate of Central Venous Catheter					
Spec.	Name	Main lumen ml/min	Branch lumen ml/min		
			Branch lumen 1	Branch lumen 2	Branch lumen 3
		White distal	Blue medial	Yellow medial	Brown proximal
1-22Ga-08cm		10	/	/	/
1-22Ga-13cm		8	/	/	/
1-22Ga-16cm		7	/	/	/
1-22Ga-20cm		6	/	/	/
1-22Ga-30cm		5	/	/	/
1-20Ga-13cm		20	/	/	/
1-20Ga-16cm		15	/	/	/
1-20Ga-20cm		12	/	/	/
1-18Ga-13cm		35	/	/	/
1-18Ga-16cm		29	/	/	/
1-18Ga-20cm		25	/	/	/
1-16Ga-16cm		60	/	/	/
1-16Ga-20cm		58	/	/	/
1-16Ga-30cm		50	/	/	/
1-14Ga-16cm		110	/	/	/
1-14Ga-20cm		100	/	/	/
1-14Ga-30cm		90	/	/	/
1-12Ga-16cm		130	/	/	/
1-12Ga-20cm		125	/	/	/
1-12Ga-30cm		118	/	/	/
2-4Fr-05cm		19	18		
2-4Fr-08cm		18	17	/	/
2-4Fr-13cm		16	14	/	/
2-4Fr-16cm		14	12	/	/
2-4Fr-20cm		12	9	/	/
2-5Fr-13cm		30	15	/	/
2-5Fr-16cm		25	14	/	/
2-5Fr-20cm		22	13	/	/
2-6Fr-13cm		63	19	/	/
2-6Fr-16cm		58	17	/	/
2-6Fr-20cm		52	14	/	/
2-7Fr-16cm		77	29	/	/
2-7Fr-20cm		70	20	/	/
2-7Fr-30cm		55	12	/	/
2-8Fr-16cm		105	54	/	/
2-8Fr-20cm		96	44	/	/
2-8Fr-30cm		86	31	/	/
3-4Fr-08cm		18	5	5	/

Annex 2: List of flow rate of Central Venous Catheter					
Spec.	Name	Main lumen ml/min	Branch lumen ml/min		
			Branch lumen 1	Branch lumen 2	Branch lumen 3
		White distal	Blue medial	Yellow medial	Brown proximal
3-4Fr-13cm		17	4	4	/
3-4Fr-16cm		15	3	3	/
3-4Fr-20cm		12	2.5	2.5	/
3-5.5Fr-08cm		34	15	15	/
3-5.5Fr-13cm		30	13	13	/
3-5.5Fr-16cm		23	11	11	/
3-5.5Fr-20cm		20	9	9	/
3-7Fr-16cm		60	28	28	/
3-7Fr-20cm		55	16	16	/
3-7Fr-30cm		40	11	11	/
3-8.5Fr-16cm		94	43	43	/
3-8.5Fr-20cm		89	34	34	/
3-8.5Fr-30cm		77	29	29	/
4-8.5Fr-20cm		56	85	29	29
4-8.5Fr-30cm		45	59	19	19

Remark: For catheters with nominal OD less than 1.0mm, the flow rate of each lumen shall not be less than 80% of the nominal value; For catheters with nominal OD greater than or equal to 1.0mm, the flow rate in each lumen shall not be less than 90% of the nominal value.

Specification Confirmation Letter

**Center for Public Procurement in Healthcare
Republic of Moldova**

WHEREAS,

We, **Henan Tuoren Medical Device Co., Ltd.**, having factories at ***Middle of Weft 7 Road, Nanpu District, Changyuan Henan 453400 China***, confirm that the ***Central Venous Catheter Sets*** offered by company FPC SOGNO meet the following requirements:

Product Code	Product name	Auction requirement	TUOREN Confirmation
CVMY25013 CVMY25016 CVMY25020	121. Double Lumen Central Venous Catheter Set 16G	Insertion needle length: ≥ 30 mm	Introducer Needle 20Ga - YModel, length 38mm
ICVM11616 ICVM11620	122. Single Lumen Central Venous Catheter Set 16G	Insertion needle length: ≥ 30 mm	Introducer Needle 18Ga - I Model, length 50mm
CVMY24013 CVMY24016 CVMY24020	123. Double Lumen Central Venous Catheter Set 18G	Insertion needle length: 38-50mm	Introducer Needle 20Ga - YModel, length 38mm
CVMY12013 CVMY12016 CVMY12020	125. Single Lumen Central Venous Catheter Set 20G	Insertion needle length: 38-50mm	Introducer Needle 20Ga - YModel, length 38mm
CVMY12208 CVMY12213 CVMY12216 CVMY12220	127. Single Lumen Central Venous Catheter Set 22G	Insertion needle length: 38-50mm	Introducer Needle 20Ga – Y Model, length 38mm

CVMY37016	129. Triple Lumen Central Venous Catheter Set	Triple lumen: 16G/18G/18G	Lumen sizes: 16G (flow rate 60ml/min) / 18G (flow rate 28ml/min) / 18G (flow rate 28ml/min)
		Guide wire	Guide Wire 0.035"*60cm
		Catheter fixator	Fastener 7Fr
		Heparinized cap	Heparinized Injection site cap
		Lateral clamp	Extension line clamp
		Stopper	Extension line clamp, Heparinized Injection site cap
		Latex & DEHP Free	Made of Latex & DEHP Free Polyurethane
CVMY48520	130. Quadruple Lumen Central Venous Catheter Set	Triple lumen: 14/16/18/18G	Lumen sizes: 14G (flow rate 85ml/min) / 16G (flow rate 56ml/min) / 18G (flow rate 29ml/min) / 18G (flow rate 29ml/min)
		Guide wire	Guide Wire 0.035"*60cm
		Catheter fixator	Fastener 7Fr
		Heparinized cap	Heparinized Injection site cap
		Lateral clamp	Extension line clamp
		Stopper	Extension line clamp, Heparinized site cap
		Latex & DEHP Free	Made of Latex & DEHP Free Polyurethane
CVMY28020 CVMY28030	157. Double Lumen Central Venous Catheter Set 12G	With safety valve that protects against gas embolism or hemorrhage in case of accidental disconnection of the infusion	With needle free positive pressure connector, that protects against gas embolism or hemorrhage in case of accidental disconnection of the infusion
CVMY27020 CVMY27030	158. Double Lumen Central Venous Catheter Set 13G	With safety valve that protects against gas embolism or hemorrhage in case of accidental disconnection of the infusion	With needle free positive pressure connector, that protects against gas embolism or hemorrhage in case of accidental disconnection of the infusion

CVMY37016 CVMY37020 CVMY38516 CVMY38520	171. Triple Lumen Central Venous Catheter Set 16-18G	Ergonomically shaped needle for easy handling, with mark indicating the bevel position	Ergonomically shaped needle for easy handling, with mark indicating the bevel position
		Metal guide: Markings every 10cm on the guide	Stainless Steel Guide Wire: Markings every 10cm on the guide
		Has gauge (G) size markings on the catheter	Yes
		The tip of the catheter is atraumatic (soft type), without the risk of bending during manipulation, the hardness of the catheter ensures correct administration of the drug without traumatizing the tissue	Yes
		Soft, movable and removable catheter fixation wings for correct fixation of the catheter on the skin surface	Yes
		Cap with safety valve that protects against gas embolism or hemorrhage in case of accidental disconnection of the infusion	With needle free positive pressure connector, that protects against gas embolism or hemorrhage in case of accidental disconnection of the infusion
		3-piece syringe, latex free	Hipodermic Syringe 5mL – I Model, 3-piece, latex free. Hipodermic Syringe 5mL – II Model, 3-piece, latex free.

Name: Kofi Guo

Title: Sales manager



Signature & stamp:

Central Venous Catheter (CVC) Kit

Instructions for Use



1 Product Name

Central Venous Catheter (CVC) Kit

2 Product Description

Central Venous Catheter (CVC) Kit contains central venous catheter, guide wire, I model/Y model introducer needle, blue introducer syringe, dilator, injection site cap, extension line clamp, fastener, scalpel, hypodermic syringe, and hypodermic needle. The central venous catheter (CVC) is a long, soft, thin, hollow tube placed into a large vein in the neck (internal jugular vein), chest (subclavian vein) or groin (femoral vein). CVC is intended to provide short-term (≤ 30 days) central venous channel for treatment of patients who requiring central venous channel.

For REF, see also Annex 1.

3 Intended Purpose

Central Venous Catheter (CVC) Kit is used to insert into the central venous system for the infusions of fluids or medications, blood sampling or pressure measurement.

4 Indications

There is no specific disease or medical condition. It is only intended for the diseases or conditions requiring central venous access or lacking usable peripheral IV sites, including but not limited to the following:

- Central venous pressure measurement;
- Blood sampling;
- Total parenteral nutrition (TPN);
- Multiple infusions of fluids, medications or chemotherapy.

5 Clinical benefits

CVC has the ability to gain access to the central circulation system through a single puncture site for applications that include infusion of fluids or medications, blood sampling or pressure measurement.

- CVC could be used for short-term infusions, compared to transient devices.
- Multi-lumen catheters could establish multiple channels compared to single-lumen catheters.

6 Contra-indications

There is no absolute contraindication, but it should be used with caution in patients with known hemorrhagic diseases.

7 Complications

Including but not limited to the following:

- Cardiac tamponade secondary to vascular, atrial or ventricular perforation;
- Pleural (i.e., pneumothorax, hemothorax) and mediastinal injuries;
- Air embolism;

- Catheter embolism;
- Catheter occlusion;
- Thoracic duct lacerations;
- Bacteraemia;
- Septicemia;
- Thrombosis;
- Inadvertent arterial puncture;
- Nerve injury;
- Hematoma;
- Hemorrhage;
- Fibrin sheath formation;
- Exit site infection;
- Vessel erosion;
- Catheter tip malposition;
- Dysrhythmias;
- Extravasation

8 Patient target groups

The device is intended for patients fitting the intended uses. It can be used for the patients over 1yr old (i.e. excluding Neonates and Infants).

9 Intended users


This product can be used by professional doctors, who has been properly trained for this operation.


10 Instructions for use

a) Read the instructions for use carefully before use and select the appropriate specification according to the *list of model specifications in Annex 1* and the *list of flow rates of CVC in Annex 2*.

b) Prepare skin for the puncture site and pave sterile sheet;

c) Before surgery, inject normal saline or heparin saline into the lumen to check whether the CVC is abnormal. Then, clamp the extension tube or put a Injection site cap on the extension tube hub, and the extension tube connecting the distal lumen should be kept open to pass through the guide wire.

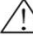 Refer to the following figure for the identification of hub and lumen:

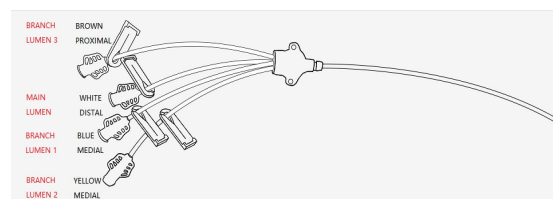


d) Insert the introducer needle into the vein and draw back to confirm that there is good venous blood reflux, but the

color of the reflux blood is not always the only dependable sign to check whether the introducer needle enters into the vein;

e) Advance the guide wire into vein through introducer needle with the aid of control handle until the required depth is met, and the "J" shaped head may need a light twirl to enter. **Do not have any angle between the guide wire and the introducer needle to deliver or withdraw the guide wire;**


⚠ Before using the control handle, gently press the guide wire at the exposed part of the handle (as shown in figure ① below), and then pull out the protective cap (as shown in figure ② below) to prevent the guide wire from slipping out of the handle; The control handle should be handled with care, do not remove by force.

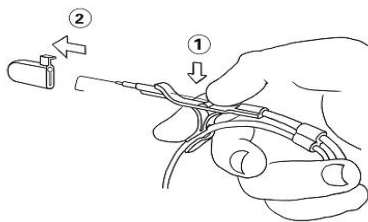

f) Hold the guide wire tightly and take out the introducer needle, and hold the guide wire tightly at any time;

g) Insert the dilator into the blood vessel along the guide wire, and turn the dilator, and gently dilate the channel in order to insert the catheter. For withdrawal of the dilator (while inserting and removing the dilator, keep the guide wire fixed), do not leave the dilator in its original position as the indwelling catheter, so as to avoid blood vessel wall perforation;

h) Hold the distal end of the catheter pass through the guide wire, when it is close to the skin, hold tightly the distal end of the catheter and then turn gently it into vein;

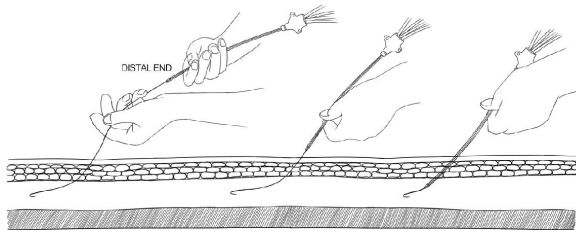

i) Advance the catheter into the indwelling position using the length mark on the catheter as a positioning reference, then record the catheter length and check the catheter position regularly;

⚠ Marking system of the catheter: Indicate distance from the distal end. From the first mark, the distance between marks is 1cm.

j) Pull out the guide wire after the catheter reaches a predetermined depth. When resistance encounters, take it slowly, and if necessary, turn the guide wire at a certain angle to move it away from a certain distance; if it still cannot be taken out, the catheter and the guide wire should be taken out together. Do not take them out forcibly since the guide wire may be broken in the human body. Check

whether the guide wire is intact after it is taken out;

k) Connect a syringe to the extension tube hub, and check the catheter position, observe whether the blood flow is unobstructed when taking out, then connect a qualified connecting tube to the extension tube hub as required. Block off the unused catheter lumen with a injection site cap. Stop the blood flow with a Extension line clamp when replacing the connecting tube and Injection site cap;

l) Check the catheter distal tip position with X-ray. Fix the catheter to the patient with sutures through fasteners. Do not suture directly outside the catheter to prevent the catheter from being punctured. Cover the puncture site with medical gauze, and record the length of catheter insertion, and check frequently whether the catheter is moved;

m) Nursing during the catheter indwelling:

1) If the three-way valve or injection site cap used at the connector is found loosening or falling off, remove it immediately, and replace it with a new three-way valve or Injection site cap. All operations at the connector such as infusion, drug administration, blood drawing and connection of infusion pump should strictly comply with aseptic operations to prevent the iatrogenic infection. For patients receiving intravenous high nutrient solution, strengthen the inspection during the period of infusion; after the infusion, flush the lumen with normal saline or replace the infusion tube with a new one before infusing other fluids;

2) When monitoring CVP for a long period of time, the catheter should be flushed with 5 ml (30 u) of dilute heparin solution prepared with saline every 24 h to keep the pressure measurement system unobstructed and reduce the incidence of infection and prevent the formation of fibrin sheath around the catheter distal end;

⚠ The residual catheter track remains the air entry point. When the catheter is successfully placed, the medical dressing used to close the wound should be kept in place until the wound is epithelialized. After the wound is epithelialized, the medical dressing used during catheter indwelling is recommended to be replaced every 24 hours to keep the wound dry and clean.

⚠ When removing the Injection site cap, holding the flanking part of the catheter hub can effectively reduce the disassembly resistance.

n) Pay attention to the patient's condition when removing the catheter. If there is a serious infection, use other methods to remove the catheter, such as surgery. Do not forcibly pull out the catheter when removing the catheter in a regular way to prevent the catheter from breaking in the human body.

⚠ Ask patient to take a breath and hold it if removing internal jugular or subclavian catheter.

11 Performance Characteristics

Central Venous Catheter is accessory devices used during interventional procedures. Clinical performance is determined by whether the catheter successfully placed into the central venous circulatory system.

For the performance characteristics of catheter:

- The catheter is made TPU material with good physical properties, smooth surface, peak tensile force, easy for

puncture and placement.

- The catheter is visible under X-ray, accurate positioning.

12 Precautions
















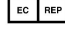




- a) Avoid repeating insertion, or it might cause hematoma and infection;
- b) This product is supplied sterile and sterilized by ethylene oxide. The shelf life of this product is 3 years. Please use it within the shelf life;
- c) Please check whether the package is intact before use. Do not use if it is damaged;
- d) This product is for single use. Do not reuse, reprocess or re-sterilize. Reuse of the device poses a potential risk of serious injury and/or infection. Dispose of the product waste as medical waste according to local regulations or hospital management regulations after use.
- e) The product can only be used by professional doctors, who has been properly trained for this operation.
- f) The product is allowed to be indwelled in human body for 30 days at most, and 14 days the best;
- g) The device used together with this product should have a standard luer lock.
- h) Adverse side effects: improper care of the puncture site may lead to infection; prolonged indwelling time may lead to phlebitis or thrombophlebitis.
- i) Do not press harshly to avoid the device damage;
- j) Fragile, handle with care;
- k) Store in a ventilated and dry place to avoid moisture as much as possible, and keep away from sunlight to avoid fire or heat deformation of the device;
- l) Intended condition of use: this product is use for operating room or catheter room environment. (Temperature range: 21 to 25 °C , relative humidity range : 30% to 60% RH, Atmospheric pressure range: 70 to 104 kPa). Please do not use it in the environment beyond the above conditions;
- m) Please note that fluid leakage during catheter indwelling may result in blood loss or delay in treatment. Please pay attention to check whether the catheter leakage during indwelling. Replace the catheter immediately once leakage is found.
- n) Do not place catheter into or allow it to remain in the right atrium or right ventricle. X-ray exam or other method in compliance with institutional policies and procedures must show catheter tip located in lower 1/3 of the Superior Vena Cava (SVC), in accordance with institutional guidelines.
- o) Do not use too much force when clamping. Also avoid repeated clamping in the same location which may weaken the catheter.
- p) This product is not for high pressure injection for such application can result in inter-lumen crossover or rupture with potential for injury.

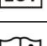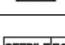13 Warnings

When selecting nitroglycerin, vitamin K1, paclitaxel, shenmai, fluorouracil and fat emulsion for drug compatibility studies, this product has been proven with good compatibility and no adsorption of drugs, which can ensure the drug efficacy.

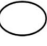The cobalt is above 0.10% w/w in 304 stainless steel of the Guide wire, Blue introducer syringe and Hypodermic needle.

Explanation of graphics, symbols and abbreviations of medical devices used in this product:

	CE mark		Manufacturer
	Authorized representative in the European Community/European Union		Catalogue number
	Use-by date		Keep away from sunlight
	Keep dry		Do not re-use
	Date of manufacture		Batch code
	Caution		Consult instructions for use or electronic instructions for use
	Sterilized using ethylene oxide		Don't use if package is damaged, and consult instructions for use
	Medical device		Do not re-sterilize
	Fragile, handle with care		Single sterile barrier system
	Unique device identifier		Contains hazardous substances: Cobalt (CAS NO.7440-48-4)

14 Additional Information

The SSCP of this product will be available in the European database on medical devices (Eudamed), where it is linked to the Basic UDI-DI (694242711359FE). The URL to the Eudamed public website is: <https://ec.europa.eu/tools/eudamed>

For a patient/user/third party in the European Union and in countries with identical regulatory regime (Regulation 2017/745/EU on Medical Devices); if during the use of this device or as a result of its use, a serious incident has occurred, please report it to the manufacturer and/or its authorized representative and to your national authority. The contacts of national competent authorities (Vigilance Contact Points) and further information can be found on the following European Commission website: https://ec.europa.eu/growth/sectors/medical-devices/contact_s_en

Manufacturer: Henan Tuoren Best Medical Device Co., Ltd.

Address of manufacture: Middle of West 7 Road, Nanpu District, Changyuan, Henan, 453400, China

Tel: +86 0373-8814000 **Fax:** +86 0373-8816222

Email: info@etuoren.com **Website:** www.tuoren.com

EU Representative: Lotus NL B.V.

Address: Koningin Julianaplein 10, 1e Verd, 2595AA, The Hague, Netherlands.

Tel: +31644168999 **Fax:** 008621-58857798

Email: peter@lotusnl.com

Document No.: BST-CE-06-04

Version: A2

Date: 2024-10-25

Website for eIFU: en.tuoren.com/ifu