

RELIABLE HEARING IN ANY SITUATION

PRODUCT CATALOG



REXTON



RELY ON REXTON

Life becomes a challenge when hearing starts to go. At Rexton, we won't let hearing loss hold people back. That's why we work hard to deliver proven hearing technology in the most reliable way.

Rexton hearing aids get the job done every single time. They will perform reliably at home, at work, during exercise, and users can wear them in all weather. That's because they are built tough, using the highest quality materials, and tested rigorously to make sure they can answer all the challenges that life throws at you.

Noisy rooms? Drop them accidentally? Caught in the rain? Working in a harsh environment? Don't worry about it. And they can handle much more. Simply put, they'll keep users prepared for every situation. They are made Lifeproof.

RELY ON
LIFEPROOF
REXTON



A ROBUST PORTFOLIO FOR RELIABLE HEARING



Reach technology

REACH R-Li and R-Li T 9

Every Voice Matters

REACH SR 10

Sleek and robust

REACH inoX-CIC Li 12

The first invisible hearing aid with rechargeability

REACH Custom IIC/CIC 14

A next-level custom made invisibility



ReCharge technology

RECHARGE R-Li 18

New Energy for every occasion

RECHARGE B-Li P 20

New Energy for every occasion

RECHARGE B-Li HP 22

New Energy for every occasion



BiCore technology

BICORE B-Li M Rugged 26

Durability and a smooth listening experience

BICORE B M, P and HP 28

A full range of BTEs to rely on

BICORE R-Li, R-Li T and R312 32

A full range of RICs

BICORE SR 36

Smaller and lighter

BICORE Custom-Li 38

Our first rechargeable custom with BiCore technology and connectivity





MotionCore technology

M-CORE B-Li M, P, HP 42

Versatile, powerful and rechargeable

M-CORE R-Li and R312 45

Trusted RIC solutions

M-CORE SR 48

Discreet design with exchangeable receivers

M-CORE iX-CIC 49

With the hidden power of MotionCore



MyCore technology

MYCORE inoX CIC 8C 52

Discreet and ready to fit

MYCORE STERLING 8C 53

Proven design – reliable technology

MYCORE EMERALD S 8C 55

Variety of accessories offers full control and usability



Accessories & Apps 59



Portfolio overview 86

Glossary 88



SOUNDPRO 3.0

SoundPro 3.0 automatically steers all sound processing features, including the new Multi-Voice Focus* based on the analysis of the acoustic scene. When changing the strategy to Natural Focus or Speech Focus, the recommended balance of features is modified. The Balanced Focus is recommended for enhanced automatic performance.



DIRECT STREAMING

Reach hearing aids can stream calls and audio from iPhone via MFi and Android via ASHA and BT LE AUDIO**.



TELECARE

Full live remote fitting sessions via Connexx resulting in higher user convenience and satisfaction.



HEADSET MODE

Rexton Reach microphones get the job done so iOS and BT LE Audio users can enjoy hands free call.



MY VOICE 2.0

General sound and own voice processed in parallel for a greater sound impression and user satisfaction.



AUTO ECHO REDUCER

Automated environment scan to reduce effects of disturbing reverberation on the perceived sound signal.



ATEX

Rexton Reach have ATEX*** certification which verifies that they are tested and approved to be intrinsically safe for use in potentially explosive environments.

* Only available for hearing aids with two microphones in a binaural fitting.

** Phones with Bluetooth 5.0+ and Android 10 (or higher) are capable of providing wireless streaming via ASHA, but ASHA must be enabled by each phone manufacturer and service carrier.

*** They can be safely used in areas where there may be flammable substances in the air such as paint and chemical manufacturing, petrol stations and mills. The ATEX group/class is printed on the ATEX card, which must be available at hand when the hearing aids are used in explosive atmospheres.



REACH TECHNOLOGY



Every voice matters, and Rexton Reach makes sure every voice is heard. It is the rechargeable hearing aid with Multi-Voice Focus that enables users to keep track of multiple voices simultaneously, even when there is background noise. And when the user needs to reach out to someone much further away, no problem. Reach connects seamlessly with smartphones.

MULTI-VOICE FOCUS

Multi-Voice Focus scans the surrounding environment 1000 times a second with four independent focus beams. This enables it to constantly adapt to the changing dynamics of the conversation and make the experience fluid and comfortable. It will boost the prominence of those in the group who are speaking while keeping track of those who are silent, and all without switching between modes, and with no delays.

VOICE STABILIZER

Voice Stabilizer controls the compression to secure that all voices within the focus beam are amplified in the best linear ways when the sound goes too loud. It keeps the signal-to-noise ratio positive to help users focusing on the signal with most contrast. Benefit: voices remain untouched while surrounding noise is managed.

CONNECTIVITY – LE AUDIO

Reach has MFi and ASHA protocols, so it connects seamlessly with smartphones. It is also future proof and will require only a simple firmware update to connect via Bluetooth® LE Audio.

LIFEPROOF

Rexton hearing aids are built using robust components and tested extensively to ensure they can withstand whatever life throws at them. You can rely on them to withstand everyday challenges and harsh environments.

REACH R-LI AND R-LI T


Rexton Reach is the rechargeable hearing aid that enables users to keep track of multiple voices simultaneously while connecting seamlessly to smart devices.

- › Features Multi-Voice Focus for greater benefit when conversing in a group
- › Rechargeable Li-ion battery for a full days use.
R-Li has 24h runtime (incl. 5h of streaming), and R-Li T has 34h runtime (incl. 5h streaming)
- › Charging on-the-go with built-in power bank in the travel charger
- › Programable via Noahlink Wireless
- › ATEX
- › Lifeproof
- › Water resistant, IP68
- › R-Li and R-Li T offers Call control when connected with Headset mode for Mfi and LE Audio

Colors



Fitting range is based on 2 ccm coupler:

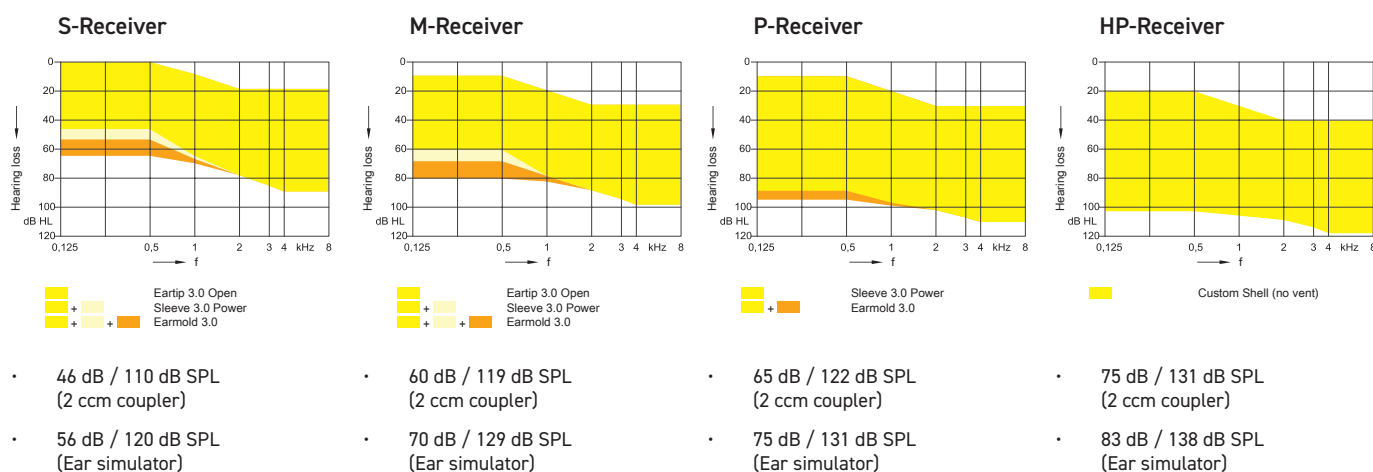

Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	48/20/6	36/18/6	32/16/6	24/12/4	16/8/4
Battery	Lithium-ion	Lithium-ion	Lithium-ion	Lithium-ion	Lithium-ion
T-coil	R-Li T	R-Li T	R-Li T	R-Li T	R-Li T
Bluetooth	•	•	•	•	•
Tinnitus function (Sound therapy and Notch)	•	•	•	•	-
Wireless sync	•	•	•	•	•
IP68	•	•	•	•	•
ATEX	•	•	•	•	•
EarWear	3.0	3.0	3.0	3.0	3.0
Technology Highlights					
Soundpro 3.0	•	•	•	•	•
Multi-Voice Focus	•	•	•	-	-
Focus Streams	4	2	1	1	1
Voice Stabilizer	•	•	•	•	•
Speech Preservation Technology	•	•	•	•	•
My Voice 2.0	•	•	•	-	-
Motion Sensor for scene analysis	•	•	•	-	-
Directionality					
Automatic adaptive	•	•	•	•	•
iOmni Sound Locator	•	•	•	•	•
iFocus 360 (Front/Back)	Automatic & Manual	Automatic & Manual	-	-	-
iFocus 360 (Left /Right)	Automatic & Manual	Manual	-	-	-
Auto Stereo iLock (Narrow)	•	•	•	-	-
Noise Reduction					
Noise Management	•	•	•	•	•
Sound Smoothing	•	•	•	•	•
Wind Noise Reduction	•	•	•	•	•
Auto Echo Reducer	•	-	-	-	-
Reverb Reducer	Manual	Manual	-	-	-
Sound Quality & Control					
Bandwidth extension	•	-	-	-	-
Bandwidth compression	•	•	•	•	•
Feedback preventer	•	•	•	•	•
Music Enhancer (presets)	3	3	1	1	-
Smart Automatic Equalizer	•	•	•	•	•
Smart Automatic Acclimatization	•	•	•	•	•
Automatic Classifier	•	•	•	•	•
Data Logging	•	•	•	•	•
Audio Streaming					
BT LE Audio	•	•	•	•	•
Made for iPhone	•	•	•	•	•
ASHA (Android 10 or higher)	•	•	•	•	•
Headset Mode					
iOS	•	•	•	•	•
BT LE Audio	•	•	•	•	•
App					
Rexton App	•	•	•	•	•
Spatial configurator	•	•	-	-	-
Rexton Assist	•	•	•	•	•
TeleCare					
Full live remote tuning	•	•	•	•	•
Basic portal remote tuning	•	•	•	•	•
Accessories		Chargers		CROS	
Smart Key	•	Charging Station	•	Reach CROS R-Li	•
Smart Transmitter 2,4	•	Travel Charger	•	Reach CROS SR	-
Smart Mic	•	Charging+ Station	•	Reach CROS iX-CIC Li	-
		IF CIC Travel Charger	-		
		Universal Cable Charger	•		
Fitting					
Noahlink Wireless	•				
ConnexxAir	-				

• available - not available

REACH SR

Sleek and yet robust, this is what you get with Rexton Reach SR. Combining the Reach technology with the Lifeproof of Rexton in an expressive design.



- › Features Multi-Voice Focus for greater benefit when conversing in a group
- › Rechargeable Li-ion battery for a full days use with 17h runtime (incl. 5h of streaming)
- › Charging on-the-go with built-in power bank in the charger
- › Programable via NoahLink Wireless
- › ATEX
- › Lifeproof
- › Water resistant, IP68
- › SR offers Tap control when connected with Headset mode for Mfi and LE Audio

Colors

Black | Black Gloss



Cosmic Blue | Rose Gold



Snow White | Rose Gold






Snow White | Snow White Gloss



Black | Silver



Fitting range is based on 2 ccm coupler:

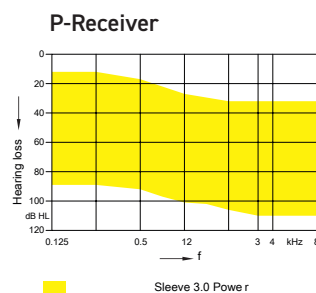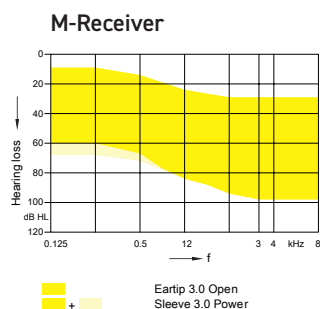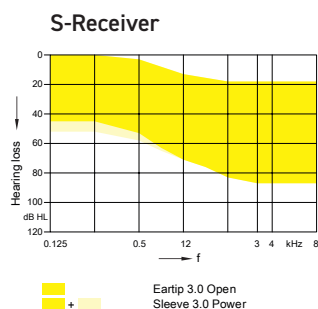

- 46 dB / 110 dB SPL (2 ccm coupler)
- 56 dB / 120 dB SPL (Ear simulator)

- 60 dB / 119 dB SPL (2 ccm coupler)
- 70 dB / 129 dB SPL (Ear simulator)

- 65 dB / 122 dB SPL (2 ccm coupler)
- 75 dB / 131 dB SPL (Ear simulator)

Features / Performance Level		80	60	40	
Channels/Controls/Programs		48/20/6	36/18/6	32/16/6	
Battery		Lithium-ion	Lithium-ion	Lithium-ion	
T-coil		-	-	-	
Bluetooth		•	•	•	
Tinnitus function (Sound therapy and Notch)		•	•	•	
Wireless sync		•	•	•	
IP68		•	•	•	
ATEX		•	•	•	
EarWear		3.0	3.0	3.0	
Technology Highlights					
Soundpro 3.0		•	•	•	
Multi-Voice Focus		•	•	•	
Focus Streams		4	2	1	
Voice Stabilizer		•	•	•	
Speech Preservation Technology		•	•	•	
My Voice 2.0		•	•	•	
Motion Sensor for scene analyzis		•	•	•	
Directionality					
Automatic adaptive		•	•	•	
iOmni Sound Locator		•	•	•	
iFocus 360 (Front/Back)		Automatic & Manual	Automatic & Manual	-	
iFocus 360 (Left /Right)		Automatic & Manual	Manual	-	
Auto Stereo iLock (Narrow)		•	•	•	
Noise Reduction					
Noise Management		•	•	•	
Sound Smoothing		•	•	•	
Wind Noise Reduction		•	•	•	
Auto Echo Reducer		•	-	-	
Reverb Reducer		Manual	Manual	-	
Sound Quality & Control					
Bandwidth extension		•	-	-	
Bandwidth compression		•	•	•	
Feedback preventer		•	•	•	
Music Enhancer (presets)		3	3	1	
Smart Automatic Equalizer		•	•	•	
Smart Automatic Acclimatization		•	•	•	
Automatic Classifier		•	•	•	
Data Logging		•	•	•	
Audio Streaming					
BT LE Audio		•	•	•	
Made for iPhone		•	•	•	
ASHA (Android 10 or higher)		•	•	•	
Headset Mode					
iOS		•	•	•	
BT LE Audio		•	•	•	
App					
Rexton App		•	•	•	
Spatial configurator		•	•	-	
Rexton Assist		•	•	•	
TeleCare					
Full live remote tuning		•	•	•	
Basic portal remote tuning		•	•	•	
Accessories		Chargers		CROS	
Smart Key	•	Charging Station	-	Reach CROS R-Li	-
Smart Transmitter 2,4	•	Travel Charger	•	Reach CROS SR	•
Smart Mic	•	Charging+ Station	•	Reach CROS iX-CIC Li	-
		IF CIC Travel Charger	-		
		Universal Cable Charger	-		
Fitting					
Noahlink Wireless	•				
ConnexxAir	-				

• available - not available

REACH inoX-CIC Li

THE FIRST INVISIBLE HEARING AID WITH RECHARGEABILITY

The first invisible hearing aid with rechargeability, and still reliable!
Reduced complexity and rechargeability in a nearly invisible solution.



- > Virtually invisible design for utmost discretion
- > Rechargeable Li-Ion battery for a full days use with up to 28h runtime
- > Rexton App compatible
- > Charging "on-the-go" with build-in power bank
- > Programmable via ConnexxAir
- > ATEX
- > Water resistant, IP68
- > Lifeproof
- > Compatible with earmold 3.0

Faceplate colors


Black



Mocha



Fitting range is based on 2 ccm coupler:

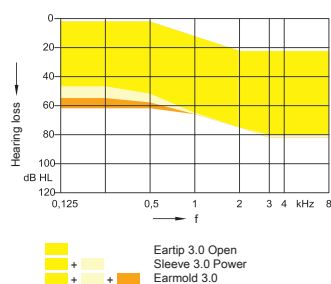

- 50 dB / 114 dB SPL
(2 ccm coupler)
- 61 dB / 126 dB SPL
(Ear simulator)



Features / Performance Level		80	60	40	30	20	
Channels/Controls/Programs		48/20/6	36/18/6	32/16/6	24/12/4	16/8/4	
Battery		Lithium-ion	Lithium-ion	Lithium-ion	Lithium-ion	Lithium-ion	
T-coil		-	-	-	-	-	
Bluetooth		-	-	-	-	-	
Tinnitus function (Sound therapy and Notch)		•	•	•	•	-	
Wireless sync		•	•	•	•	•	
IP68		•	•	•	•	•	
ATEX		•	•	•	•	•	
EarWear		3.0	3.0	3.0	3.0	3.0	
Technology Highlights							
Soundpro 3.0		•	•	•	•	•	
Multi-Voice Focus		-	-	-	-	-	
Focus Streams		-	-	-	-	-	
Voice Stabilizer		•	•	•	•	•	
Speech Preservation Technology		-	-	-	-	-	
My Voice 2.0		-	-	-	-	-	
Motion Sensor for scene analyzis		-	-	-	-	-	
Directionality							
OneMic directionality		•	•	•	-	-	
Automatic adaptive		•	•	•	•	•	
iOmni Sound Locator		-	-	-	-	-	
iFocus 360 (Front/Back)		-	-	-	-	-	
iFocus 360 (Left /Right)		-	-	-	-	-	
Auto Stereo iLock (Narrow)		-	-	-	-	-	
Noise Reduction							
Noise Management		•	•	•	•	•	
Sound Smoothing		•	•	•	•	•	
Wind Noise Reduction		•	•	•	•	•	
Auto Echo Reducer		•	-	-	-	-	
Reverb Reducer		Manual	Manual	-	-	-	
Sound Quality & Control							
Bandwidth extension		•	-	-	-	-	
Bandwidth compression		•	•	•	•	•	
Feedback preventer		•	•	•	•	•	
Music Enhancer (presets)		3	3	1	1	-	
Smart Automatic Equalizer		•	•	•	•	•	
Smart Automatic Acclimatization		•	•	•	•	•	
Automatic Classifier		•	•	•	•	•	
Data Logging		•	•	•	•	•	
Audio Streaming							
BT LE Audio		-	-	-	-	-	
Made for iPhone		-	-	-	-	-	
ASHA (Android 10 or higher)		-	-	-	-	-	
Headset Mode							
iOS		-	-	-	-	-	
BT LE Audio		-	-	-	-	-	
App							
Rexton App		•	•	•	•	•	
Spatial configurator		-	-	-	-	-	
Rexton Assist		-	-	-	-	-	
TeleCare							
Full live remote tuning		-	-	-	-	-	
Basic portal remote tuning		•	•	•	•	•	
Accessories		Chargers		CROS			
Smart Key	•	Charging Station		-		Reach CROS R-Li	-
Smart Transmitter 2,4	-	Travel Charger		-		Reach CROS SR	-
Smart Mic	-	Charging+ Station		-		Reach CROS iX-CIC Li	•
		IF CIC Travel Charger		•			
		Universal Cable Charger		-			
Fitting							
Noahlink Wireless	-						
ConnexxAir	•						

• available - not available

REACH Custom IIC AND CIC

A NEXT-LEVEL CUSTOM MADE INVISIBLE

The Reach Custom IIC and CIC are designed to optimize the perception of small size in several ways and with our new I3 faceplate the invisibility of these devices are next level.



- > Our smallest customized solution
- > 36% better battery lifetime compared to predecessor (~75h)
- > Ultra HD e2e binaural link enables OneMic Directionality
- > Rexton App compatible
- > Programmable via ConnexxAir and/or programming adaptor
- > Lifeproof

Colors

Black



Brown




Mocha




Fitting range is based on 2 ccm coupler:

S-Receiver*

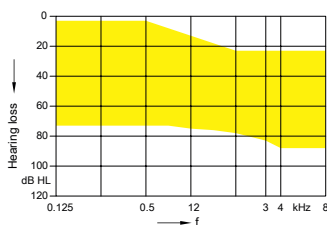


- 50 dB / 113 dB SPL (2 ccm coupler)
- 60 dB / 124 dB SPL (Ear simulator)

M-Receiver**

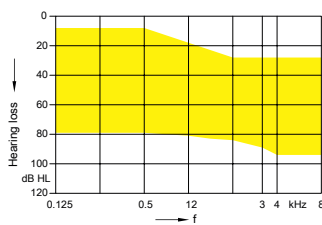

- 55 dB / 118 dB SPL (2 ccm coupler)
- 65 dB / 127 dB SPL (Ear simulator)

P-Receiver**

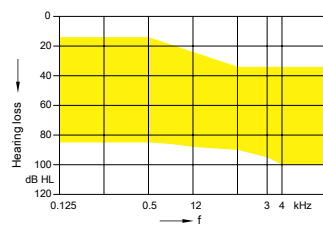

- 60 dB / 124 dB SPL (2 ccm coupler)
- 70 dB / 133 dB SPL (Ear simulator)

* Available for both IIC and CIC

** Available for CIC

Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	48/20/6	36/18/6	32/16/6	24/12/4	16/8/4
Battery	10	10	10	10	10
T-coil	-	-	-	-	-
Bluetooth	-	-	-	-	-
Tinnitus function (Sound therapy and Notch)	•	•	•	•	-
Wireless sync	•	•	•	•	•
IP68	-	-	-	-	-
ATEX	-	-	-	-	-
EarWear	3.0	3.0	3.0	3.0	3.0
Technology Highlights					
Soundpro 3.0	•	•	•	•	•
Multi-Voice Focus	-	-	-	-	-
Focus Streams	-	-	-	-	-
Voice Stabilizer	•	•	•	•	•
Speech Preservation Technology	-	-	-	-	-
My Voice 2.0	-	-	-	-	-
Motion Sensor for scene analysis	-	-	-	-	-
Directionality					
OneMic directionality	•	•	•	-	-
Automatic adaptive	•	•	•	•	•
iOmni Sound Locator	-	-	-	-	-
iFocus 360 (Front/Back)	-	-	-	-	-
iFocus 360 (Left /Right)	-	-	-	-	-
Auto Stereo iLock (Narrow)	-	-	-	-	-
Noise Reduction					
Noise Management	•	•	•	•	•
Sound Smoothing	•	•	•	•	•
Wind Noise Reduction	•	•	•	•	•
Auto Echo Reducer	•	-	-	-	-
Reverb Reducer	Manual	Manual	-	-	-
Sound Quality & Control					
Bandwidth extension	•	-	-	-	-
Bandwidth compression	•	•	•	•	•
Feedback preventer	•	•	•	•	•
Music Enhancer (presets)	3	3	1	1	-
Smart Automatic Equalizer	•	•	•	•	•
Smart Automatic Acclimatization	•	•	•	•	•
Automatic Classifier	•	•	•	•	•
Data Logging	•	•	•	•	•
Audio Streaming					
BT LE Audio	-	-	-	-	-
Made for iPhone	-	-	-	-	-
ASHA (Android 10 or higher)	-	-	-	-	-
Headset Mode					
iOS	-	-	-	-	-
BT LE Audio	-	-	-	-	-
App					
Rexton App	•	•	•	•	•
Spatial configurator	-	-	-	-	-
Rexton Assist	-	-	-	-	-
TeleCare					
Full live remote tuning	-	-	-	-	-
Basic portal remote tuning	•	•	•	•	•
Accessories		CROS			
Smart Key	•	Reach CROS R-Li	-		
Smart Transmitter 2.4	-	Reach CROS SR	-		
Smart Mic	-	Reach CROS iX-CIC Li	•		
Fitting					
Noahlink Wireless	-				
ConnexxAir	•				
Programming adaptor	•				

• available - not available

RECHARGE TECHNOLOGY



Rexton ReCharge will bring new energy to hearing situations everywhere. Bringing proven hearing aid technology and featuring a handy cable charger for recharging at home or on the move. No more hassle with tiny batteries. Just plug in, power up and away you go.

EASY TO USE

No more hassle with tiny batteries since ReCharge is rechargeable. They are easy to fit and programmable with Noahlink Wireless.

AFFORDABILITY

The first rechargeable hearing aid with proven hearing technology in the essential portfolio. The charger is designed for optimal usability and affordability.

RECHARGEABILITY

With a unique performance level offering, you can recharge at home or on the go. Plug in and power up with the handy cable-charger which can recharge one or two hearing aids.

LIFEPROOF

Rexton hearing aids are built using robust components and tested extensively to ensure they can withstand whatever life throws at them. They are durable enough to withstand everyday challenges and harsh environments. All ReCharge products are IP68 and ATEX certified.

TECHNOLOGY

- › Localization in quiet environments is improved by mimicking the ear acoustics.
- › Extended Dynamic Range
- › High input sounds with improved sound quality for less distortion.
- › Adaptive compression for more natural speech.
- › 48-channel first fit for a more precise and easy fitting from scratch.



RECHARGE R-LI

JUST PLUG IN AND POWER UP

The simplicity and affordability of Rexton ReCharge will have a major impact. New energy for work, hobbies, travel, nights out, for time with the family, and for conversations everywhere. Rexton ReCharge is powered by proven technology. No more hassles with handling and changing tiny batteries. Just plug in and power up, and everything is ready to go.




- > Programmable via Noahlink Wireless
- > Takes you through the day and night with 35h runtime
- > Unique rechargeability offering in the 30, 20, 10 and 5 performance levels
- > Addresses mild to severe hearing loss
- > Water-resistant, IP68 rating
- > Lifeproof
- > ATEX

Colors



Fitting range is based on 2 ccm coupler:


S-Receiver

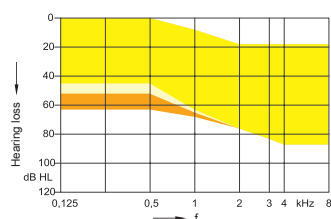

Eartip 2.0 open
 + Sleeve 2.0
 + Earmold 2.0

- 45 dB / 108 dB SPL (2 ccm coupler)
- 56 dB / 119 dB SPL (Ear simulator)


M-Receiver

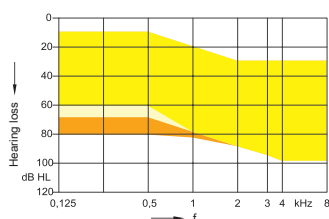

Eartip 2.0 open
 + Sleeve 2.0
 + Earmold 2.0

- 60 dB / 119 dB SPL (2 ccm coupler)
- 70 dB / 129 dB SPL (Ear simulator)


P-Receiver

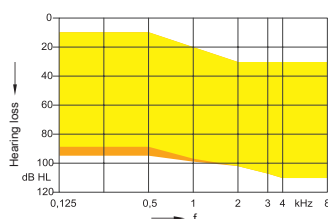

Sleeve 2.0
 + Earmold 2.0

- 70 dB / 124 dB SPL (2 ccm coupler)
- 80 dB / 134 dB SPL (Ear simulator)

HP-Receiver

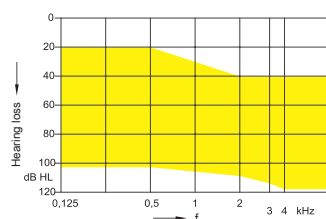

Custom Shell (no vent)

- 75 dB / 130 dB SPL (2 ccm coupler)
- 82 dB / 138 dB SPL (Ear simulator)

Features / Performance Level	30	20	10	5
Channels/Controls/Programs	16/8/4	16/8/4	12/6/4	8/4/4
Battery	Lithium-ion	Lithium-ion	Lithium-ion	Lithium-ion
T-coil	-	-	-	-
Tinnitus Sound therapy	•	-	-	-
Wireless sync	•	•	-	-
IP68	•	•	•	•
ATEX	•	•	•	•
EarWear	miniReceiver 2.0 + click dome/sleeve/Earmold			
Technology Highlights				
First Fit 48 channels	•	•	•	•
Voice Ranger	•	•	-	-
Adaptive compression	•	•	•	•
Extended Dynamic Range	•	•	•	•
Directionality				
Multichannel Adaptive Directional Microphone	•	-	-	-
Automatic Directional Microphone	•	•	•	-
Fixed Directional Microphone	•	•	•	•
iOmni Sound Locator	•	•	-	-
Noise Reduction				
Noise Management	•	•	•	•
Sound Smoothing	•	-	-	-
Wind Noise Reduction	•	-	-	-
Sound Quality & Control				
Bandwidth compression	•	•	-	-
Feedback preventer	•	•	•	•
Smart Automatic Equalizer	•	•	-	-
Smart Automatic Acclimatization	•	•	-	-
Automatic Classifier	-	-	-	-
Data Logging	•	•	•	•
App				
Rexton App	•	•	•	•
TeleCare				
Full live remote tuning	-	-	-	-
Basic portal remote tuning	•	•	•	•

Accessories	
Smart Key	•
Smart Transmitter 2,4	-
Smart Mic	-

Fitting	
Noahlink Wireless	•
ConnexxAir	-

Chargers	
Charging Station	-
Travel Charger	-
Charging+ Station	-
IF CIC Travel Charger	-
Universal Cable Charger	•

• available - not available

RECHARGE B-LI P

JUST PLUG IN AND POWER UP

The simplicity and affordability of Rexton ReCharge B-Li P will secure new energy for an entire day. Moderate to severe hearing losses can benefit from proven technology for clear and comfortable sounds.




- > Programmable via Noahlink Wireless
- > Takes you through the day and night with 39h runtime
- > Unique rechargeability offering in the 30, 20, 10 and 5 performance levels
- > Addresses mild to severe hearing loss
- > Water-resistant, IP68 rating
- > Lifeproof
- > ATEX

Colors



Fitting range is based on 2 ccm coupler:

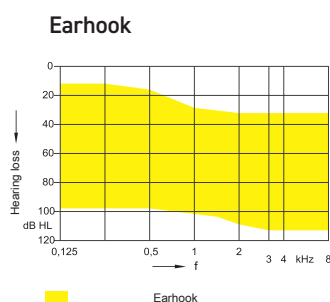

- 75 dB / 134 dB SPL
(2 ccm coupler)
- 79 dB / 139 dB SPL
(Ear simulator)

Features / Performance Level	30	20	10	5
Channels/Controls/Programs	16/8/4	16/8/4	12/6/4	8/4/4
Battery	Lithium-ion	Lithium-ion	Lithium-ion	Lithium-ion
T-coil	ReCharge B-Li P T	ReCharge B-Li P T	-	-
Tinnitus Sound therapy	•	-	-	-
Wireless sync	•	•	-	-
IP68	•	•	•	•
ATEX	•	•	•	•
EarWear	Earmold only (no thin tube available)			
Technology Highlights				
First Fit 48 channels	•	•	•	•
Voice Ranger	•	•	-	-
Adaptive compression	•	•	•	•
Extended Dynamic Range	•	•	•	•
Directionality				
Multichannel Adaptive Directional Microphone	•	-	-	-
Automatic Directional Microphone	•	•	•	-
Fixed Directional Microphone	•	•	•	•
iOmni Sound Locator	•	•	-	-
Noise Reduction				
Noise Management	•	•	•	•
Sound Smoothing	•	-	-	-
Wind Noise Reduction	•	-	-	-
Sound Quality & Control				
Bandwidth compression	•	•	-	-
Feedback preventer	•	•	•	•
Smart Automatic Equalizer	•	•	-	-
Smart Automatic Acclimatization	•	•	-	-
Automatic Classifier	-	-	-	-
Data Logging	•	•	•	•
App				
Rexton App	•	•	•	•
TeleCare				
Full live remote tuning	-	-	-	-
Basic portal remote tuning	•	•	•	•

Accessories	
Smart Key	•
Smart Transmitter 2,4	-
Smart Mic	-

Fitting	
Noahlink Wireless	•
ConnexxAir	-

Chargers	
Charging Station	-
Travel Charger	-
Charging+ Station	-
IF CIC Travel Charger	-
Universal Cable Charger	•

• available - not available

RECHARGE B-LI HP

JUST PLUG IN AND POWER UP

The powerful long-lasting Rexton ReCharge B-Li HP will secure new energy for up to 88 hours. Severe to profound hearing losses can benefit from proven technology for clear and comfortable sounds.




- > Takes you through the day and night with 88h runtime
- > Unique rechargeability offering in the 30, 20, 10 and 5 performance levels
- > Addresses severe to profound hearing loss
- > Programmable via Noahlink Wireless
- > Water-resistant, IP68 rating
- > Lifeproof
- > ATEX

Colors



Fitting range is based on 2 ccm coupler:

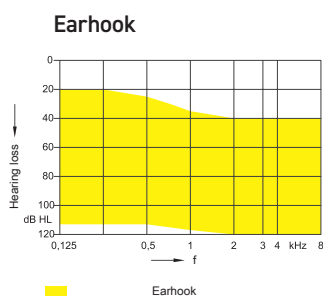

- 80 dB / 140 dB SPL
(2 ccm coupler)
- 83 dB / 144 dB SPL
(Ear simulator)

Features / Performance Level	30	20	10	5
Channels/Controls/Programs	16/8/4	16/8/4	12/6/4	8/4/4
Battery	Lithium-ion	Lithium-ion	Lithium-ion	Lithium-ion
T-coil	ReCharge B-Li HP T	ReCharge B-Li HP T	-	-
Tinnitus Sound therapy	•	-	-	-
Wireless sync	•	•	-	-
IP68	•	•	•	•
ATEX	•	•	•	•
EarWear	Earmold only (no thin tube available)			
Technology Highlights				
First Fit 48 channels	•	•	•	•
Voice Ranger	•	•	-	-
Adaptive compression	•	•	•	•
Extended Dynamic Range	•	•	•	•
Directionality				
Multichannel Adaptive Directional Microphone	•	-	-	-
Automatic Directional Microphone	•	•	•	-
Fixed Directional Microphone	•	•	•	•
iOmni Sound Locator	•	•	-	-
Noise Reduction				
Noise Management	•	•	•	•
Sound Smoothing	•	-	-	-
Wind Noise Reduction	•	-	-	-
Sound Quality & Control				
Bandwidth compression	•	•	•	•
Feedback preventer	•	•	•	•
Smart Automatic Equalizer	•	•	-	-
Smart Automatic Acclimatization	•	•	-	-
Automatic Classifier	-	-	-	-
Data Logging	•	•	•	•
App				
Rexton App	•	•	•	•
TeleCare				
Full live remote tuning	-	-	-	-
Basic portal remote tuning	•	•	•	•

Accessories	
Smart Key	•
Smart Transmitter 2,4	-
Smart Mic	-

Fitting	
Noahlink Wireless	•
ConnexxAir	-

Chargers	
Charging Station	-
Travel Charger	-
Charging+ Station	-
IF CIC Travel Charger	-
Universal Cable Charger	•

• available - not available

BICORE TECHNOLOGY



Catching the important word of a colleague, registering the sound of the machinery, the people or the environment in any situation, is at the center of the Rexton BiCore technology.

Rexton is reliable and solid listening, even in demanding situations thanks to the breakthrough BiCore Technology – getting the job done is what is most important.

There are three special reasons why you can rely on every hearing aid in the range:

CLEAN SOUND

BiCore uses innovative sound processing to faithfully reconstruct the environmental soundscape for the hearing aid user.

The robust input stage and new A/D converter reliably handle sounds of all levels without distortion, ensuring that what comes in also goes out the same way, clean and smooth.

RELIABLE TECHNOLOGY

Patented Speech Preservation Technology (SPT) minimizes unwanted noise and preserves speech sound.

SPT takes advantage of the spatial separation between target speech sounds and competing noise. This directional analysis allows for separate processing of the speech and noise before they are mixed together.

An ultra-quick wireless binaural link ensures left and right hearing aids are always synchronized and react instantly and precisely to changes in the sound environment.

LIFEPROOF

All BiCore hearing aids are Lifeproof with robust construction and proven design.

Rexton BiCore hearing aids are tested extensively to ensure they can withstand whatever life throws at them. They are durable enough to withstand everyday challenges and harsh environments.

All Rexton products are IP68 certified.



MADE FOR IPHONE AND ANDROID (ASHA)

BiCore hearing aids can stream calls and audio from iPhones and Android phones.*



STREAMING

You may stream directly from your TV and other audio devices using the Smart Transmitter 2.4 or the Smart Mic.



TELECARE

Full live remote fitting sessions via Connexx resulting in higher user convenience and satisfaction.



HEADSET MODE TAP CONTROL

Rexton BiCore microphones get the job done so iOS users can enjoy hands free call.



MY VOICE 2.0

General sound and own voice processed in parallel for a greater sound impression and user satisfaction.



AUTO ECHO REDUCER

Automated environment scan to reduce effects of disturbing reverberation on the perceived sound signal.



ATEX

Rexton BiCore hearing aids has ATEX certification which verifies that they are tested and approved to be intrinsically safe for use in potentially explosive environments. They can be safely used in areas where there may be flammable substances in the air such as paint and chemical manufacturing, petrol stations and mills. The ATEX group/class is printed on the ATEX card, which must be available at hand when the hearing aids are used in explosive atmospheres.

* Phones with Bluetooth 5.0+ and Android 10 (or higher) are capable of providing wireless streaming via ASHA, but ASHA must be enabled by each phone manufacturer and service carrier.

BICORE B-LI M RUGGED

BICORE B-LI M RUGGED IS OUR TOUGHEST HEARING AID YET AND CAN HANDLE ALL THE KNOCKS AND SPILLS THAT CAN OCCUR IN EVERYDAY LIFE.

Based on BiCore technology for a smooth listening experience.



It has proven BiCore hearing technology that delivers a range of special features:

- > Ultra-fast adaption to changing environments
- > Patented speech preservation technology
- > Advanced binaural signal processing to reduce noise efficiently
- > Stream audio directly from your iOS or Android smartphone
- > Stream directly from your TV and other audio devices using the Smart Transmitter 2.4 or the Smart Mic

Unique tough building offering:

- > Water resistant construction
- > Sweat proof
- > Drop proof
- > Scratch resistant
- > IP68 certified
- > ATEX certified

Colors

Dark Champagne | Granite

Granite | Black

Grey | Granite

Rose Gold | Grey

Sandy Brown | Granite


Silver | Grey

Beige | Beige

Black | Black



Fitting range is based on 2 ccm coupler:

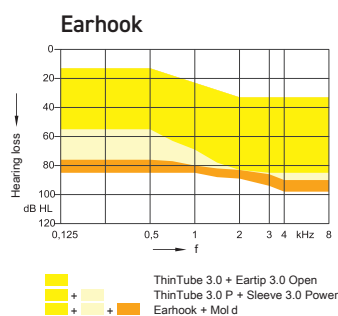

- 60 dB / 133 dB SPL (2 ccm coupler)
- 67 dB / 138 dB SPL (Ear simulator)

TOUGH DESIGN IN A SUPER-COMPACT SIZE

To optimize toughness, a long list of enhancements have been made to the hearing aid structure. These include making the receiver airtight, sealing charging contacts, reinforcing receiver suspension, thickening the housing, and sealing all its parting lines. Rugged is also far smaller than conventional behind-the-ear hearing aids. Its dimensions are more comparable to those of a RIC.



Features / Performance Level	80	60	40
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6
Soundpro 2.0	•	•	•
My Voice 2.0	•	•	•
Direct Streaming	Made for iPhone / Android (ASHA)	Made for iPhone / Android (ASHA)	Made for iPhone / Android (ASHA)
Auto Volume	•	•	•
Wireless sync	•	•	•
Directionality	Automatic Adaptive, iOmni, Front/Back Left/Right automatic & manual, Narrow	Automatic Adaptive, iOmni, Front/Back automatic & manual, Left/Right manual, Narrow	Automatic Adaptive, iOmni, Narrow
Feedback Preventer	•	•	•
Noise Reduction	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing
Wind Noise Reduction	•	•	•
Reverb Reducer	Manual	Manual	-
Auto Echo Reducer	•	-	-
Bandwidth: Extension/Compression	• / •	- / •	- / •
Music Enhancer (presets)	3	3	1
Tinnitus function (Sound therapy and Notch)	•	•	•
Xphone	•	•	•
Acclimatization/Data Logging	• / •	• / •	• / •
Small earhook	•	•	•
IP68	•	•	•
ATEX	•	•	•
Additional tests	•	•	•

Accessories & Apps

Charging Station B-M	•	•	•
Smart Key	•	•	•
Smart Transmitter 2.4	•	•	•
Smart Mic	•	•	•
Rexton App	•	•	•
BiCore CROS R312	•	•	•
BiCore CROS R-Li	•	•	•

Battery

Battery runtime (without streaming)	up to 39 h
Battery runtime (incl. 5 h streaming)	up to 36 h

• available - not available


BICORE B M

RELY ON THE BICORE BTEs FOR SOUND THAT'S CLEAN AND TRUE EVEN IN NOISY ENVIRONMENTS

The BiCore BTEs combines new enhanced premium microphone and receiver technologies with the power of proven Speech Preservation Technology®.



- › BiCore processing with Speech Preservation Technology™
- › Direct audio streaming via MFI or ASHA protocols (except PL10)
- › Battery operated (Zn-air size 13)
- › Programmable via Noahlink Wireless
- › Tamper-proof battery door (optional feature), tamper-proof rocker switch, tamper-proof earhook.

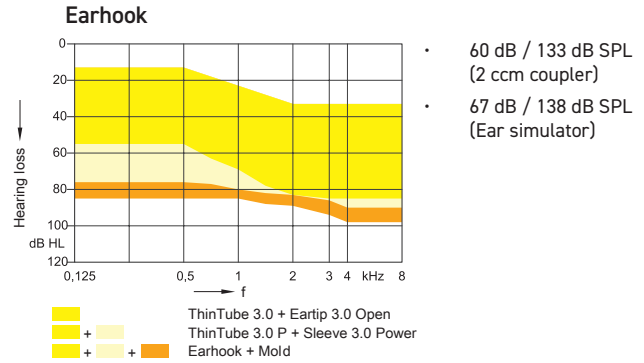

Colors



Features / Performance Level	30	20	10
Channels/Controls/Programs	16/8/4	16/8/4	12/6/4
Soundpro 2.0	•	•	•
My Voice 2.0	-	-	-
Direct Streaming	Made for iPhone / Android (ASHA)	Made for iPhone / Android (ASHA)	-
Auto Volume	•	•	-
Wireless sync	•	•	•
Directionality	Automatic adaptive, iOmni	Automatic adaptive, iOmni	Automatic fixed omni
Feedback Preventer	•	•	•
Noise Reduction	Noise Management Sound Smoothing	Noise Management	Noise Management
Wind Noise Reduction	•	-	-
Reverb Reducer	-	-	-
Auto Echo Reducer	-	-	-
Bandwidth: Extension/Compression	- / •	- / •	- / -
Music Enhancer (presets)	1	-	-
Tinnitus function (Sound therapy and Notch)	•	-	-
Xphone	•	•	-
Acclimatization/Data Logging	• / •	• / •	- / •
IP68	•	•	•
ATEX	•	•	•
T-coil	•	•	•

Accessories & Apps

Smart Key	•	•	•
Smart Transmitter 2.4	•	•	-
Smart Mic	•	•	-
Rexton App	•	•	•
BiCore CROS R-Li	•	•	-
BiCore CROS R312	•	•	-
BiCore CROS SR	-	-	-

Battery

Battery runtime (without streaming)	up to 142 h
Battery runtime (incl. 20 h streaming)	up to 129 h


BICORE B P

RELY ON THE BICORE BTEs FOR SOUND THAT'S CLEAN AND TRUE EVEN IN NOISY ENVIRONMENTS

The BiCore BTEs combines new enhanced premium microphone and receiver technologies with the power of proven Speech Preservation Technology®.



- › BiCore processing with Speech Preservation Technology™
- › Direct audio streaming via MFI or ASHA protocols (except PL10)
- › Battery operated (Zn-air size 13)
- › Programmable via Noahlink Wireless
- › Tamper-proof battery door (optional feature), tamper-proof rocker switch, tamper-proof earhook.

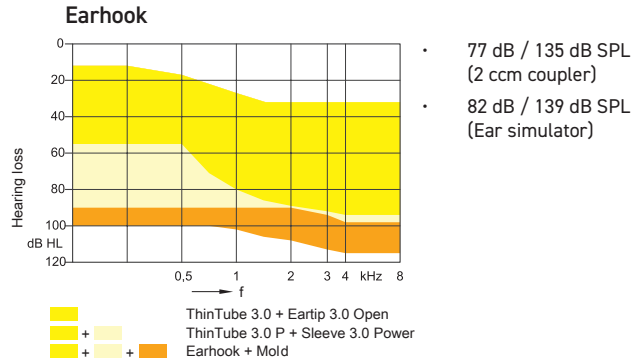

Colors



Features / Performance Level	30	20	10
Channels/Controls/Programs	16/8/4	16/8/4	12/6/4
Soundpro 2.0	•	•	•
My Voice 2.0	-	-	-
Direct Streaming	Made for iPhone / Android (ASHA)	Made for iPhone / Android (ASHA)	-
Auto Volume	•	•	-
Wireless sync	•	•	•
Directionality	Automatic adaptive, iOmni	Automatic adaptive, iOmni	Automatic fixed omni
Noise Reduction	Noise Management Sound Smoothing	Noise Management	Noise Management
Feedback Preventer	•	•	•
Wind Noise Reduction	•	-	-
Reverb Reducer	-	-	-
Auto Echo Reducer	-	-	-
Bandwidth: Extension/Compression	- / •	- / •	- / -
Music Enhancer (presets)	1	-	-
Tinnitus function (Sound therapy and Notch)	•	-	-
Xphone	•	•	-
Acclimatization/Data Logging	• / •	• / •	- / •
IP68	•	•	•
ATEX	•	•	•
T-coil	•	•	•

Accessories & Apps

Smart Key	•	•	•
Smart Transmitter 2.4	•	•	-
Smart Mic	•	•	-
Rexton App	•	•	•
BiCore CROS R-Li	•	•	-
BiCore CROS R312	•	•	-
BiCore CROS SR	-	-	-

Battery

Battery runtime (without streaming)	up to 148 h
Battery runtime (incl. 20 h streaming)	up to 115 h

BICORE B HP


RELY ON THE BICORE BTEs FOR SOUND THAT'S CLEAN AND TRUE EVEN IN NOISY ENVIRONMENTS

The BiCore B HP combines new enhanced premium microphone and receiver technologies with the power of proven Speech Preservation Technology®. And with 277 hours of performance per battery, your customers can rely on it in every situation.



- › BiCore processing with Speech Preservation Technology™
- › Direct audio streaming via MFI or ASHA protocols (except PL10)
- › Battery operated (Zn-air size 675)
- › Programmable via Noahlink Wireless
- › Tamper-proof battery door (optional feature), tamper-proof rocker switch, tamper-proof earhook.

Earhook

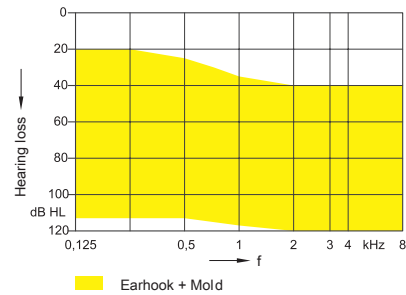

- 82 dB / 140 dB SPL (2 ccm coupler)
- 85 dB / 143 dB SPL (Ear simulator)

Colors



Features / Performance Level	30	20	10
Channels/Controls/Programs	16/8/4	16/8/4	12/6/4
Soundpro 2.0	•	•	•
My Voice 2.0	-	-	-
Direct Streaming	Made for iPhone / Android (ASHA)	Made for iPhone / Android (ASHA)	-
Auto Volume	•	•	-
Wireless sync	•	•	•
Directionality	Automatic adaptive, iOmni	Automatic adaptive, iOmni	Automatic fixed omni
Noise Reduction	Noise Management Sound Smoothing	Noise Management	Noise Management
Feedback Preventer	•	•	•
Wind Noise Reduction	•	-	-
Reverb Reducer	-	-	-
Auto Echo Reducer	-	-	-
Bandwidth: Extension/Compression	- / •	- / •	- / -
Music Enhancer (presets)	1	-	-
Tinnitus function (Sound therapy and Notch)	•	-	-
Xphone	•	•	-
Acclimatization/Data Logging	• / •	• / •	- / •
IP68	•	•	•
ATEX	•	•	•
T-coil	•	•	•

Accessories & Apps

Smart Key	•	•	•
Smart Transmitter 2.4	•	•	-
Smart Mic	•	•	-
Rexton App	•	•	•
BiCore CROS R-Li	•	•	-
BiCore CROS R312	•	•	-
BiCore CROS SR	-	-	-

Battery

Battery runtime (without streaming)	up to 277 h
Battery runtime (incl. 60 h streaming)	up to 217 h

• available - not available



BICORE R-LI AND R-LI T

DURABILITY WITH A SMOOTH LISTENING EXPERIENCE

The BiCore R-Li and R-Li T are rechargeable, they have Lifeproof microphones, 60% lower binaural latency rate (vs. the previous generation) and one of our lowest power consumption.




- › Integrated motion sensor
- › Lithium-ion battery
- › Programmable rocker switch
- › Color Conversion Kits
- › Smart Demo device available
- › Programming via:
Noahlink Wireless

Colors



Fitting range is based on 2 ccm coupler:

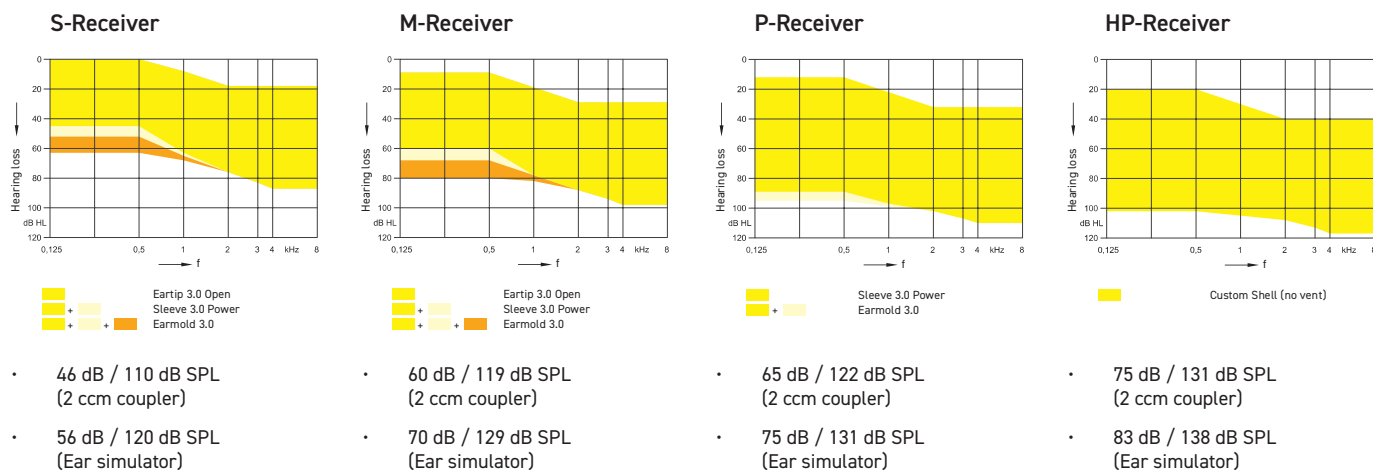

Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6	16/8/4	16/8/4
Soundpro 2.0	•	•	•	•	•
My Voice 2.0	•	•	•	-	-
Direct Streaming	Made for iPhone / Android (ASHA)	Made for iPhone / Android (ASHA)	Made for iPhone / Android (ASHA)	Made for iPhone / Android (ASHA)	Made for iPhone / Android (ASHA)
Auto Volume	•	•	•	•	•
Wireless sync	•	•	•	•	•
Directionality	Automatic Adaptive, iOmni, Front/Back Left/Right automatic & manual, Narrow	Automatic Adaptive, iOmni, Front/Back automatic & manual, Left/Right manual, Narrow	Automatic Adaptive, iOmni, Narrow	Automatic Adaptive, iOmni	Automatic Adaptive, iOmni
Noise Reduction	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing	Noise Management Sound Smoothing	Noise Management
Wind Noise Reduction	•	•	•	•	-
Reverb Reducer	Manual	Manual	-	-	-
Auto Echo Reducer	•	-	-	-	-
Bandwidth: Extension/Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer (presets)	3	3	1	1	-
Feedback Preventer	•	•	•	•	•
Tinnitus function (Sound therapy and Notch)	•	•	•	•	-
Xphone	•	•	•	•	•
Acclimatization/Data Logging	• / •	• / •	• / •	• / •	• / •
T-coil	BiCORE R-Li T	BiCORE R-Li T	BiCORE R-Li T	BiCORE R-Li T	BiCORE R-Li T
IP68	•	•	•	•	•
ATEX	•	•	•	•	•

Accessories & Apps					
Smart Key	•	•	•	•	•
Smart Transmitter 2.4	•	•	•	•	•
Smart Mic	•	•	•	•	•
Rexton App	•	•	•	•	•
Travel Charger RIC/ Charging Station R/ Charging+ Station R	•	•	•	•	•
BiCore CROS R-Li	•	•	•	•	•
BiCore CROS R312	•	•	•	•	•
BiCore CROS SR	-	-	-	-	-
Battery*					
Battery runtime (without streaming)	up to 28 h (R-Li) up to 39 h (R-Li T)				
Battery runtime (incl. 5 h streaming)	up to 24 h (R-Li) up to 36 h (R-Li T)				

• available - not available * Battery runtime value based on S-Receiver.

BICORE R312

DURABILITY WITH A SMOOTH LISTENING EXPERIENCE

Our battery RIC delivers sound that's true and clean. All have Lifeproof microphones, a 60% lower binaural latency rate (vs. the previous generation) and user-friendly rocker switch bumps for better usability and tactile feedback.




- › Integrated motion sensor
- › Battery size 312
- › Programmable rocker switch
- › Color Conversion Kits
- › Smart Demo device available
- › Programming via:
Noahlink Wireless

Colors



Fitting range is based on 2 ccm coupler:

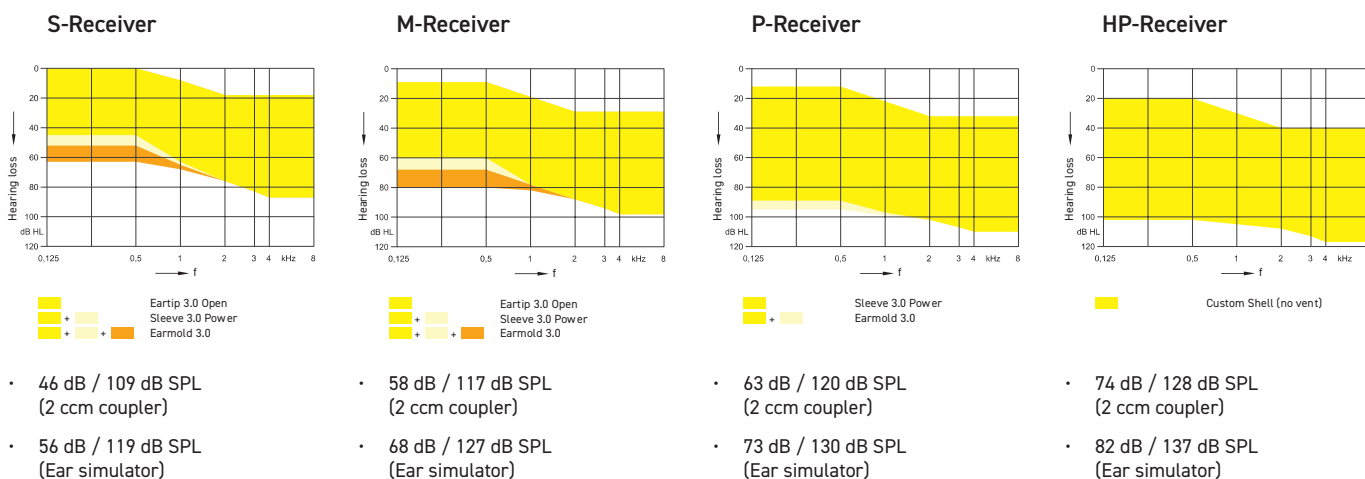

Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6	16/8/4	16/8/4
Soundpro 2.0	•	•	•	•	•
My Voice 2.0	•	•	•	-	-
Direct Streaming	Made for iPhone / Android (ASHA)	Made for iPhone / Android (ASHA)	Made for iPhone / Android (ASHA)	Made for iPhone / Android (ASHA)	Made for iPhone / Android (ASHA)
Auto Volume	•	•	•	•	•
Wireless sync	•	•	•	•	•
Directionality	Automatic Adaptive, iOmni, Front/Back Left/Right automatic & manual, Narrow	Automatic Adaptive, iOmni, Front/Back automatic & manual, Left/Right manual, Narrow	Automatic Adaptive, iOmni, Narrow	Automatic Adaptive, iOmni	Automatic Adaptive, iOmni
Feedback Preventer	•	•	•	•	•
Noise Reduction	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing	Noise Management Sound Smoothing	Noise Management
Wind Noise Reduction	•	•	•	•	-
Reverb Reducer	Manual	Manual	-	-	-
Auto Echo Reducer	•	-	-	-	-
Bandwidth: Extension/Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer (presets)	3	3	1	1	-
Tinnitus function (Sound therapy and Notch)	•	•	•	•	-
Xphone	•	•	•	•	•
Acclimatization/Data Logging	• / •	• / •	• / •	• / •	• / •
T-coil	-	-	-	-	-
IP68	•	•	•	•	•
ATEX	•	•	•	•	•
Accessories & Apps					
Smart Key	•	•	•	•	•
Smart Transmitter 2.4	•	•	•	•	•
Smart Mic	•	•	•	•	•
Rexton App	•	•	•	•	•
BiCore CROS R-Li	•	•	•	•	•
BiCore CROS R312	•	•	•	•	•
BiCore CROS SR	-	-	-	-	-
Battery*					
Battery runtime (without streaming)	up to 89 h				
Battery runtime (incl. 20 h streaming)	up to 72 h				

• available - not available * Battery runtime value based on S-Receiver.

BICORE SR

GET THE PERFECT BLEND OF WEARING COMFORT AND STYLE WITH OUR SMALL AND LIGHT SLIM-RIC

Our BiCore SR is smaller and lighter than our previous models, has Lifeproof microphones, and a 10% longer run time. Together with the Travel Charger it can provide up to 4 days of autonomy, longer than any previous generation.



- › Integrated motion sensor
- › Lithium-ion battery
- › Smart Demo device available
- › Programming via:
Noahlink Wireless

Colors

Black | Silver



Black | Black Gloss



Cosmic Blue | Rose Gold




Snow White | Rose Gold




Snow White | Snow White Gloss




Fitting range is based on 2 ccm coupler:

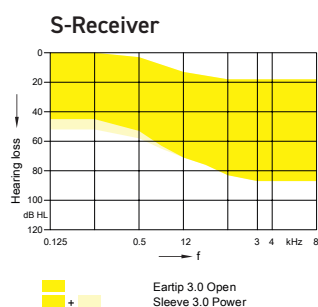

- 46 dB / 110 dB SPL (2 ccm coupler)
- 56 dB / 120 dB SPL (Ear simulator)

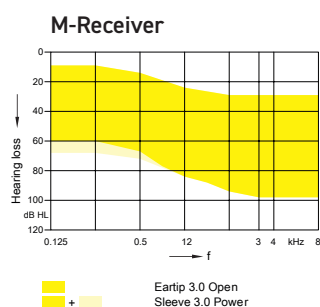

- 60 dB / 119 dB SPL (2 ccm coupler)
- 70 dB / 129 dB SPL (Ear simulator)

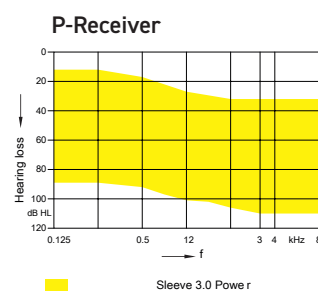

- 65 dB / 122 dB SPL (2 ccm coupler)
- 75 dB / 131 dB SPL (Ear simulator)

Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6	16/8/4	16/8/4
Soundpro 2.0	•	•	•	•	•
My Voice 2.0	•	•	•	-	-
Direct Streaming	Made for iPhone / Android (ASHA)	Made for iPhone / Android (ASHA)	Made for iPhone / Android (ASHA)	Made for iPhone / Android (ASHA)	Made for iPhone / Android (ASHA)
Auto Volume	•	•	•	•	•
Wireless sync	•	•	•	•	•
Directionality	Automatic Adaptive, iOmni, Front/Back Left/Right automatic & manual, Narrow	Automatic Adaptive, iOmni, Front/Back automatic & manual, Left/Right manual, Narrow	Automatic Adaptive, iOmni, Narrow	Automatic Adaptive, iOmni	Automatic Adaptive, iOmni
Noise Reduction	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing	Noise Management Sound Smoothing	Noise Management
Feedback Preventer	•	•	•	•	•
Wind Noise Reduction	•	•	•	•	-
Reverb Reducer	Manual	Manual	-	-	-
Auto Echo Reducer	•	-	-	-	-
Bandwidth: Extension/Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer (presets)	3	3	1	1	-
Tinnitus function (Sound therapy and Notch)	•	•	•	•	-
Xphone	•	•	•	•	•
Acclimatization/Data Logging	• / •	• / •	• / •	• / •	• / •
T-coil	-	-	-	-	-
IP68	•	•	•	•	•
ATEX	•	•	•	•	•

Accessories & Apps					
Smart Key	•	•	•	•	•
Smart Transmitter 2.4	•	•	•	•	•
Smart Mic	•	•	•	•	•
Rexton App	•	•	•	•	•
Slim-RIC Travel Charger 2 Charging+ Station SR	•	•	•	•	•
BiCore CROS R-Li	-	-	-	-	-
BiCore CROS R312	-	-	-	-	-
BiCore CROS SR	•	•	•	•	•

Battery*	
Battery runtime (without streaming)	up to 20 h
Battery runtime (incl. 5 h streaming)	up to 17 h

• available - not available * Battery runtime value based on S-Receiver.

BICORE CUSTOM LI

OUR FIRST CUSTOM LI WITH BICORE TECHNOLOGY GIVES OPTIMAL DISCRETION AND HEARING PERFORMANCE

You can rely on our very first rechargeable custom model to deliver sound that's true and clean with Bluetooth streaming, contactless charging, iOS and Android compatibility, and a sleek design.




- › Integrated motion sensor
- › Lithium-ion battery
- › Programmable push button
- › Programming via:
Noahlink Wireless

Colors




Fitting range is based on 2 ccm coupler:

S-Receiver*

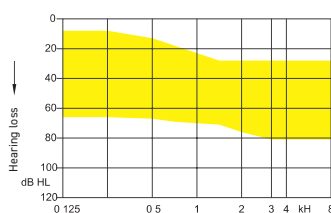


- 50 dB / 113 dB SPL (2 ccm coupler)
- 61 dB / 125 dB SPL (Ear simulator)

M-Receiver**

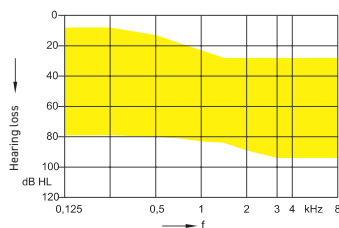

- 55 dB / 118 dB SPL (2 ccm coupler)
- 66 dB / 128 dB SPL (Ear simulator)

P-Receiver**

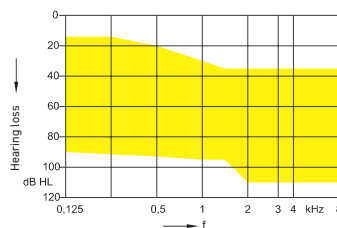

- 65 dB / 124 dB SPL (2 ccm coupler)
- 75 dB / 135 dB SPL (Ear simulator)

* Available for ITC

** Available for both ITC and ITE

Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6	16/8/4	16/8/4
Soundpro 2.0	•	•	•	•	•
My Voice 2.0	-	-	-	-	-
Direct Streaming	Made for iPhone / Android (ASHA)	Made for iPhone / Android (ASHA)	Made for iPhone / Android (ASHA)	Made for iPhone / Android (ASHA)	Made for iPhone / Android (ASHA)
Auto Volume	•	•	•	•	•
Wireless sync	•	•	•	•	•
Directionality	Automatic Adaptive, iOmni, Front/Back Left/Right automatic & manual, Narrow	Automatic Adaptive, iOmni, Front/Back automatic & manual, Left/Right manual, Narrow	Automatic Adaptive, iOmni, Narrow	Automatic Adaptive, iOmni	Automatic Adaptive, iOmni
Noise Reduction	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing	Noise Management Sound Smoothing	Noise Management
Wind Noise Reduction	•	•	•	•	-
Reverb Reducer	Manual	Manual	-	•	-
Auto Echo Reducer	•	-	-	-	-
Bandwidth: Extension/Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer (presets)	3	3	1	1	-
Feedback Preventer	•	•	•	•	•
Tinnitus function (Sound therapy and Notch)	•	•	•	•	-
Xphone	•	•	•	•	•
Acclimatization/Data Logging	• / •	• / •	• / •	• / •	• / •
T-coil	-	-	-	-	-
IP68	•	•	•	•	•

Accessories & Apps					
Smart Key	•	•	•	•	•
Smart Transmitter 2.4	•	•	•	•	•
Smart Mic	•	•	•	•	•
Rexton App	•	•	•	•	•
Custom Charger Station	•	•	•	•	•
BiCore CROS R-Li	•	•	•	•	•
BiCore CROS R312	•	•	•	•	•
BiCore CROS SR	-	-	-	-	-

Battery	
Battery runtime (without streaming)	up to 28 h
Battery runtime (incl. 5 h streaming)	up to 24 h

• available - not available

MOTIONCORE TECHNOLOGY

RX
M-CORE

Most hearing aids analyze hearing situations solely on acoustic criteria, then classify them into a small range of predefined acoustic environments. But when people rely on their senses to keep them safe and sound, hearing clearly in every situation is a must – especially when they're on the job and moving through the day, where environments can change quickly. Our M-Core hearing aids use proven MotionCore technology that automatically recognizes and adapts to different hearing situations. This ensures users can hear speech clearly and get the most reliable impression of their environment.



UP TO 20 TIMES MORE HEARING SITUATIONS THAN BEFORE AUTOMATICALLY ADAPTED



HIGH RESOLUTION SOUNDPRO

Using a combination of inputs including real-time motion data, signal to noise ratio, and sound directionality, High Resolution Soundpro can continuously adjust the sound settings to provide the most reliable impression of the user's environment, automatically.



MOTION SENSOR

MotionCore technology features an integrated real-time motion sensor that detects a user's movement, and uses the data to even more accurately determine their sound situation. With this, M-Core hearing aids can automatically detect up to 127 different sound situations.



MY VOICE

The hearing aid recognizes the user's specific voice pattern, and integrates it naturally into what they hear. In turn, this provides the user with a natural-sounding voice that leads to high acceptance rates. And for that the user will communicate more and be more socially active, greatly improving their quality of life



TELECARE

Full live remote fitting sessions via Connexx resulting in higher user convenience and satisfaction.



STREAMING

You may direct stream from your TV and other audio devices using the Smart Transmitter 2.4 or the Smart Mic.



REXTON ASSIST

Rexton Assist works on a live, learning network that continuously evolves to provide solutions to any challenge faced by the hearing aid user. Rexton Assist brings convenience for both the user and the dispenser, improving the hearing aid performance based on the needs of the user, at any time, in any place.


M-CORE B-LI M

VERSATILE AND RECHARGEABLE

Our popular BTE features MotionCore technology and a rechargeable lithium-ion battery, for an accurate impression of your environment in every situation.



- > Integrated motion sensor
- > Lithium-ion battery
- > Programmable rocker switch
- > ThinTube 3.0
- > Color Conversion Kits
- > Smart Demo device available
- > Programming via: Noahlink Wireless
- > Rexton Assist

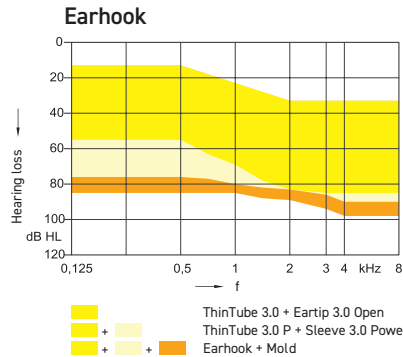

- 60 dB / 133 dB SPL (2 ccm coupler)
- 67 dB / 138 dB SPL (Ear simulator)

Colors



Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6	16/8/4	16/8/4
Soundpro	•	•	•	•	•
My Voice	•	•	•	-	-
Direct Streaming	Made for iPhone	Made for iPhone	Made for iPhone	Made for iPhone	Made for iPhone
Auto Volume	•	•	•	•	•
Wireless sync	•	•	•	•	•
Directionality	Automatic/Adaptive iOmni Front & Back Left & Right Narrow	Automatic/Adaptive iOmni Front & Back	Automatic/Adaptive iOmni	Automatic/Adaptive	Automatic/Adaptive
Noise Reduction	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing	Noise Management Sound Smoothing	Noise Management
Wind Noise Reduction	•	•	•	•	-
Feedback Preventer	•	•	•	•	•
Reverb Reducer	•	-	-	-	-
Bandwidth: Extension/Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer (presets)	3	3	1	1	-
Tinnitus function (Sound therapy and Notch)	•	•	•	•	-
Xphone	•	•	•	-	-
Acclimatization/Data Logging	• / •	• / •	• / •	• / •	• / •
T-coil	-	-	-	-	-
Small ear hook	•	•	•	•	•

Accessories & Apps					
Charging Station	•	•	•	•	•
Smart Key	•	•	•	•	•
Smart Transmitter 2.4	•	•	•	•	•
Smart Mic	•	•	•	•	•
Rexton App	•	•	•	•	•
M-Core CROS R	•	•	•	-	-
M-Core CROS R-Li	•	•	•	-	-
M-Core CROS iX-CIC	-	-	-	-	-

Battery

Battery runtime (without streaming)	up to 24 h
Battery runtime (incl. 5 h streaming)	up to 21 h


M-CORE B-LI P

POWERFUL AND RECHARGEABLE

M-Core B-Li P remains a powerful device that also offers the MotionCore technology.



- > Integrated motion sensor
- > Lithium-ion battery
- > Programmable rocker switch
- > ThinTube 3.0
- > Color Conversion Kits
- > Smart Demo device available
- > Programming via: Noahlink Wireless
- > Rexton Assist

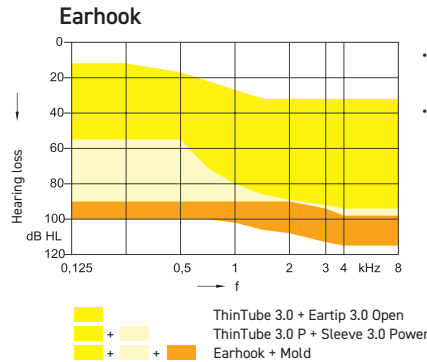

- 77 dB / 135 dB SPL (2 ccm coupler)
- 82 dB / 140 dB SPL (Ear simulator)

Colors



Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6	16/8/4	16/8/4
Soundpro	•	•	•	•	•
My Voice	•	•	•	-	-
Direct Streamin	Made for iPhone	Made for iPhone	Made for iPhone	Made for iPhone	Made for iPhone
Auto Volume	•	•	•	•	•
Wireless sync	•	•	•	•	•
Directionality	Automatic / Adaptive iOmni Front & Back Left & Right Narrow	Automatic / Adaptive iOmni Front & Back	Automatic / Adaptive iOmni	Automatic / Adaptive	Automatic / Adaptive
Noise Reduction	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing	Noise Management Sound Smoothing	Noise Management
Wind Noise Reduction	•	•	•	•	-
Feedback Preventer	•	•	•	•	•
Reverb Reducer	•	-	-	-	-
Bandwidth: Extension/Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer (presets)	3	3	1	1	-
Tinnitus function (Sound therapy and Notch)	•	•	•	•	-
Xphone	•	•	•	-	-
Acclimatization/Data Logging	• / •	• / •	• / •	• / •	• / •
T-coil	•	•	•	•	•
Small ear hook	•	•	•	•	•

Accessories & Apps					
Charging+ Station B-P / Charging Station B-P	•	•	•	•	•
Smart Key	•	•	•	•	•
Smart Transmitter 2.4	•	•	•	•	•
Smart Mic	•	•	•	•	•
Rexton App	•	•	•	•	•
M-Core CROS R	•	•	•	-	-
M-Core CROS R-Li	•	•	•	-	-
M-Core CROS iX-CIC	-	-	-	-	-

Battery	
Battery runtime (without streaming)	up to 30 h
Battery runtime (incl. 5 h streaming)	up to 27 h


M-CORE B-LI HP

EXTRA POWERFUL AND RECHARGEABLE

M-Core B-Li HP is a powerful device, suitable for those with severe to profound hearing loss, and combines the MotionCore technology with rechargeability.



- > Integrated motion sensor
- > Lithium-ion battery
- > Programmable rocker switch
- > Color Conversion Kits
- > Smart Demo device available
- > Programming via:
Noahlink Wireless
- > Rexton Assist

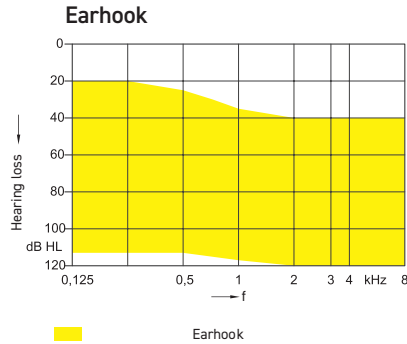

- 82 dB / 140 dB SPL
(2 ccm coupler)
- 86 dB / 144 dB SPL
(Ear simulator)

Colors



Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	34/20/6	32/16/6	24/12/6	16/8/4	16/8/4
Soundpro	•	•	•	•	•
My Voice	•	•	•	-	-
Direct Streaming	Made for iPhone	Made for iPhone	Made for iPhone	Made for iPhone	Made for iPhone
Auto Volume	•	•	•	•	•
Wireless sync	•	•	•	•	•
Directionality	Automatic/Adaptive iOmni Front & Back Left & Right Narrow	Automatic/Adaptive iOmni Front & Back	Automatic/Adaptive iOmni	Automatic/Adaptive	Automatic/Adaptive
Noise Reduction	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing	Noise Management Sound Smoothing	Noise Management
Wind Noise Reduction	•	•	•	•	-
Feedback Preventer	•	•	•	•	•
Reverb Reducer	•	-	-	-	-
Bandwidth: Extension/Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer (presets)	3	3	1	1	-
Tinnitus function (Sound therapy)	•	•	•	•	-
Xphone	•	•	•	-	-
Acclimatization/Data Logging	• / •	• / •	• / •	• / •	• / •
T-coil	•	•	•	•	•
Small ear hook	•	•	•	•	•

Accessories & Apps					
Charging+ Station B-HP / Charging Station B-HP	•	•	•	•	•
Smart Key	•	•	•	•	•
Smart Transmitter 2.4	•	•	•	•	•
Smart Mic	•	•	•	•	•
Rexton App	•	•	•	•	•
M-Core CROS R	•	•	•	-	-
M-Core CROS R-Li	•	•	•	-	-
M-Core CROS iX-CIC	-	-	-	-	-

Battery		
Battery runtime (without streaming)		up to 61 h
Battery runtime (incl. 5 h streaming)		up to 57 h


M-CORE R-LI

THE SMALL AND RECHARGEABLE RIC

Combining hearing comfort and rechargeability in a 15% smaller housing than the predecessor.



- › Integrated motion sensor
- › Lithium-ion battery
- › Programmable rocker switch
- › Color Conversion Kits
- › Smart Demo device available
- › Programming via:
Noahlink Wireless
- › Rexton Assist

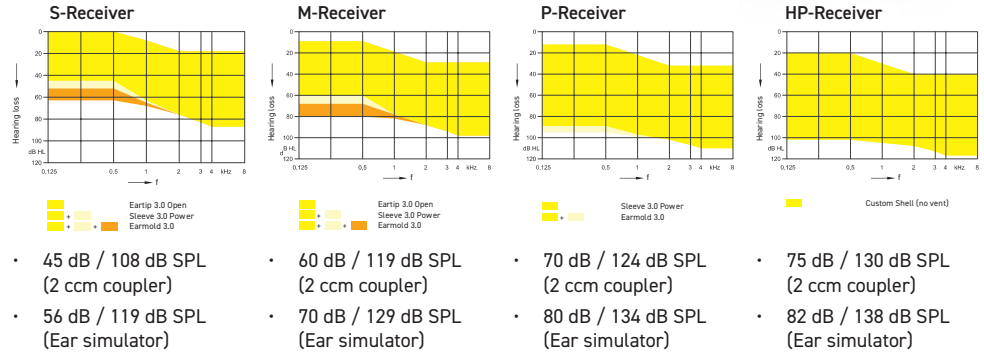

Colors



Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6	16/8/4	16/8/4
Soundpro	•	•	•	•	•
My Voice	•	•	•	-	-
Direct Streaming	Made for iPhone	Made for iPhone	Made for iPhone	Made for iPhone	Made for iPhone
Auto Volume	•	•	•	•	•
Wireless sync	•	•	•	•	•
Directionality	Automatic / Adaptive iOmni Front & Back Left & Right Narrow	Automatic / Adaptive iOmni Front & Back	Automatic / Adaptive iOmni	Automatic / Adaptive	Automatic / Adaptive
Noise Reduction	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing	Noise Management Sound Smoothing	Noise Management
Wind Noise Reduction	•	•	•	•	-
Feedback Preventer	•	•	•	•	•
Reverb Reducer	•	-	-	-	-
Bandwidth: Extension/Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer (presets)	3	3	1	1	-
Tinnitus function (Sound therapy and Notch)	•	•	•	•	-
Xphone	•	•	•	-	-
Acclimatization/Data Logging	• / •	• / •	• / •	• / •	• / •
T-coil	-	-	-	-	-

Accessories & Apps					
Charging Station	•	•	•	•	•
Smart Key	•	•	•	•	•
Smart Transmitter 2.4	•	•	•	•	•
Smart Mic	•	•	•	•	•
Rexton App	•	•	•	•	•
M-Core CROS R	•	•	•	-	-
M-Core CROS R-Li	•	•	•	-	-
M-Core CROS iX-CIC	-	-	-	-	-

Battery*

Battery runtime (without streaming)	up to 23 h
Battery runtime (incl. 5 h streaming)	up to 21 h

• available - not available

* Battery runtime value based on S-Receiver.


M-CORE R312

BATTERY DRIVEN WITH OPTIONAL T-COIL

Automatically offers the right setting in every situation, thanks to the integrated motion sensor.



- › Integrated motion sensor
- › Battery size 312
- › Programmable rocker switch
- › Smart Demo device available
- › Color Conversion Kits with completely integrated telecoil
- › Programming via: Noahlink Wireless
- › Rexton Assist

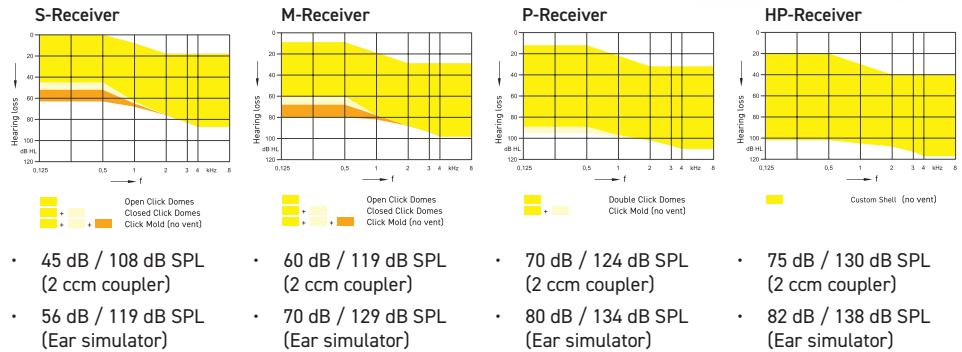

Colors



Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6	16/8/4	16/8/4
Soundpro	•	•	•	•	•
My Voice	•	•	•	-	-
Direct Streaming	Made for iPhone	Made for iPhone	Made for iPhone	Made for iPhone	Made for iPhone
Auto Volume	•	•	•	•	•
Wireless sync	•	•	•	•	•
Directionality	Automatic / Adaptive iOmni Front & Back Left & Right Narrow	Automatic / Adaptive iOmni Front & Back	Automatic / Adaptive iOmni	Automatic / Adaptive	Automatic / Adaptive
Noise Reduction	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing	Noise Management Sound Smoothing	Noise Management
Wind Noise Reduction	•	•	•	•	-
Feedback Preventer	•	•	•	•	•
Reverb Reducer	•	-	-	-	-
Bandwidth: Extension/Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer (presets)	3	3	1	1	-
Tinnitus function (Sound therapy and Notch)	•	•	•	•	-
Xphone	•	•	•	-	-
Acclimatization/Data Logging	• / •	• / •	• / •	• / •	• / •
T-coil	•	•	•	•	•

Accessories & Apps					
Smart Key	•	•	•	•	•
Smart Transmitter 2.4	•	•	•	•	•
Smart Mic	•	•	•	•	•
Rexton App	•	•	•	•	•
M-Core CROS R	•	•	•	-	-
M-Core CROS R-Li	•	•	•	-	-
M-Core CROS iX-CIC	-	-	-	-	-

Battery*

Battery runtime (without streaming)

up to 78 h



M-CORE SR


WITH EXCHANGEABLE RECEIVERS

Our discreet Slim-RIC design features proven MotionCore technology that automatically recognizes and adapts to up to 20 times more hearing situations than before.




- > Integrated motion sensor
- > Lithium-ion battery
- > Smart Demo device available
- > Programming via:
Noahlink Wireless
- > Rexton Assist

S-Receiver

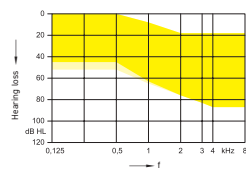


- 45 dB / 108 dB SPL (2 ccm coupler)
- 56 dB / 119 dB SPL (Ear simulator)

M-Receiver

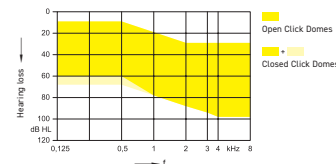

- 60 dB / 119 dB SPL (2 ccm coupler)
- 70 dB / 129 dB SPL (Ear simulator)

P-Receiver

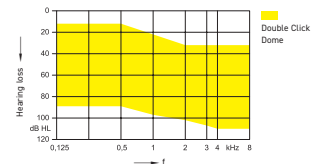

- 70 dB / 124 dB SPL (2 ccm coupler)
- 80 dB / 134 dB SPL (Ear simulator)

Colors

Black | Silver



Black | Black Gloss



Cosmic Blue | Rose Gold



Snow White | Rose Gold



Snow White | Snow White Gloss



* Colors not available for all performance levels

Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6	16/8/4	16/8/4
Soundpro	•	•	•	•	•
My Voice	•	•	•	-	-
Direct Streaming	Made for iPhone	Made for iPhone	Made for iPhone	Made for iPhone	Made for iPhone
Auto Volume	•	•	•	•	•
Wireless sync	•	•	•	•	•
Directionality	Automatic / Adaptive iOmni Front & Back Left & Right Narrow	Automatic / Adaptive iOmni Front & Back	Automatic / Adaptive iOmni	Automatic / Adaptive	Automatic / Adaptive
Noise Reduction	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing	Noise Management Sound Smoothing	Noise Management
Wind Noise Reduction	•	•	•	•	-
Feedback Preventer	•	•	•	•	•
Reverb Reducer	•	-	-	-	-
Bandwidth: Extension/Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer (presets)	3	3	1	1	-
Tinnitus function (Sound therapy and Notch)	•	•	•	•	-
Xphone	•	•	•	-	-
Acclimatization/Data Logging	• / •	• / •	• / •	• / •	• / •
T-coil	-	-	-	-	-

Accessories & Apps					
Slim-RIC Travel Charger	•	•	•	•	•
Smart Key	•	•	•	•	•
Smart Transmitter 2.4	•	•	•	•	•
Smart Mic	•	•	•	•	•
Rexton App	•	•	•	•	•
M-Core CROS R	-	-	-	-	-
M-Core CROS R-Li	-	-	-	-	-
M-Core CROS iX-CIC	-	-	-	-	-

Battery*

Battery runtime (without streaming)	up to 19 h
Battery runtime (incl. 5 h streaming)	up to 15 h


M-CORE IX-CIC

WITH THE HIDDEN POWER OF MOTIONCORE



Almost invisible design featuring the proven MotionCore technology that automatically recognizes and adapts to up to 10 times more hearing situations than before.

- > 10 size battery
- > Smart Demo device available
- > Programming via: ConnexxAir Flex connector

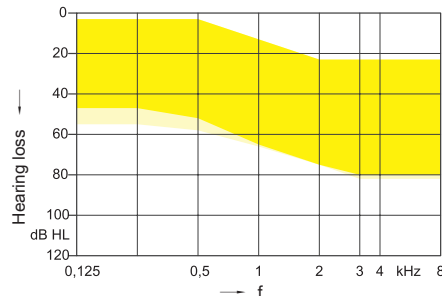

- 50 dB / 114 dB SPL (2 ccm coupler)
- 60 dB / 125 dB SPL (Ear simulator)

Colors



Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6	16/8/4	16/8/4
Soundpro	•	•	•	•	•
My Voice	-	-	-	-	-
Direct Streaming	-	-	-	-	-
Auto Volume	-	-	-	-	-
Wireless sync	•	•	•	•	•
Directionality	Automatic/Adaptive iOmni Front	Automatic/Adaptive iOmni Front	Automatic/Adaptive iOmni	-	-
Noise Reduction	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing	Noise Management Sound Smoothing	Noise Management
Wind Noise Reduction	•	•	•	•	-
Feedback Preventer	•	•	•	•	•
Reverb Reducer	•	-	-	-	-
Bandwidth: Extension/Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer (presets)	3	3	1	1	-
Tinnitus function (Sound therapy and Notch)	•	•	•	•	-
Xphone	•	•	•	-	-
Acclimatization/Data Logging	• / •	• / •	• / •	• / •	• / •
T-coil	-	-	-	-	-

Accessories & Apps					
Smart Key	•	•	•	•	•
Smart Transmitter 2.4	-	-	-	-	-
Smart Mic	-	-	-	-	-
Rexton App	•	•	•	•	•
M-Core CROS R	-	-	-	-	-
M-Core CROS R-Li	-	-	-	-	-
M-Core CROS iX-CIC	•	•	•	-	-

Battery

Battery runtime

up to 70 h



MYCORE

A PERFECT COMBINATION OF NATURAL SOUND AND DIRECT STREAMING



DIRECT STREAMING

Optimally connected to smartphones, TV's and more



MY VOICE

Natural perception of the user's own voice



STEREO ILOCK

Improves the signal-to-noise ratio in difficult listening situations



VOICE RANGER

Easier one-on-one listening in noisy situations



DYNAMIC EXTENDER

Automatically optimizes the dynamic bandwidth of sounds



TELECARE

Full live remote fitting sessions via Connexx resulting in higher user convenience and satisfaction.


MYCORE inoX CIC 8C

DISCREET AND READY TO FIT

inoX combines the discretion of a custom model without having to take an ear impression.



- › Battery size 10
- › Smart Demo device available
- › Programming via:
 - › ConnexxAir
 - › Flex Connector

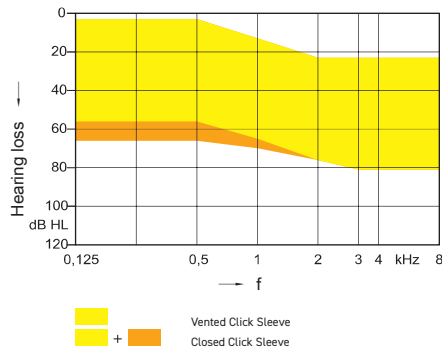

- 50 dB / 113 dB SPL (2 ccm coupler)
- 60 dB / 125 dB SPL (Ear simulator)

Colors

Shell	Red	Blue	Mocha	Black
Faceplate				

Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6	16/8/4	16/8/4
My Voice	-	-	-	-	-
Voice Ranger	•	•	•	•	•
Direct Streaming/ Auto Volume	- / -	- / -	- / -	- / -	- / -
Wireless sync	•	•	•	•	•
Directionality	Automatic / Adaptive OneMic Directionality	Automatic / Adaptive OneMic Directionality	Automatic / Adaptive OneMic Directionality	-	-
Noise Reduction	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing	Noise Management Sound Smoothing	Noise Management
Feedback Preventer	•	•	•	•	•
Wind Noise Reduction	•	•	•	•	-
Reverb Reducer	•	-	-	-	-
Bandwidth: Extension/Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer: Live/recorded/playing	• / • / •	• / - / -	-	-	-
Tinnitus	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy	-
Xphone	•	•	•	-	-
Acclimatization/Data Logging	• / •	• / •	• / •	• / •	• / •

Accessories & Apps					
Smart Key	•	•	•	•	•
CROS inoX 8C	•	•	•	-	-
CROS RIC 8C	-	-	-	-	-
Smart Transmitter 2.4	-	-	-	-	-
Smart Mic	-	-	-	-	-
Rexton App	•	•	•	•	•

Battery	
Battery runtime	up to 60 h


MYCORE STERLING IIC AND CIC 8C

ULTRA-DISCREET WITH MYCORE FEATURES


Discretion ensured by small components and custom-made housing. Discreet Smart Key or Rexton App allows users to change volume, programs and more with discretion.




- › Battery size 10
- › Programming via:
 - › ConnexxAir
 - › FlexConnector (CIC only)

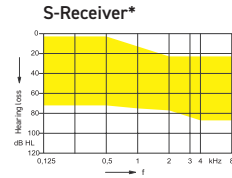

- 50 dB / 113 dB SPL (2 ccm coupler)
- 60 dB / 124 dB SPL (Ear simulator)

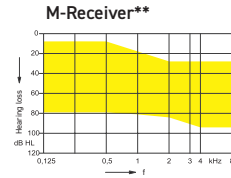

- 55 dB / 118 dB SPL (2 ccm coupler)
- 65 dB / 128 dB SPL (Ear simulator)

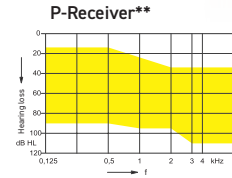

- 65 dB / 124 dB SPL (2 ccm coupler)
- 75 dB / 135 dB SPL (Ear simulator)

* Available for both IIC and CIC

** Available for CIC

Colors

	Mocha	Black	Brown
IIC			
CIC			

Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6	16/8/4	16/8/4
My Voice	-	-	-	-	-
Voice Ranger	•	•	•	•	•
Direct Streaming/ Auto Volume	- / -	- / -	- / -	- / -	- / -
Wireless sync	•	•	•	•	•
Directionality	Automatic / Adaptive OneMic Directionality	Automatic / Adaptive OneMic Directionality	Automatic / Adaptive OneMic Directionality	-	-
Noise Reduction	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing	Noise Management Sound Smoothing	Noise Management
Feedback Preventer	•	•	•	•	•
Wind Noise Reduction	•	•	•	•	-
Reverb Reducer	•	-	-	-	-
Bandwidth: Extension/Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer: Live/recorded/playing	• / • / •	• / - / -	-	-	-
Tinnitus	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy	-
Xphone	•	•	•	-	-
Acclimatization/Data Logging	• / •	• / •	• / •	• / •	• / •

Accessories & Apps					
Smart Key	•	•	•	•	•
CROS inoX 8C	•	•	•	-	-
CROS RIC 8C	-	-	-	-	-
Smart Transmitter 2.4	-	-	-	-	-
Smart Mic	-	-	-	-	-
Rexton App	•	•	•	•	•

Battery	
Battery runtime	up to 55 h


MYCORE STERLING ITC AND ITE 8C

PROVEN DESIGN – RELIABLE TECHNOLOGY


This hearing system provides a pleasant sound experience and the binaural signal processing ensures an easier listening in loud environments. The small-sized remote control or Rexton App allows users to change volume, programs, and more with discretion.




- › Battery size 312 or 13
- › Bluetooth or T-coil option available
- › Programming via:
 - › Noahlink Wireless (BT)
 - › ConnexxAir
 - › Flex Connector

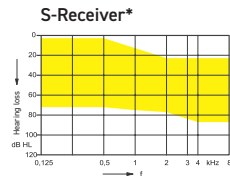

- 50 dB / 113 dB SPL (2 ccm coupler)
- 61 dB / 124 dB SPL (Ear simulator)

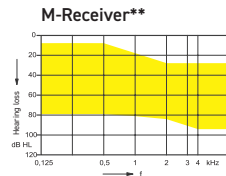

- 55 dB / 118 dB SPL (2 ccm coupler)
- 65 dB / 128 dB SPL (Ear simulator)

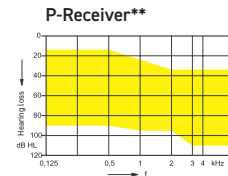


- 65 dB / 124 dB SPL (2 ccm coupler)
- 75 dB / 133 dB SPL (Ear simulator)

* Available for ITC

** Available for both ITC and ITE

Colors

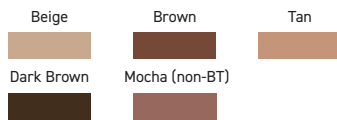

Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	48/20/6	32/16/6	24/12/6	16/8/4	16/8/4
My Voice	-	-	-	-	-
Voice Ranger	•	•	•	•	•
Direct Streaming (BT version only)	Made for iPhone	Made for iPhone	Made for iPhone	Made for iPhone	Made for iPhone
Auto Volume	•	•	•	•	•
Wireless sync	•	•	•	•	•
Directionality	Automatic / Adaptive Front & Back Left & Right Narrow - Stereo iLock	Automatic / Adaptive Front & Back Narrow - iLock	Automatic / Adaptive Narrow - iLock	Automatic / Adaptive	Automatic / Adaptive
Noise Reduction	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing	Noise Management Sound Smoothing	Noise Management
Feedback Preventer	•	•	•	•	•
Wind Noise Reduction	•	•	•	•	-
Reverb Reducer	•	-	-	-	-
Bandwidth: Extension/Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer: Live/recorded/playing	• / • / •	• / - / -	-	-	-
Tinnitus	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy	-
Xphone	•	•	•	-	-
Acclimatization/Data Logging	• / •	• / •	• / •	• / •	• / •

Accessories & Apps					
Smart Key	•	•	•	•	•
CROS inoX 8C (non BT version only)	•	•	•	-	-
CROS RIC 8C (BT version only)	•	•	•	-	-
Smart Transmitter 2.4 (BT version only)	•	•	•	•	•
Smart Mic (BT version only)	•	•	•	•	•
Rexton App	•	•	•	•	•

Battery		
Battery runtime (without streaming)		up to 100 h
Battery runtime (incl. 5 h streaming)		up to 70 h


MYCORE EMERALD S 8C

A PERFECT COMBINATION OF SMALL SIZE AND PERFORMANCE

Supports direct streaming from iPhone via *Bluetooth®* and easier streaming to most *Bluetooth®*-enabled phones via Smart Mic.



- › Battery size 312
- › Smart Demo device available
- › Programming via:
 - › Noahlink Wireless
 - › HiPro II

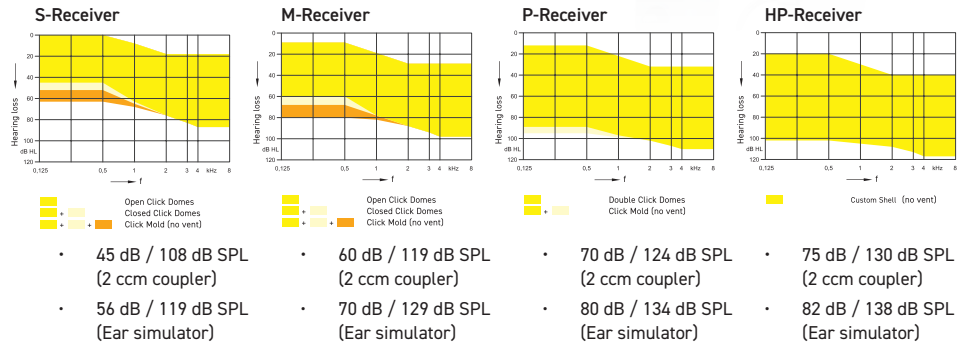

Colors: Beige, Dark Brown, Dark Granite, Granite, Grey, Pearl White, Sandy Brown, Silver

Features / Performance Level	80	60	40	30	20
Channels/Controls/Programs	48 / 20 / 6	32 / 16 / 6	24 / 12 / 6	16 / 8 / 4	16 / 8 / 4
My Voice	•	•	•	-	-
Voice Ranger	•	•	•	•	•
Direct Streaming	Made for iPhone	Made for iPhone	Made for iPhone	Made for iPhone	Made for iPhone
Auto Volume	•	•	•	•	•
Wireless Sync	•	•	•	•	•
Directionality	Automatic / Adaptive iOmni Front & Back Left & Right Narrow - iLock	Automatic / Adaptive iOmni Front & Back Narrow - iLock	Automatic / Adaptive iOmni Narrow - iLock	Automatic / Adaptive iOmni	Automatic iOmni
Noise Reduction	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing Directional	Noise Management Sound Smoothing	Noise Management Sound Smoothing	Noise Management
Feedback Preventer	•	•	•	•	•
Wind Noise Reduction	•	•	•	•	-
Reverb Reducer	•	-	-	-	-
Bandwidth: Extension/Compression	• / •	- / •	- / •	- / •	- / •
Music Enhancer: Live/recorded/playing	• / • / •	• / - / -	• / - / -	-	-
Tinnitus	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy	Sound Therapy Notch Therapy
Xphone	•	•	•	-	-
Acclimatization/Data Logging	• / •	• / •	• / •	• / •	• / •
Accessories & Apps					
Smart Key	•	•	•	•	•
CROS inoX 8C	-	-	-	-	-
CROS RIC 8C	•	•	•	•	-
Smart Transmitter 2,4	•	•	•	•	-
Smart Mic	•	•	•	•	•
Rexton App	•	•	•	•	•
Battery*					
Battery runtime	up to 70 h				

• available - not available

* Battery runtime value based on S-Receiver.







ACCESSORIES & APPS





MORE THAN A HEARING AID, A COMPLETE EXPERIENCE

Rextons commitment is to provide a reliable hearing experience and that goes beyond a hearing device.

With our App and accessories, users can get the most out of their experience by controlling and fine tuning their hearing aid at anytime and anywhere by themselves.

Through a range of accessories users can also stream directly to their devices and just enjoy their preferred streaming applications. Our charger range offers flexibility and power options that suits different tastes and needs.

In one click users can connect to Healthcare Professional for assistance or get help through an automated artificial intelligent system with 100% of efficiency.

Rely on Rexton for a complete and seamless hearing experience all day.

CHARGERS

Rexton offers different charger types so that you can choose the one that best suits your needs.

STANDARD

Get the job done with Rexton stationary charger. Practical and easy-to use, the stationary charger solution needs to be plugged into an outlet.



DRY & CLEAN

Go beyond and offer extra care to the hearing device with an integrated drying function and UV-cleaning as it charges.



ON THE GO

The travel charger gives the freedom of built-in power banks that holds several extra charges. This means that the user don't have to stay close to a power outlet to charge their hearing aids.



CABLE CHARGER

Recharge at home or on the go. Plug in and power up with the handy cable charger. Can recharge one or two hearing aids.



CHARGING STATION B-M

ROBUST AND FUNCTIONAL



A highly powerful hearing aid needs a powerful and long lasting charger. Its' galvanic charging system ensures a lossless charging process.

The housing has a robust and scratch resistant design and is very easy to use. Just place the hearing aids into the intended slots to charge it overnight or make use of the 30 minute fast charging functionality for a quick boost of extra wearing time if you're in a hurry.



COMPATIBILITY

- B-Li M Rugged

SPECIAL FEATURES

- Robust and scratch resistant design
- 30-minute fast charge functionality
- Automatic on/off
- Fixed USB charging cable
- Included power adapter

CUSTOM CHARGER STATION

CONTACTLESS CHARGING ENABLED BY MAGNETIC RESONANCE TECHNOLOGY



Featuring a new technology, the new Custom Charger Station offers optimal charging performance and ease of use with its intuitive placement slots.



COMPATIBILITY

- BiCore Custom-Li

SPECIAL FEATURES

- Intuitive placement in charger
- 30-minute fast charge functionality
- Automatic on/off
- USB-C and power adapter included

CHARGING STATION R

ROBUST AND FUNCTIONAL



A highly powerful hearing aid needs a powerful and long lasting charger. It's galvanic charging system ensures a lossless charging process.

The housing has a robust and scratch resistant design and is very easy to use. Just place the hearing aids into the intended slots to charge it overnight or make use of the 30 minute fast charging functionality for a quick boost of extra wearing time if you're in a hurry.



COMPATIBILITY

- Reach R-Li/R-Li T
- BiCore R-Li/R-Li T

SPECIAL FEATURES

- Robust and scratch resistant design
- 30-minute fast charge functionality
- Automatic on/off
- Fixed USB charging cable
- Included power adaptor

CHARGING STATION B-P

ROBUST AND FUNCTIONAL



The galvanic charging system ensures a lossless charging process.

The housing has a robust and scratch resistant design and is very easy to use. Just place the hearing aids into the intended slots to charge it overnight or make use of the 30 minute fast charging functionality.



COMPATIBILITY

- M-Core B Li P

SPECIAL FEATURES

- Robust and scratch resistant design
- 30-minute fast charge functionality
- Automatic on/off
- Fixed USB charging cable
- Included power adapter

CHARGING STATION B-HP

ROBUST AND FUNCTIONAL



The galvanic charging system ensures a lossless charging process.

The housing has a robust and scratch resistant design and is very easy to use. Just place the hearing aids into the intended slots to charge it overnight or make use of the 30 minute fast charging functionality for a quick boost of extra wearing time if you're in a hurry.



COMPATIBILITY

- M-Core B Li HP

SPECIAL FEATURES

- Robust and scratch resistant design
- 30-minute fast charge functionality
- Automatic on/off
- Fixed USB charging cable
- Included power adapter

CHARGING STATION



Inductive charger with a lid covering the hearing aids. Simply place the hearing aids in the slots to start the charging. The LEDs on the charger indicates how much battery the hearing aids have.

An automatic switch-off guarantees overcharge protection. The charger also has a drying function while charging.



COMPATIBILITY

- M-Core B-Li M
- M-Core R-Li

SPECIAL FEATURES

- Design with case cover
- 30-minute fast charge functionality
- Automatic on/off
- Inductive charging
- Micro-USB and power adaptor included

CHARGING+ STATION B-P / B-HP

ADVANCED PRACTICALITY



The premium version compared to the standard charger with additional features.

It is equipped with an integrated drying functionality, an integrated UV-cleaning light and also offers a 30 minute fast charging functionality.

It's robust and scratch resistant finish is ideal for the daily use and for a long lasting operation.



COMPATIBILITY

- M-Core B-Li P, M-Core B-Li HP

SPECIAL FEATURES

- Robust and scratch resistant design
- 30-minute fast charge functionality
- Automatic on/off
- Integrated drying functionality
- Integrated UV-cleaning functionality
- Fixed USB charging cable
- Included power adapter

CHARGING+ STATION R

ADVANCED PRACTICALITY



The premium version compared to the standard charger with additional features.

It is equipped with an integrated drying functionality, an integrated UV-cleaning light and also offers a 30 minute fast charging functionality.

It's robust and scratch resistant finish is ideal for the daily use and for a long lasting operation.



COMPATIBILITY

- Reach R-Li/R-Li T
- BiCore R-Li/R-Li T

SPECIAL FEATURES

- New robust and scratch resistant design
- 30-minute fast charge functionality
- Automatic on/off
- Integrated drying functionality
- Integrated UV-cleaning functionality
- USB-C and power adapter included

CHARGING+ STATION SR

ADVANCED PRACTICALITY



The premium charging version equipped with an integrated drying functionality and UV-cleaning light.

It's robust and scratch resistant finish is ideal for the daily use and for a long lasting operation.



COMPATIBILITY

- Reach SR
- BiCore SR

SPECIAL FEATURES

- 30-minute fast charge functionality
- Charges in 4 hours
- Integrated drying functionality
- Integrated UV-cleaning functionality
- Automatic on/off
- USB-C and power adapter included

TRAVEL CHARGER RIC

PORTABLE AND LIGHT



This portable charger provides up to 3 charges on the go – without connecting to a power outlet.



COMPATIBILITY

- Reach R-Li/R-Li T
- BiCore R-Li/R-Li T

SPECIAL FEATURES

- Holds 3 charges
- 30-minute fast charge functionality
- Charges in 4 hours
- Automatic on/off
- USB-C and power adapter included

SLIM-RIC TRAVEL CHARGER

WITH Qi CHARGING



Featuring an ultra-convenient Qi wireless charging system for the internal power bank and advanced lithium-ion technology, the sleek and compact design provides up to four days hearing on-the-go.

The fast charge feature provides up to five hours wearing time in just 30 minutes, while the case can be fully recharged in just three hours.



COMPATIBILITY

- M-Core SR

SPECIAL FEATURES

- Holds 3 charges
- 30-minute fast charge functionality
- The size is optimized to accommodate:
 - Devices with miniReceivers 2.0 length 0-3
 - Devices with miniReceivers 2.0 types S, M and P
- Automatic on/off

SLIM-RIC TRAVEL CHARGER 2

WITH Qi CHARGING



Featuring an ultra-convenient Qi wireless charging system, this portable charger provides up to 4 days hearing on the go – without connecting to the power outlet.

A 30-minute fast charge gives an extra 5.5 hours of runtime.



COMPATIBILITY

- Reach SR
- BiCore SR

SPECIAL FEATURES

- Holds 4 charges (1 more vs. the previous generation)
- 30-minute fast charge functionality
- LEDs to show battery status of HI and Charger
- Automatic on/off
- USB-C and power adapter included

IF CIC TRAVEL CHARGER

WITH Qi CHARGING



Small and portable charger with a built-in power bank for 4 extra charges on the go.

With Qi charging you are not dependent on power outlets and the right and left indicators inside the charger makes it easy to use for anyone.

COMPATIBILITY

- Reach inoX-CIC Li

SPECIAL FEATURES

- Holds 4 charges
- 30-minute fast charge functionality
- Charges in 4 hours
- Automatic on/off
- USB-C and power adaptor included

CABLE CHARGER

A HANDY AND LIGHT-WEIGHT CHARGER



Rexton ReCharge is the new rechargeable hearing aid with a handy cable charger, and it makes the convenience of rechargeability more affordable for all. With ReCharge you get new energy for work, hobbies, travel, nights out, for time with the family, and for conversations everywhere. Just plug in, power up and off you go!

COMPATIBILITY

- Mandatory for ReCharge products
- Reach R-Li/R-Li T
- BiCore R-Li/R-Li T
- BiCore B-Li M Rugged
- M-Core B Li P & HP


SPECIAL FEATURES

- Use at home or on the move
- Simple to use
- Plug in, power up and away you go
- One-charger-fits-all
- Optional stand available for optimized stability



THE COMPREHENSIVE REXTON APP

ONE FOR ALL

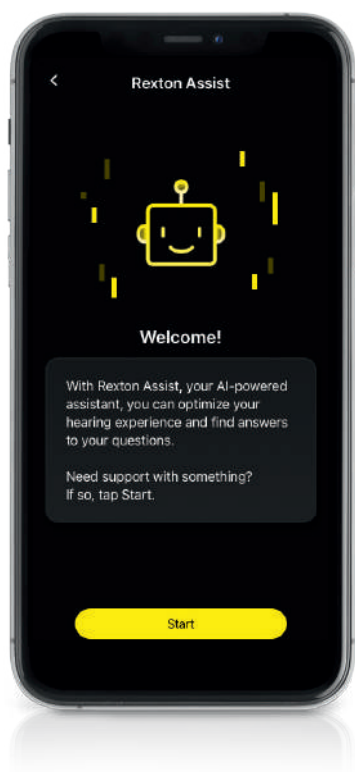

The Rexton App allows users to control their hearing system conveniently as it turns a smartphone into a discreet remote control for the hearing aids.

Telephone calls are forwarded directly from the smartphone and take priority over other audio signals. Streamed music can be listened to in stereo, directly and without any additional device*.

SPECIAL FEATURES

- Hearing programs and volumes can be adjusted discreetly
- The volume can be set individually for the right and left hearing system
- Directionality control
- Remote support
- Real time battery status

Download the Rexton App now!




iOS



Android


*Check the compatibility of your hearing aid and smartphone device

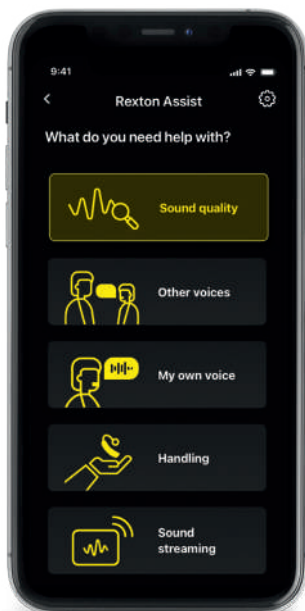

REXTON ASSIST

Our cloud-based Artificial Intelligence neural network allows users to make real-time programming changes when they are on the go. Users can personalize their hearing aids' performance when and where challenges are occurring without having to wait for an in-person or virtual appointment.

After receiving solutions curated for that specific situation, users provide feedback to the living network so it can evolve based on users satisfaction level.

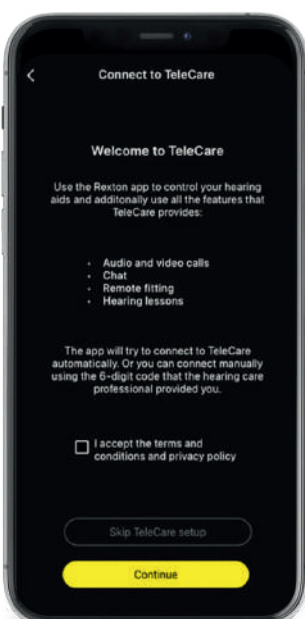

REXTON TELECARE

With Telecare**, hearing care professionals can increase their service journey and provide remote support for their patients through:

- Instant messages
- Audio & Video calls
- Fine tuning

**Telecare interface can be different across technologies and form factors

AUDIO STREAMING DEVICES

Rexton offers audio streaming devices as interfaces to keep users connected for a pleasant hearing experience, no matter what. Compatible with Bluetooth devices only.

THE PRACTICAL SMART TRANSMITTER 2.4

THE BLUETOOTH® TRANSMITTER FOR TELEVISIONS AND STEREO SYSTEMS



The Smart Transmitter 2.4 allows direct transfer of stereo sound from audio sources such as television sets and stereo systems to the hearing systems. Simply place the hearing systems on the transmitter and the pairing process begins. The transmitter can be paired with multiple hearing systems at the same time for shared enjoyment with friends.

SPECIAL FEATURES

- No separate remote control needed
- Connecting/pairing is particularly easy
- Several hearing systems can connect simultaneously
- Analog and digital audio inputs (Cinch and TOSLINK)

THE COMPACT SMART MIC

DIRECT AUDIO STREAMING FROM PRACTICALLY ANY BLUETOOTH-ENABLED ANDROID SMARTPHONE



The Smart Mic enables hands-free telephone calls without your customer needing to pick up his/her smartphone. In challenging hearing environments it can easily be used as an additional microphone.

SPECIAL FEATURES

- Direct streaming from Android smartphones
- Listen to music in stereo quality
- Can be used as a conference microphone

REMOTE CONTROL DEVICE

Control the volume of your hearing aid easily and quickly with a touch.

THE HANDY SMART KEY

SMALL AND PRACTICAL



With Smart Key, programs and volume can be easily adjusted. A keylock prevents buttons from being inadvertently pressed. A status LED indicates whether the remote control is enabled.

SPECIAL FEATURES


- Small and handy
- Control via high-frequency sounds

EARWEAR 2.0

EarWear is Rexton solution for an optimal fitting. With a complete range of receivers, thin tubes and ear tips, we offer different options that can fit users unique ears with comfort for a great acoustical experience.

MINIRECEIVER 2.0

The miniReceiver 2.0 is a smart RIC receiver that uses SIM technology. It securely and reliably stores and transmits data between the receiver and the hearing aid. The SIM card can also detect an incorrect receiver and automatically reduce gain to protect the user. And the design enables more ergonomic placement of the receiver for better microphone placement.

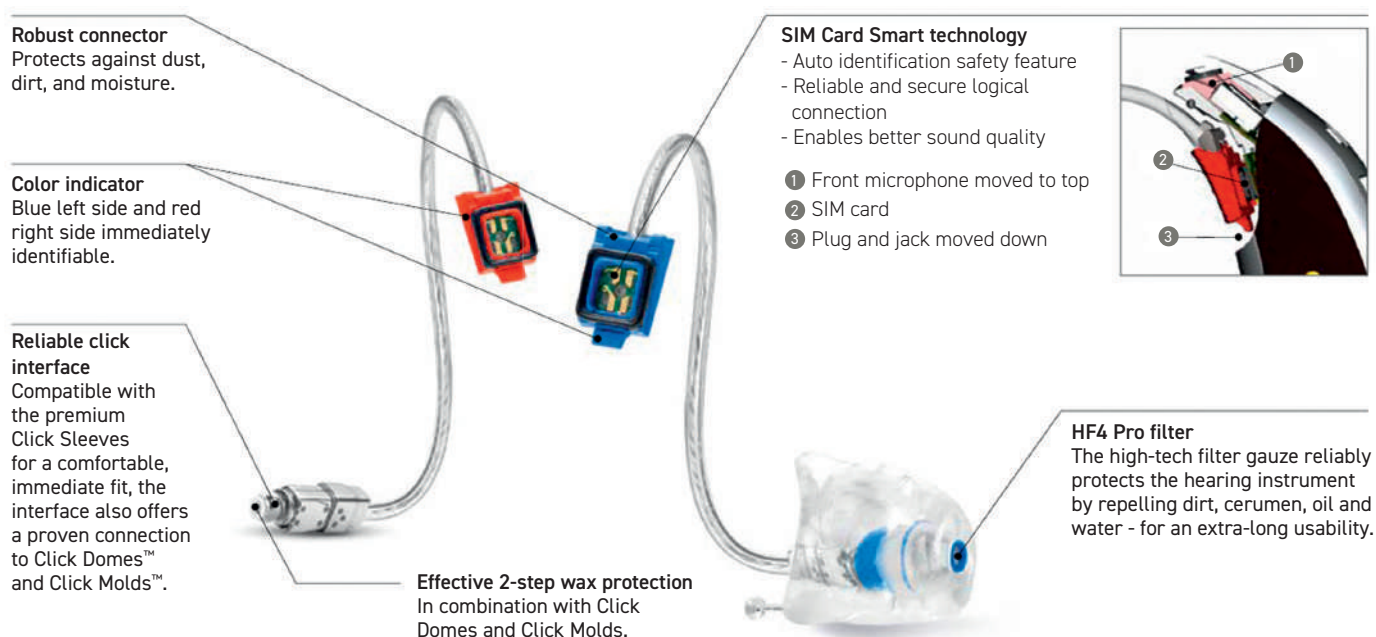


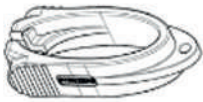




COMPATIBLE PLATFORMS





POWER	LENGTH
S, M, P, HP	0, 1, 2, 3, 4* *only available for M-Receiver

SUPPORTING TOOLS

ITEM		PART NUMBER
Removal tool		10674370
Measuring gauge for S and M		10174571
Measuring gauge for P and HP		10674512
Click Dome removal tool		10605854
Click Mold service set		10597930
Wax Guard		10962332

TYPES OF MINIRECEIVER 2.0

	SIZE	UNITS	ITEM NUMBER
 S-RECEIVER	miniReceiver 2.0 0L S	1	10669907
	miniReceiver 2.0 0R S	1	10669908
	miniReceiver 2.0 1L S	1	10669909
	miniReceiver 2.0 1R S	1	10669910
	miniReceiver 2.0 2L S	1	10669911
	miniReceiver 2.0 2R S	1	10669912
	miniReceiver 2.0 3L S	1	10669913
	miniReceiver 2.0 3R S	1	10669914
 M-RECEIVER	miniReceiver 2.0 0L M	1	10824820
	miniReceiver 2.0 0R M	1	10824821
	miniReceiver 2.0 1L M	1	10824822
	miniReceiver 2.0 1R M	1	10824823
	miniReceiver 2.0 2L M	1	10824824
	miniReceiver 2.0 2R M	1	10824825
	miniReceiver 2.0 3L M	1	10824826
	miniReceiver 2.0 3R M	1	10824827
 P-RECEIVER	miniReceiver 2.0 4L M	1	10933054
	miniReceiver 2.0 4R M	1	10933053
	miniReceiver 2.0 0L P	1	10669923
	miniReceiver 2.0 0R P	1	10669924
	miniReceiver 2.0 1L P	1	10669925
	miniReceiver 2.0 1R P	1	10669926
	miniReceiver 2.0 2L P	1	10669927
	miniReceiver 2.0 2R P	1	10669928
	miniReceiver 2.0 3L P	1	10669929
	miniReceiver 2.0 3R P	1	10669930







HP

Please order the miniReceiver size HP with a custom shell in the same way as a custom hearing aid using the order form.

EAR TIPS FOR MINIRECEIVER 2.0



CLICK DOMES

	TYPE	SIZE
	The open dome is for users who require maximum comfort and no occlusion	4 mm
		8 mm
		10 mm
	The semi-open dome combines the benefits of an open and a closed dome for users who desire the comfort of an open fitting but need more power.	8-12 mm
	The closed dome is for users who require more loudness than the open dome can offer.	6 mm
		8 mm
		10 mm
	The double dome is for users who require broadband amplification and more retention. Can also be used with a miniReceiver size P.	8/10 mm
		10/12 mm

CLICK SLEEVES


TYPE	SIZE
Vented 	XS (6pcs)
	S (6pcs)
	M (6pcs)
	L (6pcs)
Closed 	XS (6pcs)
	S (6pcs)
	M (6pcs)
	L (6pcs)

CLICK MOLD TYPES

TYPE

Short mold Open mold, short

Long mold Open mold, long

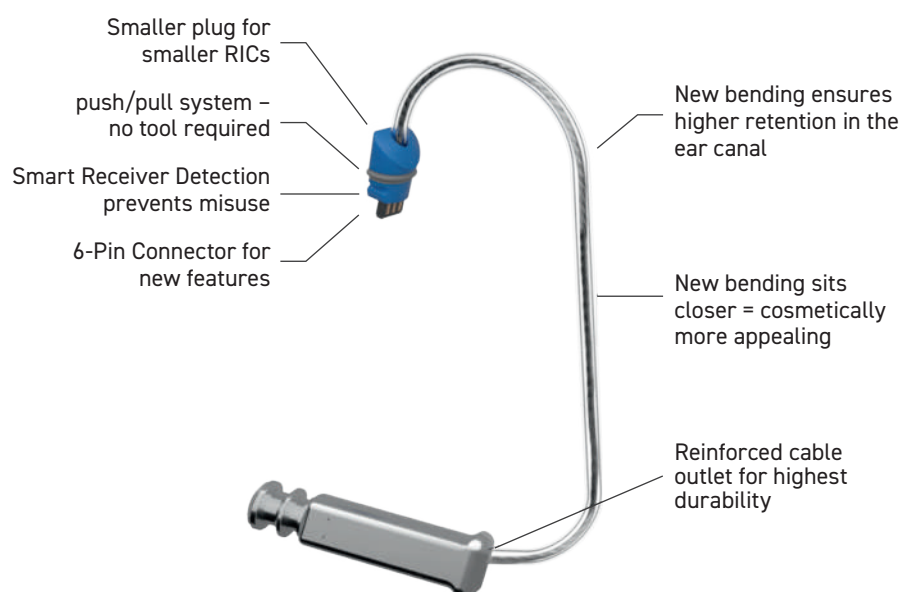EARWEAR 3.0

EarWear 3.0 provides a comprehensive solution which is highly reliable in its acoustic performance and durability. It improves the cosmetic appeal as well as physical fit and retention in the ear canal.



COMPATIBLE HEARING AIDS

- Reach R-Li/R-Li T
- Reach SR
- BiCore R-Li/R-Li T
- BiCore SR
- BiCore R312

POWER


- S, M, P and HP (HP is not recommended for SR due to size).

LENGTHS

- 0, 1, 2, 3, 4 and 5 (5 is not recommended for SR due to size).

Improved design for even better protection against earwax

WaxGuard 3.0 NanoCare reliably protects the receiver from moisture and earwax

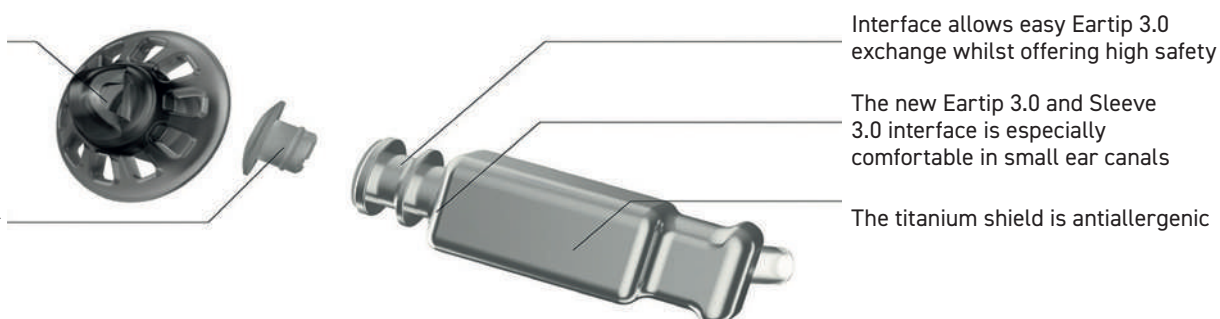

COMPATIBLE PLATFORMS



S



M



P



HP



	UNITS	ITEM NUMBER
miniReceiver 3.0 0R S	1	10992548
miniReceiver 3.0 0L S	1	10992549
miniReceiver 3.0 1R S	1	10992550
miniReceiver 3.0 1L S	1	10992552
miniReceiver 3.0 2R S	1	10992553
miniReceiver 3.0 2L S	1	10992554
miniReceiver 3.0 3R S	1	10992555
miniReceiver 3.0 3L S	1	10992556
miniReceiver 3.0 4R S	1	10992557
miniReceiver 3.0 4L S	1	10992558
miniReceiver 3.0 5R S	1	21023207
miniReceiver 3.0 5L S	1	21023211
miniReceiver 3.0 0R M	1	10992559
miniReceiver 3.0 0L M	1	10992560
miniReceiver 3.0 1R M	1	10992561
miniReceiver 3.0 1L M	1	10992562
miniReceiver 3.0 2R M	1	10992563
miniReceiver 3.0 2L M	1	10992564
miniReceiver 3.0 3R M	1	10992565
miniReceiver 3.0 3L M	1	10992566
miniReceiver 3.0 4R M	1	10992567
miniReceiver 3.0 4L M	1	10992568
miniReceiver 3.0 5R M	1	21023215
miniReceiver 3.0 5L M	1	21023216
miniReceiver 3.0 0R P	1	10992569
miniReceiver 3.0 0L P	1	10992570
miniReceiver 3.0 1R P	1	10992571
miniReceiver 3.0 1L P	1	10992572
miniReceiver 3.0 2R P	1	10992573
miniReceiver 3.0 2L P	1	10992574
miniReceiver 3.0 3R P	1	10992575
miniReceiver 3.0 3L P	1	10992576
miniReceiver 3.0 4R P	1	10992577
miniReceiver 3.0 4L P	1	10992578
miniReceiver 3.0 5R P	1	21023217
miniReceiver 3.0 5L P	1	21023218

Please order the miniReceiver size HP with a custom shell in the same way as a custom hearing aid using the order form.

THINTUBE 3.0

New cosmetically highly appealing options – also for severe hearing losses.



COMPATIBLE HEARING AIDS

- BiCORE B-Li M Rugged
- BiCORE B M
- BiCORE B P
- M-Core B-Li M
- M-Core B-Li P

LENGTHS

- 0, 1, 2, 3 and 4

miniReceiver 3.0 and ThinTube 3.0 use the same Eartip 3.0 and Sleeve 3.0 for maximum simplicity and the same cable/tube geometry for best cosmetics.



ThinTube 3.0 (0.9mm)



	UNITS	ITEM NUMBER
Set ThinTube 3.0 S0	1R + 1L	10991609
Set ThinTube 3.0 S1	1R + 1L	10991610
Set ThinTube 3.0 S2	1R + 1L	10991611
Set ThinTube 3.0 S3	1R + 1L	10991612
Set ThinTube 3.0 S4	1R + 1L	10991613

ThinTube 3.0 P (1.4mm)







Set ThinTube 3.0 S0 P	1R + 1L	10992492
Set ThinTube 3.0 S1 P	1R + 1L	10992493
Set ThinTube 3.0 S2 P	1R + 1L	10992494
Set ThinTube 3.0 S3 P	1R + 1L	10992495
Set ThinTube 3.0 S4 P	1R + 1L	10992496

EAR TIP SOLUTIONS FOR MINIRECEIVER 3.0 AND THINTUBE 3.0





SLEEVE 3.0 VENTED

Open option, still highly trustable retention and comfort in the ear canal

	SIZE/UNITS	ITEM NUMBER
	XS (6 pcs)	10992863
	S (6 pcs)	10992864
	M (6 pcs)	10992865
	L (6 pcs)	10992866





SLEEVE 3.0 CLOSED

Similar to tulips, with one vent for delivering the right amount of gain


	XS (6 pcs)	10992867
	S (6 pcs)	10992868
	M (6 pcs)	10992869
	L (6 pcs)	10992870

SLEEVE 3.0 POWER


Fully closed option, replaces the double domes, for a reliable acoustical performance

	XS (6 pcs)	10997192
	S (6 pcs)	10997193
	M (6 pcs)	10997194
	L (6 pcs)	10997195

EARTIP 3.0 OPEN

	SIZE/UNITS	ITEM NUMBER
	5mm (6 pcs)	10993678
	7mm (6 pcs)	10993680
	10mm (6 pcs)	10993681

EARTIP 3.0 TULIP

	8mm (6 pcs)	10993682
	12mm (6 pcs)	10993683

EARMOLD 3.0



Short



Long

SUPPORTING TOOLS

ITEM		PART NUMBER
Measuring Gauge 3.0		10991403
Waxguard 3.0		10993649
Concha Lock 3.0 (10 pcs)		10993677
Optional tool to exchange eartips		10984033
Short custom mold exchange tool		10967340



PORTFOLIO OVERVIEW

RIC













Slim-RIC

BTE



Product	Reach R-Li NEW	Reach R-Li T NEW	ReCharge R-Li NEW	BiCore R-Li	BiCore R-Li T	BiCore R312	M-Core R-Li	M-Core R312	Reach SR NEW	BiCore SR	M-Core SR	BiCore B-Li M Rugged
Battery	Li-ion	Li-ion	Li-ion	Li-ion	Li-ion	312	Li-ion	312	Li-ion	Li-ion	Li-ion	Li-ion
Programming Interface	Noahlink Wireless	Noahlink Wireless	Noahlink Wireless	Noahlink Wireless	Noahlink Wireless	Noahlink Wireless	Noahlink Wireless	Noahlink Wireless	Noahlink Wireless	Noahlink Wireless	Noahlink Wireless	Noahlink Wireless
TeleCare	Full Live Tuning	Full Live Tuning	Basic Portal Tuning	Full Live Tuning	Full Live Tuning	Full Live Tuning	Full Live Tuning	Full Live Tuning	Full Live Tuning	Full Live Tuning	Full Live Tuning	Full Live Tuning
Bluetooth Enabled	Yes	Yes	No	Yes	Yes	Yes	Yes	Yes	Yes	Yes	Yes	Yes
T-coil	No	Yes	No	No	Yes	No	No	Yes	No	No	No	Yes
Rexton Assist	Yes	Yes	No	Yes	Yes	Yes	Yes	Yes	Yes	Yes	Yes	Yes

* not available for M products

Instant Fit						Custom					
											
BiCore B M,P,HP	M-Core B-Li M,P,HP	MyCore Emerald S 8C	ReCharge B-Li P NEW	ReCharge B-Li HP NEW	Reach inoX-CIC Li NEW	M-Core iX-CIC	MyCore inoX CIC 8C	BiCore Custom Li	Reach Custom IIC/CIC NEW	MyCore Sterling IIC/CIC 8C	MyCore Sterling ITC/ITE 8C
13/13/675	Li-ion	312	Li-ion	Li-ion	Li-ion	10	10/312	Li-ion	10	10	312/13
Noahlink Wireless	Noahlink Wireless	Noahlink Wireless, HiPro II	Noahlink Wireless	Noahlink Wireless	ConnexxAir	Connexx- Air, HiPro II	Connexx- Air, HiPro II	Noahlink Wireless	Connexx- Air, HiPro II	Connexx- Air, HiPro II	Noahlink Wireless, Connexx- Air, HiPro II
Full Live Tuning	Full Live Tuning	Basic Portal Tuning	Basic Portal Tuning	Basic Portal Tuning	Basic Portal Tuning	Basic Portal Tuning	Basic Portal Tuning	Full Live Tuning	Basic Portal Tuning	Basic Portal Tuning	Basic Portal Tuning
Yes	Yes	Yes	No	No	No	No	No	Yes	No	No	Yes
Yes*	Yes*	No	Yes	Yes	No	No	No	No	No	No	Yes (for non BT)
Yes	Yes	No	No	No	No	No	No	Yes	No	No	No

GLOSSARY

A

ADAPTIVE COMPRESSION

Adaptive compression enables automatic adjustments in sound amplification according to the volume of incoming signals. This function ensures that our users can hear both soft and loud sounds comfortably.

ASHA PORTOCOLS

Phones with Bluetooth 5.0+ and Android 10 (or higher) are capable of providing wireless streaming via ASHA (Audio Streaming for Hearing Aids)

ATEX

ATEX certification verifies that our hearing aids are tested and approved to be intrinsically safe for use in potentially explosive environments. They can be safely used in areas where there may be flammable substances in the air such as paint and chemical manufacturing, petrol stations and mills. The ATEX group/class is printed on the ATEX card, which must be available at hand when the hearing aids are used in explosive atmospheres.

AUTO ECHO REDUCER

Automated environment scan to reduce effects of disturbing reverberation on the perceived sound signal.

AUTO STEREO ILOCK

With a narrow focus beam, the hearing aids can efficiently improve the signal-to-noise ratio in difficult listening situations. Giving a clear focus to the direct front of the user.

AUTO VOLUME

Adapts amplification levels based on changes in the surrounding sound environment, ensuring optimal comfort and clarity for users in various listening conditions.

AUTOMATIC ADAPTIVE

It detects sounds from a specific direction, automatically change its focus direction, and adjust noise filters and amplification settings accordingly.

B

BANDWIDTH EXTENSION/COMPRESSION

Adjusts the frequency range of audio signals by either expanding or compressing it. This ensures sound quality by adding higher frequencies or narrowing the range while preserving clarity and fidelity.

BINAURAL LATENCY RATE

Binaural latency rate refers to the time lapse between a sound reaching one ear and the other in a two-ear auditory setup. It's essential for spatial awareness and pinpointing sound sources. In hearing aids, reducing this delay ensures both ears process sound simultaneously, enhancing the perception of natural sound.

BT LE AUDIO

BT LE Audio facilitates wireless audio streaming with minimal power usage, ensuring high-quality sound transmission. This capability enables easy connection to compatible devices, enhancing communication and audio streaming in hearing aids.

C

CONNEXXAIR

ConnexxAir is our programming interface between our wireless Rexton hearing aids and our Rexton fitting software.

CROS

CROS (Contralateral Routing of Signals): This hearing aid system is designed for one-sided deafness, transferring sounds from the impaired ear to the better-hearing ear, enhancing overall auditory perception.

D

DATA LOGGING

Automatically records usage details like environmental sounds and user preferences. This information helps professionals optimize performance based on actual usage patterns.

DIRECTIONAL MICROPHONE

A microphone that concentrates on capturing sounds from particular directions, effectively diminishing background noise from other directions.

DYNAMIC EXTENDER

Increases the dynamic range of sounds. Enhances sound quality, especially important when listening to music.

E

EXTENDED DYNAMIC RANGE

An extended dynamic range of 117 dB SPL is important to secure stable and high-quality sound when listening to music or while in a very noisy environment, keeping the speech clarity.

F

FEEDBACK PREVENTER

Identifies and handles feedback at a high speed. Lets users rely on and enjoy a whistle-free hearing experience.

FOCUS STREAMS

Virtual beams created to support independent directional processing of the listening environment.

I

IFOCUS 360

Tracks the direction of the loudest voice signal and steers the focus beam to the front or back, left or right depending on the direction of speech.

IOMNI SOUND LOCATOR

iOmni helps locating the sound source by replicating the natural pinna effect of an outer ear.

IP68

An IP68 rating means strong defense against dust and water, guaranteeing durability and resilience, and ensuring the lifeproof of our hearing aids.

M

MFI

Rexton Mfi (made for iPhone, iPad or iPod touch) hearing aids connect directly to your iPhone, iPad or iPod so that you can stream your phone calls and music directly into your hearing aids.

MOTION SENSOR FOR SCENE ANALYSIS

Integrated real-time motion sensor that detects a user's movement, and uses the data to even more accurately determine their sound scene.

MULTI-VOICE FOCUS

Developed of and based on Speech Preservation Technology. It scans the surrounding environment 1000 times a second with 4 independent focus beams. This enables it to constantly adapt to the changing dynamics of the conversation and make the experience fluid and comfortable. It will boost the prominence of those in the group who are speaking while keeping track of those who are silent, and all without switching between modes, and with no delays.

MUSIC ENHANCER

Provides exceptional sound quality created for music. Whether attending a concert, listening to music at home, or even performing live on stage, Music Enhancer delivers a listening experience like no other.

MY VOICE/MY VOICE 2.0

The hearing aid recognizes the user's specific voice pattern and integrates it naturally into what they hear. In turn, this provides the user with a natural-sounding voice that leads to high acceptance rates. And for that the user will communicate more and be more socially active, greatly improving their quality of life.

N

NOAHLINK WIRELESS

Facilitates wireless programming and instant fitting for hearing aids and compatible devices.

O

ONEMIC DIRECTIONALITY

Thanks to Rextons wireless connection between the hearing aids, they can exchange the input data with one another. And although they only have one microphone each, in a pair, these hearing aids can then create a directionality which will improve the speech intelligibility and sound experience of the user.

R

REXTON ASSIST

Our cloud-based Artificial Intelligence assistant, is built right into the Rexton App. By clicking the icon you can answer a few short questions and Rexton Assist will help you make tailored, real-time sound adjustments on the go and support with basic troubleshooting through videos.

REVERB REDUCER

Dedicated program for pleasant communication in an echoing room, the reverb reducer can efficiently minimize intrusive echoes so that the user can focus on speech sources.

REXTON TELECARE

Enables remote fitting sessions via Connexx resulting in higher user convenience and satisfaction.

S

SMART AUTOMATIC ACCLIMATIZATION

Smart automatic acclimatization adapts amplification settings over time, aiding users in comfortably adjusting to new hearing aid technology.

SMART AUTOMATIC EQUALIZER

The smart automatic equalizer adapts frequency settings to suit the user's requirements and the surrounding environment.

SOUND SMOOTHING

Letting the hearing aids take care of any high pitched sounds, without distorting it or amplifying it.

SOUND THERAPY AND NOTCH

Rexton tinnitus management systems, available in several different setups to accomodate the user that need extra tinnitus management throughout the day.

SOUNDPRO

Soundpro can automatically and continuously adjust the sound settings to provide the most reliable impression of the user's environment.

SOUNDPRO 2.0

Automatically steers all sound processing features, including Speech Preservation Technology and My Voice 2.0. Additionally, uses an advanced classification system for enhanced steering of features.

SOUNDPRO 3.0

Automatically steers all sound processing features, including the new Multi-Voice Focus*, based on the analyzis of the acoustic scene. When changing the strategy to Natural Focus or Speech Focus, the recommended balance of features is modified. The Balanced Focus is recommended for enhanced automatic performance.

* Only available for hearing aids with two microphones in a binaural fitting.

SPEECH PRESERVATION TECHNOLOGY

Removes unwanted noise and emphasizes speech signals by offering a dedicated pathway to process speech independently from the one to process noise, so that you can focus on the conversation, even in busy environments.

STEREO ILOCK

A technology that allows for a very narrow directional listening beam, focusing on the immediate front. Improves the signal-to-noise ratio in difficult listening situations.

T

T-COIL

Captures sounds from phones and public venues equipped with loop systems, improving phone calls and listening in noisy environments for individuals with hearing loss.

V

VOICE RANGER

Coordinates all noise and speech management tools to adjust to each sound environment for a more relaxed listening experience.

VOICE STABILIZER

Controls and adjust the compression to secure that all voices within the focus beam are amplified in the best linear ways when the sound gets too loud. It keeps the signal-to-noise ratio positive to help users focusing on the signal with most contrast.

W

WIND NOISE REDUCTION

Smoothens wind noise fluctuations automatically to deliver a comfortable and pleasant listening experience even in windy environments.

WIRELESS SYNC

Wireless sync allows our hearing aids to connect to compatible devices without the need for cables.

X

XPHONE

Dedicated program that simultaneously sends phone signals to both hearing aids when the phone is held to one ear.



RELY ON REXTON

It's a tough world out there, and even tougher when your hearing starts to go. At Rexton, we won't let hearing loss hold you back. Good hearing is critical: From getting the job done, to getting home safely and being there for the people who count on you, we know what's at stake. That's why we work hard to deliver proven hearing technology in the most reliable way. We understand what reliability means in real life, and make practical, easy-to-use products you can count on. Our hearing aids have been getting the job done since 1955, so whatever life throws at you, you can rely on Rexton.

We believe that nothing hearing related should ever be left to chance, both for the people with hearing loss, and those who depend upon them. We understand their challenges and we meet them with proven hearing technology that is made Lifeproof so it can be counted upon to perform reliably in even the harshest environments. At work, at home, during leisure time, exercise and in all weather, you can always RELY ON REXTON.

**GETTING
THE JOB DONE
SINCE 1955**

[rexton.com](https://www.rexton.com)

The Bluetooth® word mark and logos are owned by the Bluetooth SIG, Inc. and any use of such marks by Sivantos GmbH is under license. Other trademarks and trade names are those of their respective owners. Apple and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc. Google Play is a trademark of Google Inc. "Made for iPhone" and "Made for iPad" mean that an electronic accessory has been designed to connect specifically to iPhone or iPad, respectively, and has been certified by the developer to meet Apple performance standards. Apple is not responsible for the operation of this device or its compliance with safety and regulatory standards. Please note that the use of this accessory with iPhone or iPad may affect wireless performance. iPad, iPhone are trademarks of Apple Inc., registered in the U.S. and other countries.

Order No. 04538-99C1-7600 © 2024 Sivantos GmbH