




## NOVIOMAGUS FEMORAL HEAD REAMER SET - SPLITTING

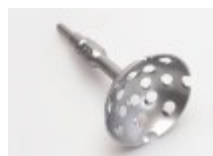

The Noviomagus Femoral Head Reamer Set has been developed to fixate bone parts, remove cartilage and split the femoral head in a safe, easy and efficient way. Literature suggests that cartilage and sclerotic bone material osseointegrate and thus removal is necessary. The inverted reamers are designed to fit standard power tools. They can easily be assembled and disassembled which is ideal for cleaning and sterilization and it has a low centre of gravity, making it a stable device. When using the inverted reamers in combination with a power tool, the splash shield prevents the removal of tissue. The slot in the bone vice makes it easy to cut the bone or split a femoral head in two halves, so the bone can be instantaneously milled in the Noviomagus Bone Mill.

- Strong and solid fixation of bone material
- Tissue removal in less than 30 seconds
- Less risk in personal harm or glove tears
- Long lifespan without wear
- Preparation for high quality bone chips



00:00

ADDITIONAL INFORMATION

PRODUCT DETAILS

The Noviomagus Femoral Head Reamer Set is developed and manufactured by Spierings Orthopaedics E the Netherlands. The Noviomagus Femoral Head Reamer Set consists of the following parts:

Sterilisation Tray (NFHRS)SBM-1109-45

STANDARD OPERATING KIT NOVIOMAGUS FEMORAL HEAD REAMER SET

Article	Article Number
Base Plate	SBM-1109-14
Detachable Clamp	SBM-0503-15
Screw Handle	SBM-0503-16
Splash Shield	SBM-1109-17
At least one Inverted Reamer selected below	SBM-XXXX-XX

The Inverted Reamers for the Noviomagus Femoral Head Reamer Set are available with a standard AO-C as the Zimmer-Stryker-Jacobs Hall-Coupling. The sterilization tray can hold up to 10 inverted reamers.

PRODUCTS	DEVELOPMENT	COMPANY	BLOG/NEWS	CONTACT
Inverted	Reamer AO ø42	SBM-0503-18	Inverted	Reamer Z-S-J Hall ø42
Inverted	Reamer AO ø44	SBM-0904-21	Inverted	Reamer Z-S-J Hall ø44
Inverted	Reamer AO ø46	SBM-0904-22	Inverted	Reamer Z-S-J Hall ø46
Inverted	Reamer AO ø48	SBM-0503-19	Inverted	Reamer Z-S-J Hall ø48
Inverted	Reamer AO ø50	SBM-0904-23	Inverted	Reamer Z-S-J Hall ø50
Inverted	Reamer AO ø52	SBM-0904-24	Inverted	Reamer Z-S-J Hall ø52
Inverted	Reamer AO ø54	SBM-0503-20	Inverted	Reamer Z-S-J Hall ø54
Inverted	Reamer AO ø56	SBM-0904-25	Inverted	Reamer Z-S-J Hall ø56
Inverted	Reamer AO ø58	SBM-0904-26	Inverted	Reamer Z-S-J Hall ø58

Our account managers and distributors are available to advise the best solution and reamer selection based on your specific needs. Feel free to reach out to us using the [contact form](#). In case you are looking for [orthopedic instruments](#) for surgical use, please visit [www.spierings.biz](http://www.spierings.biz)

FAQ