

Protective Medical Coverall TYPE 5-B/6-B

Product description:

Disposable medical laminate protective coverall, made of laminated spunbond, density 60 gr / m². It is provided with a hood; zipper with protective flap and adhesive tape; and protective shoes - an integral part of the coverall; elastic at the hood, sleeves, ankles and waist. Provides an effective barrier against fluids (waterproof).

Instructions for use:

The equipment is worn over the work coat. The zipper opens, and the coverall is dressed starting with the feet, then the hands are inserted into the sleeves of the coverall, the elastic ends are arranged correctly, the hood is positioned on the head and, at the end, the zipper is closed and the protective flap is fixed with adhesive tape.

Undressing instructions:

The coverall will be undressed in the designated area, away from pollution and dangers; the help of another person is recommended, avoiding contamination of skin and clothing. Carefully remove the self-adhesive flap for the zipper, unzip the coverall and remove the hood making sure the gloves do not touch the inside of the coverall; then the coverall is totally removed. The gloves should be removed according to the manufacturer's instructions and with safety procedures. Dispose of the coverall in a suitable waste container.

Storage instructions:

Storage is done in the original packaging. Store at a temperature between 0 ° C and + 20 ° C; in a dry and clean place, away from direct sunlight, any flame, water droplets and solvent vapors.



AGAINST
INFECTIVE AGENTS
EN 14126



ANTI-
STATIC
EN 1149-5



DRY
PARTICLES
EN ISO 13928-1

TURKEY

ISO

CE

2163

Made in Turkey

REF

GNC-2213




GENC TAHHÜT GİYİM TEKSTİL
İNŞAAT SAN. VE TİC. LTD. ŞTİ.

Address: General At İsmi güran cadde'si merkezi 9
numaralı blok 2/21 35098 merkezi İstanbul Turkey
Tel: +90212627150
e-mail: info@gencgiyim.com.tr


Authorized representative:

"ARTBUTIC" LLC
Uzinelor 8/1 str. Chişinău MD-
2023 Republic of Moldova,
E-mail: vlvartoffice@gmail.com
Tel: +37369500084
+37369500057

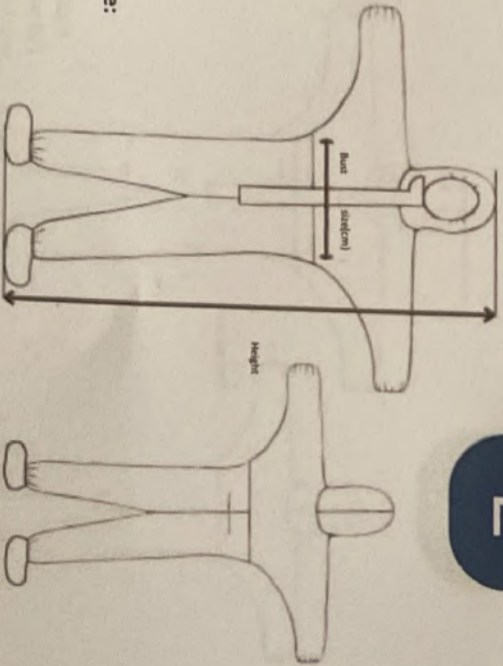

Registration number: DM000305279

Size	Height (cm)	Bust size (cm)
L	172-176	130-136
XL	176-182	136-140
2XL	182-188	140-144
3XL	188-194	144-148

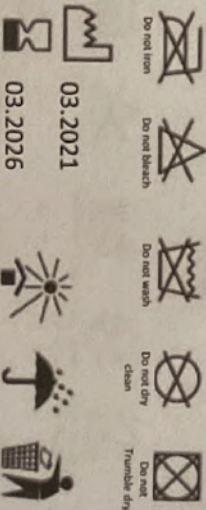

03.2021
03.2026

GENC TAHHÜT GİYİM TEKSTİL
İNŞAAT SAN. VE TİC. LTD. ŞTİ.
Merkez: Refahiye Mh.28 Sk. Solim Apt. Zemin Kat No:9/2
Tel: 0414 215 52 52 - Faks: 0414 215 52 52
Şube: General At İsmi güran Cad. No: 9/1 35098 merkezi İstanbul
Tel: 0212 482 71 50 - 51 - Faks: 0212 482 71 53
Tic. Sic. No: 273183
Mail: muhasebe@gencgiyim.com.tr - genclik@gencgiyim.com.tr - info@gencgiyim.com.tr
Mersis No: 0353050322500011-Ş. Urfa Tic. Odası: 11430

Combinezon Medical de Protecție TIP 5-B/6-B

Descrierea produsului:

Combinezon de protecție laminat medical, de unică folosință, fabricat din spunbond laminat, densitatea 60gr/m2. Este prevăzut cu glugă; fermoar cu clapetă protectoare și bandă adezivă; și botoșei - parte integrală a combinezonului; elastic la glugă, mâneci, glezne și în talie. Oferă o barieră eficientă împotriva fluidelor (impermeabil).

Instrucțiuni de utilizare:

Echipamentul se îmbracă peste haina de lucru. Se deschide fermoarul, se îmbracă începând cu picioarele, apoi se introduc mâinile în mânecile combinezonului și se aranjează corect terminațiile elastice, se poziționează gluga pe cap și, la final, se închide fermoarul și se fixează clapeta protectoare cu bandă adezivă.

Instrucțiuni de dezechipare:

Dezechiparea combinezonului se va efectua în zona desemnată, departe de poluare și pericole; se recomandă ajutorul unei alte persoane, evitând contaminarea pielii și a hainelor. Se îndepărtează cu grijă clapeta autoadezivă a fermoarului, se desface fermoarul și se scoate gluga costumului, asigurându-vă că mânușile nu ating interiorul costumului; ulterior, se scoate costumul definitiv. Se vor îndepărta mânușile în conformitate cu instrucțiunile producătorului și cu procedurile de siguranță. Se va arunca combinezonul într-un recipient adecvat pentru deșuri.

Instrucțiuni de depozitare:

Depozitarea se realizează în ambalajul original. A se păstra la temperatura cuprinsă între 0 °C și +20 °C; într-un loc uscat și curat, ferit de umiditate; departe de lumina directă a soarelui, de orice flacără, picături de apă și vapori de solvenți.



AGAINST
INFECTIVE AGENTS
EN 14126

ANTI-
STATIC
EN 1149-5

DRY
PARTICLES
EN ISO 15925-1

TURKEY

ISO

CE

2163

REF

GNC-2213




GENC TAHHÜT GİYİM TEKSTİL
İNŞAAT SAN. VE TİC. LTD. ŞTİ.

Adres: General Ali Rıza güran caddeşi merkezi 15
mahallesi no: 2/21 3609 merkezi istanbul Turkey
Tel: +902124627159
E-mail: info@gencgiyim.com.tr

Reprezentant autorizat:


"Artbute" SRL
Republica Moldova, MD-2023
mun. Chișinău, str. Uzinelor 8/1
E-mail: viviatoffice@gmail.com
Tel.: +37369500084
+37369500057

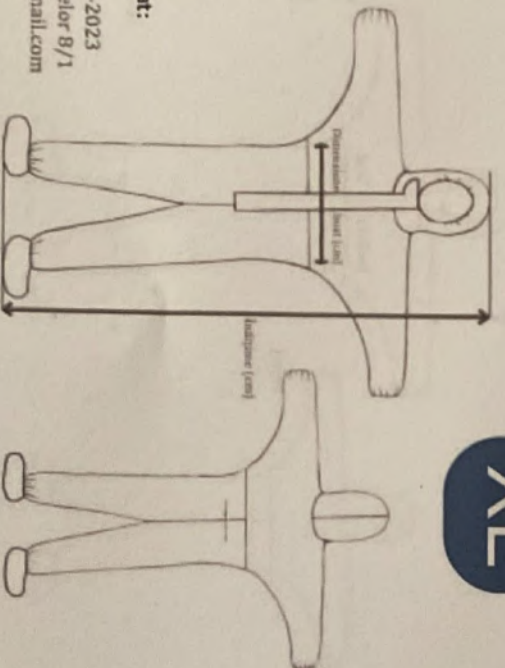

XL

Număr de înregistrare: DM000305280

Mărime	Înălțimea (cm)	Dimensiune bust (cm)
L	172-176	130-136
XL	176-182	136-140
2XL	182-188	140-144
3XL	188-194	144-148

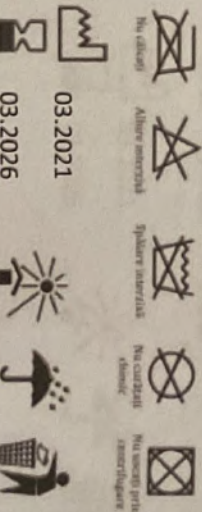

03.2021
03.2026

GENC TAHHÜT GİYİM TEKSTİL
İNŞAAT SAN. VE TİC. LTD. ŞTİ.
Merkez: Rafahya Mh.28 Sk. Salm Apt. Zemin Kat No:0/2
Tel: 0414 215 52 25 - Faks: 0414 217 16 35 Halkiye S. URFA
Şube: General Ali Rıza Güran Cad. Merkez İş Merkezi No:2/21
Tel: 0212 462 11 50 - 51 Faks: 0212 462 71 53 Merkez İSTANBUL
Tic. Sic. No: 27 33 856 3225
Mail: muhasebe@gencgiyim.com.tr info@gencgiyim.com.tr www.gencgiyim.com.tr
Mersis No: 0330058322500011 Ş. Urfa Tic. Odası: 11430

Protective Medical Coverall TYPE 5-B/6-B

Product description:

Disposable medical laminate protective coverall, made of laminated spunbond, density 60 gr / m². It is provided with a hood, zipper with protective flap and adhesive tape; and protective shoes - an integral part of the coverall; elastic at the hood, sleeves, ankles and waist. Provides an effective barrier against fluids (waterproof).

Instructions for use:

The equipment is worn over the work coat. The zipper opens, and the coverall is dressed starting with the feet, then the hands are inserted into the sleeves of the coverall, the elastic ends are arranged correctly, the hood is positioned on the head and, at the end, the zipper is closed and the protective flap is fixed with adhesive tape.

Undressing instructions:

The coverall will be undressed in the designated area, away from pollution and dangers; the help of another person is recommended, avoiding contamination of skin and clothing. Carefully remove the self-adhesive flap for the zipper, unzip the coverall and remove the hood making sure the gloves do not touch the inside of the coverall; then the coverall is totally removed. The gloves should be removed according to the manufacturer's instructions and with safety procedures. Dispose of the coverall in a suitable waste container.

Storage instructions:

Storage is done in the original packaging. Store at a temperature between 0 ° C and + 20 ° C; in a dry and clean place, away from moisture; away from direct sunlight, any flame, water droplets and solvent vapors.



AGAINST
INFECTIVE AGENTS
EN 14129



ANTI-
STATIC
EN 11485



DRY
PARTICLES
EN ISO 13982-1

Turkey

Disposal
Procedures

ISO
2163

CE
2163

Made in Turkey

REF

GNC-2213

GENC Giyim

GENC TAHRÜT GİYİM TEKSTİL
İNŞAAT SAN. VE TİC. LTD. ŞTİ.

Address: General Ali Rıza paşa caddeşi merkezi
mühürli no: 2/21 34098 merkezi İstanbul Turkey
Tel: +902124827150
e-mail: info@gencgiyim.com.tr

Authorized representative:

"ARTBUTIC" LLC


Uzineleor 8/1 str. Chişinău MD.

2023 Republic of Moldova,

E-mail: viviaroffice@gmail.com

Tel.: +37369500084


+37369500057

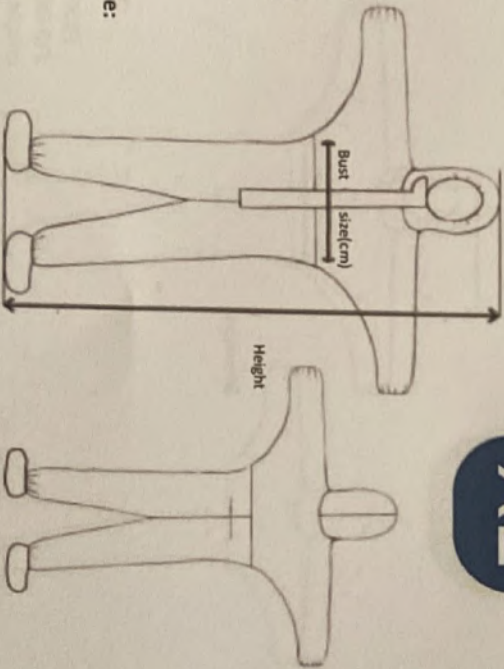

XL

Registration number: DM000305280

Size	Height (cm)	Bust size (cm)
L	172-176	130-136
XL	176-182	136-140
2XL	182-188	140-144
3XL	188-194	144-148

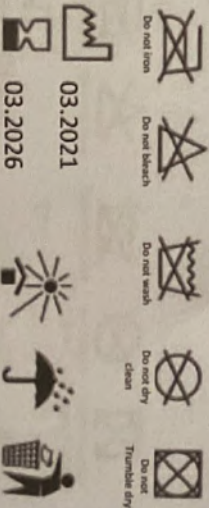

03.2021

03.2026

GENC TAHRÜT GİYİM TEKSTİL
İNŞAAT SAN. VE TİC. LTD. ŞTİ.
Merkez: Refahiye Mh 881. Sk. Solim Apt. Zemin Kat No:9/2
Tel: 0414 212 48 22 Fax: 0414 212 16 36 Halkıyo/Ş.URFA
Şube: General Ali Rıza paşa caddeşi merkezi
Tel: 0212 482 71 50 Fax: 0212 482 71 53 Merkez/İSTANBUL
Mühürli no: 2/21 34098
Mail: muhurat@gencgiyim.com.tr - info@gencgiyim.com.tr
Mersis No: 036305832500011-Ş. Urfa Tic. Odası: 11430

Combinezon Medical de Protecție TIP 5-B/6-B

Descrierea produsului:

Combinezon de protecție laminat medical, de unică folosință, fabricat din spunbond laminat, densitatea 60gr/m2. Este prevăzut cu glugă, fermoar cu clapetă protectoare și bandă adezivă; și botoșei - parte integrată a combinezonului; elastic la glugă, mâneci, glezne și în talie. Oferă o barieră eficientă împotriva fluidelor (impermeabil).

Instrucțiuni de utilizare:

Echipamentul se îmbracă peste haina de lucru. Se deschide fermoarul, se îmbracă începând cu picioarele, apoi se introduc mâinile în mânecile combinezonului și se aranjează corect terminațiile elastice, se poziționează gluga pe cap și, la final, se închide fermoarul și se fixează clapeta protectoare cu bandă adezivă.

Instrucțiuni de dezechipare:

Dezechiparea combinezonului se va efectua în zona desemnată, departe de poluare și pericole; se recomandă ajutorul unei alte persoane, evitând contaminarea pielii și a hainelor. Se îndepărtează cu grijă clapeta autoadezivă a fermoarului, se desface fermoarul și se scoate gluga costumului, asigurându-vă că mânușile nu ating interiorul costumului; ulterior, se scoate costumul definitiv. Se vor îndepărta mânușile în conformitate cu instrucțiunile producătorului și cu procedurile de siguranță. Se va arunca combinezonul într-un recipient adecvat pentru deșuri.

Instrucțiuni de depozitare:

Depozitarea se realizează în ambalajul original. A se păstra la temperatura cuprinsă între 0 °C și +20 °C, într-un loc uscat și curat, ferit de umiditate; departe de lumina directă a soarelui, de orice flacără, picături de apă și vapori de solvenți.



AGAINST
INFECTIVE AGENTS
EN 14126

ANTI-
STATIC
EN 1149-5

DRY
PARTICLES
EN ISO 15982-1

Turkey

ISO

CE

2163

Made in Turkey

REF


GNC-2213



GENC TAHHÜT GİYİM TEKSTİL
İNŞAAT SAN. VE TİC. LTD. ŞTİ.

Adres: General At Rıza Hasan Cad.İ. merkezi 11
merkezi no: 2/21 3409 merkezi İstanbul Turkey
Tel: +902124827150
E-mail: info@gencgiyim.com.tr


Reprezentant autorizat:
"Arbutic" SRL
Republica Moldova, MD-2023
mun. Chișinău, str. Uzinelor 8/1
E-mail: viviaroffice@gmail.com
Tel.: +37369500084
+37369500057

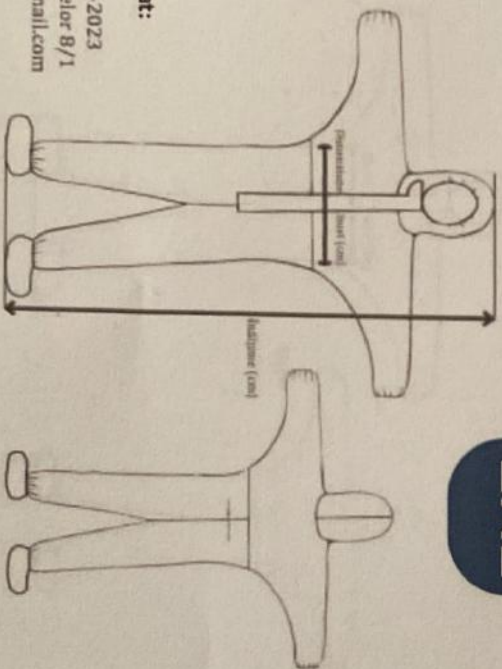

2XL

Număr de înregistrare: DM000305281

Mărime	Înălțimea (cm)	Dimensiune bust (cm)
L	172-176	130-136
XL	176-182	136-140
2XL	182-188	140-144
3XL	188-194	144-148

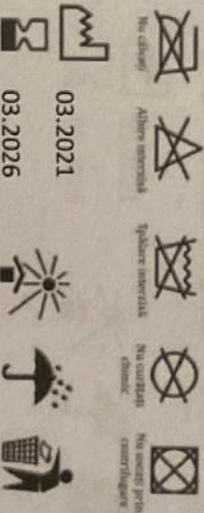

03.2021
03.2026

GENC TAHHÜT GİYİM TEKSTİL
İNŞAAT SAN. VE TİC. LTD. ŞTİ.
Merkez: Râzîyîyî, No: 2/21, Sk. Şahin, Apt. Emin Kat No:9/2
Tel: 0414 2 5 12 22 Fax: 0414 2 5 12 23
Sube: General At Rıza Hasan Cad.İ. merkezi 11
Tel: 0212 482 7150 - 51 Fax: 0212 482 7153 Merkez/İSTANBUL
Nöbetçi: 0212 482 7150 - 51
Mail: muhasibe@gencgiyim.com.tr - genc.giyim@hotmail.com - www.gencgiyim.com.tr
Mersis No: 383058322500011-Ş. Ürfâ Tic. Odası: 11430

Protective Medical Coverall TYPE 5-B/6-B

Product description:

Disposable medical laminate protective coverall, made of laminated spunbond, density 60 gr / m². It is provided with a hood; zipper with protective flap and adhesive tape; and protective shoes - an integral part of the coverall; elastic at the hood, sleeves, ankles and waist. Provides an effective barrier against fluids (waterproof).

Instructions for use:

The equipment is worn over the work coat. The zipper opens, and the coverall is dressed starting with the feet, then the hands are inserted into the sleeves of the coverall, the elastic ends are arranged correctly, the hood is positioned on the head and, at the end, the zipper is closed and the protective flap is fixed with adhesive tape.

Undressing instructions:

The coverall will be undressed in the designated area, away from pollution and dangers; the help of another person is recommended, avoiding contamination of skin and clothing. Carefully remove the self-adhesive flap for the zipper, unzip the coverall and remove the hood making sure the gloves do not touch the inside of the coverall; then the coverall is totally removed. The gloves should be removed according to the manufacturer's instructions and with safety procedures. Dispose of the coverall in a suitable waste container.

Storage instructions:

Storage is done in the original packaging. Store at a temperature between 0 ° C and + 20 ° C; in a dry and clean place, away from moisture; away from direct sunlight, any flame, water droplets and solvent vapors.



AGAINST
INFECTIVE AGENTS
EN 14126



ANTI-
STATIC
EN 1149-5



DRY
PARTICLES
EN ISO 13882-1

TURKEY

ISO

2163

CE

Made in Turkey


REF

GNC-2213

GENC Giyim

Address: General Ali Rıza Güran cadde'si merkezi 15
mimarlar no: 1/1 34098 merkezi (Halkalı) Turkey
Tel: +902124627150
e-mail: info@gencgiyim.com.tr


Authorized representative:
"ARTBUTIC" LLC
Uzineleor 8/1 str. Chişinău MD-
2023 Republic of Moldova,
E-mail: vivartoffice@gmail.com
Tel.: +37369500084
+37369500057

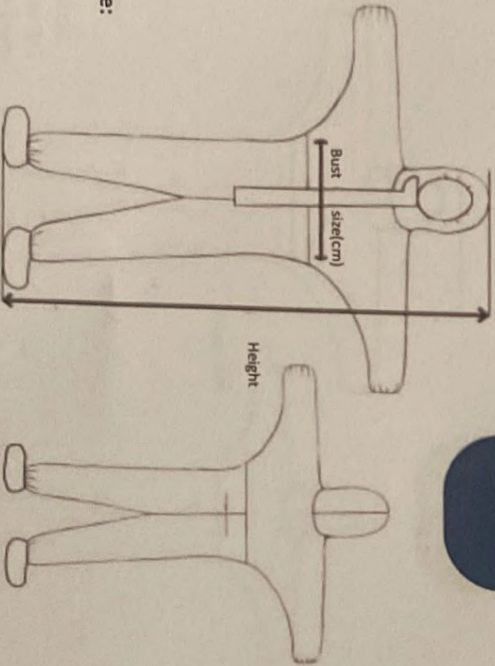

2XL

Registration number: DM000305281

Size	Height (cm)	Bust size (cm)
L	172-176	130-136
XL	176-182	136-140
2XL	182-188	140-144
3XL	188-194	144-148

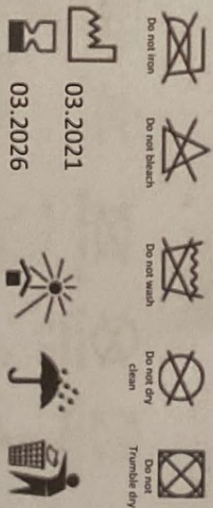

03.2021

03.2026

GENC TAHHÜT GİYİM TEKSTİL
İNŞAAT SAN. VE TİC. LTD. ŞTİ.
Merkez: Refahiye Mh. 281. Sk. Selim Apt. Zemin Kat No:9/2
Tel: 0414 218 52 22 - Fax: 0414 218 52 36 Halkılıy/Ş.ÜRFA
Şube: General Ali Rıza Güran Cad. Merkez İş Merkezi No:2/21
Tel: 0212 463 71 10 - Ş.Ü. Fax: 012 463 71 63 Merkezi/İSTANBUL
Tic. Sic. No: 306 093 1225
Müş. Muhasebe/Giyim San. ve Tic. Ltd. Şti. www.gencgiyim.com.tr
Mersis No: 0803050322000011-Ş.Ürfa Tic. Odası: 11430

CamScanner ile tarandı

Protective Medical Coverall TYPE 5-B/6-B

Product description:

Disposable medical laminate protective coverall, made of laminated spunbond, density 60 gr / m². It is provided with a hood, zipper with protective flap and adhesive tape; and protective shoes - an integral part of the coverall; elastic at the hood, sleeves, ankles and waist. Provides an effective barrier against fluids (waterproof).

Instructions for use:

The equipment is worn over the work coat. The zipper opens, and the coverall is dressed starting with the feet, then the hands are inserted into the sleeves of the coverall, the elastic ends are arranged correctly, the hood is positioned on the head and, at the end, the zipper is closed and the protective flap is fixed with adhesive tape.

Undressing instructions:

The coverall will be undressed in the designated area, away from pollution and dangers; the help of another person is recommended, avoiding contamination of skin and clothing. Carefully remove the self-adhesive flap for the zipper, unzip the coverall and remove the hood making sure the gloves do not touch the inside of the coverall; then the coverall is totally removed. The gloves should be removed according to the manufacturer's instructions and with safety procedures. Dispose of the coverall in a suitable waste container.

Storage instructions:

Storage is done in the original packaging. Store at a temperature between 0 ° C and + 20 ° C; in a dry and clean place, away from moisture; away from direct sunlight, any flame, water droplets and solvent vapors.




AGAINST
INFECTIVE AGENTS
EN 14126



ANTI-
STATIC
EN 1149-5



DRY
PARTICLES
EN ISO 15828-1

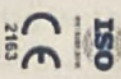

Made in Turkey

REF

GNC-2213

GENC Giyim

GENC TAHHÜT GİYİM TEKSTİL

İNŞAAT SAN. VE TİC. LTD. ŞTİ.

Address: General Ali Rıza gıyım caddeisi mermer 5

mekan no: 2/21 34098 mermer ilçesi İstanbul Turkey

Tel: +902124627150

e-mail: info@gencgiyim.com.tr

Authorized representative:

"ARTBUTIC" LLC


Uzinelor 8/1 str. Chişinău MD-

2023 Republic of Moldova,

E-mail: viviantoffice@gmail.com

Tel: +37369500084


+37369500057

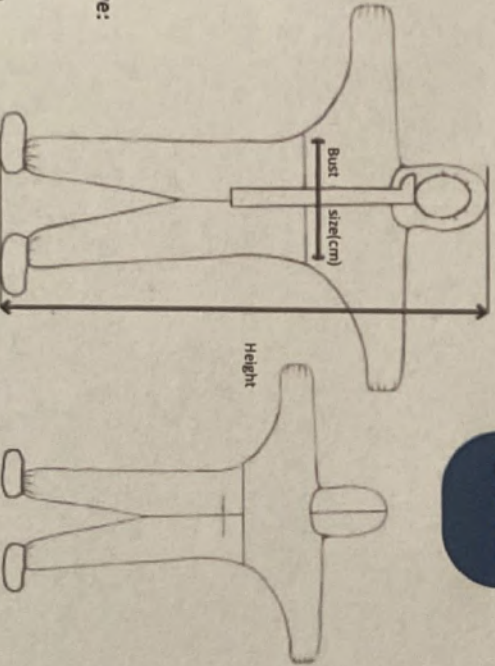


3XL

Registration number: DM000305282


Size	Height (cm)	Bust size (cm)
L	172-176	130-136
XL	176-182	136-140
2XL	182-188	140-144
3XL	188-194	144-148

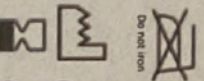


Do not iron

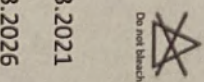


Do not bleach

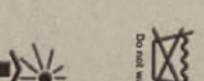


Do not wash

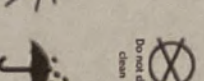


Do not dry

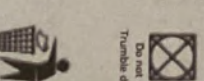

Do not tumble dry

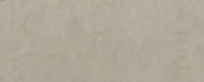

Do not clean

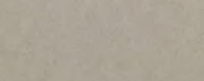

Do not

GENC TAHHÜT GİYİM TEKSTİL
İNŞAAT SAN. VE TİC. LTD. ŞTİ.
Merkez: Refahiyah 2/1 Sk. Solim Apt. Zemin Kat No:9/2
Tel: 0414 216 5272 - Fax: 0414 217 16 36 Halilive@URFA
Şube: General Ali Rıza gıyım cad. mermer 5 Mecidiyeköy/İSTANBUL
Tel: 0212 462 7150 - Fax: 0212 462 7153 Mermer@ISTANBUL
E-posta: info@gencgiyim.com.tr - 0303 103 3225
Mail: muhasebe@gencgiyim.com.tr - gencgiyim@hotmail.com - www.gencgiyim.com.tr
Mersis No: 0633056322000011-Ş. Urfa Tic. Odası: 11430