


Home » Supplies

Supplies

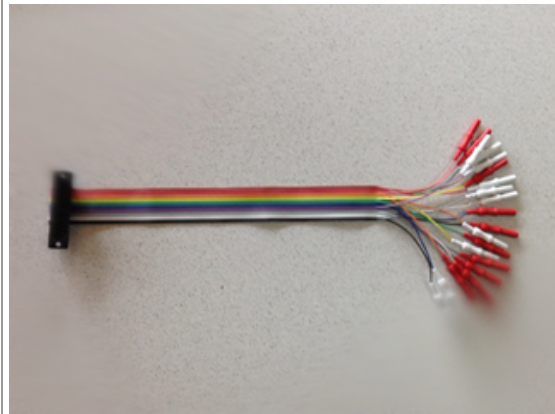

E2S ELECTRODE BOARD ADAPTER*

The Electrode Board Adapter converts the EEG instrument from individual electrodes to the Electro-Cap. Each wire is color coded for easy identification. Once the EBA is plugged into the machine, a single plug connects the appropriate cap. Special and custom EBA's are available.



E3-L BODY HARNESS - 44" - 52" E3-M BODY HARNESS - 30" - 44" E3-S BODY HARNESS - 21" - 30"

An adjustable elastic band that fastens around the patient's chest with Velcro. Cap straps attach the cap to the body harness. Proper use of the body harness will hold the cap securely in place reducing the possibility of movement artifact and high impedances.



E4S QUICK INSERT ELECTRODE*

A temporary electrode which will override the cap wiring in the event of a faulty electrode during a test. A Quick Insert Electrode may be placed into the mount of the faulty electrode and connected directly to the EEG machine until the cap can be properly repaired.