

EU Declaration of Conformity

Unit type Minicentrifuges-vortexes

Models FV-2400, FVL-2400N, MSC-3000, MSC-6000, CVP-2

Serial number 14 digits styled XXXXXYYMMZZZZ, where XXXXXX is model code, YY and MM – year and month of production, ZZZZ – unit number.

Manufacturer SIA BIOSAN
Latvia, LV-1067, Riga, Ratsupites str. 7/2

The objects of the declaration described above is in conformity with the following relevant Union harmonization legislations:

LVD 2014/35/EU	LVS EN 61010-1:2011 Safety requirements for electrical equipment for measurement, control, and laboratory use. General requirements. LVS EN 61010-2-020:2016 Particular requirements for laboratory centrifuges.
EMC 2014/30/EU	LVS EN 61326-1:2013 Electrical equipment for measurement, control and laboratory use. EMC requirements. General requirements.
RoHS3 2015/863/EU	Directive on the restriction of the use of certain hazardous substances in electrical and electronic equipment.
WEEE 2012/19/EU	Directive on waste electrical and electronic equipment.

I declare that the Declaration of Conformity is issued under sole responsibility of the manufacturer and belongs to the above-mentioned objects of the declaration.

Svetlana Bankovska
Managing director



Signature

07.02.2020.
Date

MSC-6000, Centrifuge/Vortex Multispin

DESCRIPTION

Centrifuge/vortex Multi-Spin **MSC-6000** is product of many years evolution of Spin-Mix-Spin technology that is intended for collecting micro volumes of reagents on the microtube's bottom (first centrifugation spin), following mixing (mix) and collecting the reagents again from the walls and cap of the microtube (second spin). We named this repetitive algorithm of operation that is aimed at reducing the mistakes during sample preparation for PCR analysis a "sms-algorithm".

Multi-Spin is a fully automatic device for reproducing sms-algorithm for 12 tubes at one time, thus saving time considerably. A must-have instrument for PCR and DNA analyses laboratory.

Multi Spin is four devices combined in one:

1. Centrifuge — Maximum RCF: up to 2350 × g
2. Vortex (3 mixing modes — soft, medium, hard; regulated time; Vortexing regulation timer 1–20 s)
3. Centrifuge/Vortex;
4. SMS-cycler for realization of the "sms-algorithm".

Table on the right compares FVL-2400N, MSC-3000 and **MSC-6000**

Multi-Spin allows considerable time saving compared to Combi-Spin by automatically performing cycling program of sample mixing and spinning according to the set spin-mix-spin cycle for 12 microtubes simultaneously.

SPECIFICATIONS

Speed control range	1000–6000 rpm (increment 100 rpm)
RCF max.	2350 × g
Spin timer	1 s - 30 min
Vortexing intensity	Soft, medium, hard
Vortexing time	0–20 s (increment 1 s)
SMS-cycle regulation	1 - 999 cycles
Timer sound signal	+
Display	LCD, 2 x 16 signs
Safety	Lid lock
Overall dimensions (W×D×H)	190x235x125 mm
Weight	2.5 kg
Input current/power consumption	24 V, 24 W (1 A)
External power supply	Input AC 100–240 V; 50/60 Hz; Output DC 24 V



CAT. NUMBER

Including rotors R-1.5, R-0.5/0.2	Including rotors R-1.5, R-0.5/0.2
BS-010211-AAL	230VAC 50/60Hz Euro plug
BS-010211-AAQ	230VAC 50/60Hz UK plug
BS-010211-AA4	230VAC 50/60Hz AU plug
BS-010211-AAU	100VAC 50/60Hz US plug
BS-010211-AAM	120VAC 60Hz US plug
BS-010211-CK	IQ OQ document
BS-010211-DK	PQ document



R-0.5/0.2
BS-010205-BK
rotor

Rotor for 12 x 0.5 ml and 12 x 0.2 ml microtest tubes



R-1.5
BS-010205-AK
rotor

Rotor for 12 x 1.5/2 ml microtest tubes




R-2/0.5
BS-010205-CK
rotor

Rotor for 8 x 2/1.5 ml and 8 x 0.5 ml microtest tubes



R-2/0.5/0.2
BS-010205-DK
rotor

Rotor for 6x2/1.5 ml + 6x0.5 ml + 6x0.2 ml microtest tubes

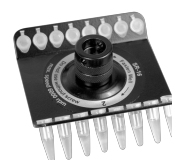

SR-16
BS-010202-AK
rotor

Rotor for 2 x 8-section 0.2 ml microtube strips



SR-32
BS-010205-FK
rotor

Rotor for 4 x 8-section 0.2 ml microtube strips

* not compatible with Combi-Spins produced before 2015

Anexa 6. Mini-centrifugă pentru eprubete 0.2 - 1.5 ml

Biosan MSC-6000

Parametri solicitati	Parametri oferiti
Anul de productie: 2021 Viteza reglabila 2000-6000RPM, pas 100 RPM Timer 10 sec-60min Rotor inclus pentru microtuburi 12x1.5ml, 12x0.5ml, 12x0.2ml Afisare digitală Sistem de blocare a capacului in timpul functionarii. Protecție la pornire cu capacul deschis. Capacul sa fie transparent sau dotat cu fereastră sau alta cale de acces, necesar pentru efectuarea procedurii de verificare periodica a turatiilor, conform normelor si standardelor în vigoare Tensiune alimentare: 210-230 V/ 50 Hz Cerințe de certificare: Certificat CE sau declarație de conformitate CE cu anexele corespunzătoare pentru produsele oferite, valabil, copie confirmată prin semnatura și ștampila participantului. Declarație de la Ofertant – confirmată prin semnatura și ștampila, în care să certifice termenul de garanție pentru echipament și accesorii nu mai mic de 36 luni din momentul instalării/darii în exploatare a bunului. Instalare, darea in exploatare de către participantul câștigător-obligatoriu Ghid rapid al utilizatorului (max. 4 pagini A4), in limba de stat-obligatoriu.	Anul de productie: 2021 Viteza reglabila 1000-6000RPM, pas 100 RPM Timer 1 sec-30min Rotor inclus pentru microtuburi 12x1.5ml, 12x0.5ml, 12x0.2ml Afisare digitală Sistem de blocare a capacului in timpul functionarii. Protecție la pornire cu capacul deschis. Capacul sa fie transparent sau dotat cu fereastră sau alta cale de acces, necesar pentru efectuarea procedurii de verificare periodica a turatiilor, conform normelor si standardelor în vigoare Tensiune alimentare: 210-230 V/ 50 Hz Cerințe de certificare: Certificat CE sau declarație de conformitate CE cu anexele corespunzătoare pentru produsele oferite, valabil, copie confirmată prin semnatura și ștampila participantului. Declarație de la Ofertant – confirmată prin semnatura și ștampila, în care să certifice termenul de garanție pentru echipament și accesorii nu mai mic de 36 luni din momentul instalării/darii în exploatare a bunului. Instalare, darea in exploatare de către participantul câștigător-obligatoriu Ghid rapid al utilizatorului (max. 4 pagini A4), in limba de stat-obligatoriu.