



## TECHNICAL DATA SHEET

### Disposable Sterile Reinforced Surgical Gown

1. Product name  
Disposable Sterile Reinforced Surgical Gown.

2. Product code  
Model 101-501.

3. Product description

Reinforced sterile surgical gown made of SMS (Spunbond + Meltblown + Spunbond) three-layer non-woven, according to EN 20811:1981 the resistance to liquid penetration is greater than 50 cm H<sub>2</sub>O. Reinforced with waterproof material on sleeves and front of gown to provide maximum protection against liquids, dust particles, alcohol and blood. With long sleeves and elasticated knitted cuffs, neck fastening system with „velcro” ( hedgehog) tape and back overlap fastening with ties. Product complies with: ISO 13485, EN 13795, EN 20811:1981, TS EN ISO 811(ISO 811:2018), TS 393 EN ISO 13938-1:2002.

Length of gown according to table. The gown is sterile, disposable and cannot be resterilized. It retains its sterility as long as the package is unopened and undamaged.


4. Material

Fibre composition - 100% polypropylene

Material density - 35 gr/m<sup>2</sup>

Sterilized - Ethylene Oxide.

5. Construction features

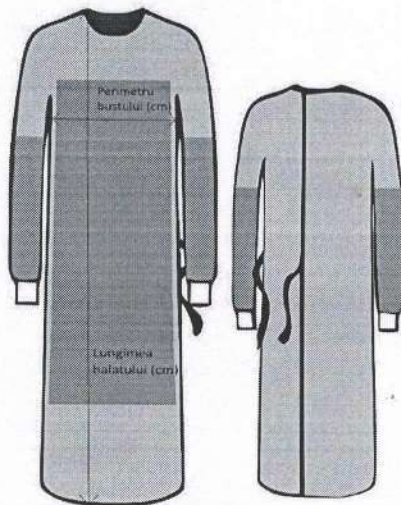

3TEKSTİ  
Gıda İlaç Sanayi Tic. Ltd. Şti  
Z. Organize Sanayi Bölgesi No:  
Sanayiye Giriş / Çıkış / T.C. (TÜRKİYE)  
Tel: +90 342 342 337 21 Fax: +90 342 337 21  
Sakarya V.D. 09\* 351 6111

Size	Model number	Gown length (cm)	Breast perimeter (cm)
S	101-501-01	117	116
M	101-501-02	125	124
L	101-501-03	132	132
XL	101-501-04	140	140
2XL	101-501-05	150	148

6. Instructions for use

The equipment is worn over the work coat. It is fastened at the back with straps and at the neck with „velcro” fasteners (hedgehog).

7. Packing method

The product is individually packed in a polythene bag with the size and design of the dressing gown.

8. Storage instructions

The storage is in the original packaging. Store in a dry and clean place away from moisture, from 0°C to +20°C, away from direct sunlight, water droplets and solvent vapors.

9. Warranty term

3 years from date of production.

**FACTORY CENTER**

2. Organized industrial Zone No:1 Başpınar/Gaziantep

Tel.: +90 (342) 337 24 16; +90 (342) 337 21 11

E-mail: info@3teks.com.tr

