
# Danis Procedure Pack

# **SX-ELLA Stent Danis**







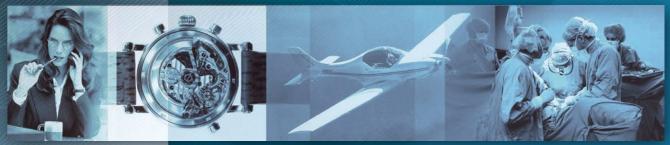

Removable stent system for treatment of bleeding esophageal varices Easy implantation without the need for endoscopy / X-Ray control Secure hemostasis through direct compression of varices Oral intake maintained from time of implantation



# **Procedure Pack with SX-ELLA Stent Danis**

Self-expandable nitinol stent preloaded in a ready-to-use delivery system. Complete procedure pack including all items required for implantation.

#### Indications

SX-ELLA Stent Danis is intended to stop acute and/or refractory bleeding from esophageal varices.

#### **Contraindications**

- Diagnosed malignancy or stricture affecting the esophagus
- ☑ Case history of irradiation therapy of the chest or the esophagus
- Malignancy of the throat or larynx
- Malignancy of bronchus or fistula
- ☑ Stomach malignancy
- ☑ Suspicion of bleeding due to foreign body injuring the upper gastrointestinal tract

# **Features / Benefits**

- ☑ Variable stent braiding is adjusted to esophageal peristalsis reduced migration rates.
- Standardized compression of the varices results in effective hemostasis.
- $\ensuremath{\boxtimes}$  Atraumatic stent ends.
- ☑ Radiopaque markers at both stent ends and at the midpoint for possible X-ray control of the stent position.
- Retrieval loops made of medical grade alloy at both stent ends high mechanical durability and acid resistance.
- ☑ Easy endoscopic stent removal within 7 days.

## **Benefits of the method**

- Readily implantable without endoscopy / X-Ray control, also in emergency situations.
- Delta Patient can continue an oral diet after stent implantation.

## **Additional description**

Danis Procedure Pack (019-08S-25-135) consists of a plastic carry case containing all items required for stent implantation:

☑ Delivery system preloaded with stent, guide wire, plastic syringes, mouthpiece, kidney dish, slobber cloth, waste bag with clip, gloves.

#### Danis Procedure Pack Basic (019-08S-25-135-B) consists of:

 $\blacksquare$  Delivery system preloaded with stent, guide wire, plastic syringe 50 ml.


The delivery system 9.4 mm (28 F) / 6.6 mm (20 F) is available in a standard length of 60 cm.




#### Available sizes

REF No.	SX-ELLA Stent Danis			Delivery system	
	Stent flares diameter [mm]	Stent body diameter [mm]	Nominal length [mm]	Active length [cm]	Outer diameter [F]
019-08S-25-135	30	25	135	60	28/20
019-08S-25-135-B	30	25	135	60	28/20



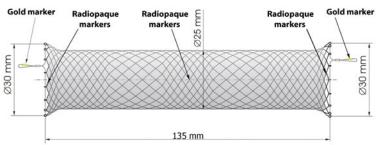

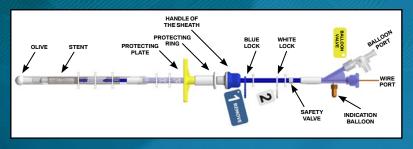

Delivery system - Gastric balloon deployed



Distributor:



Tel.: +420 495 279 111 E-mail: info@ellacs.eu URL: www.ellacs.eu

## **Danis Procedure Pack**