

# PCECG-500 ECG workstation

## PC operating environment :

### Hardware:

CPU : Pentium CPU 1G Hz or better

Memory: 1G or better

Communication port : USB or RS232

Display : 1024\*768 or better

### Software:

Operation system : Microsoft windows XP/ Vista/ win7

Licensing type: USB dongle software license.

## Modularized design (The relationship between four systems and modules)

Name of system	ECG Receiver	ECG Manager	ECG Vet Manager & Acquirer	ECG Manager& Acquirer
Type	PCECG-500(R)	PCECG-500(M)	PCECG-500(VM)	PCECG-500(A)
Brief introduction to the functions	Can receive ECG files from several ECGs and support a variety file format conversion, storage and output	Can receive ECG files from several ECGs, and manage files, such as register, edit, waveform review, analysis and measurement	Except for all the function of receiver and manager systems, it can also support real-time lead acquisition, display and analysis.	Except for all the function of receiver and manager systems, it can also support real-time lead acquisition, display and analysis.
Module composition	File receiver module	•	•	•
	Data management module for human		•	•
	Data management module for veterinary			•
	Acquirer module			•
Connected device	ECG (human / vet)	ECG (human)	ECG (Vet) & collecting box	ECG (human) & collecting box
Standard configuration	Software CD, RS232 cable, Dongle(R)	Software CD, RS232 cable, Dongle(M)	Software CD, RS232 cable, USB cable, Dongle (VA)	Software CD, RS232 cable, USB cable, Dongle (A)
Options	Network converter	Network converter	Network converter	Network converter



# PCECG-500 ECG workstation



## Main features

- Can be used in hospital, clinique, physical examination center, community health center and pet center
- Application of CSE and AHA database, reliable automatic measurement and analysis result
- One key operation to install USB drive software, database and PCECG-500 system
- 12-lead simultaneous acquisition and display
- High-precision electronic ruler, supporting adjustment and analysis of any measuring point.
- Up to 3600 seconds of ECG data can be stored, any 10 seconds can be picked out to analyze and print.
- Multiple, practical and user-definable ECG diagnosis dictionaries and modules
- Professional management database of 4G capacity, faster, safer and more stable
- 9 formats for option: XML, SCP, PDF, JPEG, DICOM3.0, HL7 aECG, WORD, BITMAP, MFER
- Multi-format printout of 12\*1, 6\*2+1 or 4\*3+1 on A4 paper vertically or horizontally.
- Quick search function: the doctor can set all kinds of searching terms, convenient for patient data recording and statistics.
- Can be connected to PC and hospital network

Shenzhen Carewell Electronics Co., Ltd.

Add: 4F, BLD 9, Baiwangxin High-Tech Industrial Park, Songbai RD, Xili Street, Nanshan District, Shenzhen, China

Tel: +86-755-86952278 Fax: +86-755-86952287

E-mail: info@carewell.com.cn Website: www.carewell.com.cn

※ Specification subject to change without prior notice



# PCECG-500

## ECG workstation




### ● Network management

#### Time-saving :


- Anytime, anywhere, fast collection for patient data
- One key installation

#### Labor-saving :

- Up to 9 kinds of output formats
- For patients: stay in the wards and take ECG examination
- For nurses: upload files by PCECG-500
- For doctors: to review, analyze and print out the report in the office
- Mass data storage capacity and long-time preservation

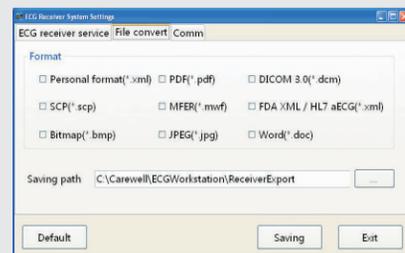

File convert setting

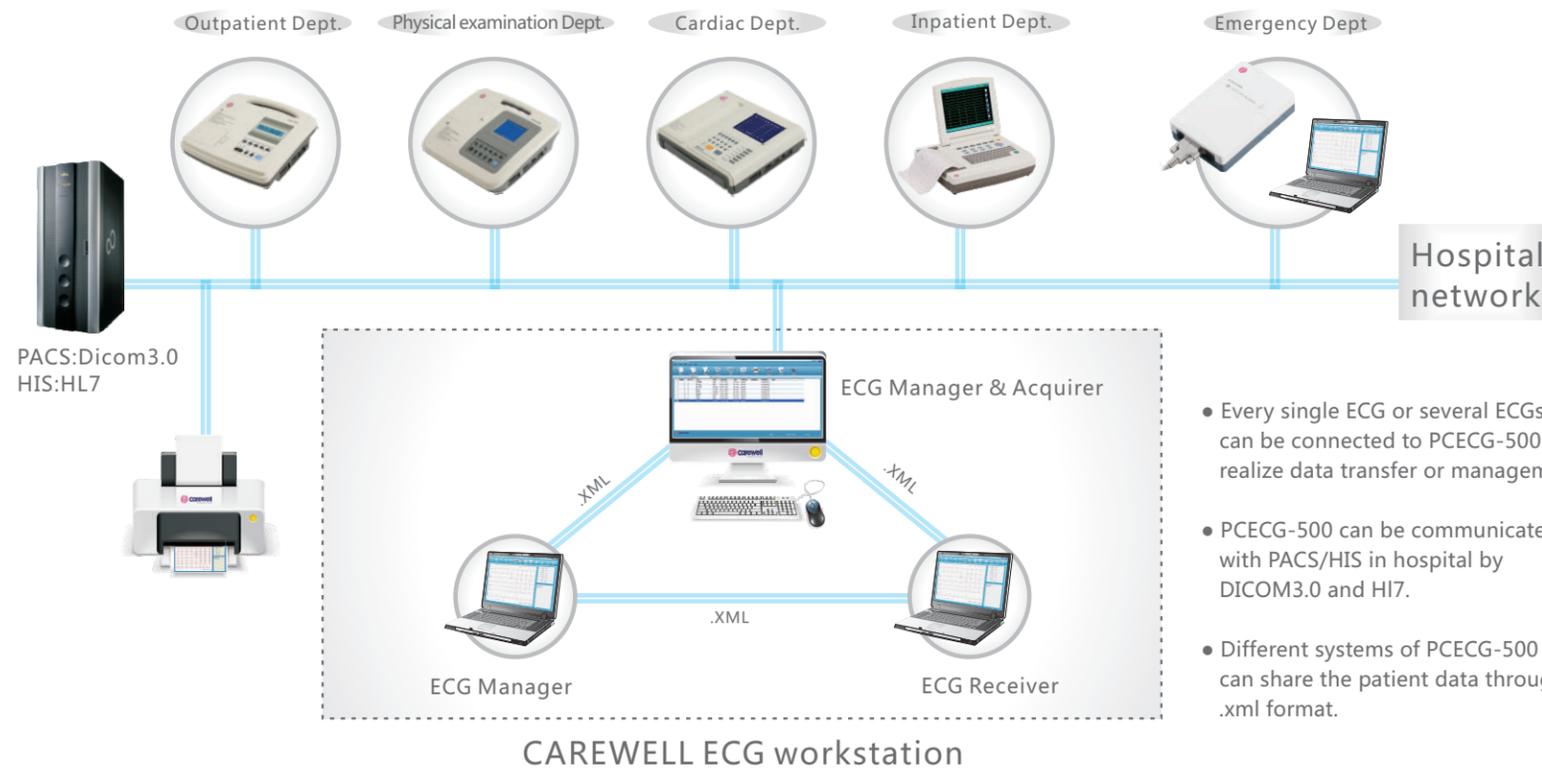


### ● Network management


#### Money-saving :

- A4 paper printout
- Fast and convenient second analysis on the network without printout
- Modulized design, four workstation systems for different ECG rooms
- Flexible upgrading solutions by changing a dongle


- Every single ECG or several ECGs can be connected to PCECG-500 to realize data transfer or management.
- PCECG-500 can be communicated with PACS/HIS in hospital by DICOM3.0 and HL7.
- Different systems of PCECG-500 can share the patient data through .xml format.

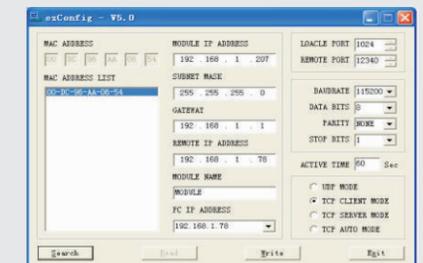


Network convert setting

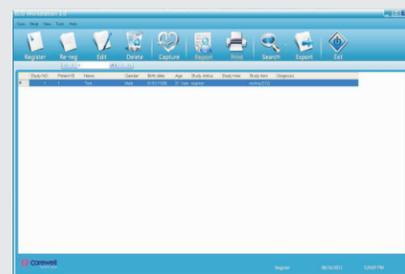


Main interface

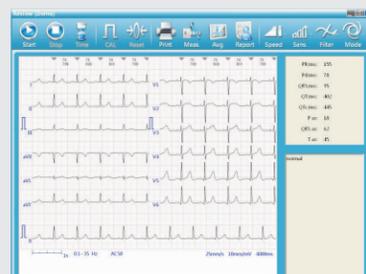


Acquisition interface

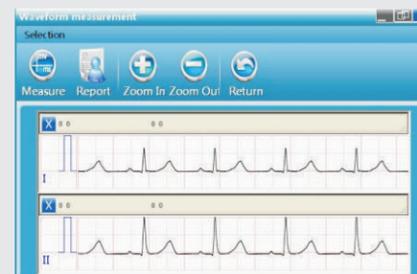


Measurement interface

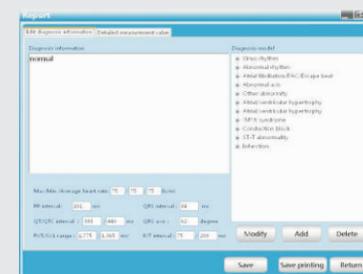

Analysis & Diagnosis interface

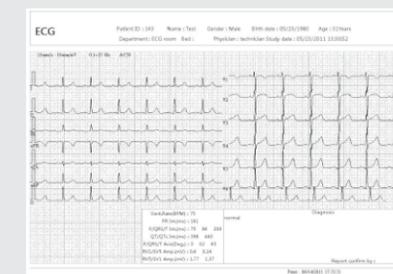

Report interface

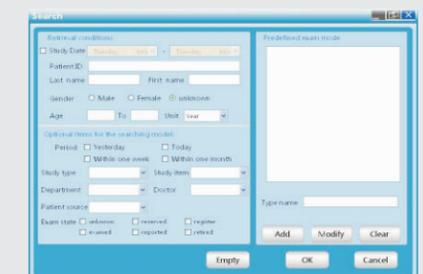

Search interface