



NeuroSlider®
Microcatheter DLC



Superior torqueability and pushability
Smooth and safe device delivery
Long-lasting tip shape retention


ORDERING INFORMATION | NeuroSlider® DLC

Product Name	Reference Number	ID (Inch)	OD dist. / prox. (French)	Usable Length (cm)	Tip Shape	Tip Marker
NeuroSlider® 17 DLC	01-000282	0.0165	1.9 / 2.3	155	Straight (shapeable)	2
	01-000283	0.0165	1.9 / 2.3	160	Straight (shapeable)	2
	01-000284	0.0165	1.9 / 2.3	167	Straight (shapeable)	2
NeuroSlider® 21 DLC	01-000292	0.021	2.2 / 2.6	155	Straight (shapeable)	2
	01-000293	0.021	2.2 / 2.6	160	Straight (shapeable)	2
	01-000294	0.021	2.2 / 2.6	167	Straight (shapeable)	2
NeuroSlider® 27 DLC	01-000276	0.027	3.0 / 3.1	155	Straight (shapeable)	1
NeuroSlider® 27 DLC pro*	01-000277	0.027	3.0 / 3.1	155	Straight (shapeable)	1
NeuroSlider® 39 DLC	01-000262	0.039	4.0 / 4.1	125	Straight (shapeable)	1
	01-000263	0.039	4.0 / 4.1	135	Straight (shapeable)	1
	01-000264	0.039	4.0 / 4.1	145	Straight (shapeable)	1

* For availability and device compatibility please contact your local representative from Acandis®

www.3we.de 2107

NeuroSlider® Microcatheter DLC

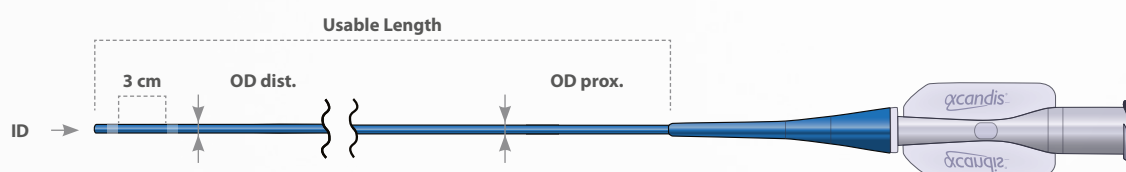

All changes or modifications, may they be technical or other, or changes in the availability of products are expressively reserved.

CE 0297

ADVANCE.

- Dual layer hydrophilic coating ensures outstanding lubricity and durability.
- Braiding / coiling reinforcement induces superior torqueability and significant reduction of ovalisation and elongation.


NAVIGATE.

- Shapeable tip with lasting shape retention allows excellent distal navigation even in tortuous anatomies.
- Multi polymer construction consisting of different flexibility zones with smooth transitions from maximum stability at the hub to maximum flexibility at the tip permits precise and effective navigation.

DELIVER.

- Inner PTFE liner minimises friction and allows a controlled and safe delivery of therapeutic and diagnostic agents.
- Advanced hub design with a transparent window results in a precise device transfer into the hub.

NeuroSlider® Microcatheter DLC

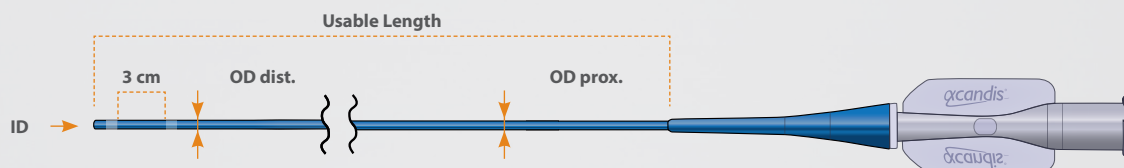

ORDERING INFORMATION

Product Name	Reference Number	ID (Inch)	OD dist. / prox. (French)	Usable Length (cm)	Tip Shape	Tip Marker
NeuroSlider® 17	01-000272	0.0165	1.9 / 2.1	155	Straight (shapeable)	2
NeuroSlider® 17 DLC	01-000282	0.0165	1.9 / 2.3	155	Straight (shapeable)	2
	01-000283	0.0165	1.9 / 2.3	160	Straight (shapeable)	2
	01-000284	0.0165	1.9 / 2.3	167	Straight (shapeable)	2
NeuroSlider® 21	01-000273	0.021	2.4 / 2.5	155	Straight (shapeable)	2
NeuroSlider® 21 DLC	01-000292	0.021	2.2 / 2.6	155	Straight (shapeable)	2
	01-000293	0.021	2.2 / 2.6	160	Straight (shapeable)	2
	01-000294	0.021	2.2 / 2.6	167	Straight (shapeable)	2
NeuroSlider® 27 (DLC)	01-000274	0.027	3.0 / 3.6	155	Straight (shapeable)	1

All changes or modifications, may they be technical or other, or changes in the availability of products are expressly reserved.

Distributed by:

acandis®
ENGINEERING STROKE SOLUTIONS

CE 0297

ACANDIS GmbH
Theodor-Fahrner-Str. 6
75177 Pforzheim
Germany

Tel: +49 7231 155 00 0
Fax: +49 7231 155 00 129
E-Mail: info@acandis.com
www.acandis.com