


BALTON[®]

www.balton.pl



Урология

NEPHROSTOMY KIT

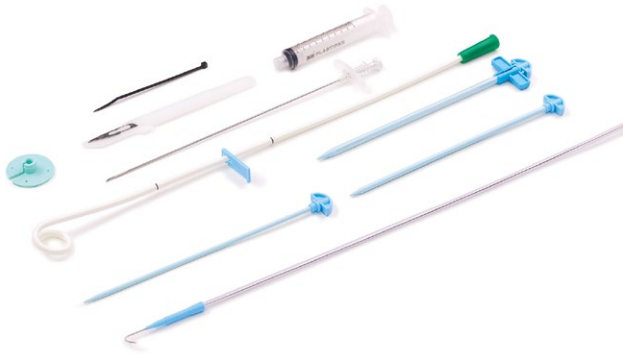

Kit components:

- Pigtail type catheter
- guidewire "J"
- two-part needle or three-part needle
- dilators
- dilator with sleeve to peel away
- fixing collar
- syringe 10 ml LL
- scalpel
- band clip

The kit for nephrostomy is used for making the nephrostomy fistula with the puncture method. We draw your attention to our modification of the superstiff guidewire. In this modification is a very thin tube of stainless steel with a softened part at the end. This provides a better possibility of safe manoeuvres with the dilator with a sleeve pushed on the superstiff guidewire.

The best quality polyurethane is used for the catheter. This plastic material ensures an optimally long time of catheter maintenance in nephrostomy.


Reference:	Catheter size	Catheter length [cm]	Dilators size	Needle type	Needle size/ length	Guidewire diameter [inch]
ZNEF5F	5F	45	5F, 6F*	Three-part needle	20Gx19 cm	.035"
ZNEF6F	6F		7F, 8F*			
ZNEF8F	8F		7F, 10F*	Two-part needle	18Gx20 cm	.038"
ZNEF9F	9F		7F, 10F*			
ZNEF10F	10F		7F, 10F, 13F*			
ZNEF11F	11F		7F, 10F, 13F*			
ZNEF12F	12F		7F, 10F, 13F*			
ZNEF14F	14F		10F, 13F, 15F*			

When ordering use Ref. Number, e.g.: ZNEF14F.

* Dilator with sleeve to peel away size



NEPHROSTOMY KIT FOR CATHETERS EXCHANGE

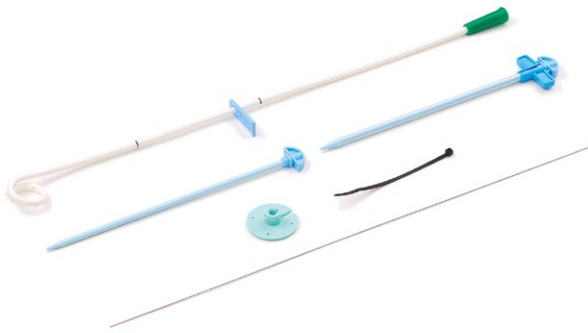

Kit components:

- Pigtail type catheter
- leading wire .038" x 80 cm
- dilators
- dilator with sleeve to peel away
- fixing collar
- band clip

Reference:	Catheter size	Catheter length [cm]	Dilator with sleeve to peel away size	Dilators size
ZWK9F	9F	45	10F	7F
ZWK12F	12F		13F	10F
ZWK14F	14F		15F	13F

When ordering use Ref. Number, e.g.: ZWK14F.

