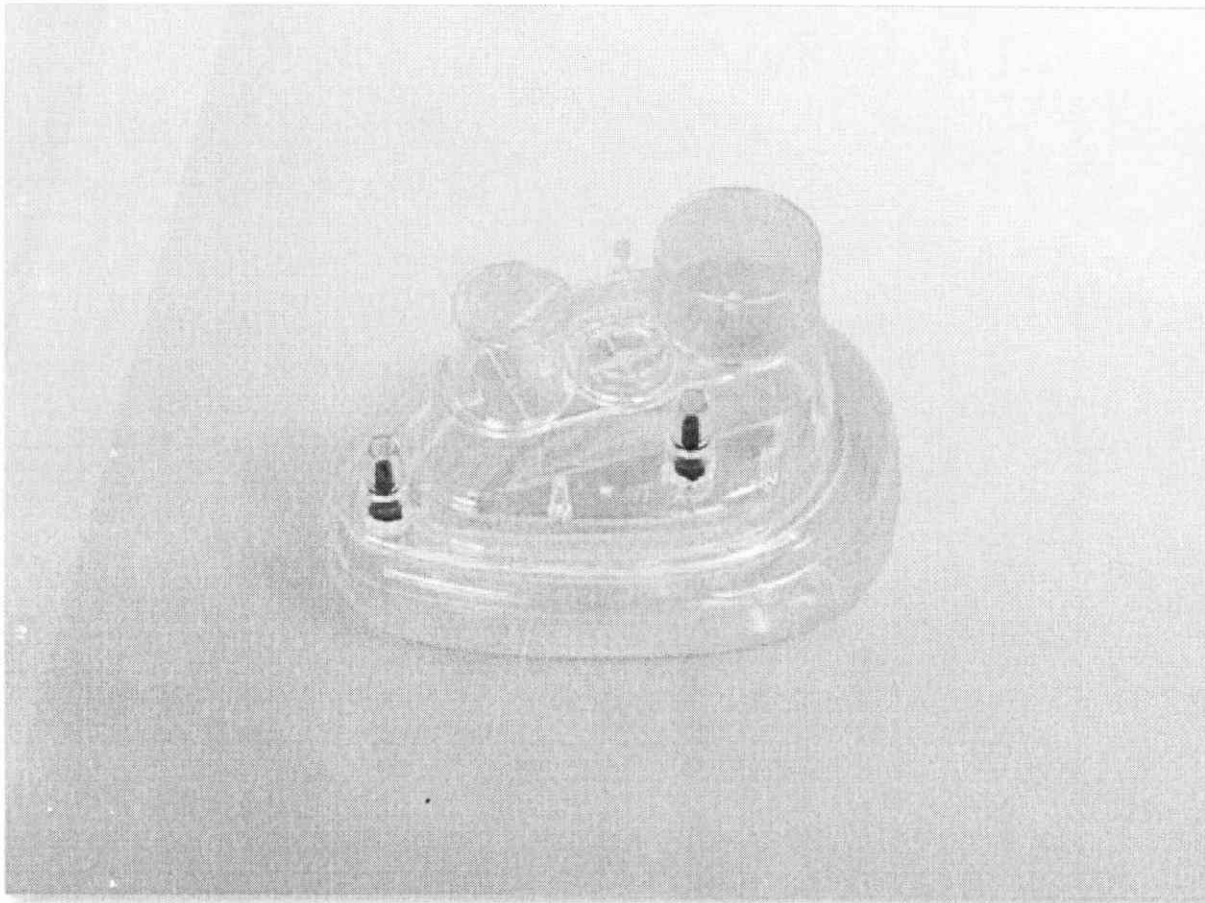


**DAHLHAUSEN®**

Lo-ful 84

**CPAP Masks**

CPAP mask with cuff. The mask cuff can be adapted to the shape of the head by the amount of air administered into the mask.


The mask has a double port. The smaller port with a 22 mm diameter is used to attach it to a tube for connection to the CPAP generator. The larger port with 30 mm is used to attach it the PEEP-valve required.

Art.no.	Product type
5160005521	Smal/medium
5160005520	Large

P. J. Dahlhausen & Co. GmbH  
Emil-Hoffmann-Str. 53  
50996 Köln

Tel.: (02236) 39 13 - 0  
Fax: (02236) 39 13 - 109

A handwritten signature in black ink, likely belonging to a representative of P. J. Dahlhausen &amp; Co. GmbH, located in the bottom right corner of the page.

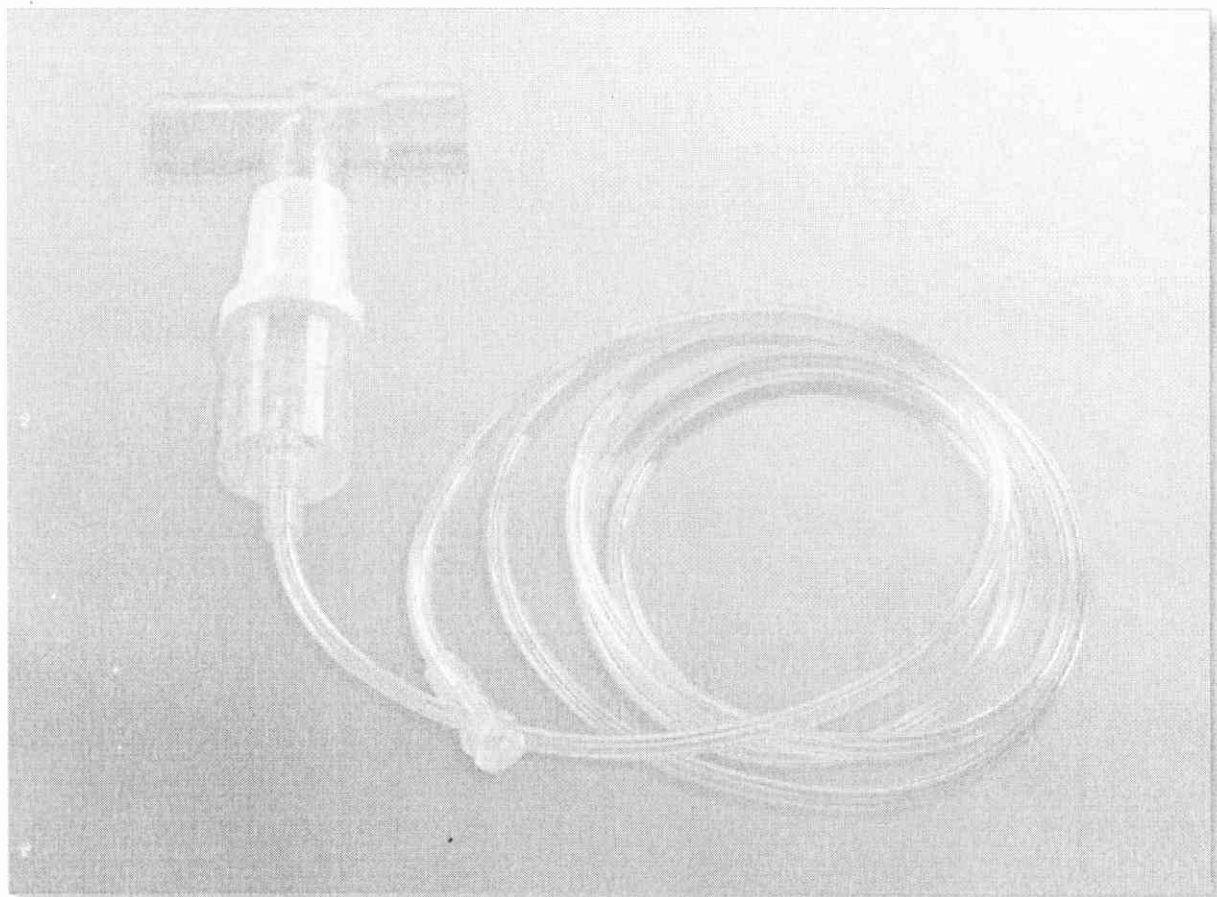**DAHLHAUSEN®**

### **Nebulizer Kit with standard connector**

Nebulizer Kit for respiration systems is used for nebulization of medicines and their inhalation.

The nebulizer has the following outputs:

Flow rate at 3,5 bar (LPM):

4 (minimum)

**6 (recommended)**

10 (maximum)

Air pressure (bar):

1 (minimum)

**3,5 (recommended)**

6 (maximum)

MMAD (NaCl 9 g/l) ( $\mu\text{m}$ ):

2,53 (at 4 LPM)

**1,18 (at 6 LPM)**

0,88 (at 10 LPM)

Nebulized quantity (NaCl 9g/l) (mg/min)

100 (at 4 LPM)

**200 (at 6 LPM)**

340 (at 10 LPM)

Volumes of solution 2-6 ml.

**Allergenic potential:** DEHP-free

**Application area:** Adult

<b>Length (oxygen tube):</b>	210 cm
<b>Material (connector):</b>	Polypropylene
<b>Material (oxygen tube):</b>	Polyvinyl Chloride
<b>Material (valve):</b>	Polystyrene, Silicone
<b>Material (connecting piece):</b>	T-piece: Polystyrene, Polypropylene, Polyethylene
<b>Material (nebulizer):</b>	Polyethylene
<b>Set composition:</b>	Non-return valve NRV), nebulizer chamber, safety connecting tube with two standard connectors

Art.no.

0100008905

P. J. Dahlhausen & Co. GmbH  
Emil-Hoffmann-Str. 53  
50996 Köln

Tel.: (02236) 39 13 - 0  
Fax: (02236) 39 13 - 109

