




magic
over 20 years of
flow dependence

ordering information

Reference	Total length (cm)	Guidewire max	Distal tip O.D.	Turquoise pursil O.D. (F) x L (cm)	White pursil O.D. (F) x L (cm)	Pink pursil O.D. (F) x L (cm)	Torver O.D. (F) x L (cm)
MAGIC1,2F	165	.008"	1.2F	1.2 x 12	1.8 x 3	2.4 x 30	2.7 x 114
MAGIC1,2FM	165		1.2F	1.2 x 3	1.8 x 12	2.4 x 30	
MAGIC1,2FM180	180		1.2F	1.2 x 3	1.8 x 17	2.4 x 40	
MA1,2FOLIVE	165		1.2F	1.2 x 12	1.8 x 3	2.4 x 30	
MAGIC1,5F	160	.009"	1.5F	/	1.5 x 15	2.4 x 25	2.7 x 114
MA1,5FMP	165		1.5F	/	1.5 x 20		
MA1,5FOLIVEMP	165		1.5F	/	1.5 x 20		
MAGICSTD	155	.010"	1.8F	/	1.8 x 10	2.4 x 25	
MAGICMP	165		1.8F	/	1.8 x 20		

magic compatibilities

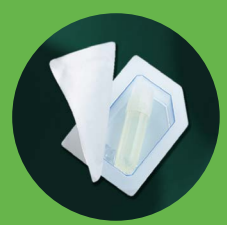

Cyanoacrylate Glue*



SPIF Flow Coil
(e.g. pressure cooker)

*MAGIC is glue compatible.

MAGIC, micro-catheters are intended for selective and hyperselctive vascular catheterization for diagnostic and therapeutic purposes. Class III CE0297 in compliance with Medical Device Regulation (MDR 2017/745). Manufactured by BALT EXTRUSION SAS. Carefully read the instructions for use before use. Not reimbursed. First CE marking:1998 The occlusion devices SPIF Flow Coil are intended for the embolization of vascular malformations of the vasculature. SPIF is class III CE0297 in compliance with Medical Device Directive (MDD 93/42/EEC amended by 2007/47/EC). Manufactured by BALT Extrusion SAS. Carefully read the instructions for use before use. First CE marking: 1998 (Magic); 2008 (SPIF). The content of this document, in particular data, information, trademarks and logos is BALT SAS's sole property. © 2022 BALT SAS and affiliates, all rights reserved. All representation and/or reproduction, whether in part or in full, is forbidden and would be considered a violation of BALT SAS and its affiliates' copyrights and other intellectual proprietary rights. This document with associated pictures is non-contractual and is solely dedicated to healthcare professionals and BALT's distributors (BALT's supplier's distributors). The products commercialized by BALT shall exclusively be used in accordance with the instructions for use included in the boxes. DC034GB (05/22)



Balt
10, rue de la Croix Vigneron,
95160 Montmorency, France
Tél. : +33 (0)1 39 89 46 41
Fax : +33 (0)1 34 17 03 46
www.baltgroup.com




magic

Flow dependent microcatheter

Intended for selective and hyperselective vascular catheterization for diagnostic or therapeutic purposes.

unique alloy pursil and torver

avm* treatment

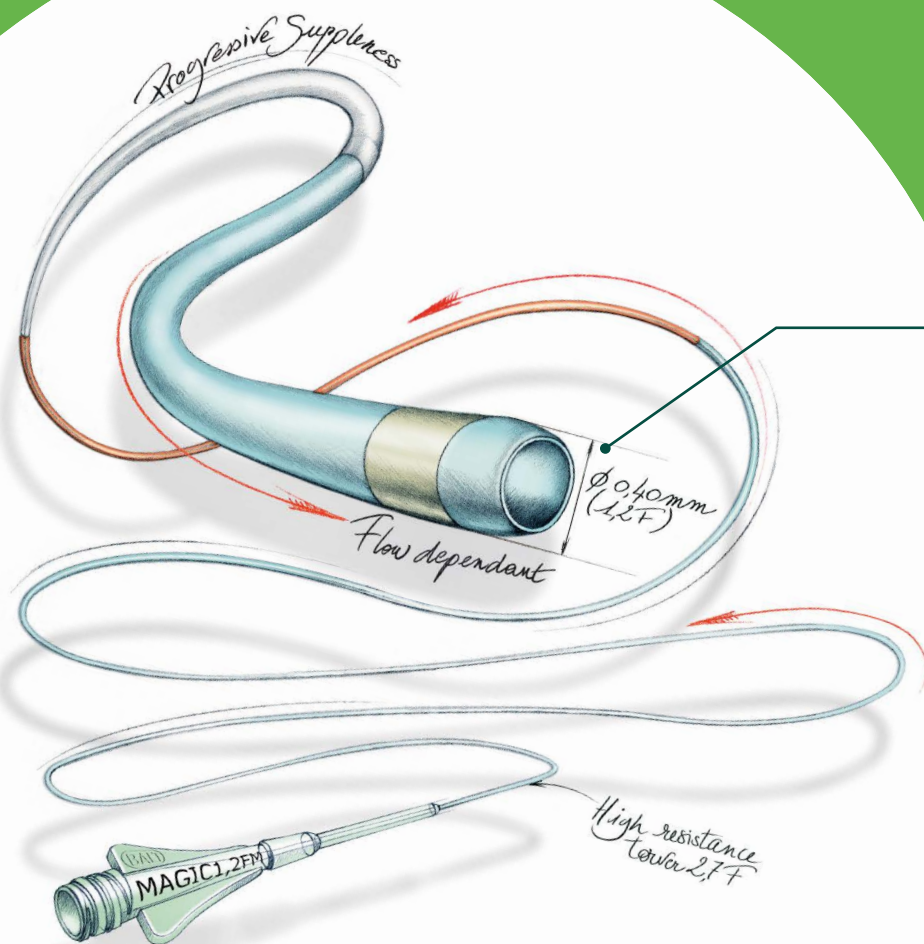

Flow dependent

MAGIC progression is facilitated by blood flow

Gentle navigation of the microcatheter inside distal and tortuous vasculature


● Distal access

The small size of the tip allows distal intracerebral hyperselective catheterism of small vessels

Unique alloy

In 1982, BALT developed an alloy that combines the mechanical resistance of polyurethane (PUR) and the flexibility of silicone (SIL): **PURSIL**

In 1987, another alloy called **TORVER** was introduced for the proximal part of the catheter

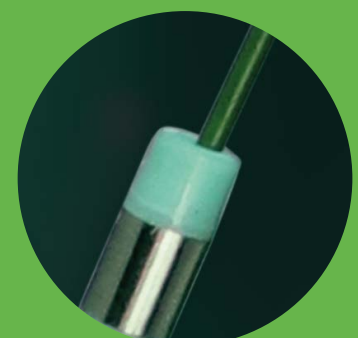

Up to **4** progressive suppleness transition zones

magic progression is facilitated by
blood flow


* Arteriovenous Malformation

Visibility


Technical features

Partea proximală
construită din polimer
transparent pentru
vizualizarea refluxului
sanguin și bulelor

Partea distală conține un
inel din Platinum
radioopac

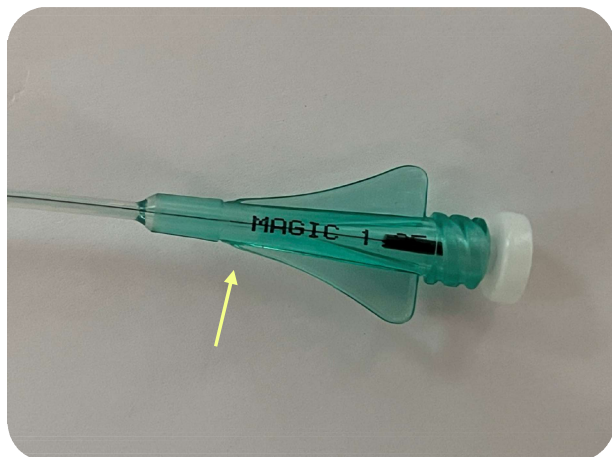

Transparent proximal part made of
polymer
→ visualization of blood reflux &
bubbles

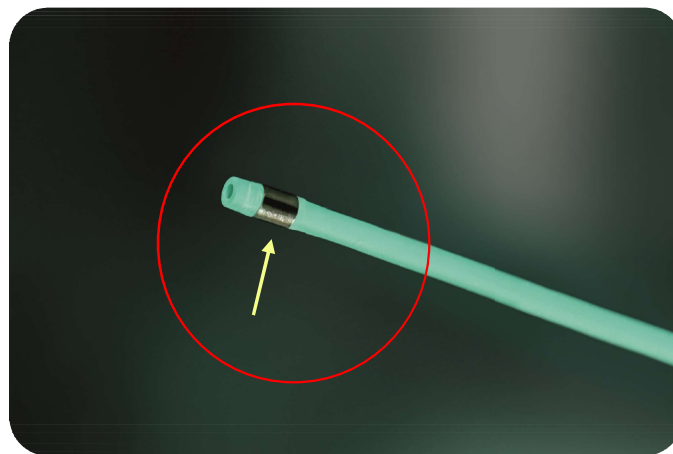

Radiopaque ring made of Platinum at
the distal tip



Radiopaque
microcatheter

The content of this document, in particular data, information, trademarks and logos is BALT SAS and affiliates' sole property. ©2018 BALT and affiliates, all rights reserved. Consequently, all representation and/or reproduction, whether in part or in full, is forbidden and would be considered a violation of BALT SAS and affiliates' copyrights and other intellectual proprietary rights. This document with associated pictures is non-contractual and is solely dedicated to BALT's distributors (BALT's supplier's distributors). The products commercialized by BALT shall exclusively be used in accordance with the instructions for use included in the boxes. The distribution, the sale and the use of BALT's products is under distributor's liability.

INTERNAL USE –CONFIDENTIAL



Indications

Technical features

N°1 flow-dependant microcatheter for distal navigation


AVF treatments

AVM treatments

For GLUE injections

Indicat pentru tratamentul AVF, AVM și injectarea agentului adeziv

Conține mandren

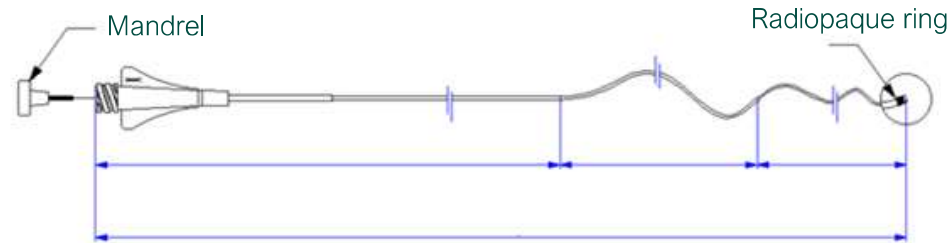

The content of this document, in particular data, information, trademarks and logos is BALT SAS and affiliates' sole property. ©2018 BALT and affiliates, all rights reserved. Consequently, all representation and/or reproduction, whether in part or in full, is forbidden and would be considered a violation of BALT SAS and affiliates' copyrights and other intellectual proprietary rights. This document with associated pictures is non-contractual and is solely dedicated to BALT's distributors (BALT's supplier's distributors). The products commercialized by BALT shall exclusively be used in accordance with the instructions for use included in the boxes. The distribution, the sale and the use of BALT's products is under distributor's liability.

INTERNAL USE –CONFIDENTIAL



Compatibilities

Technical features

Magic is compatible with the following Balt devices:



Can be inserted in BALT 6F guiding catheter and BALT 6F sheath introducer



As well as 6F guiding catheter (→ minimum of 1.20 mm / .047" internal lumen)

Magic *cannot be used* with DMSO-based embolization product

DMSO incompatibil

The content of this document, in particular data, information, trademarks and logos is BALT SAS and affiliates' sole property. ©2018 BALT and affiliates, all rights reserved. Consequently, all representation and/or reproduction, whether in part or in full, is forbidden and would be considered a violation of BALT SAS and affiliates' copyrights and other intellectual proprietary rights. This document with associated pictures is non-contractual and is solely dedicated to BALT's distributors (BALT's supplier's distributors). The products commercialized by BALT shall exclusively be used in accordance with the instructions for use included in the boxes. The distribution, the sale and the use of BALT's products is under distributor's liability.

INTERNAL USE –CONFIDENTIAL



Range details

Technical features

Format din 4 zone de tranziție (green, pink, white, turquoise) cu suplețe progresivă și rigiditate sporită de la distal la proximal, nu necesită ghid pentru navigare

Diametre distal:

- 1.2F (cu Ø min. a lumenului intern de 0.23mm, .009", GW max. .007")
- 1.5F (cu Ø min. a lumenului intern de 0.28mm, .011", GW max. .008")
- 1.8F (cu Ø min. a lumenului intern de 0.32mm, .013", GW max. .009")

Lungimi: 155, 160, 165, 180 cm

	REFERENCE	Green Torver Ø (F) x L (cm)	Pink PURSIL Ø (F) x L (cm)	White PURSIL Ø (F) x L (cm)	Turquoise PURSIL Ø (F) x L (cm)	Total length (cm)	Tip Ø	Dead space (mL)	Min Ø internal lumen	Max Ø olive tip
MAGIC1,2F	MAGIC1,2F	2,7 x 114	2,4 x 30	1,8 x 3	1,2 x 12	165	1,2F	0,38	0,23mm .009" GW max .007"	0,80mm
	MAGIC1,2FM	2,7 x 114	2,4 x 30	1,8 x 12	1,2 x 3	165	1,2F	0,39		
	MAGIC1,2FM180	2,7 x 114	2,4 x 40	1,8 x 17	1,2 x 3	180	/	0,40		
	MA1,2FOLIVE	2,7 x 114	2,4 x 40	1,8 x 3	1,2 x 12	165	Olive tip 1,2F	0,38		
MAGIC1,5F	MAGIC1,5F	2,7 x 114	2,4 x 25	1,5 x 15	/	160	1,5F	0,38	0,28mm .011" GW max .008"	1,02mm
	MA1,5FMP	2,7 x 114	2,4 x 25	1,5 x 20	/	165	1,5F	0,38		
	MA1,5FOLIVEMP	2,7 x 114	2,4 x 25	1,5 x 20	/	165	Olive tip 1,5F	0,38		
MAGIC1,8F	MAGICSTD	2,7 x 114	2,4 x 25	1,8 x 10	/	155	1,8F	0,38	0,32mm .013" GW max .009"	1,12mm
	MAGICMP	2,7 x 114	2,4 x 25	1,8 x 20	/	165	1,8F	0,38		

The content of this document, in particular data, information, trademarks and logos is BALT SAS and affiliates' sole property. ©2018 BALT and affiliates, all rights reserved. Consequently, all representation and/or reproduction, whether in part or in full, is forbidden and would be considered a violation of BALT SAS and affiliates' copyrights and other intellectual proprietary rights. This document with associated pictures is non-contractual and is solely dedicated to BALT's distributors (BALT's supplier's distributors). The products commercialized by BALT shall exclusively be used in accordance with the instructions for use included in the boxes. The distribution, the sale and the use of BALT's products is under distributor's liability.

