

FemoSeal™ Vascular Closure System



Description

Characteristics

Preparation

Item specifications

References

Description

Active vascular closure

The FemoSeal vascular closure system uses

3 bioabsorbable

components to actively seal the arteriotomy:

Two bioabsorbable polymer discs

- An inner disc deployed intra-arterially
- An outer locking disc deployed on the outer wall of the artery

Suture

Hemostasis is achieved by sealing the arteriotomy between the two discs, which are held together by bioabsorbable sutures complete with friction lock



Description

Characteristics

Preparation

Item specifications

References

Characteristics

Fast, effective, easy to use

- Rapid hemostasis with an easy to use deployment system

Fully bioabsorbable

- All components of FemoSeal vascular closure system are fully bioabsorbable: two polymer discs held together by a suture

Rapid hemostasis and early ambulation

- Designed to achieve effective and fast hemostasis. Sealing femoral artery punctures following catheterization procedures promotes rapid ambulation with low vascular complication rates.¹

Characteristics

Preparation

Item specifications


References



Characteristics

Preparation

Item specifications

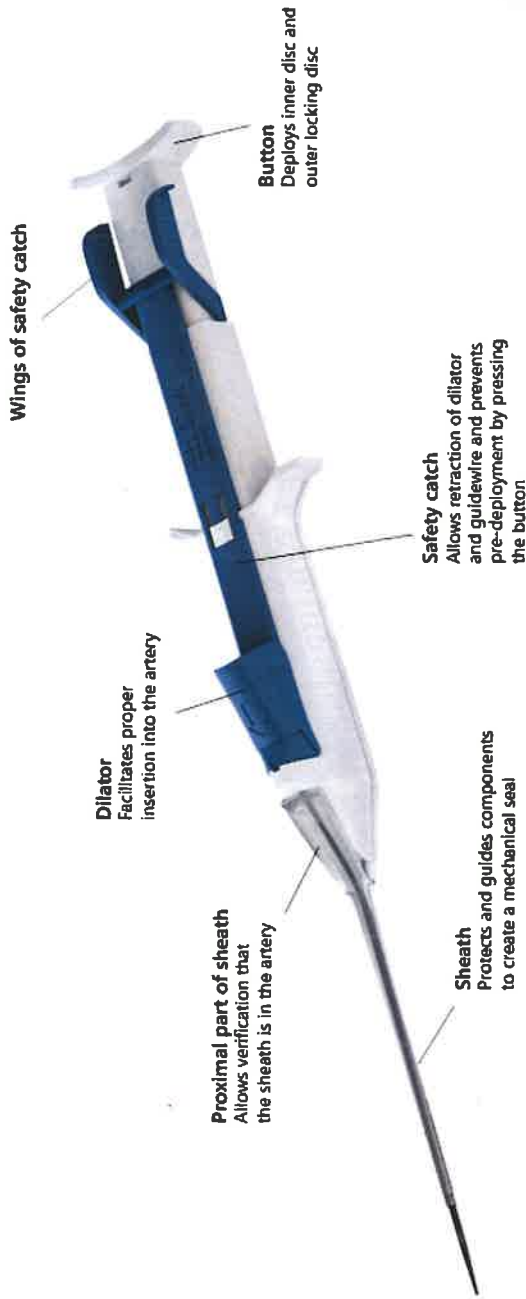

Privacy - Terms

Item specifications

Item reference	French size	Guidewire diameter (in)
C.11202	7	0.038