



The excellent choice of the professionals guarantees the improved accuracy and convenience about the optometry measurement.

POTEC

CE



Auto Ref-keratometer
PRK-7000

SPECIFICATIONS

MEASUREMENT MODE	
K/R Mode	Continuous Keratometry and Refractometry
REF Mode	Refractometry
KER Mode	Keratometry
CLBC Mode	Contact Lens Base Curve Measurement
K(P) Mode	Peripheral Keratometry
REFRACTOMETRY	
Vertex Distance(VD)	0.0, 12.0, 13.5, 15.0
Sphere(SPH)	-25.00 ~ +22.00 D (when VD=12 mm) (Increments : 0.12 and 0.25 D)
Cylinder(CYL)	0.00 ~ ±10.00 D (Increments : 0.12 and 0.25 D)
Axis(AX)	1 ~ 180° (Increments : 1°)
Cylinder Form	+, MIX
Pupil Distance(PD)	10 ~ 88 mm
Minimum Pupil Diameter	Ø 2.0 mm
KERATOMETRY	
Radius of Curvature	5.0 ~ 10.2 mm (Increments : 0.01 mm)
Corneal astigmatism	33.00 ~ 67.50 D (when cornea equivalent refractive index is 1.3375) (Increments : 0.05/0.12/0.25 D)
Corneal astigmatism Axis	0.00 ~ ±15.00 D (Increments : 0.05/0.12/0.25 D) 1 ~ 180° (Increments : 1°)
OTHERS	
Corneal Diameter	2.0 ~ 14.0 mm (Increments: 0.1 mm)
Memory of Data	10 measured value for each eye
Internal Printer	Thermal line printer with Auto-Cutter function
Monitor	14.48cm(5.7 inch) TFT LCD (640x480 pixels, Touch-Screen function)
Power supply	AC 100-240V, 50/60Hz
Dimensions	Approximately 260(W)x500(D)x450(H) mm
Weight	Approximately 20kg

* Designs and details can be changed without prior notice when its improvement.

POTEC POTEC Co., Ltd.

40-4, Techno 2-ro, Yuseong-gu, Daejeon
305-509, Korea
Tel +82 42 632 3536 / Fax +82 42 632 3537
webmaster@potec.biz

Distributed by

www.potec.biz

www.potec.biz

Make smart technology yours

Auto Ref-keratometer
PRK-7000

- Stable measurement result
- Convenient measurement
- Various functions
- The best design for user






Auto Ref-keratometer PRK-7000

The excellent choice of the professionals guarantees the improved accuracy and convenience about the optometry measurement.

optometry

1 Refraction & Keratometry Measurement

PRK-7000 provides an extensive dioptric measurement range (-25D to +22D), and the radius of curvature for keratometry is 5.0mm to 10.2mm. Easy alignment and friendly operation allows you to measure the refraction and keratometry in sequence, and the results can be validated immediately. Owing to its easy alignment and friendly operation which you can measure the refraction and keratometry in a sequence, the results can be checked all at once.

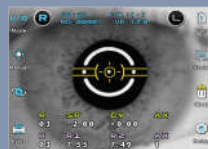

| Refraction & Keratometry Measurement

2 Peripheral Keratometry

Peripheral corneal curvatures can be measured by having the examinee look at the peripheral eye fixation lamps. Measuring the corneal periphery will help you examine irregular astigmatism, and also determine a better fitting for a contact lens



| Peripheral Keratometry

Auto Ref-keratometer PRK-7000



optometry

3 Intuitive Diameter Measurement

Using the freeze function, measurement of the diameter of the cornea, pupil or hard contact lenses worn by the patient can be performed. By simply touching and dragging the screen with your finger measurements are exactly calculated.



| Intuitive Diameter Measurement

4 IOL Measurement

When refraction results in an error reading due to an intraocular lens or cataract the measurement can be performed with the IOL icon switched 'ON'



| TFT-LCD with Touch Screen

1 TFT-LCD with Touch Screen and Convenient control

By adopting high resolution VGA TFT-LCD with Touch Screen function, we virtually removed all input keys except the measurement button on the Joystick and Chinrest control switches. Tests are performed rapidly and conveniently by simply pressing a button or the icon indicated on the screen. Also, user can use the Key button as convenient function.

2 Motorized chinrest

The chinrest is motorized with the use of a switch conveniently positioned for the operator. This makes for simple adjustment of the chinrest from patient to patient. There is the motorized chinrest which makes user easy to control by the simple switch pressing up and down.



interface

1 Instant Display Data

By simply touching the icon on the screen results of up to 10 measurements stored in memory can be viewed or printed by the built in thermal auto cut-off printer.

2 Interactive SETUP Change

Simply touch the icon on the screen and changed settings can be seen. The touch screen supports interactive setting changes quickly and conveniently saving user time innovatively.



| Instant Display Data


| Interactive SETUP Change

3 Assistant measurement function

PRK-7000 provides the assistant measurement function for the eyes which have different measuring between center of miring and center of REF focus.

4 The guidance function for measuring

PRK-7000 indicates the directions for focusing the center of the measuring eye and shows the conditions of focusing for the measurement convenience

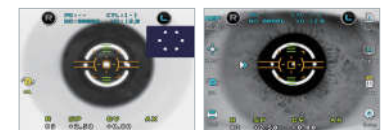

| Assistant measurement function

| The guidance function for measuring


convenience

3 Printer with Auto Cutter Function

PRK-7000 immediately provide a complete printout of the measurement results and an auto cutter is provided for convenience. And there is the economy mode of printer which makes printer paper saved by being well arranged and amended letter size.

4 Data transfer

Data can be transferred to external devices (personal computer, etc.) via an RS-232 interface. A USB interface is also provided for future compatibility. And user can send the images of patient's eyes through USB port and have the faster S/W upgrade by the PC program which is provided as option.

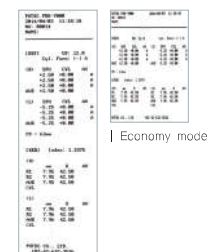

| Economy mode

| Normal mode

5 The tilting LCD for user

The 90° tilting LCD provides easy control by user standing or sitting.

6 The convenience stage-lock

User can fasten the main body by simple control.



| The tilting LCD for user



CT 1707/CT 1708

MOTORIZED TABLE

Specifications:

1. Table size: 72*49(CM)
2. Table elevation range: 73-94(CM)
3. Power: 100V-230V 50HZ/60HZ
4. Weight: 37KGS



CT 1707



Specifications:

1. Table size: 91*48(CM)
2. Table elevation range: 69-90(CM)
3. Power: 100V-230V 50HZ/60HZ
4. Weight: 35KGS



CT 1708



CT 1716/CT 1717

MOTORIZED TABLE

Specifications:

1. Table size: 56*48CM
2. Table elevation height: 63~81CM
3. Input voltage: 220V/50HZ, 110V/60HZ
4. Net weight: 18KGS

CT 1716



Specifications:

1. Table Size: 580 X 400mm(REGULAR) or 480 X 400mm
2. Range of elevation: 630mm - 820mm
3. Max Loading: 75kg
4. Electric Voltage : 220V/50Hz or 110V/60Hz
5. Power : 120W
6. Weight : 19kg

CT 1717



CT 1720/CT 1723

MOTORIZED TABLE

Specifications:

1. Table Size: 540 X 470mm
2. Range of elevation: 630mm - 820mm
3. Max Loading: 75kg
4. Electric Voltage : 220V/50Hz or 110V/60Hz
5. Power : 120W
6. Weight : 19kg



CT 1720



Specifications:

1. Table Size: 58*40CM
2. Size of drawer: 54*34.5*5(CM)
(Trial lens set is available)
3. Range of elevation: 63-82(CM)
4. Max Loading: 75KGS
5. Electric Voltage: 220V/50Hz or 110V/60Hz
6. Power: 120W
7. Weight: 19KGS



CT 1723



CT 1713/CT 1714

MOTORIZED TABLE

Specifications:

1. Table size: 96*48(CM)
2. Table elevation range: 77-98(CM)
3. Input: 100V-230V 50HZ/60HZ
4. Weight: 27KGS



CT 1713



Specifications:

1. Table size: 93*53(CM)
2. Table elevation range: 69-90(CM)
3. Input: 100V-230V 50HZ/60HZ
4. Weight: 40KGS



CT 1714



CT 1730/CT1731B

MOTORIZED TABLE

Specifications:

1. Table size: 92*55CM
2. Table elevation height: 69~90CM
3. Input voltage: 110V~240V
4. Net weight: 47KGS

CT 1730



Feature :

- Suitable for three instruments.
- Low noise, strong load-bearing.
- Large amplitude movements: 35cm.

Specifications:

1. Table Size: 140*68cm
2. Table elevation height: 69~104cm.
3. Input voltage: 110V/60HZ, 220V/50HZ
4. Net weight: 61KG
5. Max. loading: 90KG

CT 1731B



CT 1736/CT 1737

MOTORIZED TABLE

Specifications:

1. Table Size: 100*55CM
2. Range of elevation: 73-94(CM)
3. Electric Voltage: 110V/220V 50Hz-60Hz
4. Weight: 61KGS

CT 1736



Specifications:

1. Table Size: 98*60CM
2. Range of elevation: 76-97(CM)
3. Electric Voltage: 110V-230V 50Hz/60Hz
4. Weight: 57KGS

CT 1737

