

Cervix Examination Brush



Cervix Examination Brush

Use this brush for cell collection from the endocervical t-zone in Pap smears. It also can be used with thin prep technology and glass slides. The brush has a polyethylene shaft and convenient removable tip. Overall length 200mm.

Cat. No.	Description	Qty/Case (Pcs)
110203	Non Sterile, Bulk Pack	3000
110204	Individual Sterile Pack	3000

Boen Healthcare Co., Ltd.


Technical data sheet of cervix examination brush



1. Product Description

Disposable cervix examination brush is intended for the gynecological examination. It is composed of brush head and handle. The brush head is made of PE and the handle is made of PP. The brush can be with push tube or without push tube, the push tube is made of PP.

2. Product Construction

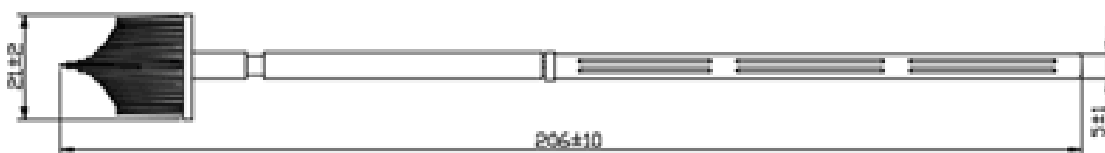

Type B

3. Manufacturing Process

Raw material → Injection → Assembly → Packing and Sealing → Sterilization → Resolving → Inspection → Storage

Injection, assembly, packing and sealing is key process, it should be made in one hundred thousand clean room. Sterilization is a special process.

4. Technical Requirement

Cervix examination brush should be smooth, no burr, no flying edges.

The brush head should be flexible.

The brush head and the handle can be connected firmly.

After being sterilized by EO, the residue should be less than 10mg/kg after 7 days.

5. Packing Details

It is individually packed in blister paper bag, 100pcs/box, 2000pcs/ctn.