



Specifications

Power Supply	100~240VAC, 50/60Hz, 200VA
Flow Rate	0.1~10ml/s
Pressure Limit	325PSI
Auto Filling Rate	8ml/s
Air Purge Rate	A: 1.0~9.9ml/s B: 2.5~9.9ml/s
Piston Rod Speed	9.9ml/s
Injection Delay	1~300s
Scan Delay	1~300s
Pause Phase	1~900s
Phases	6
Injection Programs	100
Injection Records	100
Volume	A: 1~65ml B: 1~115ml
Multiple Language User Interface	

High Quality Disposables



· DC Power



- 65ml Syringe for Contrast
- 115ml Syringe for Saline



imaStar
Up to 3 Tesla
MR Environments



MR Injector

Leading Contrast Media Injectors Forward

antmed

Shenzhen Antmed Co., Ltd.

18 Jinhui Ave., Pingshan New District, Shenzhen 518122, China
Tel: +86 755 8606 0992 Fax: +86 755 8606 0991
E-mail: info@antmed.com http://www.antmed.com

CE 0123 
P/N: EN-MR-201912

antmed
Partner For Life

MR Injector

Introducing an all-new ImaStar MR Injector, offering an extended list of product features, reliable performance, and a safe patient experience.



Ease of Operation and Clinical Efficiency

- Fill contrast and saline with a single press of a button
- Optional IV pole

Independent Dual Motor Drive Systems

- Dual motor drive capability for time saving loading and filling of syringes
- Allows independent injection of contrast and saline



Color Touch-Screen

- 12.1" color touch-screens for easy operation
- Clear and concise information display to ensure user operational confidence

Injection Record Storage


- Stores up to 100 injection history records



Bluetooth communication

Real-time Pressure Curve

- Display provides real-time pressure graph for injection monitoring

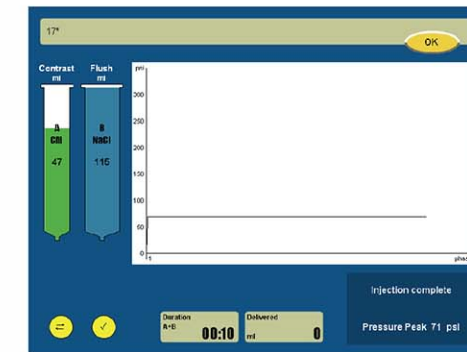


A Variety of Configuration Options

- Stores up to 100 customized injection programs
- Program storage capacity can be increased based on user needs



KVO Function

- Independent KVO (Keep Vein Open) ensures efficient examination procedure
- Flexible KVO interval time settings
- 115ml saline syringe capacity allows for longer KVO

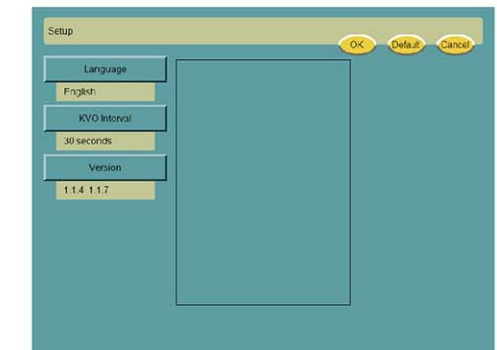

Safety

- To ensure the safety of the injection, there is a double confirmation feature with separate prompts to confirm air purging and injection ready status before any injection
- Automatically slows down the injection when approaching to pressure limit
- Injection stops when pressure limit is reached. Alerts with both light and sound
- Waterproof design inside the head to minimize damage from contrast media and saline leakage