





GENÇ

FDA

**SURGICAL
ISOLATION GOWNS
AND COVERALLS CATALOGUE**

2020

 www.gencgiyim.com.tr

CE

ISO 9001:2015

ISO 14001:2015

ISO 113485:2016

OHSAS 18001:2007

ISO 26000:2010

ISO 22301:2014

ISO 27001:2013

GUIDE TO PROTECTIVE MATERIALS

To help choose the right product, here is a brief description of the features and benefits of the main materials available at Genc Giyim;

Product Range and Materials Descriptions

GNC - 2207 Spunbond Polypropylene

Standard, cost-effective and comfortable
For basic infection control, this non-woven fabric bonds fibers together to form a single layer that is appropriate only for very minimal fluid exposure.

Medium Weight
Multi-Ply Protection
and Comfort Levels

Protection Good
Comfort Best




GNC - 2209 SMS Material

Balanced mix of protection and comfort
Strong and breathable, Spunbond/Meltblown/Spunbond (SMS) is a multi-layer fabric composed of inner layers of meltblown polypropylene between outer layers of spunbond polypropylene that is ideal for extended wear. Light/medium weight SMS is appropriate for low amounts of fluid and heavy weight SMS maybe appropriate for moderate contact with fluids.

Medium Weight
Multi-Ply Protection
and Comfort Levels

Protection Better
Comfort Best

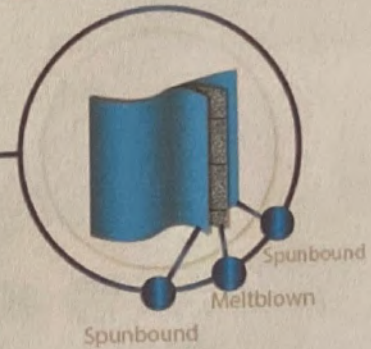


GNC-2211 Coated Polypropylene ("Poly-coated")

Soft, spunbond polypropylene is coated with a layer of polyethylene (plastic) film.

Medium Weight
Multi-Ply Protection
and Comfort Levels

Protection Best
Comfort Good

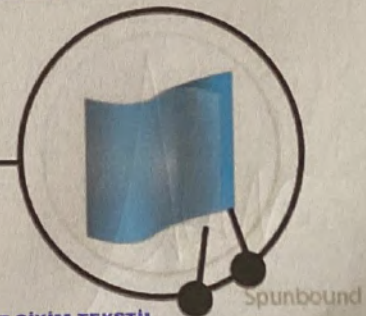

GENÇ TAHHÜT GİYİM TEKSTİL
İNŞAAT SAN. VE TİC. LTD. ŞTİ.

Merkez: Refahiye Mh.281. Sk. Erim Apt. Zemin Kat No:9/2
Tel: 0414 215 52 22 - Fax: 0414 217 16 36 H. Aliye/S. Ümrüpa
Şube: Generali Ali Rza, Gürcan, Sk. Merkez No:2/21
Tel: 0212 462 71 50 - Fax: 0212 172 11 63 Merter/İSTANBUL
Topçumutluları, D. No: 33 056 3225

Mail: muhasabe@gencgiyim.com.tr - gencgiyim@hotmail.com - www.gencgiyim.com.tr
Mersis No: 039305422500011-Ş. Ürfa Tic. Odası: 11430

GNC - 2207 Isolation Gown - Light Weight PP

Personal Protective Equipment (PPE) is a two-way street. Originally developed to protect the healthcare worker, PPE when used properly can also represent the first line of defense against contact transfer of pathogens like MRSA and viruses like Covid-19.

Lightweight spunbond polypropylene provides comfort and protection in very minimal fluid settings. For basic infection control, this non-woven fabric bonds fibers together to form a single layer that is appropriate only for very minimal fluid exposure.

Our products provide high protection and great comfort.

Single layer spunbond isolation gown designed for dry particle protection, featuring full coverage. Also these gowns are able to repel light amounts of fluids.

Meets AAMI Level 1 requirements

Light PP fabric

Natural rubber latex-free


Low lint generation

Anti-static treated fabric

Regulatory Compliance

Compliance with
EN ISO 13485:2016,
EN 14971:2012 and
EN ISO 15223-1:2016

 **GENÇTAHHÜT DİJİN TENSTİL
İNŞAAT SAN. VE TİC. LTD. ŞTİ.**
Merkez: Refahiye Mh.241. Sk. Suluç. Kat: 2. Kat. No: 9/2
Tel: 0414 215 52 21 Fax: 0414 215 52 21 Haliççe/İSTANBUL
Şube: General Ali Rıza Gürbüz Cad. No: 2/3
Tel: 0212 482 71 50 - 3 Fax: 0212 182 71 53 Merkez/İSTANBUL
Tölgelievleri V.D. 393 058 3225
Mail: tr@usex@gncgyt.com @ gncgyt@holmail.com - www.gncgyt.com
Mersis No: 0353003322500011-Ş.Urla Tic. Odası: 11430

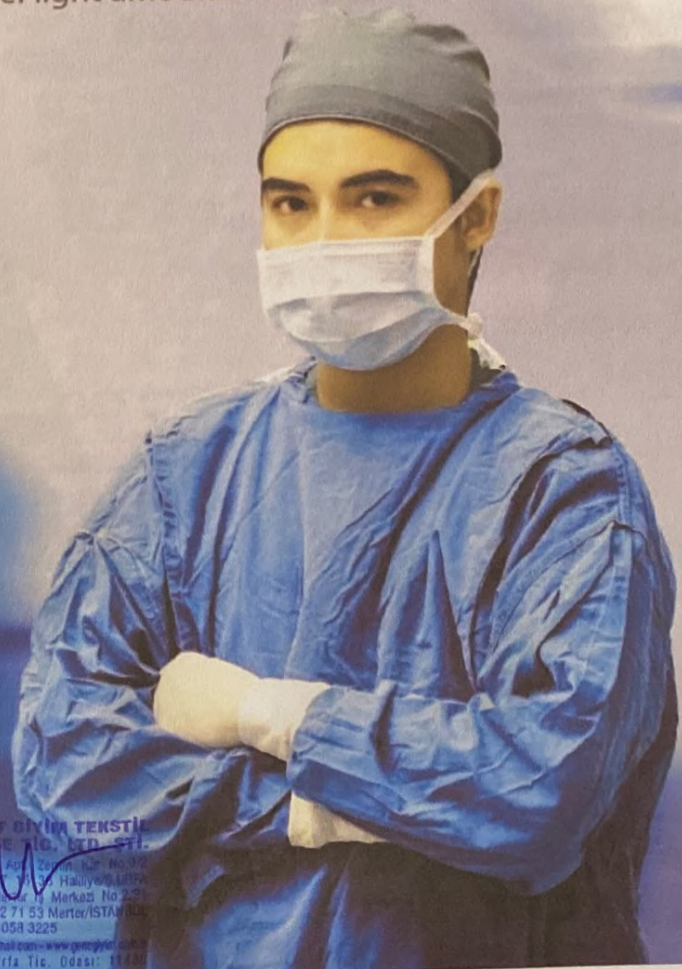

GNC - 2207 Fluid Resistant, Isolation Gown Light Weight PP



Size
S to XL


GENÇ

Basic spunbond polypropylene gown Comfort and protection for use in very minimal fluid settings. Generously sized for greater coverage and flexibility.

Suitable for health care workers and visitors, helps to prevent contamination with particulates

Knit cuffs, ties at the neck and waist, breathable back

Full coverage, breathable and comfortable - Waist and neck ties, Elastic or Knit Cuffs

 **GENÇ TAHHÜT GİYİM TEKSTİL
İNŞAAT SAN. VE TİC. LTD. ŞTİ.**
Adres: Refahiye Mh.2611 Ş. Salim Apt. Zemin Kat No:9/2
Tel: 0414 216 52 22 - Fax: 0414 217 10 36 Haliliye/Ş.URFA
Sakarya General Ali Rıza Gürün Cad. Merkez İş. Merkezi No:2/21
Tel: 0212 482 71 60 - 61 Fax: 0212 482 71 52 Merter/İSTANBUL
Töpcümevleri V. D. 93 058 0217
E-mail: info@gencgiyim.com gençgiyim@gmail.com - www.gençgiyim.com.tr
Tic. Sic. No: 039305832500011 Ş. Urfa Tic. Odası: 11430



GNC-2207

GNC ISOLATION GOWN

GNC-2209 Isolation Gown- Medium Weight SMS

We know you rely on disposable apparel to keep your staff and patients protected and comfortable. For this reason, we use the finest materials available for protection, offering the greatest level of comfort. We apply the highest quality standards to our products and in our manufacturing facilities.

Using the latest advancements in SMS technology, GNC gowns offer outstanding softness and enhanced breathability. This light weight line has great moisture evaporation, leaving the wearer cool and dry. Those who spend much of the day in a surgical gown appreciate the comfort of GNC. Count on three fabric layers, meltblown layer sandwiched between two strong spunbound outer layers, for excellent barrier protection against fluids and pathogens.

Our products provides high protection and great comfort.

For use when expected risk of exposure to fluid is between low and moderate

- ✓ Meets AAMI Level 2 requirements
- ✓ Medium-weight SMS fabric
- ✓ Recommended for
- ✓ Su turing, blood draw, inserting I.V. lines, specimen handling, drawing arterial blood
- ✓ IC U, Med/Surg Unit, Hyperbaric and Dialysis Units, Labs and Pathology, Nursery

Regulatory Compliance

Compliance with
EN 13795.,
EN ISO 13485:2016,
EN 14971:2012
and EN ISO 15223-1:2016

Latex-free protection

The "Latex-Free" symbol on all our gowns provides additional peace of mind for patients and clinicians.



GENÇ TAHHÜT GİYİM TEKSTİL
İNŞAAT SAN. VE TİC. LTD. ŞTİ.

Merkez: Refahiye Mh. 2811. Sk. Solim Apt. Zemin Kat No:9/2
Tel: 0414 215 52 22 - Fax: 0414 217 76 36 H. Urfa
Şube: General Ali Rıza Çarçan Cd. - Merkez 5. Market No:2/21
Tel: 0212 482 71 50 - Fax: 0212 482 71 50 Martir/İSTANBUL
Topçumudahi V.1: 303 95 32
Mail: mulasebe@gencyim.com.tr - gnc@icloud.com - www.gencyim.com.tr
Mersis No: 039305322500011-Ş. Urfa Tic. Odası: 11430

GNC-2209 Fluid Resistant, Isolation Gown Medium Weight SMS




Fluid Resistant, Isolation Gown - Medium Weight ,
SMS Suitable for procedures with moderate risks
of exposure to nonbiological fluids.

Knit cuffs, ties at the neck and waist, breathable
back

This Gown offers increased protection from
non-biological fluids without sacrificing the
level of comfort and breathability.

- Class 6 for Resistance to wet bacterial
penetration (EN ISO 22610 per EN 14126)
- Class 3 for Resistance to penetration
by contaminated solid particles
(EN ISO 22612 per EN 14126)
- Resistance to liquid penetration
>200cmH2O



 **GENÇ YAHHUT GİYİM TEKSTİL
İNŞAAT SAN. VE TİC. LTD. ŞTİ.**
Merkez: Rafahiye, M...
Tel: 0414 215 52...
Şube: General Al...
Tel: 0212 482 71...
Mail: muhasebe@...
Mersis No: ...

GNC - 2209
GNC ISOLATION GOWN

GNC - 2211 SURGICAL LAMINATED GOWN AGAINST LIQUIDS -PP + PE Microporous

We know you rely on disposable apparel to keep your staff and patients protected and comfortable. For this reason, we use the finest materials available for protection, offering the greatest level of comfort. We apply the highest quality standards to our products and in our manufacturing facilities.

GNC-2213 gown is a high performance gown. With microporous film technology, this surgical gown offers a superb combination of protection and comfort. Moisture vapor is 1,000 times smaller than the smallest viruses; hence 2211 gown blocks fluids while still allowing moisture vapour to escape, keeping the wearer dry and comfortable.

It provides great fluid protection. It is ideal for any situation demanding a high level of fluid protection. The Single Use Bilaminare Isolation Gowns are impervious and suitable for procedures with high risks of exposure to fluids.

- ✓ Meets AAMI Level 3 requirements
- ✓ High-weight PP fabric+PE Lamination
- ✓ Recommended for

Ideal for highly infectious situations and all operating room and ward applications. To prevent the spread of germs on emergency and isolation wards, in care homes, and in hospitals, for personnel, patients and visitors



Regulatory Compliance

Compliance with
EN 13795.,
EN ISO 13485:2016,
EN 14971:2012
and EN ISO 15223-1:2016


 **GENC TAHHÜT GIYİM TEKSTİL
İNŞAAT SAN. VE TİC. LTD. ŞTİ.**
Merkez: Refahiye Mh. 2811. Sk. Sağlık Apt. Katman No:9/2
Tel: 0414 215 54 21 - Fax: 0414 211 10 35 Haliçliysu, ÜRFA
Şube: General Ali Rıza Gökalp Cad. İnşaat İş Merkezi No:2/21
Tel: 0212 482 71 51 - Fax: 0212 340 11 53 Merter/İSTANBUL
Toplantı yeri: 093 053 3225
Mail: muhasebe@gncgiyim.com.tr | gncgiyim@hotmail.com | www.gncgiyim.com.tr
Mersis No: 03930010022500011-Ş. Urfa Tic. Odası: 11430

Latex-free protection

The "Latex-Free" symbol on all our gowns provides additional peace of mind for patients and clinicians.

GNC - 2211

SURGICAL LAMINATED GOWN AGAINST LIQUIDS - PP + PE Microporous

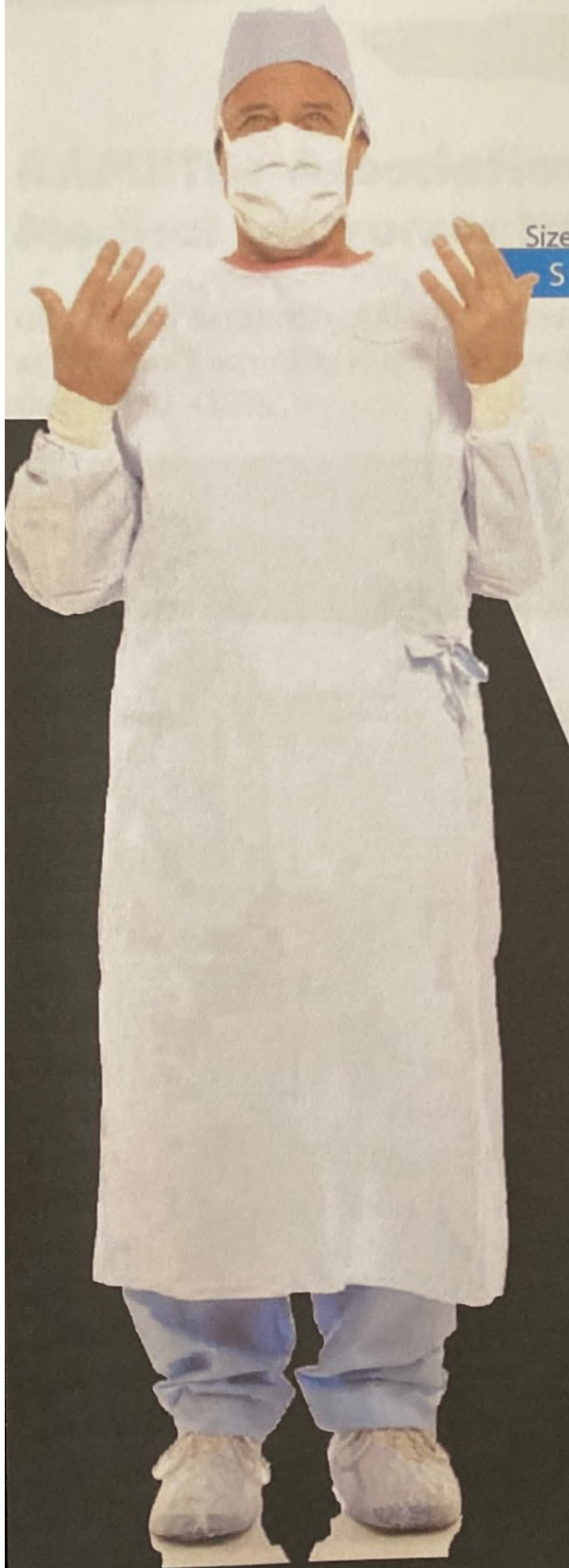

Size
S to XL


GENÇ

Surgical Laminated Isolation Gown - PP + PE
Microporous Suitable for procedures with high
risks of exposure to fluids.

Knit cuffs, ties at the neck and waist,
breathable back

- Class 6 for Resistance to wet bacterial penetration (EN ISO 22610 per EN 14126)
- Class 3 for Resistance to penetration by contaminated solid particles (EN ISO 22612 per EN 14126)
- Resistance to liquid penetration >500cmH2O



GNC-2211

 **GENÇ SAĞLIK GİYİM TEKSTİL İNŞAAT SAN. VE TİC. LTD. ŞTİ.**
Merkez: Baranlıye Mh.2811 Şişli, Salim İpot. Zeytinli, Kat: No:9/2
Tel: 0414 215 52 20 Faks: 0414 215 52 21 Halkiye/ŞURFA
Şube: General Ali Rıza Göktenel Cd. No: 2/21 Merkez: İSTANBUL
Tel: 0212 482 71 50 - 51 - Faks: 0212 482 71 53
Topçum Yılları: 03 058 3225
Mail: muhasebe@gencgiyim.com.tr - info@gencgiyim.com.tr - www.gencgiyim.com.tr
Mersis No: 039305322500011-Ş. Urfa Tic. Odası: 11430

USA Levels of Barrier Protection Standards

AAMI(The Association for the Advancement of Medical Instrumentation®) PB70 guidelines

GNC-2209 is classified by AAMI PB70 guidelines as PPE Level 2 Level 2 while GNC-2211 is classified as PPE Level 3 according to the Directive 89/686/EEC or MDD Class I non sterile according to the Directive 93/42/EEC.

ANSI/AAMI PB70 BARRIER PERFORMANCE	TEST METHOD	TEST DEFINITIO	REQUIREMENT	ANTICIPATED FLUID XPOSURE
AAMI Level 1	Water resistance: Impact penetration AATCC 42	AATCC42 Measures the resistance of fabrics to the liquid penetration of water by impact	Water impact < 4.5 g	Minimal fluid levels
AAMI Level 2	Water resistance: Impact penetration AATCC 42 Water resistance: Hydrostatic pressure AATCC 27	AATCC 42 Measures the resistance of fabrics to the liquid penetration of water by impact AATCC 127 Measures the resistance of fabrics to the liquid penetration of water by impact under constant and increasing hydrostatic pressure	Spray impact < 1.0 g Hydrostatic Pressure > 20 cm	Low fluid levels
AAMI Level 3	Water resistance: Impact penetration AATCC 42 Water resistance: Hydrostatic pressure AATCC 27	AATCC 42 Measures the resistance of fabrics to the liquid penetration of water by impact AATCC 127 Measures the resistance of fabrics to the liquid penetration of water by impact under constant and increasing hydrostatic pressure	Spray impact < 1.0 g Hydrostatic Pressure > 50 cm	Moderate fluid levels
AAMI Level 4	ASTM F1671, Standard Test Method for Resistance of Materials Used in Protective Clothing to Penetration by Blood-Borne Pathogens Using Phi-X174 Bacteriophage Penetration as a Test System	ASTM F1671 Measures the resistance of materials used in protective clothing to penetration by blood borne pathogens using a surrogate microbe under conditions of continuous liquid contact.	Pass	High fluid levels

GENÇ TAHHÜT GİYİM TEKSTİL İNŞAAT SAN. VE TİC. LTD. ŞTİ.
 Merkez: Refahiye, Mh. 51. Sk. Selim Ar. Zemin Kat. No:9/2
 Tel: 0414 215 42 22 Faks: 0414 117 89 93 Haliçliye/SURFA
 Şube: General Atik, Rıhtım Çarşı, Mh. 13. Merkez No:2/21
 Tel: 0212 482 71 50-51-52 C21/187713 Merkez/İSTANBUL
 Tıpçimeydanı V.D. 3309543225
 Mersis No: 0959200322500011900119000111430

Association for the Advancement of Medical Instrumentation. Liquid Barrier Performance Level Classification of Protective Apparel and Drapes Intended for use in Health Care Facilities. ANSI/AAMI PB70:2012. Arlington, VA: AAMI.

EU PERFORMANCE REQUIREMENTS AND PERFORMANCE LEVELS FOR SURGICAL GOWNS EN 13795

GNC-2209 AND GNC-2211 COMPLY WITH PERFORMANCE REQUIREMENTS FOR SURGICAL GOWNS BELOW;

PERFORMANCE REQUIREMENTS FOR SURGICAL GOWNS

CHARACTERISTIC	UNIT	STANDARD PERFORMANCE		HIGH PERFORMANCE	
		CRITICAL AREA	LESS CRITICAL AREA	CRITICAL AREA	LESS CRITICAL AREA
RESISTANCE TO MICROBIAL PENETRATION - DRY	LOC ₁₀₀ (CFU)	N/A	≤2 ^{a,c}	N/A	≤2 ^{a,c}
RESISTANCE TO MICROBIAL PENETRATION - WET	BI	≥2,8 ^b	N/A	6,0 ^{b,d}	N/A
CLEANLINESS - MICROBIAL	LOC ₁₀₀ (CFU/cm ²)	≤2 ^c	≤2 ^c	≤2 ^c	≤2 ^c
CLEANLINESS - PARTICULATE MATTER	IPM	≤3,5	≤3,5	≤3,5	≤3,5
LINTING	LOC ₁₀₀ (Lint Counts)	≤4,0	≤4,0	≤4,0	≤4,0
RESISTANCE TO LIQUID PENETRATION	Cm H ₂ O	≥20	≥10	≥100	≥10
BURSTING STRENGTH - DRY	kPa	≥40	≥40	≥40	≥40
BURSTING STRENGTH - WET	kPa	≥40	N/A	≥40	N/A
TENSILE STRENGTH - DRY	N	≥20	≥20	≥20	≥20
TENSILE STRENGTH - WET	N	≥20	N/A	≥20	N/A




CRITICAL AREA
LESS CRITICAL AREA

- a** Test conditions: challenge concentration 10⁸ CFU/g talc. and 30 minutes vibration time.
- b** The Least Significant Difference (LSD) for BI when estimated using EN ISO 22610, was found to be 0,98 at the 95% confidence level. This is the minimum difference needed to distinguish between two materials thought to be different. This means materials varying by up to 0,98 BI are probably not different; materials varying by more than 0,98 BI probably are different. (The 95% confidence level means that an observer would be correct 19 times out of 20 to accept these alternatives).
- c** For the purpose of this standard, log₁₀ CFU² 2 means maximum 300 CFU.
- d** BI = 6,0 for the purpose of this standard means: no penetration. BI = 6,0 is the maximum achievable value.

GNC TAARHUT GYIM TEKSTİL
İNŞAAT SAN. VE TİC. LTD. ŞTİ.
Merkez: Refahiye Anı 211, Sk. / Soğanlı, Kat: Zemin Kat No:9/2
Tel: 0414 215 60 22 / Faks: 0414 215 16 36 Haliççe/Ş. URFA
Şube: General / Rıza Çiftçiyer Cad. / Mimarlar İş. Merkez No:2/21
Tel: 0212 482 71 50 - 51 / Kat: C-12 / 482 71 53 Merkez/İSTANBUL
Tic. Sic. No: 27000 / Ş. No: 393 058 3225
Mail: muhasebe@gncplyim.com.tr / gnc.plyim@hotmail.com / www.gncplyim.com.tr
Mersis No: 0693058322500011-Ş. Urfa Tic. Odası: 11430


GENÇ

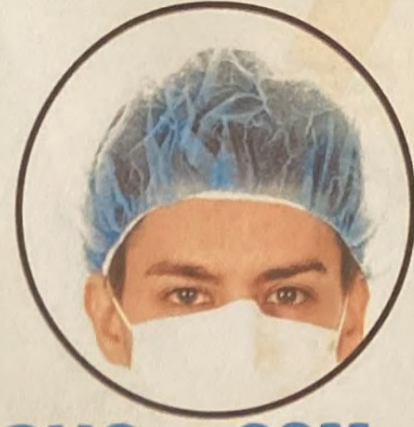


GNC - 2311 Head Cover

Nonwovens spunband fabric
12 gr/m2
Elastic
500 pairs



GNC - 2411 Shoe Covers

%100 Polythene
Elastic (with mouth rubber)
500 pairs with package

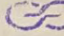 **GENÇ TAHHÜT GİYİM TEKSTİL
İNŞAAT SAN. VE TİC. LTD. ŞTİ.**

Merkez: Refahiye Mh/28. Sk. Selim Apt. Zemin Kat No:0/2
Tel: 0414 215 52 22 - Faks: 0414 517 116 30 - Hattıy/Ş. URFA
Şube: General Ali Rıza Güneşli Cad. No:2/21
Tel: 0212 482 71 51 - 51 - Faks: 0212 482 1153 Merter/İSTANBUL
Tic. Sic. No: 393 050 3225
E-posta: rahesebe@gencgiyim.com.tr - gencgiyim@hotmail.com - www.gencgiyim.com.tr
Tic. Sic. No: 0333066322500011-Ş. Urfa Tic. Odası: 11430



DISPOSABLE CHEMICAL PROTECTION


Protection against infectious diseases such as COVID-19

**GENÇ GÜVENLİK VE KORUMA ÜRÜNLERİ
İNŞAAT SAN. VE TİC. LTD. ŞTİ.**
Merkez: Refahiye Mh.2801. Sk. Sağlık Apt. Zemin Kat No:9/2
Tel: 0414 215 52 22 - Faks: 0414 217 15 55 - HATAY-Ş.URFA
Şube: General At. Plaza Güvenlik İş Merkezi No:2/21
Tel: 0212 482 71 50 - 51 - Faks: 0212 414 71 53 Merter/İSTANBUL
Tic. Sic. No: 270907 / D. No: 593 058 3225
Mail: muhasebe@genc.com.tr - genclik@genc.com.tr - genclik@hotmail.com - www.genclik.com.tr
Tic. Sic. No: 270907 / D. No: 593 058 3225 / T. No: 0212 414 71 53 / F. No: 0212 414 71 53 / P. No: 11430

Workers may come into contact with biological substances in a wide range of situations. Appropriate personal protective equipment must be used to prevent potential infection.

Clothing designed to provide protection against infectious diseases establishes a physical barrier between the wearer's skin and the source of the infection, and prevents the disease from spreading to other people or environments.

Generally, micro-organisms such as bacteria, fungi and viruses are classed as biological substances. A detailed definition is provided in EU Directive 2000/54/EC. The crucial common denominator between all of these materials is that they can cause infections, trigger allergies, or produce a toxic effect.

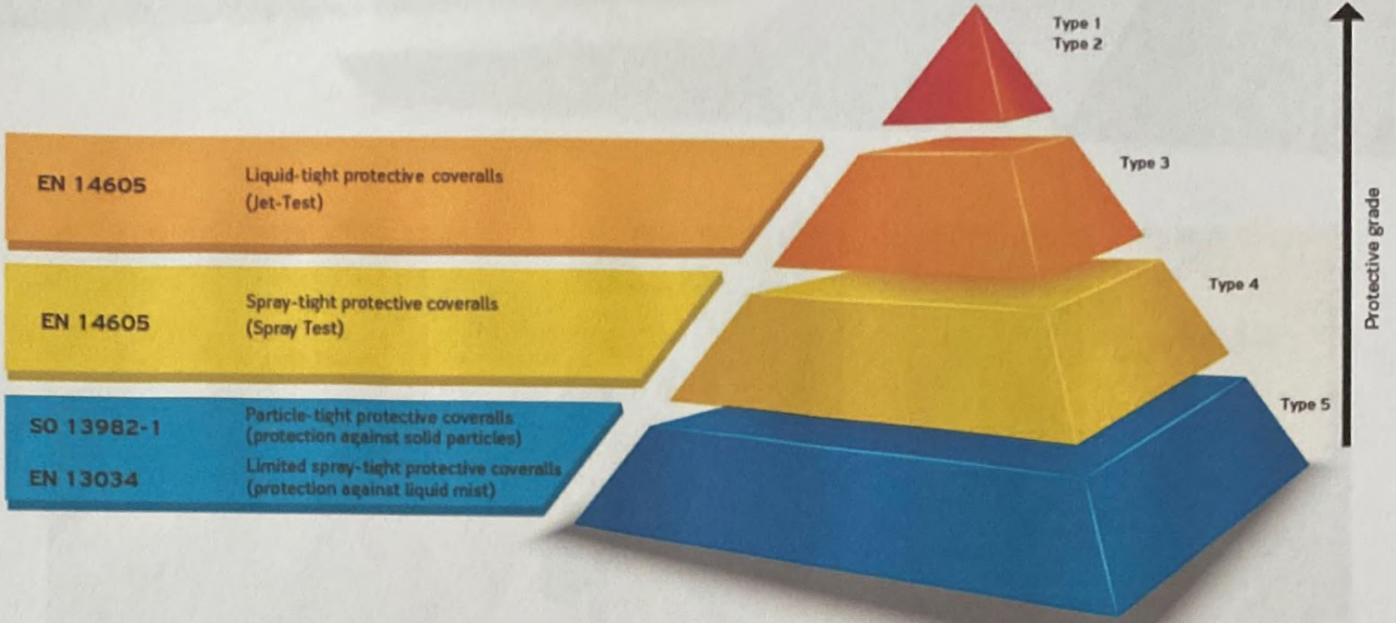

PPE Regulation EU 2016/425

The conditions for the sale of personal protective equipment (PPE) and the fundamental safety requirements that this PPE must meet are regulated in the PPE Regulation EU 2016/425.

If the requirements set out in the directive are met, the product is eligible for CE marking and can be offered for sale on the European internal market.

Depending on the level of risk that protective clothing is designed to protect against, PPE is assigned to one of three categories:

**GENC TAHHÜT GİYİM TEKSTİL
İNŞAAT SAN. VE TİC. LTD. ŞTİ.**
Merkez: Refahiye Mah. 2. Sk. Solim Apt. Zemin Kat No:9/2
Tel: 0414 214 52 22 Fax: 0414 217 11 55 HATAY / ÜRFA
Şube: General Ali Rıza İrçan Cad. Meşer İş Merkezi No:2/21
Tel: 0212 482 71 50 - 51 Fax: 0212 482 71 53 Merter / İSTANBUL
Tic. Sic. No: 290583
Mali: muhasebe@gencgiyim.com.tr genc.giyim@hotmail.com - www.gencgiyim.com.tr
Mersis No: 0393058322500011-Ş. Ürfâ Tic. Odası: 11430

1 CATEGORY I
Simple protective equipment, low risk

2 CATEGORY 2
Protective clothing that protects the wearer against hazards; PPE that does not fall into category 1 or category 3

3 CATEGORY 3
Protective clothing that protects the wearer against potentially fatal hazards or serious and irreversible damage to health, high risk

GNC Protection Clothes comply with Category III Standards. Chemical protection is always classified as Category III.

Protective clothing in accordance with EN 14126:2003



The EN 14126 standard defines the performance requirements for clothing designed to protect against infectious diseases. The defined test procedures focus on the medium that contains the micro-organisms, such as fluids, aerosols or solid dust particles. Due to the hetero-geneous nature of the micro-organisms, it is not possible to define performance criteria for them.

The tests stipulated by the standard relate only to the material; seam technologies are not taken into account. Taped seams provide a higher level of protection, as the micro-organisms are small enough to penetrate through the tiny needle holes along the seam.

The entire protective suit must also be certified as chemical protection clothing.

Wearing appropriate protective clothing does not provide complete, guaranteed protection against all chemical risks. It is also essential that you put on and remove the PPE correctly to ensure safety. Any person who assists in the removal of the clothing is also exposed to the risk of contamination.

**GENÇ TAHHÜT GİYİM TEKSTİL
İNŞAAT SAN. VE TİC. LTD. ŞTİ.**
Merkez: Refahiye Mh.2/21 Sk. Seim Apt. Zemin Kat. No:9/2
Tel: 0414 215 62 72 - Fax: 0414 211 16 76 Haliliye/Ş. URFA
Şube: General Ali Paşa Gökhan Cd. Merkez No:2/21
Tel: 0212 462 71 51 - 81 - Fax: 0212 462 71 53 Merter/İSTANBUL
Tic. Sic. No: 393 098 3225
E-mail: muhasebe@genctahhut.com.tr - info@genctahhut.com - www.genctahhut.com.tr
Mersis No: 0833006322500011-9-Urla Tic. Odası: 11430


GENC



Product Features

GNC-2215

GNC-2213

EN 14605



Liquid-tight coveralls
TYPE-3



EN 14605



Spray-tight coveralls
TYPE-4



EN 13034



Limited spray-tight coveralls (protection
against liquid mist)
TYPE-5



SO 13982-1



Particle-tight coveralls
(protection against solid particles)
TYPE-6



EN 1149-5



Protective clothing with antistatic properties
(antistatic function only guaranteed if air
humidity > 25%)



EN 14126



Protective clothing against infective agents




EN 1073-2




Protective clothing against particulate radioactive
contamination (no protection against radioactive
rays)



 **GENC TAHHÜT GİYİM TEKSTİL
İNŞAAT SAN. VE TİC. LTD. ŞTİ.**

Merkez: Refahiye M. 28. Sok. / Soliç / A. Z. Kat. No: 9/2
Tel: 0414 215 81 22 - Fax: 0414 215 81 23
Şube: General Ali Rıza Gürbüz Cad. No: 2/21
Tel: 0212 482 71 81 - 51 - Fax: 0212 482 71 83 Merter/İSTANBUL
Tic. Sic. No: 273930 / 393 059 3225

E-mail: muhasebe@gencgiyim.com.tr / gencgiyim@hotmail.com - www.gencgiyim.com.tr
Mersis No: 0390058322500011-Ş. Urfa Tic. Odası: 11430

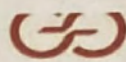

GENÇ

Certificates



**GENÇ TAHHÜT GİYİM TEKSTİL
İNŞAAT SAN. VE TİC. LTD. ŞTİ.**
Merkez: Rafahiye M.29/1. Sok. S.İ.İ.İ. Kat No:9/2
Tel: 0414 215 52 22 - Fax: 0414 215 52 16
Şube: General Ali Rıza Gürcan Cad. Merter İş Merkez No:2/21
Tel: 0212 482 71 50 - Fax: 0212 482 71 53 Merter/İSTANBUL
Tic. Sic. No: 270930 / Şirket Sic. No: 270930 / Mersis: 34090100000000000000
E-posta: muhasebe@gencyem.com.tr / gencyem@tccmail.com - www.gencyem.com.tr
No: 0303180322500011-Ş.Urlu Tic. Odası: 11430

Genc Taahhut Giyim Tekstil Ins. San.Tic. Ltd. Sti.

General Ali Riza Gurcan Cad. Merter Is Merkezi No: 2/21 Merter/Istanbul

Info@gencyem.com.tr Tel: +902124857150