F-1501 Wanda Plaza, No.18 Shijingshan Road, Beijing 100043, China Tel:+8610 88696020 Fax:+8610 88696085 A12.1502 **Laboratory Microscope** OP10 60U **JOC**




david@optoedu.com



Product Details



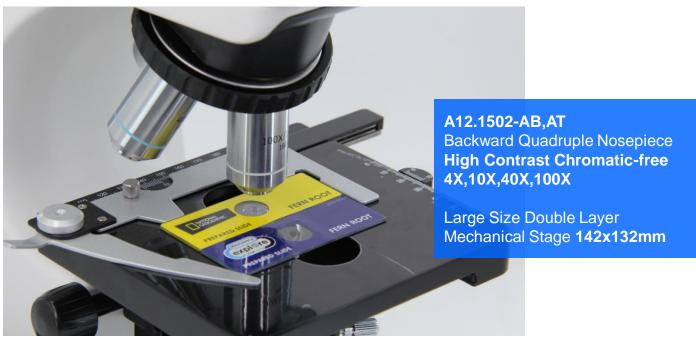







Condenser





Product Models










Product Specification

	A12.1502 Laboratory Biological Microscope	В	Т	AB	AT
Optical System	Infinity Optical System				
	Achromatic Optical System			•	•
Head	Binocular, Inclined at 30°, Interpupillary 48-75mm, 360° Rotatable				
i icau	Trinocular, Inclined at 30°, Interpupillary 48-75mm, 360° Rotatable				•
Eveniese	WF10x/22mm, Dia.30mm, High Eyepoint, Diopter Adjustment -/+5	•	•		
Eyepiece	WF10x/20mm, Dia.30mm, High Eyepoint, Diopter Adjustment -/+5			•	•
Nosepiece	Backward Quintuple Nosepiece	•	•		•
Objective	Infinite Plan High Contrast Chromatic-free 4X,10X, 20X,40X(S),100X(S,Oil)				
Objective	Plan High Contrast Chromatic-free 4X,10X,20X,40X(S),100X(S,Oil)				•
Condenser	Abbe Condenser NA1.25, with Iris Diaphragm, Center Adjustable			•	•
Focusing	Coaxial Coarse and Fine Adjustment, Fine Division 0.002mm, Coarse Stroke 37.7mm/Rotation, Moving Range 24mm, With Safety Stop Screw, Tension Adjustable	•	•	•	•
Morking Store	Large Size Double Layer Mechanical Stage 216x150mm, Graphite Surface, Moving Range 75mmx55mm, Reckless Safety Round Edge	•	•		
Working Stage	Large Size Double Layer Mechanical Stage 142x132mm, Graphite Surface, Moving Range 75mmx50mm			•	•
	Kohler Illumination, 3W/LED, Brightness Adjustable, With Blue Filter	•			
Light Source	External Illumination, 3W/LED, Brightness Adjustable, With Blue Filter Polarizing Polarizer & Analyzer			•	•
Adapter	C-Mount 0.5x,				•

System Size(mm)

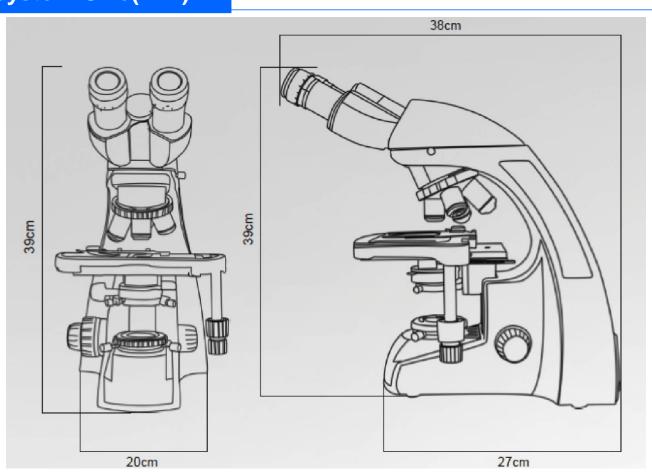

	A12.1502 Laboratory Biological Microscope Optional Accessor	ies
Eveniece	WF16X/F.N.15mm	A51.1502-16
	WF20X/F.N.12mm	A51.1502-20
Nosepiece	Backward Quintuple Nosepiece	A54.1502-5
II Jark Fidia	Dark Field Condenser, Dry, N.A.1.25	A5D.1502-D
	Dark Field Condenser, Immersion, N.A.1.25~1.36	A5D.1502-I
Polarizing	Polarizer & Analyzer	A5P.1501
Phase Contrast	Turret Disc Phase Contrast Condenser, For 10x,20x,40x,100x Infinity Plan Phase Contrast Objective 10x,20x,40x,100x Centering Telescope 10x Optional For B, T Only	A5C.1501
Fluorescent	Epi-Fluorescent Illuminator 100W Mercury Lamp House 12V100W Power Supply Box Fluorescent Filter B, G Optional For B, T Only	A5F.1501
Adapter	Eyepiece Tube Adapter	A55.1502-E
	C-Mount 1.0x	A55.1502-10



◆ Dark-field Viewing

It can be used for clearly viewing of blood or flagellum etc. fine structure.



A5D.1502-D



A5D.1502-I







◆ Polarizing Viewing It is quite suitable for viewing collagen, amyloid and crystal etc. double refracting specimen.





♦ Phase Contrast Viewing

Operators can get high contrast image of neutral background color whatever the magnifications are. It is suitable for viewing non-stained specimen.

Turret Disc Phase Contrast Condenser, For 10x,20x,40x,100x

Infinity Plan Phase Contrast Objective 10x,20x,40x,100x

Centering Telescope 10x

Optional For A12.1502-B/-T Only









