




MICROSCOPY AND
LABORATORY DIAGNOSTICS

Our global presence

“ We manufacture products for use in histopathology, hematology, cytology, and microbiology with a wide national and international distribution network present in over 70 countries.

BioGnost offers only the best quality products at competitive prices and has one of the widest ranges of medical products in South-Eastern Europe. Our experience and flexibility towards the client makes us a sought-after business partner, which has been recognized on the global market.

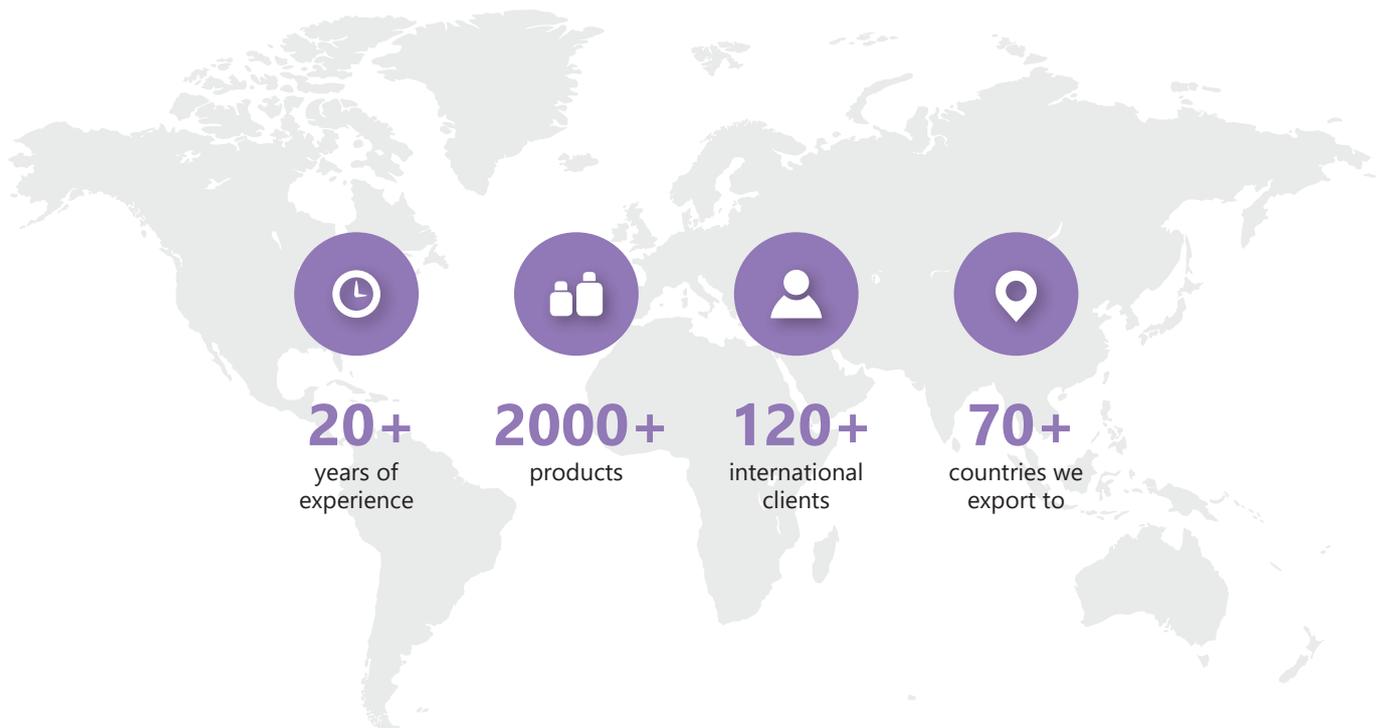

Our entrepreneurship started in 1990 with the founding of the family company. Accelerated business growth began in 1998 with production of blood typing serums followed by the Croatian Ministry of Health granting BioGnost the License for the manufacture of medical products. We are proud of one of our biggest successes - relocating to our new commercial and manufacturing building in 2012. The site enabled an enhanced work environment, including a carefully devised warehouse structure, segmented production facilities with increased manufacturing capacity, and additional R&D capabilities.

“ We continue to build on the strong momentum from last year to create high-quality products and provide excellent customer service.


Strategic priorities

Our mission is to create high-quality products and provide excellent services, as well as introduce new or expanded product lines and sales ranges in accordance with market dynamics and demands.

Product and service excellence makes BioGnost a desirable regional and international business partner. The key to success is the ever-expanding network of business partnerships, as well as recognizable social responsibility. In the coming years we plan to further our market expansion, additionally invest in research and development, and to become a competitive manufacturer in the fields of microbiological diagnostics and histopathology.

ISO 9001, ISO 13485, ISO 14001

BioGnost’s entire business is in accordance with international ISO 9001, ISO 13485 and ISO 14001 standards. The company implemented ISO standards more than 15 years ago, and is the first Croatian company to implement the ISO 13485 quality system for medical devices. Having fulfilled the requirements of these standards, the company is completely in accordance with national legal regulation and the European Union’s Medical Device Regulation (EU 2017/745).

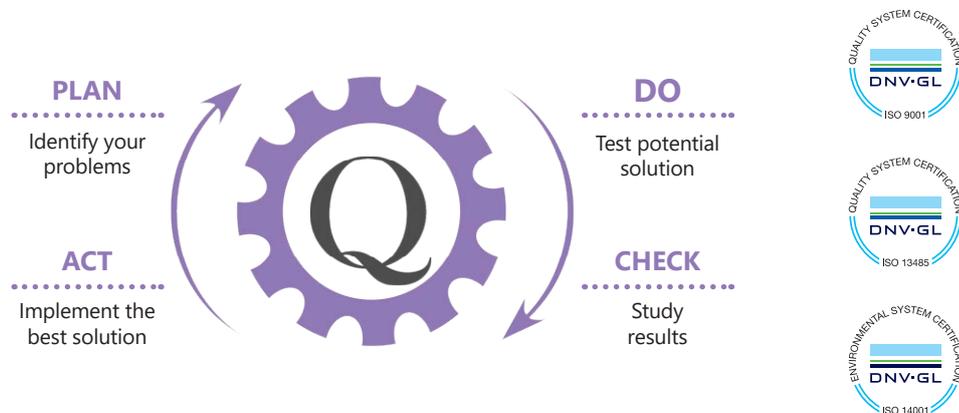

CE, IVDR

BioGnost’s products bear the CE mark, a certificate that confirms that the products are in accordance with the European Union’s standards. They are declared as *in vitro* diagnostic medical devices (IVD) classified as “Class A” products (according to the Regulation (EU) 2017/746). No control of a Notified Body is required for this class; entire company business is in accordance with applicable regulations and standards.

OEM manufacture

Besides products from our own sales range, we also offer OEM manufacture for your company. It includes a wide range of services, such as label design, product packaging, instructions for use and certificates of analysis, and products manufactured according to your specifications.

Contents



FIXATION AND PROCESSING 6

Marking systems 6

Formalin fixatives 8

Other fixatives 10

Decalcifying solutions 11

Dehydrating and rehydrating agents 12

Clearing and deparaffinization agents 13

Infiltration, embedding and demolding media 14

MOUNTING AND MICROSCOPY 16

Mounting media 16

Immersion oils 20

Cleaning solutions for microscope lens and slides 21

STANDARD HISTOLOGY STAINING REAGENTS 22

Hematoxylin reagents 22

Differentiation and bluing reagents 24

Eosin reagents 25

CYTOLOGY AND HEMATOLOGY STAINING REAGENTS 26

Papanicolaou (Pap) reagents 26

Rapid and standard hematology reagents and kit 28

LEUKOGNOST Cytochemical staining kits and ancillary reagents 32

BACTERIOLOGY AND MYCOLOGY STAINING REAGENTS 34

BioGram kits and reagents 34

TB-Stain kits and reagents 36

Solutions for staining and analysis of fungi 39

SPECIAL STAINING KITS FOR HISTOLOGY AND CYTOLOGY 40

ADDITIONAL STAINING REAGENTS AND SOLUTIONS 51

POWDER DYES AND STAINS 58

ACCESSORIES 62

Legend:

● DG cargo in transport ○ DG cargo only in air transport ● temperature controlled cargo (transport regime 2-8°C)

Tissue marking dyes, 5 color set

TMD Tissue marking dyes set 5x50 mL contains 5 dyes for tissue marking: green, yellow, black, red and blue. The set includes wooden sticks and rack.

5x50 mL 3408-SET

Tissue marking dyes, 7 color set

TMD Tissue marking dyes set 7x50 mL contains 7 dyes for tissue marking: green, yellow, black, red, blue, orange and lilac. The set includes wooden sticks and rack.

7x50 mL 3408-SET-7

Tissue marking dyes

For marking and orientation of tissue samples.

Tissue marking dye, green	50 mL	3408-1-50
	100 mL	3408-1-100
	250 mL	3408-1-250
Tissue marking dye, yellow	50 mL	3408-2-50
	100 mL	3408-2-100
	250 mL	3408-2-250
Tissue marking dye, black	50 mL	3408-3-50
	100 mL	3408-3-100
	250 mL	3408-3-250
Tissue marking dye, red	50 mL	3408-4-50
	100 mL	3408-4-100
	250 mL	3408-4-250
Tissue marking dye, blue	50 mL	3408-5-50
	100 mL	3408-5-100
	250 mL	3408-5-250
Tissue marking dye, orange	50 mL	3408-6-50
	100 mL	3408-6-100
	250 mL	3408-6-250
Tissue marking dye, lilac	50 mL	3408-7-50
	100 mL	3408-7-100
	250 mL	3408-7-250

Liquid Blocker (PAP pen)


Liquid Blocker contains hydrophobic liquid as a water repellent barrier around specimens or sections. Useful in immunohistochemistry that uses small amounts of antiserum. Multiple specimens can be separated by drawn circles or lines applied to the same slide.

1 pc ML-1

Duoline marker pen

Professional marker pen for permanent labelling of cassettes, slides and transport containers. With fine and thick marker pen tip.

1 pc MD-1

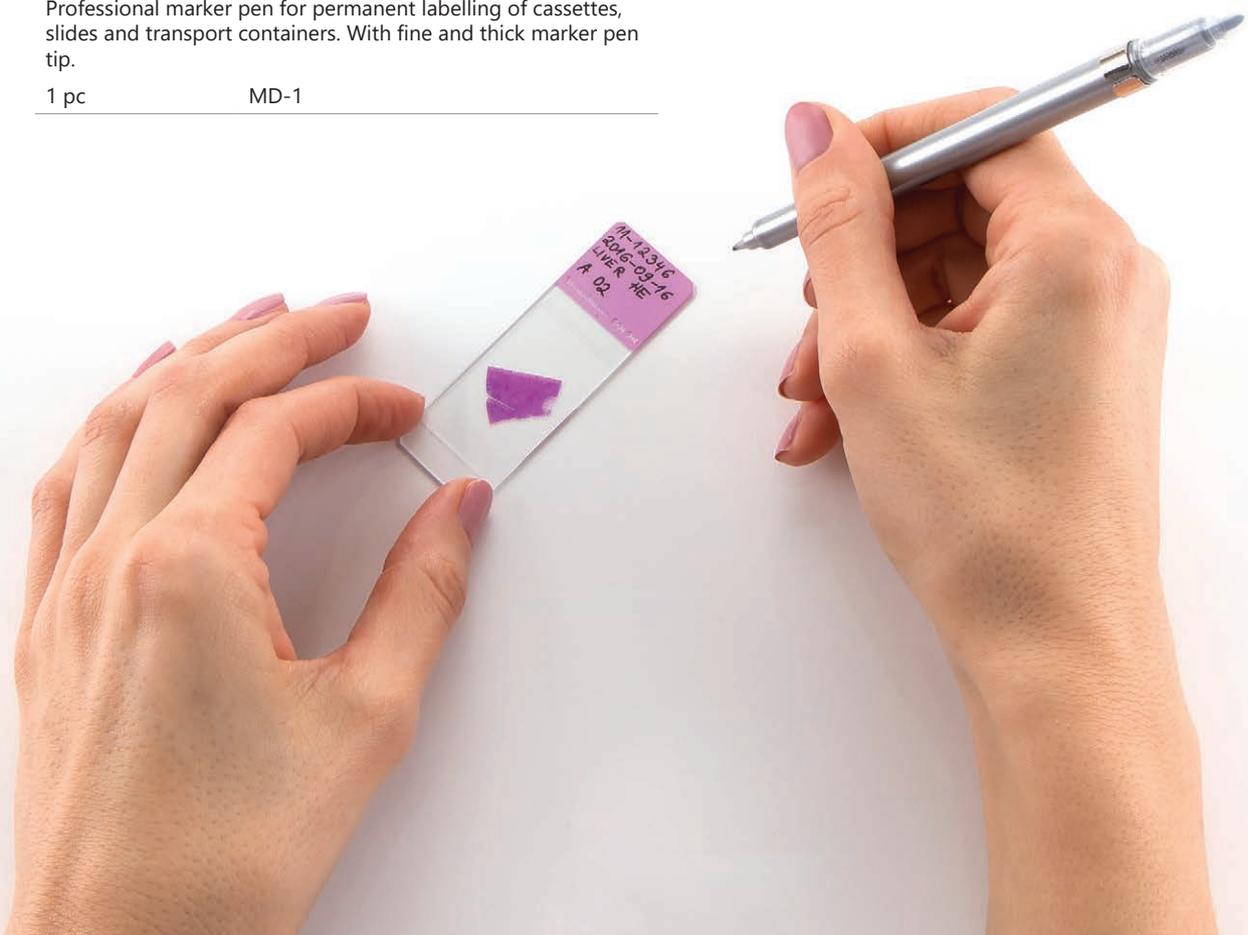

10% NB Formaldehyde

25% neutral buffered and stabilized formalin solution, pH 7.0.

1000 mL	FNB10-1L
5000 mL	FNB10-5L
10000 mL	FNB10-10L
20000 mL	FNB10-20L



ImmunoForm

Stabilized 10% formalin solution produced from paraformaldehyde without methyl alcohol. For use in histology and immunohistochemistry.

1000 mL	IMF-1L
5000 mL	IMF-5L
10000 mL	IMF-10L



10% Formaldehyde

Stabilized 25% formalin solution.

1000 mL	F10-1L
5000 mL	F10-5L
10000 mL	F10-10L
20000 mL	F10-20L



Formaldehyde 37% for use in histology

100% formalin solution stable at low temperatures.

1000 mL	F37H-1L
5000 mL	F37H-5L
10000 mL	F37H-10L



NB Formaldehyde concentrated 5x

Concentrated neutral buffered formalin for preparation of 10% buffered formalin (synonym: 4% buffered formaldehyde) solution, pH 7.0.

For preparation of 5 L of 4% buffered formaldehyde	1000 mL	FK5X-1L
--	---------	---------

For preparation of 25 L of 4% buffered formaldehyde	5000 mL	FK5X-5L
---	---------	---------

For preparation of 50 L of 4% buffered formaldehyde	10000 mL	FK5X-10L
---	----------	----------



DeHematin

Modified Kardasewitsch solution for depigmentation of tissues fixed in non-buffered formalin (after acid hematin occurs).

1000 mL	KO-OT-1L
5000 mL	KO-OT-5L



Other fixatives

BioFix

Methyl alcohol-based fixative for cryostat sections.

500 mL	BF-500
1000 mL	BF-1L



BioFix C

Modified Carnoy's fixative based on ethyl alcohol, chloroform and acetic acid.

500 mL	BFC-500
1000 mL	BFC-1L
5000 mL	BFC-5L



BioFix GL

Glyoxal-based fixative for use in histology.

500 mL	BFGL-500
1000 mL	BFGL-1L
5000 mL	BFGL-5L

BioFix HD

Modified Hartmann/Davidson's fixative that enables rapid fixation and better visualization of lymph nodes. The results provide clearer nuclear details and minimize formalin pigment occurrence.

500 mL	BFHD-500
1000 mL	BFHD-1L
5000 mL	BFHD-5L



BioForm S

Modified Schaeffer's buffered fixative solution of formaldehyde in methanol.

1000 mL	BIOS-OT-1L
5000 mL	BIOS-OT-5L



Decalcifying solutions

Decalcifiers most often contain inorganic or organic acids enabling rapid decalcification (i.e. calcium removal). Milder types of decalcifiers are formulated from chelating agents that are used for decalcification of more sensitive types of tissue.

OsteoFast 1

Light blue fixative and decalcifying solution for bone and hard tissue in histology, based on formaldehyde and hydrochloric acid.

1000 mL OF1-OT-1L
2500 mL OF1-OT-2.5L



OsteoFast 2

Light yellow decalcifying solution for bone and hard tissue in histology, based on formic acid.

1000 mL OF2-OT-1L
2500 mL OF2-OT-2.5L



OsteoSens

Decalcifying solution for sensitive calcified tissues in histology, contains EDTA.

1000 mL OS-OT-1L
2500 mL OS-OT-2.5L



CryoFix Gel

Embedding matrix for cryosectioning of tissue samples known as OCT medium (optimal cutting temperature). Efficient at temperatures of -10°C and below. Available in 5 different colors.

CryoFix Gel	100 mL	CF-100
CryoFix Gel pink	100 mL	CFR-100
CryoFix Gel blue	100 mL	CFP-100
CryoFix Gel green	100 mL	CFZ-100
CryoFix Gel yellow	100 mL	CFY-100



Cryo Spray

Medium in spray bottle for freezing histological samples.

200 mL 08-SPRAY-2



DeMold Spray

Solution in spray bottle for easy removal of paraffin tissue blocks from the mold.

100 mL DS-100



DeMold solution

Solution for refilling DeMold Spray bottle.

500 mL DS-500
1000 mL DS-1L



Tissue Soft

Solution for softening the tissue embedded in paraffin block before microtome sectioning.

500 mL TS-500
1000 mL TS-1L



BioGnost's synthetic BioMount media properties

PRODUCT NAME	SOLVENT	POLYMER	VISCOSITY @ 20°C (cSt, mPa)	REFRACTIVE INDEX @ 20°C (n _p)	USAGE IN AUTOMATED COVER SLIP MACHINES
BIOMOUNT DPX	xylene (BioClear)	S	850-950	1.5210-1.5230	-
BIOMOUNT DPX HIGH	xylene (BioClear)	S	1900-2100	1.5235-1.5255	-
BIOMOUNT DPX LOW	xylene (BioClear)	S	140-160	1.5115-1.5135	✓
BIOMOUNT DPX NEW*	xylene (BioClear)	S	140-160	1.5150-1.5180	✓
BIOMOUNT NEW	isoparaffin (BioClear New)	M	950-1050	1.4385-1.4395	-
BIOMOUNT NEW LOW	isoparaffin (BioClear New)	M	130-160	1.4345-1.4365	✓
BIOMOUNT	xylene (BioClear)	M	450-650	1.4920-1.4930	✓
BIOMOUNT HIGH	xylene (BioClear)	M	5000-5500	1.4900-1.4920	-
BIOMOUNT M	toluene	M	140-160	1.4900-1.4920	✓
BIOMOUNT C	xylene (BioClear)	M	5000-5500	1.4910-1.4925	-

Legend: M – methacrylate, S – styrene (polystyrene)
 * does not contain harmful dibutyl phthalate plasticizer
 ** contact us to find out how to decrease viscosity of a medium

Xylene and styrene-based media for permanent mounting

BioMount DPX

Moderate viscosity medium for permanent covering of microscope samples and mounting cover glass.

100 mL BM-100
 250 mL BM-250
 500 mL BM-500



BioMount DPX High

High viscosity medium for permanent covering of microscope samples and mounting cover glass.

100 mL BMH-100
 250 mL BMH-250
 500 mL BMH-500



BioMount DPX Low

Low viscosity medium for permanent covering of microscope samples and mounting cover glass. Can be used in automatic cover slip machines.

100 mL BML-100
 250 mL BML-250
 500 mL BML-500



BioMount DPX New

Low viscosity medium for permanent covering of microscope samples and mounting cover glass. Medium contains more environmentally friendly plasticizer instead of harmful dibutyl phthalate. Can be used in automatic cover slip machines.

100 mL BMLE-100
 250 mL BMLE-250
 500 mL BMLE-500





Methacrylate and xylene substitute-based media for permanent mounting

BioMount New

Moderate viscosity medium for permanent covering of microscope samples deparaffinized using xylene substitute-based medium, as well as xylene-based medium.

100 mL	BMN-100
250 mL	BMN-250
500 mL	BMN-500

BioMount New Low

Very low viscosity medium for permanent covering of microscope samples deparaffinized using xylene substitute, as well as xylene.

100 mL	BMNL-100
250 mL	BMNL-250
500 mL	BMNL-500



Methacrylate-based permanent mounting media

BioMount

Low viscosity xylene-based medium for permanent covering of microscope slides and mounting cover glass.

100 mL	BMT-100
250 mL	BMT-250
500 mL	BMT-500

BioMount High

Very high viscosity xylene-based medium for permanent covering of microscope samples and mounting cover glass.

100 mL	BMTH-100
250 mL	BMTH-250
500 mL	BMTH-500

BioMount M

Very low viscosity toluene-based medium for permanent covering of microscope samples and mounting cover glass.

100 mL	BMM-100
250 mL	BMM-250
500 mL	BMM-500



Additional permanent mounting media

BioMount C

Very high viscosity medium for permanent covering of microscope slides and mounting cover glass.

100 mL	BMC-100
250 mL	BMC-250
500 mL	BMC-500

* excellent synthetic substitute for natural covering media, such as Canada Balsam



Water-based mounting media

PRODUCT NAME	SOLVENT	POLYMER	VISCOSITY @ 20°C (cSt, mPa)	REFRACTIVE INDEX @ 20°C (n _p)	USAGE IN AUTOMATED COVER SLIP MACHINES
BIOMOUNT AQUA	water	PVP	700-1000	1.4145-1.4205	-



BioMount Aqua

Water-based medium for covering microscope slides and mounting cover slips. With application in hematology, immunohistochemistry and immunocytochemistry with Fast-Red or AEC chromogens that require an aqueous mounting medium, but also with DAB chromogen.

30 mL	BMA-30
100 mL	BMA-100
250 mL	BMA-250
500 mL	BMA-500



Immersion oil type C

Immersion oil for microscopy analysis with high viscosity of 1500-2200 cSt.

100 mL	IUC-100
250 mL	IUC-250
500 mL	IUC-500

* excellent synthetic substitute for natural immersion oils, such as Cedarwood oil

Immersion oil type FF

Immersion oil for sample analysis in fluorescent microscopy.

100 mL	IUF-100
250 mL	IUF-250
500 mL	IUF-500

Immersion oil type 37

Immersion oil for microscopy analysis with optimal refractive index and viscosity of 1250 cSt at 37°C.

100 mL	IU37-100
250 mL	IU37-250
500 mL	IU37-500

Immersion oil, tropical grade

Immersion oil for microscopy analysis in tropical climates with optimal refractive index at 27-29°C.

10x10 mL	IUT-1010
30 mL	IUT-30
100 mL	IUT-100
250 mL	IUT-250
500 mL	IUT-500



Cleaning solutions for microscope lens and slides

Clean Lens, solution 1

Solution according to Carl Zeiss formulation for cleaning microscope lens. Based on n-hexane and 2-propanol.

100 mL	CL1-OT-100
--------	------------

Clean Lens, solution 3

Solution for cleaning microscope lens and removing immersion oil from stained slide after using microscope. Based on petroleum ether and 2-propanol.

100 mL	CL3-OT-100
--------	------------

Clean Lens, solution 2

Solution for cleaning microscope lens and removing immersion oil from stained slide after using microscope. Based on diethyl ether and 2-propanol.

100 mL	CL2-OT-100
--------	------------



Standard histology staining reagents


Reagents for hematoxylin and eosin (HE) staining, the gold standard for demonstrating cellular and tissue structures, providing detailed anatomy and morphology of frozen and paraffin tissue sections, enabling clear visualization of nuclear, cytoplasmic and other cell structures.

Hematoxylin reagents

Nuclear stains for standard HE staining

Properties of BioGnost's Hematoxylin reagents

PRODUCT NAME	CYTOLOGY	HISTOLOGY	IHC	PROGRESSIVE METHOD	REGRESSIVE METHOD
HEMATOXYLIN H	-	✓	-	✓	✓
HEMATOXYLIN M	-	✓	✓	✓	-
HEMATOXYLIN ML	-	✓	✓	✓	✓
HEMATOXYLIN G1	✓	✓	✓	✓	-
HEMATOXYLIN G2	✓	✓	✓	✓	✓
HEMATOXYLIN G3	✓	✓	-	✓	✓
HEMATOXYLIN CAR	-	✓	✓	✓	-
HEMATOXYLIN COL	-	✓	✓	✓	✓
HEMATOXYLIN HP, PAP 1A	✓	-	-	✓	✓
HEMATOXYLIN HP, PAP 1B	✓	-	-	✓	✓
HEMATOXYLIN W KIT	-	✓	✓	✓	acid resistant
HEMATOXYLIN P.T.A.	-	✓	-	✓	-

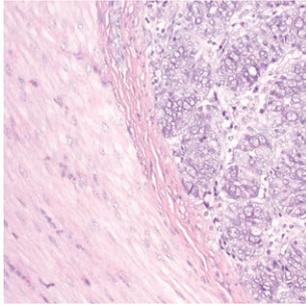

Each of BioGnost's Hematoxylin is specially stabilized in order to maintain the reagent's reactivity for prolonged periods of time, as well as to maintain the quality of the reagent. Every staining provides a crystal clear and highly detailed image of the cell's nucleus, which in turn makes abnormal or pathological changes easier to detect.

Hematoxylin H

Modified hematoxylin according to Harris, most commonly used formulation in HE staining. Strong intensity reagent for progressive and regressive staining in histopathology.

100 mL	HEMH-OT-100
500 mL	HEMH-OT-500
1000 mL	HEMH-OT-1L
2500 mL	HEMH-OT-2.5L

Hematoxylin M

Modified hematoxylin according to Mayer. Weaker intensity reagent for progressive staining in histopathology and counterstaining in immunohistochemistry.

100 mL	HEMM-OT-100
500 mL	HEMM-OT-500
1000 mL	HEMM-OT-1L
2500 mL	HEMM-OT-2.5L

Hematoxylin ML

Modified hematoxylin according to Mayer-Lillie. Strong intensity new generation reagent for progressive and regressive staining in histopathology.

100 mL	HEMML-OT-100
500 mL	HEMML-OT-500
1000 mL	HEMML-OT-1L
2500 mL	HEMML-OT-2.5L

Hematoxylin G1

Hematoxylin according to Gill, ideal for staining goblet cells. New generation reagent for delicate progressive staining in histopathology, cytology and for counterstaining in immunohistochemistry.

100 mL	HEMG1-OT-100
500 mL	HEMG1-OT-500
1000 mL	HEMG1-OT-1L
2500 mL	HEMG1-OT-2.5L

Hematoxylin G2

Hematoxylin according to Gill, contains double concentration of hematoxylin compared to Hematoxylin G1. Moderate intensity new generation reagent for progressive and regressive staining in histopathology, cytology and for counterstaining in immunohistochemistry.

100 mL	HEMG2-OT-100
500 mL	HEMG2-OT-500
1000 mL	HEMG2-OT-1L
2500 mL	HEMG2-OT-2.5L

Hematoxylin G3

Modified hematoxylin contains triple concentration of hematoxylin compared to Hematoxylin G1. Strong intensity new generation reagent for progressive and regressive staining in histopathology and cytology.

100 mL	HEMG3-OT-100
500 mL	HEMG3-OT-500
1000 mL	HEMG3-OT-1L
2500 mL	HEMG3-OT-2.5L





Hematoxylin reagents

Nuclear stains for standard HE staining

Hematoxylin W kit

Acid-resistant hematoxylin according to Weigert. Two-reagent kit that stains the nuclei blue-black, often a component of special staining kits for connective tissues.

2x100 mL	HEMW-K-100
2x500 mL	HEMW-K-500
2x1000 mL	HEMW-K-1L

Hematoxylin P.T.A.

Synonyms: Hematoxylin Phosphotungstic Acid, Mallory P.T.A.H. reagent. Simultaneously stains different types of tissue reddish-brown and blue. For differentiation between fibrin, smooth and cross-striated muscle fibers.

500 mL	HPTA-OT-500
1000 mL	HPTA-OT-1L

Hematoxylin CAR

Modified hematoxylin according to Carazzi for delicate nuclear staining without staining cytoplasmic components. Reagent for progressive staining of frozen sections in histopathology.

100 mL	HEMCA-OT-100
500 mL	HEMCA-OT-500
1000 mL	HEMCA-OT-1L

Hematoxylin COL

Modified hematoxylin according to Cole, contains iodine. Weaker intensity reagent for progressive and regressive staining in histopathology.

100 mL	HEMCO-OT-100
500 mL	HEMCO-OT-500
1000 mL	HEMCO-OT-1L

Hematoxylin Instant kit

Two-component kit that contains ready-to-use powders for preparation of hematoxylin reagent for progressive and regressive nuclear staining in histopathology and cytology.

for 1000 mL	HEMI-K-1L
for 6000 mL	HEMI-K-6L

Kit for rapid HE staining of frozen and paraffin-embedded sections

Ready-to-use eight-reagent kit (in 20 containers that can be used as staining jars) for rapid HE staining of frozen and paraffin tissue sections in histopathology. Contains xylene substitute as clearing agent and xylene substitute-based medium for permanent coverslipping.

for 100 tests	HE-RTU-100T
---------------	-------------

ErioGnost reagent

Acid fast reagent for nuclear staining, comparable to hematoxylin staining. Eco-friendly synthetic replacement for hematoxylin reagents.

100 mL	ERS-OT-100
500 mL	ERS-OT-500
1000 mL	ERS-OT-1L

Differentiation and bluing reagents for histology and cytology

Differentiation of the nucleus enables clear visualization of nuclear structures in regressive procedure, while the bluing reagent turns the nucleus blue.

Acid alcohol

Reagent for enabling crisp visualization of nuclear chromatin detail and removing of excess background dye in histological samples after regressive staining with hematoxylin.

1000 mL	KA-OT-1L
2500 mL	KA-OT-2.5L

Scott's solution

Buffered reagent with potassium carbonate and magnesium sulphate for bluing of nuclear chromatin and nuclear membranes.

1000 mL	SC-OT-1L
2500 mL	SC-OT-2.5L

HCL Pap reagent

Aqueous reagent for enabling crisp visualization of nuclear chromatin detail and removing of excess background dye in cytological after regressive staining with hematoxylin.

1000 mL	HCLPAP-OT-1L
2500 mL	HCLPAP-OT-2.5L

Bluing reagent

Blue colored buffered reagent with lithium carbonate for bluing of nuclear chromatin and nuclear membranes.

1000 mL	BR-OT-1L
2500 mL	BR-OT-2.5L

Eosin reagents

Cytoplasmic stains for standard HE staining

Eosin Y (yellowish eosin) solutions are commonly used as a counterstain to hematoxylin in the standard HE staining method. Eosin Y is an anionic dye that stains basic cell components (such as cytoplasm, collagen, muscle fibers and erythrocytes) bright red.

Eosin Y aqueous solutions

Eosin Y aqueous solutions for cytoplasmic counterstaining.

Eosin Y 0.2% aqueous	1000 mL 2500 mL	EOY-02-OT-1L EOY-02-OT-2.5L
Eosin Y 0.5% aqueous	1000 mL 2500 mL	EOY-05-OT-1L EOY-05-OT-2.5L
Eosin Y 1% aqueous	1000 mL 2500 mL	EOY-10-OT-1L EOY-10-OT-2.5L
Eosin Y 2% aqueous	1000 mL 2500 mL	EOY-20-OT-1L EOY-20-OT-2.5L

Eosin Y alcoholic solutions

Eosin Y alcoholic solutions for cytoplasmic counterstaining.

Eosin Y 0.5% alcoholic	1000 mL 2500 mL	EOYA-05-OT-1L EOYA-05-OT-2.5L
Eosin Y 1% alcoholic	1000 mL 2500 mL	EOYA-10-OT-1L EOYA-10-OT-2.5L
Eosin Y 2% alcoholic	1000 mL 2500 mL	EOYA-20-OT-1L EOYA-20-OT-2.5L

Eosin Contrast

Modified alcoholic solution for cytoplasmic counterstaining. The reagent contains Eosin Y and Phloxine B for enhanced counterstaining effect.

1000 mL 2500 mL	EOYK-OT-1L EOYK-OT-2.5L
--------------------	----------------------------

Eosin Contrast PLUS

Modified alcoholic solution for intensive cytoplasmic counterstaining. Reagent contains Eosin Y, Phloxine B, and Biebrich Scarlet dyes for additional counterstaining effect.

1000 mL 2500 mL	EOYKP-OT-1L EOYKP-OT-2.5L
--------------------	------------------------------

Eosin Y instant

Mixture of powders for preparation of eosin reagent in aqueous or alcoholic medium.

Eosin Y Instant, aqueous	for 1000 mL	EOYI-AQ-1L
Eosin Y Instant, alcoholic	for 1000 mL	EOYI-ALC-1L





Cytology and hematology staining reagents

Cytology and hematology staining kits contain reagents for fixation, staining and buffering cytology and hematology samples, as well as cytology smears of hematopoietic organs, providing precise distinction between cellular structures. Cytochemical staining kits contain reagents for cytochemical diagnosis of leukemia using bone marrow or whole blood smears, and also for staining and determining specific cell components in blood samples.

Papanicolaou (Pap) reagents

Nuclear and cytoplasmic reagents for Pap staining

Due to their simplicity of use and exceptional precision in diagnostics, staining reagents according to Papanicolaou have substantially contributed to decreasing occurrence of cervical cancer. By using a simple staining protocol, BioGnost's Pap reagents enable excellent differentiation between cellular structures.



Hematoxylin HP, Pap 1A

Modified Harris hematoxylin solution for nuclear staining according to Papanicolaou. Strong intensity reagent for progressive and regressive staining in exfoliative cytology.

100 mL	HEMHP-OT-100
500 mL	HEMHP-OT-500
1000 mL	HEMHP-OT-1L
2500 mL	HEMHP-OT-2.5L

Hematoxylin HP, Pap 1B

Hematoxylin modified according to Harris for intensive nuclear staining according to Papanicolaou. Reagent for progressive and regressive staining in exfoliative cytology, with formulation stronger than Hematoxylin HP, Pap 1A.

100 mL	HEMHP1B-OT-100
500 mL	HEMHP1B-OT-500
1000 mL	HEMHP1B-OT-1L
2500 mL	HEMHP1B-OT-2.5L

OG-6 reagent, Pap 2A

Monochromatic reagent for counterstaining mature and keratinized squamous cell. Contains Orange G dye that stains the cytoplasm orange.

100 mL	OG6-OT-100
500 mL	OG6-OT-500
1000 mL	OG6-OT-1L
2500 mL	OG6-OT-2.5L

Orange II reagent, Pap 2B

Monochromatic reagent for counterstaining mature and keratinized squamous cells. Contains Orange II dye that stains the cytoplasm intense red.

100 mL	OR2-OT-100
500 mL	OR2-OT-500
1000 mL	OR2-OT-1L
2500 mL	OR2-OT-2.5L



Rapid and standard hematology reagents and kits

Bio-Diff kits

Rapid staining solutions for hematology, cytology and histology

Rapid diagnostic kits that provide precise staining results in a short period of time, such as results provided by May-Gruenwald Giemsa (MGG) staining method. Each component of Bio-Diff and Bio-Diff RTU kits is stabilized and prepared according to the highest standards.

Bio-Diff kit

Three-reagent kit that contains fixative agent, red and blue components for fast and effective staining. Each kit contains buffer tablets for consistent staining results.

3x100 mL	BD-K-100
3x500 mL	BD-K-500
3x1000 mL	BD-K-1L
3x2500 mL	BD-K-2.5L



Bio-Diff 1 reagent

Fixative agent, component of Bio-Diff kit.

100 mL	BD1-OT-100
500 mL	BD1-OT-500
1000 mL	BD1-OT-1L
2500 mL	BD1-OT-2.5L

Bio-Diff 2 reagent

Buffered red component of Bio-Diff kit.

100 mL	BD2-OT-100
500 mL	BD2-OT-500
1000 mL	BD2-OT-1L
2500 mL	BD2-OT-2.5L

Bio-Diff 3 reagent

Buffered blue component of Bio-Diff kit.

100 mL	BD3-OT-100
500 mL	BD3-OT-500
1000 mL	BD3-OT-1L
2500 mL	BD3-OT-2.5L

STAINING DONE IN
14 SECONDS*



* Depending on the type of sample.

Bio-Diff RTU kit

Ready-to-Use three-reagent kit with reagents in containers that can be used as staining jars. Kit contains fixative agent, red and blue components for fast and effective staining and buffer tablet for consistent staining results.

3x100 mL BD-RTU-100



Bio-Diff 1 RTU reagent

Fixative component Bio-Diff RTU kit.

100 mL BD1-RTU



Bio-Diff 2 RTU reagent

Buffered red component of Bio-Diff RTU kit.

100 mL BD2-RTU

Bio-Diff 3 RTU reagent

Buffered blue component of Bio-Diff RTU kit.

100 mL BD3-RTU

STAINING DONE IN
14 SECONDS*



* depending on the type of sample

Buffer solution, pH 7.2

Phosphate buffer solution, molarity 6 mM for use in hematology and histology.

1000 mL PUF72-OT-1L
2500 mL PUF72-OT-2.5L
5000 mL PUF72-OT-5L

Buffer tablets pH 7.2

Phosphate buffer tablets for preparing buffer solutions used in hematology and histology. One buffer tablet is sufficient for preparation of 1000 mL of phosphate buffer solution.

50 pcs PT-72-50
100 pcs PT-72-100



Use buffer solution pH 6.8 for hematological and cytological diagnostics. Use buffer solution pH 7.2 for the detection of blood parasites.



May-Gruenwald and Giemsa reagents

Romanowsky staining solutions for hematology, cytology and histology

One of the most commonly used Pappenheim methods is May-Gruenwald Giemsa (MGG) for staining hematology and cytology materials, as well as cytology smears of hematopoietic organs. Giemsa solution is also used for staining histology sections.

May-Gruenwald solution

Polychromatic fixative solution containing Eosin Y, Methylene Blue and azure dyes. For staining in hematology, cytology and staining sections of hematopoietic organs in histopathology.

100 mL	MG-OT-100
10x100 mL	MG-OT-110
500 mL	MG-OT-500
1000 mL	MG-OT-1L
2500 mL	MG-OT-2.5L



Giemsa solution

Polychromatic solution of Eosin Y, Methylene Blue and azure dyes. For use in hematology, cytology, histology and for staining section samples of hematopoietic organs.

100 mL	GM-OT-100
10x100 mL	GM-OT-110
500 mL	GM-OT-500
1000 mL	GM-OT-1L
2500 mL	GM-OT-2.5L



Buffer solution, pH 6.4

Phosphate buffer solution, molarity 6 mM for use in hematology and histology.

1000 mL	PUF64-OT-1L
2500 mL	PUF64-OT-2.5L
5000 mL	PUF64-OT-5L

Buffer solution, pH 6.8

Phosphate buffer solution, molarity 6 mM for use in hematology and histology.

1000 mL	PUF68-OT-1L
2500 mL	PUF68-OT-2.5L
5000 mL	PUF68-OT-5L

Buffer tablets pH 6.8

Phosphate buffer tablets for preparing solutions used in hematology and histology. One buffer tablet is sufficient for preparation of 1000 mL of phosphate buffer solution.

50 pcs	PT-68-50
100 pcs	PT-68-100

Use buffer solution pH 6.8 for hematological and cytological diagnostics. Use buffer solution pH 7.2 for the detection of blood parasites.



LEUKOGNOST

Cytochemical staining kits and ancillary reagents

Ancillary reagents

Ancillary reagents greatly complement LeukoGnost kits - can be purchased separately or as a kit (LeukoGnost PLUS) containing a fixative, hematoxylin stain and a water-based mounting medium.

LeukoGnost PLUS

Additional set of reagents for use with LeukoGnost kits. Contains fixative, hematoxylin reagent and water-based mounting medium.

for at least 100 tests LKG-PLUS

LeukoGnost Fixative

Fixative for use in cytochemical diagnostics of leukemia

500 mL LKF-500

LeukoGnost Fixative ECO

Formalin-free fixative for use in cytochemical diagnostics of leukemia.

500 mL LFEC-500

LeukoGnost Fixative GL

Fixative for use in cytochemical diagnostics of leukemia, based on A formalin substitute.

500 mL LFGL-500

LeukoGnost HEM

Hematoxylin reagent for use with LeukoGnost kits.

100 mL	LKH-OT-100
250 mL	LKH-OT-250
500 mL	LKH-OT-500

BioMount Aqua

Aqueous covering medium for microscope sections and mounting medium for cover glasses.

30 mL	BMA-30
100 mL	BMA-100
250 mL	BMA-250
500 mL	BMA-500





Bacteriology and mycology staining reagents

Reagents for rapid and easy staining in microbiology allowing optimal differential analysis of bacteria. Phenol-free staining kits offer environment-friendly alternative without compromising quality of staining. Reagents for fungal analysis enable visualization of fungal structures (otherwise difficult to observe) and samples preservation.

BioGram kits and reagents

Kits for identification of bacteria according to Gram

Differential staining of bacteria according to Gram is widely known and used method in microbiology. BioGram 4 kit uses standard method to provide optimal staining and identification of Gram-positive and Gram-negative bacteria. BioGram ECO kit is a phenol-free, eco-friendly alternative to BioGram 4 kit with equal staining results.

BioGram 4 kit

Four-reagent kit for identification of bacteria according to Gram. Kit contains Gram Crystal Violet 1% solution, stabilized Gram Lugol solution, double amount of Gram Decolorizer solution 2 and Gram Safranin solution as counterstain.

for 100 tests	BGR4-100T
5x100 mL	BGR4-K-100
5x250 mL	BGR4-K-250
5x500 mL	BGR4-K-500

Gram Crystal Violet 1% solution

Primary stain reagent for Gram-positive bacteria. Crystal Violet is also known as Gentian Violet and Methyl Violet 10B.

100 mL	GC1-OT-100
250 mL	GC1-OT-250
500 mL	GC1-OT-500
1000 mL	GC1-OT-1L
2500 mL	GC1-OT-2.5L

Gram Lugol solution, stabilized

Solution contains polyvinylpyrrolidone iodine (PVP-I).

100 mL	GLS-OT-100
250 mL	GLS-OT-250
500 mL	GLS-OT-500
1000 mL	GLS-OT-1L
2500 mL	GLS-OT-2.5L

Gram Decolorizer solution 2

Medium speed decolorizer used after stabilized iodine solution.

100 mL	GD2-OT-100
250 mL	GD2-OT-250
500 mL	GD2-OT-500
1000 mL	GD2-OT-1L
2500 mL	GD2-OT-2.5L

Gram Safranin solution

Counterstain for Gram-negative bacteria, contains Safranin O dye.

100 mL	GSF-OT-100
250 mL	GSF-OT-250
500 mL	GSF-OT-500
1000 mL	GSF-OT-1L
2500 mL	GSF-OT-2.5L

BioGram ECO kit PHENOL-FREE

Four-reagent phenol-free kit for identification of bacteria according to Gram. Kit contains Gram Crystal violet, phenol-free reagent, Gram Sodium hydrogencarbonate, solution, stabilized Gram Lugol solution, double amount of Gram Decolorizer solution 2 and Gram Safranin solution as counterstain.

2x50 mL+4x100 mL	BGRE-K-100
2x250 mL+4x500 mL	BGRE-K-500

Gram Crystal Violet, phenol-free reagent

Primary stain reagent for Gram-positive bacteria, doesn't contain phenol. Crystal Violet is also known as Gentian Violet and Methyl Violet 10B.

100 mL	GCF-OT-100
250 mL	GCF-OT-250
500 mL	GCF-OT-500
1000 mL	GCF-OT-1L
2500 mL	GCF-OT-2.5L

Gram Sodium hydrogencarbonate, solution

Solution of sodium hydrogencarbonate for use with BioGram ECO kit.

100 mL	GNHK-OT-100
250 mL	GNHK-OT-250
500 mL	GNHK-OT-500
1000 mL	GNHK-OT-1L
2500 mL	GNHK-OT-2.5L

Additional bacteriology staining solutions according to Gram

Gram Lugol solution

Non-stabilized Lugol solution.

100 mL	GL-OT-100
250 mL	GL-OT-250
500 mL	GL-OT-500
1000 mL	GL-OT-1L

Gram Decolorizer solution 1

Slow acting decolorizer solution used after non-stabilized iodine solution.

100 mL	GD1-OT-100
250 mL	GD1-OT-250
500 mL	GD1-OT-500
1000 mL	GD1-OT-1L

Gram Decolorizer solution 3

Fast acting solution used with automated stainers.

100 mL	GD3-OT-100
250 mL	GD3-OT-250
500 mL	GD3-OT-500
1000 mL	GD3-OT-1L

Gram Safranin 5x concentrated

Concentrated counterstain solution of Safranin O dye.

For preparation of 500 mL of Safranin O counterstain solution	100 mL	SO05X-OT-100
For preparation of 2500 mL of Safranin O counterstain solution	500 mL	SO05X-OT-500
For preparation of 5000 mL of Safranin O counterstain solution	1000 mL	SO05X-OT-1L

Carbol Fuchsin Gram 10x concentrated

Concentrated Carbol Fuchsin solution for preparation of counterstain solution for modified staining. Also known as concentrated Gram Fuchsin solution.

For preparation of 1000 mL of Carbol Fuchsin Gram counterstain solution	100 mL	CFGR10-OT-100
For preparation of 5000 mL of Carbol Fuchsin Gram counterstain solution	500 mL	CFGR10-OT-500
For preparation of 10000 mL of Carbol Fuchsin Gram counterstain solution	1000 mL	CFGR10-OT-1L

Additional iodine solutions

Lugol's solution

Classic formulation of aqueous solution of iodine and potassium iodide.

100 mL	LUG-OT-100
250 mL	LUG-OT-250
500 mL	LUG-OT-500

Lugol's C solution

Aqueous solution of iodine and potassium iodide, enhanced formulation with stronger staining (compared to classic Lugol solution).

100 mL	LUGC-OT-100
250 mL	LUGC-OT-250
500 mL	LUGC-OT-500

Lugol's solution, acid

Aqueous iodine and potassium iodide solution with added acid.

100 mL	LUGK-OT-100
250 mL	LUGK-OT-250
500 mL	LUGK-OT-500



TB-Stain kits and reagents

Kits for staining acid-fast bacteria according to Ziehl-Neelsen

Most commonly known method for detecting TB bacteria according to Ziehl-Neelsen. This method enables rapid and easy staining of acid-fast bacterial strains, such as *Mycobacteria*, *Nocardia*, *Rhodococcus*, *Dietzia* and *Gordonia*. TB-Stain ECO kit is a phenol-free, environment-friendly alternative to TB-Stain Hot kit.

TB-Stain Hot kit

Three-reagent kit for staining acid-fast bacteria. Contains TB Carbol Fuchsin reagent, double amount of TB Decolorizer and Methylene Blue Loeffler's reagent as counterstain.

for 100 tests	TBH-100T
4x100 mL	TBH-K-100
4x250 mL	TBH-K-250
4x500 mL	TBH-K-500



TB Carbol Fuchsin reagent

Synonym: Ziehl-Neelsen reagent. Primary stain reagent for acid-fast bacteria used in TB-Stain Hot, TB-Stain Cold and TB-Stain Quick kits.

100 mL	TBC-OT-100
250 mL	TBC-OT-250
500 mL	TBC-OT-500
1000 mL	TBC-OT-1L
2500 mL	TBC-OT-2.5L



TB-Stain ECO kit PHENOL-FREE

Three-reagent phenol-free kit for staining acid-fast bacteria. Contains TB-Fuchsin reagent, double amount of TB Decolorizer and Methylene Blue Loeffler's reagent as counterstain.

4x100 mL	TBE-K-100
4x250 mL	TBE-K-250
4x500 mL	TBE-K-500



TB Decolorizer

For use in TB-Stain Hot and TB-Stain Cold kits.

100 mL	TBD-OT-100
250 mL	TBD-OT-250
500 mL	TBD-OT-500
1000 mL	TBD-OT-1L
2500 mL	TBD-OT-2.5L



Methylene Blue Loeffler's reagent

Blue counterstain for TB-Stain Hot kit.

100 mL	MBL-OT-100
250 mL	MBL-OT-250
500 mL	MBL-OT-500
1000 mL	MBL-OT-1L
2500 mL	MBL-OT-2.5L

TB Fuchsin reagent

Primary stain **phenol-free** reagent for acid-fast bacteria used in TB-Stain ECO kit.

100 mL	TBFR-OT-100
250 mL	TBFR-OT-250
500 mL	TBFR-OT-500
1000 mL	TBFR-OT-1L
2500 mL	TBFR-OT-2.5L

Additional TB-Stain Hot kit staining solution

Methylene Blue 10x concentrated

Concentrated counterstain solution of Methylene Blue dye for use in microbiology.

For preparation of 1000 mL of Methylene Blue counterstain solution	100 mL	MB10X-OT-100
--	--------	--------------

For preparation of 5000 mL of Methylene Blue counterstain solution	500 mL	MB10X-OT-500
--	--------	--------------

For preparation of 10000 mL of Methylene Blue counterstain solution	1000 mL	MB10X-OT-1L
---	---------	-------------





Kits for detection of acid-fast bacteria using fluorescence microscopy according to modified Truant method

TB-Stain Fluorescent kit

Three-reagent kit for detection of acid-fast bacteria using fluorescence microscopy. Contains TB Auramine-Rhodamine reagent, double amount of TB Decolorizer Fluorescent and TB Permanganate reagent as counterstain.

for 100 tests	TBF-100T
4x100 mL	TBF-K-100
4x250 mL	TBF-K-250
4x500 mL	TBF-K-500



TB-Stain Auramine O kit

Three-reagent kit for staining acid-fast bacteria using fluorescence method. Contains TB Auramine O reagent, double amount of TB Decolorizer Fluorescent and counterstain of TB Permanganate reagent.

4x100 mL	TBAO-K-100
4x250 mL	TBAO-K-250
4x500 mL	TBAO-K-500



TB Auramine O reagent

Primary stain for use in TB Stain Auramine O kit.

100 mL	TBAO-OT-100
250 mL	TBAO-OT-250
500 mL	TBAO-OT-500
1000 mL	TBAO-OT-1L
2500 mL	TBAO-OT-2.5L

TB Auramine-Rhodamine reagent

Primary stain reagent for use in TB-Stain Fluorescent kit.

100 mL	TBA-OT-100
250 mL	TBA-OT-250
500 mL	TBA-OT-500
1000 mL	TBA-OT-1L
2500 mL	TBA-OT-2.5L



TB Decolorizer Fluorescent

For use in TB-Stain Fluorescent and TB-Stain Auramine O kits.

100 mL	TBF-OT-100
250 mL	TBF-OT-250
500 mL	TBF-OT-500
1000 mL	TBF-OT-1L
2500 mL	TBF-OT-2.5L



TB Permanganate reagent

Counterstain for use in TB-Stain Fluorescent kit and TB-Stain Auramine O kits.

100 mL	TBP-OT-100
250 mL	TBP-OT-250
500 mL	TBP-OT-500
1000 mL	TBP-OT-1L
2500 mL	TBP-OT-2.5L



Solutions for staining and analysis of fungi

In order to detect and achieve better visualization of fungi and fungal structures that are otherwise difficult to observe, fungi can be stained using Lactophenol solutions. Using Myco-Mount medium they can also be stained and preserved as a permanent preparation.

Lactophenol Blue solution

For use in mycology for staining and analysis of fungi.

100 mL	LPB-OT-100
250 mL	LPB-OT-250



Lactophenol solution

Colorless solution for use in mycology and analysis of fungi.

100 mL	LP-OT-100
250 mL	LP-OT-250



For detection and fungal staining of histology sections, please see page 44 for more details on Grocott kit.

Myco-Mount Blue

Medium for simultaneous staining of fungi blue and preparation of permanent slides.

30 mL	MYMB-30
100 mL	MYMB-100



Myco-Mount Red

Medium for simultaneous staining of fungi red and preparation of permanent slides.

30 mL	MYMR-30
100 mL	MYMR-100



Potassium hydroxide, solution for fungi

Solution for use in mycology for analysis and visualization of fungi.

30 mL	KHOG-OT-30
100 mL	KHOG-OT-100





Special staining kits for histology and cytology

Special kits for histology and cytology enable clear visualization of specific structures and cellular components. All of the reagents are easy to use and precisely formulated to maximize usability and to provide optimal results for the user.

Classification of special staining kits according to use:

MUSCLE AND/OR CONNECTIVE TISSUE

- A.F.O.G kit
- Azan Trichrome kit
- Elastica-Van Gieson kit
- Gomori Trichrome kit
- H.B.F.P. kit
- Hematoxylin P.T.A. kit
- Hematoxylin W kit
- Ladewig's Trichrome kit
- Mallory Trichrome kit
- Martius Scarlet Blue (MSB) kit
- Masson Trichrome kit
- Masson-Goldner Trichrome kit
- Movat kit
- Orcein kit
- Paraldehyde Fuchsin kit
- P.A.S.M. / Jones kit
- Picro-Sirius Red kit
- Safranin O kit
- Reticulin kit
- Reticulin Contrast kit
- Van Gieson Trichrome kit
- Verhoeff kit
- Weigert-Van Gieson kit

ENZYME ACTIVITY

- LeukoGnost SPE

FUNGI, BACTERIA AND PARASITES

- Alcian Yellow Toluidine Blue kit
- BioGram Histo kit
- Field kit
- Giemsa HP kit
- Grocott kit
- M.I.F. kit
- TB-Stain Histo kit
- Warthin Starry kit

LIPIDS

- Oil Red O kit
- Sudan Black B Lipid kit

NUCLEI AND NUCLEIC ACIDS

- Feulgen kit

PIGMENTS AND MINERAL DEPOSITS

- Fouchet-Van Gieson kit
- Grimelius kit
- HemoGnost Perls kit
- Masson Fontana kit
- Orcein kit
- Perls-Van Gieson kit
- Rhodanine kit
- Von Kossa kit

CYTOLOGY

- Bio-Diff RTU kit
- Eosin and nigrosin vital kit
- Sperm-Diff RTU kit
- UriGnost S kit
- UriGnost SM kit

CARBOHYDRATES


- Alcian Blue - P.A.S. kit
- Alcian Blue pH 1.0 kit
- Alcian Blue pH 2.5 kit
- Colloidal Iron kit
- Congo Red Highman kit
- Congo Red Puchtler kit
- Mucicarmine kit
- P.A.S. kit
- P.A.S. Diastase kit

CRYOSTAT

- Gomori Trichrome kit
- Oil Red O kit
- Orcein kit
- Sudan Black B Lipid kit

CENTRAL NERVOUS SYSTEM


- Bielschowsky kit
- Luxol Fast Blue kit

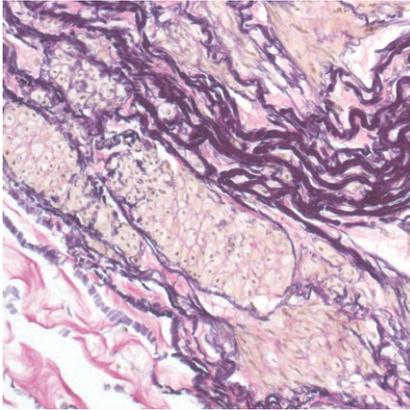

Elastica-Van Gieson kit

Four-reagent kit for staining elastic fibers and differentiation between elastic tissue, collagen and other types of connective tissue. The rapid method enables a satisfactory result with shorter section staining time.


for 100 tests WVGB-100T
4x100 mL WVGB-K-100

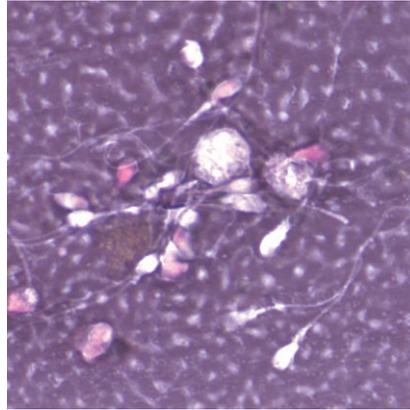

Eosin and nigrosine vital kit

Two-reagent kit containing separate dyes for rapid detection of sperm vitality and simple visualization of dead and living sperm cells. The nigrosin stain provides dark background for easier recognition of both viable and non-viable spermatozoa.


2x30 mL ENVK-K-30
2x100 mL ENVK-K-100

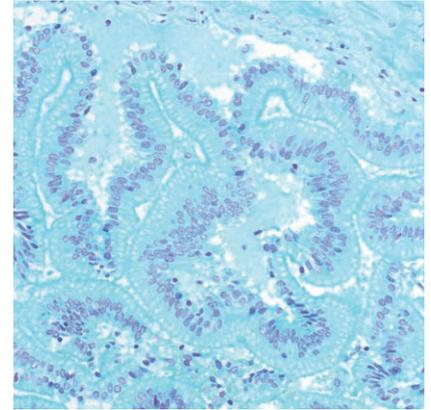

Feulgen kit

Five-reagent DNA staining kit according to Feulgen. For use in semiquantitative DNA determination in histological and cytological samples. Specimen is first treated with hydrochloric acid creating an aldehyde group of DNA that can be visualized by Schiff (BioSchiff) reagent. This reaction is specific for nuclear DNA.

for 100 tests FE-100T
6x100 mL FE-K-100

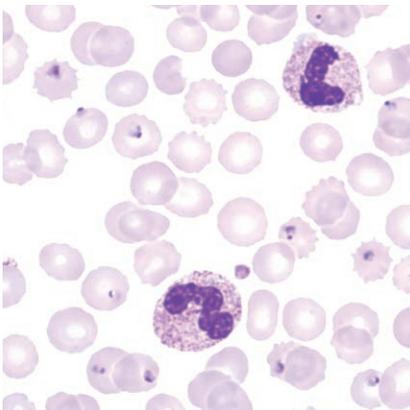

Field kit

Ready-to-use two-reagent kit for rapid and efficient staining and detection of parasites in hematology samples. Primarily used for staining thin and thick blood smears (dense drop) for purpose of diagnosing blood parasites. Reagents are stored in containers that can be used as staining jars.


for 100 tests FD-100T
2x500 mL FD-K-500



Fouchet-Van Gieson kit

Three-reagent kit for visualization of bilirubin and collagen according to Kutllick. Bilirubin is a yellow-brown pigment, but changes to green due to oxidation induced by Fouchet reaction. Green bilirubin can easily be detected on yellowish and pink colored background.

for 100 tests FVG-100T
3x100 mL FVG-K-100

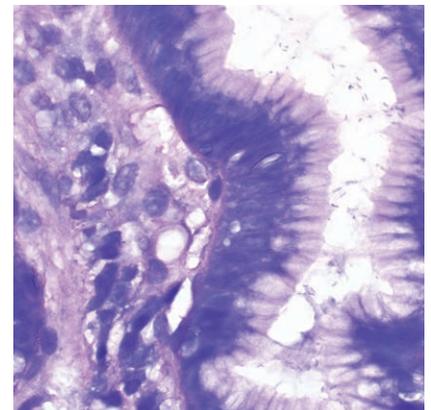


Giemsa HP kit

Four-reagent kit for staining *Helicobacter pylori* in gastroscopic sections according to Lennart. Advantages of this method for detecting *H. pylori* are sensitive and reproducible results and easy performance.

for 100 tests GMHP-100T




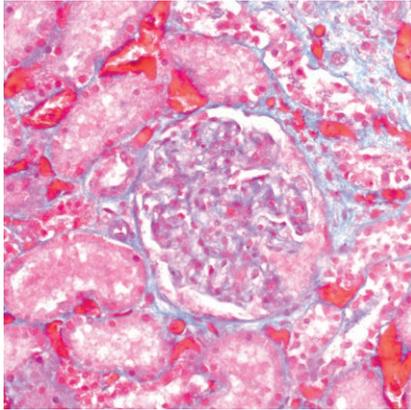

Gomori Trichrome kit

Five-reagent kit for staining muscle, collagen fibre and nuclei, contains blue counterstain. Kit can be used to contrast skeletal, cardiac or smooth muscle.


for 100 tests GT-100T
5x100 mL GT-K-100
* Green counterstain also available

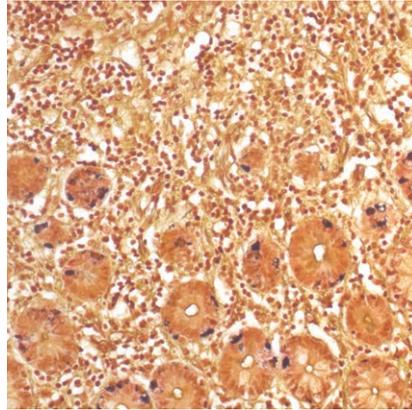

Grimelius kit

Five-reagent kit for staining argyrophilic granules. Grimelius kit can be used for detection of secretory intracytoplasmic granules specific for carcinoid tumors and for identification of neuroendocrine cells.

for 100 tests GRI-100T
5x100 mL GRI-K-100

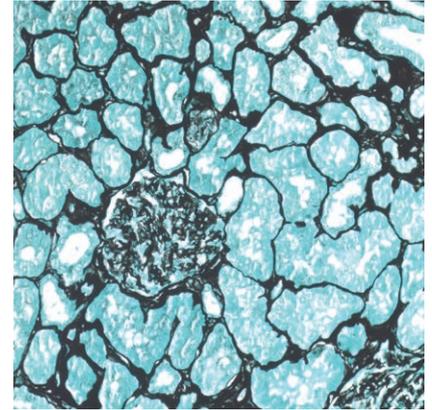


Grocott kit, stabilized

Seven-reagent kit for visualization of fungi and histological argentaaffin structures in general (such as basal membranes). Green counterstain provides clear and visually rich contrast to target structures stained black.

for 100 tests GRC-100T
300-350 tests GRC-K-100


Stable at room temperature

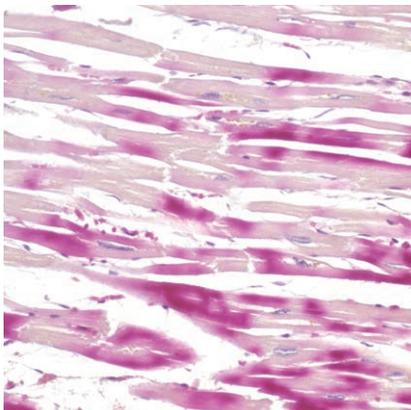

H.B.F.P. kit

Three-reagent Hematoxylin-Basic Euchsin-Picric acid staining kit for detection of cardiac muscle changes after ischemia or myocardial infarction. H.B.F.P. kit is a non-enzymatic histochemical technique used for identification of early myocardial ischemia occurring on the cardiac muscle during the first six hours of symptoms.


for 100 tests HBFP-100T
3x100 mL HBFP-K-100

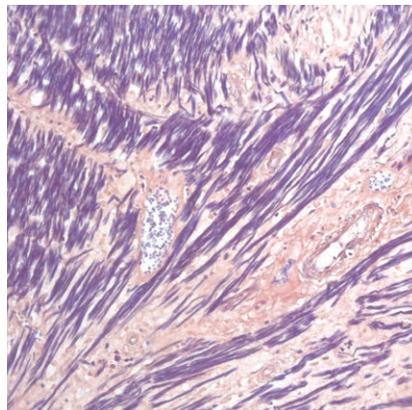

Hematoxylin P.T.A. kit

Four-reagent Hematoxylin-Phosphotungstic Acid kit for differentiation of smooth and striated muscle tissues as well as for detection of fibrin, collagen and elements of central nervous system according to Mallory.

for 100 tests HPTA-100T
4x100 mL HPTA-K-100

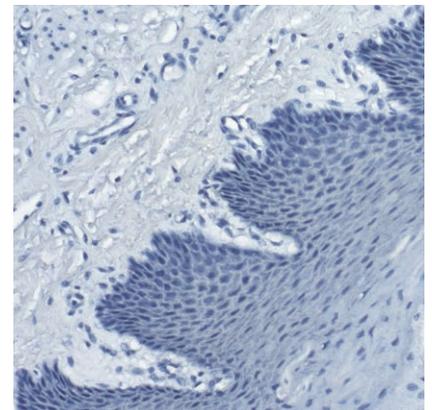


Hematoxylin W kit

Acid-resistant hematoxylin according to Weigert. Two-reagent kit that stains the nuclei blue-black, often a component of special staining kits for connective tissues.

for 100 tests HEMW-100T
2x100 mL HEMW-K-100
2x500 mL HEMW-K-500
2x1000 mL HEMW-K-1L




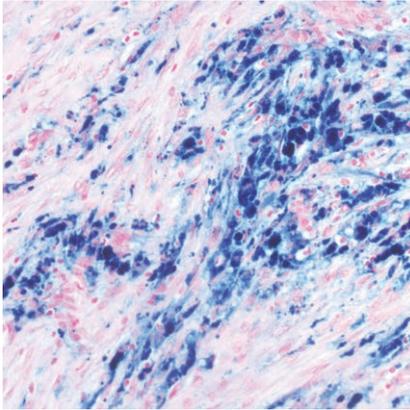

HemoGnost Perls kit

Three-reagent HemoGnost Perls (Prussian blue / Berlin blue) kit for the detection of reactive ferric (Fe³⁺) ions in tissue cells (like spleen or endometrial tissue). When staining bone marrow or other hematology tissue, use LeukoGnost PERLS kit (page 32).


for 100 tests HP-100T
3x250 mL HP-K-250

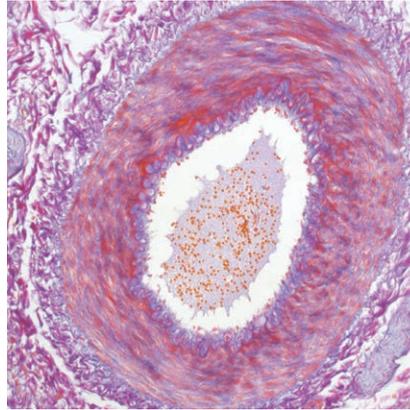

Ladewig Trichrome kit

Four-reagent kit for staining muscle (brown to red) and collagen fibers with blue counterstain. Staining procedure is easier and faster compared to Masson and Mallory Trichrome staining procedures.


for 100 tests LDW-100T
4x100 mL LDW-K-100

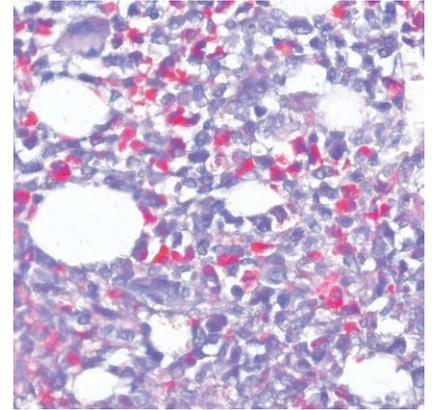

LeukoGnost SPE kit

Kit for detection of specific esterase activity (Leder stain, chloroacetate esterase) on formalin-fixed, paraffin-embedded (FFPE) samples of bone marrow tissue.

for at least 100 tests LKG-SPE

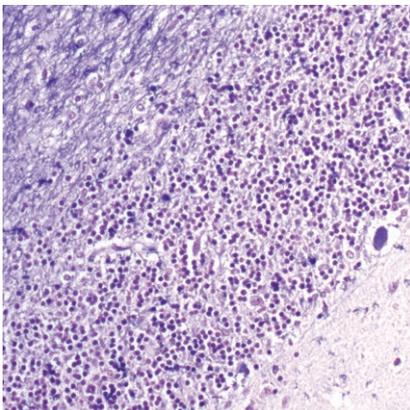

Luxol Fast Blue kit

Three-reagent kit for staining myelin and myelinated axons, Nissl bodies and phospholipids according to Kluver-Barrera. Kit is used for identifying the basic neuronal brain structure or spinal cord sections.


for 100 tests LFB-100T
3x100 mL LFB-K-100



M.I.F. kit

Two-reagent Merthiolate-Iodine-Formalin kit for fixing and staining fecal parasites (especially protozoae, cysts, helminth eggs and larvae). Stool sample is fixed with formalin and stained with two coloring agents; iodine and Eosin Y.


for 100 tests MIF-100T
for 200 tests MIF-200T

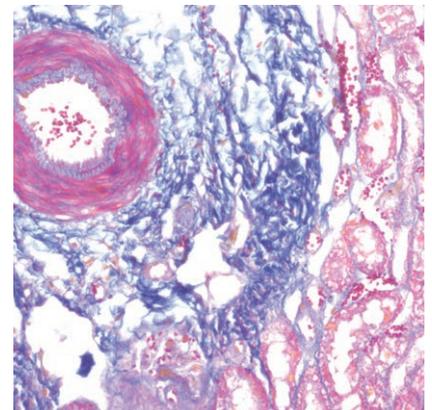

Mallory Trichrome kit

Three-reagent staining kit for connective tissue visualization and detection of collagen, cartilage, muscle, elastic fibers, mucous, pituitary cells, reticulum, bones, amyloid and erythrocytes.


for 100 tests MT-100T
3x100 mL MT-K-100

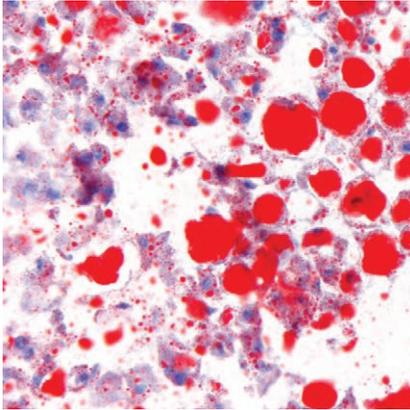

Oil Red O kit

Four-reagent kit for selective staining and detection of fat cells and neutral fats according to Johnson. Can be used with frozen sections and fresh smears to detect obesity-linked pathologies such as dyslipidemia and diabetes.


for 100 tests ORO-100T
4x250 mL ORO-K-250

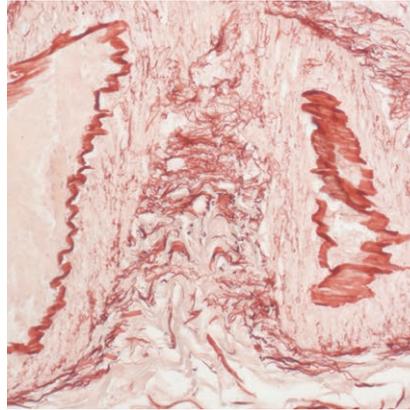

Orcein kit

Five-reagent kit for visualization of hepatitis B surface antigen (HBsAg) seen as viral inclusion bodies in hepatocytes, for elastic fibers and copper associated protein in tissue sections. Can be used with frozen sections.


for 100 tests ORC-100T
5x250 mL ORC-K-250

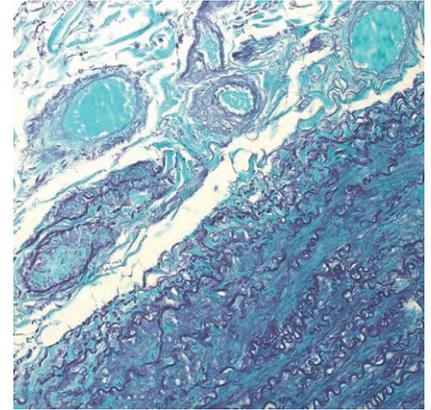

Paraldehyde Fuchsin kit

Seven-reagent kit according to Gomori for detecting pathological changes in elastic fibers. It also stains mast cell granules, beta granules in pancreatic islets, neurosecretory material, mast cell granules and beta cells in the pituitary gland.


for 100 tests PAF-100T
9x100 mL PAF-K-100

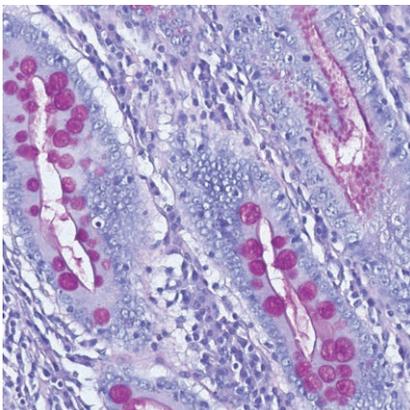

P.A.S. kit

Five-reagent Periodic Acid-Schiff kit for staining aldehydes, mucopolysaccharides, mucoproteins and lymphocytes according to Hotchkiss-McManus. P.A.S. staining may also be used for the demonstration of fungal organisms in tissue sections.


for 100 tests PAS5-100T
5x100 mL PAS5-K-100
5x500 mL PAS5-K-500

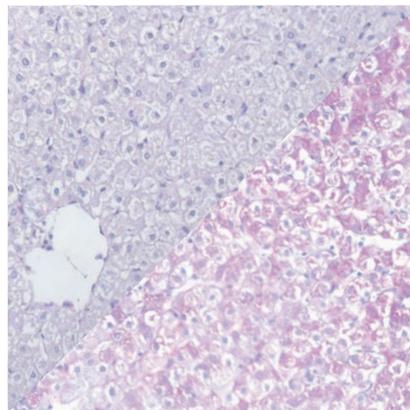

P.A.S. Diastase kit

Five-reagent kit for detecting glycogen and mucopolysaccharide structures mostly in liver. The diastase digestion procedure, followed by P.A.S. staining, is useful in the evaluation of glycogen storage disease.

for 100 tests PDIA-100T

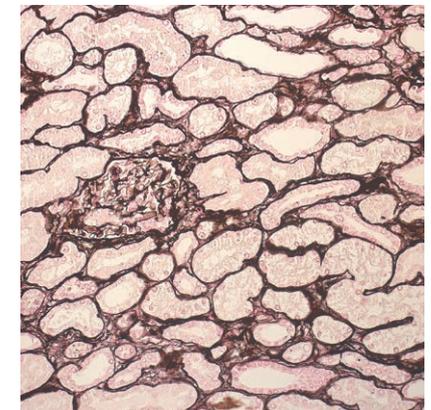


P.A.S.M. / Jones kit, stabilized

Seven-reagent Periodic Acid Silver Methenamine kit for staining kidney glomerular basement membranes. Kit includes red counterstain which provides visually rich contrast to target structures stained black.

for 100 tests PASM-100T
300-350 tests PASM-K-100


Stable at room temperature

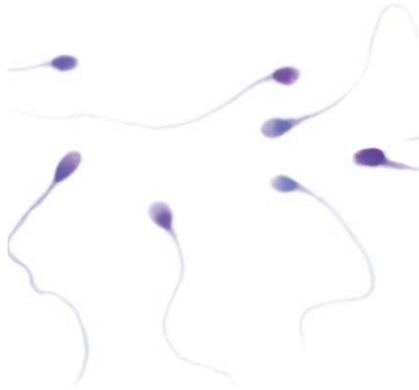

Sperm-Diff RTU kit

Ready-to-Use three-reagent kit with reagents stored in containers that can be used as staining jars. Kit contains fixative reagent, red and blue components for fast and effective staining.


3x100 mL SP-RTU-100

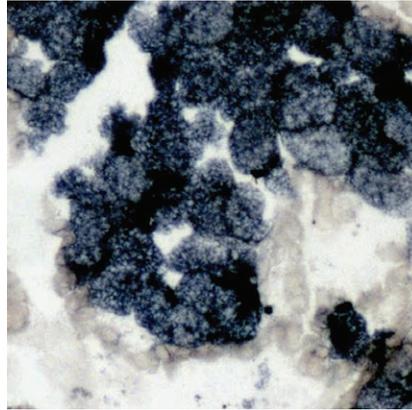

Sudan Black B Lipid kit

Four-reagent kit for specific lipid staining in cytochemistry. Sudan Black dye stains a few types of lipids, including neutral fats, phospholipids and sterols. Contains double amount of Sudan Black B, solution.


for 100 tests SBB-100T
6x100 mL SBB-K-100

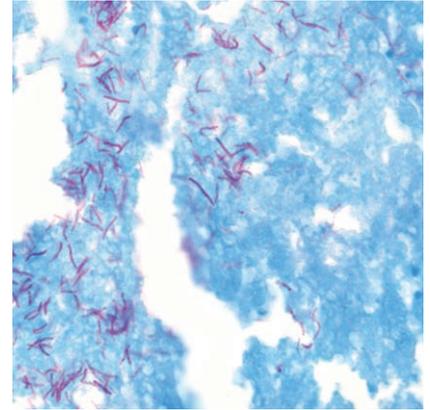

TB-Stain Histo kit

Three-reagent kit for staining acid-fast bacteria (pathogenic mycobacteria) in histology sections, sputum, smears and culture smears according to Ziehl-Neelsen. Heating of the carbol-fuchsin solution is avoided in this protocol hence the release of hazardous phenolic vapors is omitted.

for 100 tests TBHK-100T

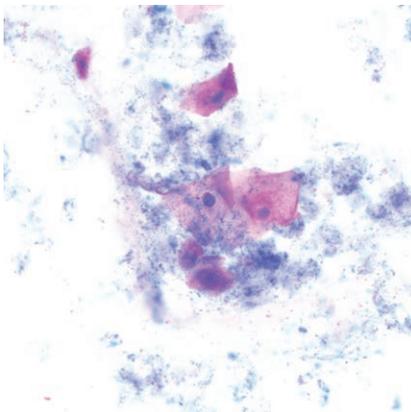

UriGnost S kit

Kit for sampling, staining, and microscopic analysis of urine sediment. Contains modified reagent according to Sternheimer (reagent manufactured according to the European Confederation of Laboratory Medicine (ECLM) and all the necessary equipment for sampling, concentrating, staining, and counting cells and kidney cylinders, as well as urine sediment analysis (test tubes with retentive bottom, Eppendorf 200 µL pipette tips and urine plates).

for 500 tests USK-500

* for more information, see page 57




UriGnost SM kit

Kit for sampling, staining and microscopic analysis of urine sediments. Contains modified reagent according to Sternheimer-Malbin and equipment required for sampling, concentration, staining, counting of kidney cells and kidney cylinders and urine sediment analysis (test tubes with retentive bottom, Eppendorf 200 µL pipette tips and urine plates).

for 500 tests USMK-500

* for more information, see page 57

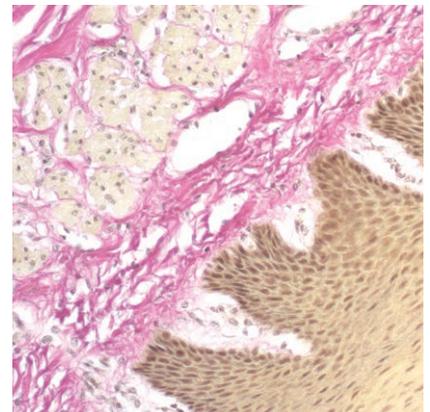


Van Gieson Trichrome kit

Three-reagent kit for staining collagen fibers, muscle tissue, keratinized epithelium, cytoplasm, glial fibers and erythrocytes. Used for differentiation between collagen and smooth fibers in tumors and various other diseases.

for 100 tests VG-100T
3x100 mL VG-K-100
3x500 mL VG-K-500




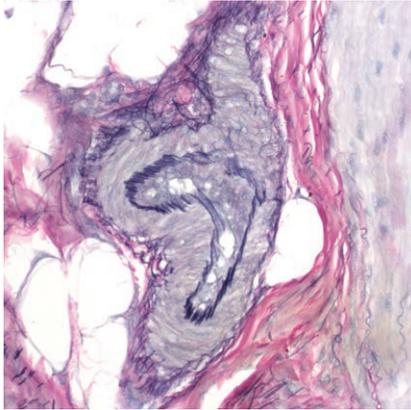

Verhoeff kit

Six-reagent kit for detecting atrophy of elastic tissue in cases of emphysema, thinning and loss of elastic fibers in arteriosclerosis and other vascular diseases, or whether blood vessels have been invaded by a tumor.


for 100 tests VER-100T
6x100 mL VER-K-100

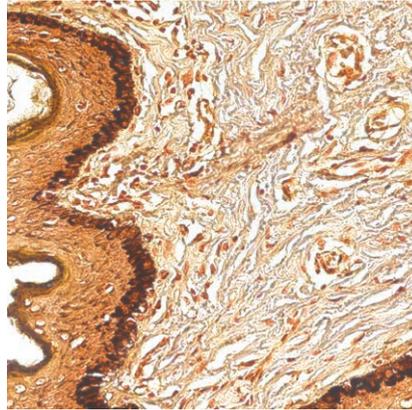

Von Kossa kit

Five-reagent kit for simple and reproducible detection of calcium deposits and calcium salt in tissue samples according to Von Kossa. Tissue calcification is associated with metabolic problems in various tissues (like bone marrow and mamma) and in tumors.

for 100 tests VK-100T
5x100 mL VK-K-100

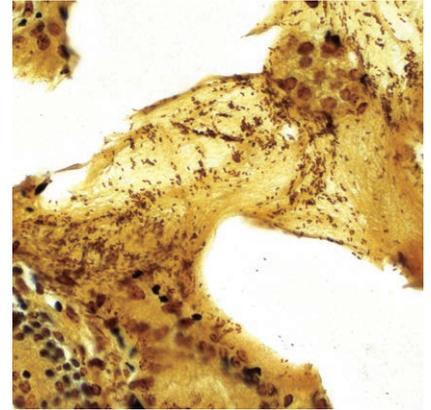

Warthin Starry kit

Five-reagent kit for staining *Spirochaeta spp.*, *Helicobacter pylori*, *Microsporidia spp.* and *Legionella pneumophila*. The kit contains 12 jars with gelatin that enables both incubation and staining of sections, as well as other reagents that enable precipitation of silver on the bacterial surface. The bacteria are found in the mucus of the surface epithelium, in the apical gastric glands and in the gastric mucosa.

up to 48 tests WS-12T



Weigert-Van Gieson kit

Six-reagent kit for staining elastic fibers with prolonged incubation period. Used for differentiation between elastic tissue, collagen and other types of connective tissue.

for 100 tests WVG-100T
6x100 mL WVG-K-100





Additional staining reagents and solutions

Ammonium water

For use with Bielschowsky kit.

100 mL AV-OT-100

Ammonium iron sulfate, solution

For use with Reticulin and Reticulin Contrast staining kits.

100 mL ASF-OT-100

Aniline, alcoholic solution

For use with Azan Trichrome kit.

100 mL ANA-OT-100
250 mL ANA-OT-250

Aniline Blue reagent

Counterstain used with Masson Trichrome kit.

100 mL ABR-OT-100
500 mL ABR-OT-500

Aniline Blue, solution

For use with MSB kit.

100 mL ABO-OT-100

Azan reagent

For use with Azan Trichrome kit.

100 mL AZR-OT-100
250 mL AZR-OT-250

Azocarmine, solution

For use with Azan Trichrome kit.

100 mL AZ-OT-100
250 mL AZ-OT-250

Basic Activation buffer

For use with Oil Red O kit.

250 mL BAP-OT-250

Benedict's reagent

For detecting glucose in urine.

100 mL BNR-OT-100
250 mL BNR-OT-250
500 mL BNR-OT-500
1000 mL BNR-OT-1L

Biebrich Scarlet-Acid Fuchsin reagent

For use with special kits.

100 mL BSAF-OT-100
500 mL BSAF-OT-500

BioSchiff reagent

Schiff's reagent for detecting aldehydes and mucous substances. A component of several special staining kits.

100 mL BS-OT-100
500 mL BS-OT-500
1000 mL BS-OT-1L

BioSchiff Forte

Reagent for use with LeukoGnost PAS kit. Enhanced formulation of Schiff reagents for staining cytoplasm of neutrophils, monocytes, lymphocytes and basophils intensive magenta.

100 mL BSF-OT-100
500 mL BSF-OT-500

Borax, solution

For use with Grocott and P.A.S.M./Jones kits.

100 mL BO-OT-100
105 mL BO-OT-105
250 mL BO-OT-250

SCHIFF'S REAGENT
STABLE AT 22°C*



* In an unopened original packaging. Keep refrigerated after first use.



Additional staining reagents and solutions

Fuchsin Acid Van Gieson reagent

For use in special staining kits.
Synonym: Picrofuchsin reagent.

100 mL FAG-OT-100

Fuchsin Basic reagent

For use in bacteriology.

100 mL FBR-OT-100
500 mL FBR-OT-500
1000 mL FBR-OT-1L

Fuchsin Basic, solution

For use with H.B.F.P. kit.

100 mL FBO-OT-100

Gelatin, gel in transport jar

For preparing developing solution with Warthin Starry staining kit.

12x8 mL GEL-128



Glutaraldehyde, 4% buffered solution

Solution used with Sudan Black B kit.

100 mL G4-OT-100



Gold chloride solutions

Aqueous solutions of gold chloride (aurichloric acid, HAuCl_4), component of many special staining kits.

Gold chloride, 0.2% solution	100 mL	ZK02-OT-100
Gold chloride, 0.6% solution	100 mL	ZK06-OT-100
Gold chloride, 1% solution	50 mL	ZK1-OT-50
Gold chloride, 2% solution	50 mL	ZK2-OT-50

Hematoxylin, Verhoeff A

For use with Verhoeff and Movat kits.

100 mL HEMV-OT-100



Hematoxylin, Weigert A

Hematoxylin solution in ethanol, component of Hematoxylin W kit.
Part of several special staining kits.

100 mL HEMA-OT-100
500 mL HEMA-OT-500
1000 mL HEMA-OT-1L



Hydroquinone, solution

For preparation of developing solution in Warthin Starry kit.

50 mL HQ-OT-50



Hydroquinone-sodium sulphite

For use with Grimelius kit

15 x 3,5 g HNS-15-35

Iodine solution, Verhoeff C

For use with Verhoeff and Movat kits.

100 mL JODV-OT-100

Iron chloride, solution

For use with Fouchet-Van Gieson kit.

100 mL FK-OT-100

Kovac indole reagent

Reagent for identification of indole-positive and indole-negative microorganisms.

100 mL KOV-OT-100



Ladewig's reagent

For use with Ladewig Trichrome kit, contains Aniline Blue, Acid Fuchsin and Orange G dyes.

100 mL LDWR-OT-100

Light Green S.F. reagent

Counterstain used with trichrome special staining kits.

100 mL LGS-OT-100

Lithium carbonate, Luxol

Lithium carbonate solution used with Luxol Fast Blue kit.

100 mL LKL-OT-100

Lithium carbonate, solution

Lithium carbonate solution used with Von-Kossa staining kit.

100 mL LK-OT-100

Lugol's solution, M.I.F.

For use with M.I.F. kit.

100 mL LUGM-OT-100

Luxol Fast Blue, solution

For use with Luxol Fast Blue kit.

100 mL LFB-OT-100



M.F. reagent

For use with M.I.F. kit.

100 mL MF-OT-100

Martius Yellow, solution

For use with MSB kit.

100 mL MAY-OT-100

Mayer's mucicarmine

For use with Mucicarmine kit.

100 mL MMU-OT-100



Metanil Yellow, solution

For use with Mucicarmine kit.

100 mL MY-OT-100



Additional staining reagents and solutions

Potassium hexacyanoferrate, solution

Potassium hexacyanoferrate (II) solution, optimally concentrated for use with HemoGnost Perls kit.

100 mL	KHC-OT-100
250 mL	KHC-OT-250

Potassium hydroxide, solution

Optimally concentrated aqueous solution of potassium hydroxide used with Congo Red Highman kit.

100 mL	KHO-OT-100
--------	------------

Potassium permanganate, solutions

For use with Hematoxylin P.T.A., Masson Fontana, Orcein, Paraldehyde Fuchsin and Reticulin staining kits.

Potassium permanganate, 0.5% solution	100 mL	KP05-OT-100
---------------------------------------	--------	-------------

Potassium permanganate, 1% solution	250 mL	KP1-OT-250
-------------------------------------	--------	------------

Propylene glycol

For use with Sudan Black B Lipid kit.

100 mL	PG-100
--------	--------

Propylene glycol, 85%

For use with Sudan Black B Lipid kit.

100 mL	PG85-100
--------	----------

P.T.A. - P.M.A. reagent

Reagent containing optimal concentration of phosphotungstic and phosphomolybdic acid. Used with Masson Trichrome and Masson-Goldner Trichrome staining kits.

100 mL	PPR-OT-100
500 mL	PPR-OT-500

Reagent for differentiation in Verhoeff kit

For use with Verhoeff kit.

100 mL	RDV-OT-100
--------	------------

Reducing buffer, solution

For use with Von Kossa kit.

100 mL	RP-OT-100
--------	-----------

Resorcine Fuchsin reagent

For use with Weigert-Van Gieson kit.

100 mL	RFR-OT-100
--------	------------

Rhodanine reagent

For use with Rhodanine kit.

100 mL	RR-OT-100
--------	-----------

Safranin O reagent

Counterstaining reagent for use in histology and cytology. Used for visualizing and staining cartilage, mucins and mastocytes in formalin-fixed and paraffin-embedded tissue samples.

100 mL	SAFO-OT-100
--------	-------------

Silver nitrate solutions

For use in special staining kits.

Silver ammonia solution (component of Reticulin and Reticulin Contrast staining kits) ●	100 mL	SA-OT-100
---	--------	-----------

Silver ammonia reagent (component of Von Kossa and Masson Fontana staining kits) ●	100 mL	SAR-OT-100
--	--------	------------

Silver ammonia reagent stock solution (component of Bielschowsky staining kit) ●	100 mL	SSAR-OT-100
--	--------	-------------

Silver nitrate, 20% solution (component of Bielschowsky staining kit) ●	100 mL	SN20-OT-100
---	--------	-------------

Silver nitrate solution, stabilized (component of Grocott and P.A.S.M. / Jones, stabilized staining kits)	100 mL	SNS-OT-100
---	--------	------------

Silver nitrate, WS solution (component of Warthin Starry staining kit) ●	100 mL	SNWS-OT-100
--	--------	-------------

Sodium acetate, solution

For use with Rhodanine kit.

100 mL	NA-OT-100
--------	-----------

Sodium hydroxide, solution

For use with Congo Red Puchtler kit.

30 mL	NH-OT-30
-------	----------

Sodium metabisulfite solutions

Component of Alcian Blue-P.A.S., P.A.S. and Feulgen staining kits.

Sodium metabisulfite	100 mL	NM-OT-100
	500 mL	NM-OT-500

Sodium metabisulfite, solution for Feulgen	100 mL	NMF-OT-100
--	--------	------------

Sodium tetraborate, solution

A component of Alcian Blue pH 2.5 and Alcian Blue-P.A.S. staining kits.

100 mL	NTB-OT-100
500 mL	NTB-OT-500

Sodium thiosulfate solution

Solutions are component of Grocott stabilized, Masson Fontana, P.A.S.M. / Jones, stabilized, Reticulin, Reticulin Contrast and Von Kossa staining kits.

Sodium thiosulfate, 2% solution	100 mL	NT2-OT-100
---------------------------------	--------	------------

Sodium thiosulfate, 5% solution	100 mL	NT5-OT-100
---------------------------------	--------	------------

Sudan Black B, solution

Solution for staining lipids, used with Sudan Black B lipid kit.

100 mL	SBO-OT-100
--------	------------

**Hematoxylin, C.I. 75290**

Natural Black 1, BSC certified stain. For preparing nuclear staining reagents.

25 g	H-P-25
100 g	H-P-100

Janus Green B, C.I. 11050

Diazine Green, BSC certified stain. For displaying mitochondria in living cells.

25 g	JG-P-25
------	---------

Jenner Stain

Jenner, BSC certified stain. Romanowsky polychromatic stain routinely used in hematology for staining blood smears.

25 g	JS-P-25
------	---------

Leishman Stain

Eosin-polychrome Methylene Blue. Romanowsky polychromatic dye. Used to distinguish and identify leukocytes and malaria parasites.

25 g	LS-P-25
------	---------

Light Green S.F., C.I. 42095

Acid Green 5, Pencil Green SF, BSC certified stain. For counterstaining according to Papanicolaou and Masson trichrome methods.

25 g	LG-P-25
------	---------

Luxol Fast Blue MBSN, C.I. 74180

Luxol Fast Blue, Solvent Blue 38. For use in histology.

25 g	LFB-P-25
------	----------

Malachite Green, C.I. 42000

Malachite Green Oxalate, Basic Green 4, Victoria Green, BSC certified stain. Used as counterstain in TB-Stain Cold kit.

10 g	MAG-P-10
25 g	MAG-P-25

Martius Yellow C.I. 10315

Golden Yellow, Martinsgelb, Acid yellow 24. For staining erythrocytes yellow in trichrome staining methods.

25 g	MAY-P-25
------	----------

May-Gruenwald

Polychromatic Romanowsky dye. For blood smears and bone marrow staining.

25 g	MG-P-25
50 g	MG-P-50

Metanil Yellow, C.I. 13065

Acid yellow 36. For staining cytoplasm, muscle connectors and identifying parietal cells of gastric glands.

25 g	MEY-P-25
50 g	MEY-P-50
100 g	MEY-P-100

Methyl Green, C.I. 42590

For staining nuclei green.

25 g	MGR-P-25
------	----------

Methyl Violet 2B, C.I. 42535

Basic Violet 1, BSC certified stain. For bacterial staining according to Gram.

25 g	MV2B-P-25
------	-----------

Methyl Violet 10B, C.I. 42555

Gentian Violet, Crystal Violet, BSC certified stain. For bacterial staining according to Gram.

25 g	MV10B-P-25
50 g	MV10B-P-50

Methylene Blue, C.I. 52015

Basic Blue 9, BSC certified stain. Used for preparing Romanowsky polychromatic stain for staining bone marrow and as a counterstain in staining of acid-fast bacteria.

25 g	MB-P-25
100 g	MB-P-100

Methylene Green, C.I. 52020

Methylene Green, BSC certified stain. For staining nuclei green.

25 g	MEG-P-25
------	----------

Neutral Red, C.I. 50040

Basic Red 5, BSC certified stain. For supravital staining of blood cells using Neutral Red-Janus Green B method.

25 g	NR-P-25
------	---------

Nigrosine WS, C.I. 50420

Nigrosine water soluble, Nigrosine B. For use in microbiology and cytology, as negative staining.

25 g	NGR-P-25
------	----------

Nuclear Fast Red, C.I. 60760

Kernechtrot. Used for staining nuclei red by creating aluminum complex.

5 g	NF-P-5
-----	--------

Oil Red O, C.I. 26125

Sudan Red 5B, BSC certified stain. For staining lipids.

25 g	ORO-P-25
------	----------

Orange II, C.I. 15510

Acid Orange 7, BSC certified stain. Used as a counterstain, as a substitute for Orange G.

25 g	O2-P-25
------	---------

Orange G, C.I. 16230

Acid Orange 10, Acid Fast Orange G, BSC certified stain. For staining keratin according to Papanicolaou, counterstain for hematoxylin.

25 g	OG-P-25
100 g	OG-P-100

Orcein, C.I. (ed.1) 1242

Natural Red 28, Orcein Synthetic. For staining elastic tissues and viral inclusion bodies of the hepatitis B surface antigen.

5 g	ORC-P-5
10 g	ORC-P-10

Phloxine B, C.I. 45410

Acid Red 92, Eosin 10B, Eosin Blue, BSC certified stain. For Methylene Blue-Phloxine staining and staining keratin and mucus according to Kreyberg.

25 g PB-P-25

Ponceau S, C.I. 27195

Acid Red 112, Fast Ponceau 2B. For use in electrophoresis.

25 g PS-P-25

Potassium permanganate

Synonym: Potassium (VII) manganate, Condy's crystals. Potassium permanganate for use in microscopy.

50 g KPH-P-50
100 g KPH-P-100

Pyronine B, C.I. 45010

Pyronine B, BSC certified stain. For bacterial, mold and RNA staining.

25 g PYB-P-25

Resorcinol, C.I. 76505

Resorcin, 1,3-benzendiol. Used for differentiation between aldoses and ketoses in Seliwanoff's test by staining ketose sugars dark red.

25 g RES-P-25

Rhodamine B, C.I. 45170

Basic Violet, Acid Red 52. For staining acid-fast bacteria.

25 g RHB-P-25

Rhodanine

DMAB-Rhodanine. For staining deposits of copper bound to proteins.

5 g RHO-P-5

Rose Bengal, C.I. 45440

Acid Red 94, BSC certified stain. For staining in microscopy and histology.

25 g RB-P-25

Safranin O, C.I. 50240

Basic Red 2, Safranin T, Phenazinium, BSC certified stain. Counterstain in the Gram method for detecting Gram-negative bacteria.

25 g SO-P-25

Sirius Red, C.I. 35780

Direct Red 80. Used for staining collagen and amyloid in histology.

25 g SR-P-25

Sudan III, C.I. 26100

Solvent Red 23, Sudan G, Fat Ponceau G, BSC certified stain. For lipid staining according to Lillie-Ashburn.

25 g S3-P-25

Sudan IV, C.I. 26105

Solvent Red 24, Scarlet Red, Biebrich Scarlet R, Fat Ponceau. For staining lipids according to Churukian's modification of Lillie-Ashburn's method.

25 g S4-P-25

Sudan Black B, C.I. 26150

Sudan Black, Solvent Black 3, Fat Black HB. For staining lipids and other hydrophobic materials.

5 g SB-P-5
10 g SB-P-10

Toluidine Blue O, C.I. 52040

Basic Blue 17, Tolonium Chloride, BSC certified stain. For staining mast cells, cartilage and Nissl bodies in the spinal cord.

25 g TB-P-25
100 g TB-P-100

Victoria Blue B, C.I. 44045

Basic Blue 26, Fat Blue B. For use in cytology and histology.

25 g VBB-P-25

Victoria Blue R, C.I. 42563

Basic Blue 8, Fat Blue 4R. For use in microbiology and histology.

25 g VBR-P-25

Water Blue, C.I. 42755

Acid Blue 22, BSC certified stain. For staining connective tissue according to Mallory and Gomori trichrome methods.

25 g WB-P-25

Wright's Stain

Colorants de Wright, BSC certified stain. Part of polychromatic Romanowsky stains commonly used in hematology.

25 g WS-P-25

Xylidine Ponceau, C.I. 16150

Ponceau 2R, Ponceau de Xylidine, BSC certified stain. For Masson trichrome staining.

25 g XP-P-25





Accessories


High quality disposable laboratory equipment including glass, plastic and metal items. The product range includes VitroGnost slides for microscopy, sample transport containers and various other accessories.

VitroGnost® Plus Ultra adhesive microscope slides

These slides are ideal for pathology, cytology, and special immunohistochemistry (IHC). VitroGnost Plus Ultra adhesive slide is a new generation slide with a denser and stronger positive charge on the surface to achieve even better adhesion of the sample. VitroGnost Plus Ultra adhesive slide can be used for many scientific purposes, especially with tissues difficult to fix, like breast tissue and samples of brain and skin. Dimensions: 26x76x1 mm.

72 pcs

103001A72

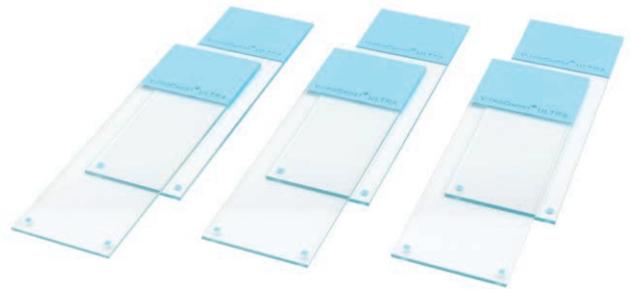

Lid for cassettes

Cassette embedding stainless steel lid, standard.

material	packaging	product code
INOX	1 pc	51050103





Bio/Cut microtome blades - high profile

Bio/Cut H35C microtome blades for cryotome, high profile (wide). Ideal for frozen sections. Suitable for both soft and hard tissues.

50 pcs Bio/Cut H35C

Histoblades Duraedge histology blades - high profile

HistoBlades Duraedge blades with Teflon coating - high profile.

50 pcs 7310

Bio/Cut microtome blades - low profile

Bio/Cut microtome blades manufactured in order to provide most precise sectioning.

Bio/Cut 22A microtome blades for ultra-thin sections, kidneys and lymph nodes, narrow blade 50 pcs Bio/Cut 22A

Bio/Cut 35A microtome blades for soft and hard tissues, low profile 50 pcs Bio/Cut 35A

Bio/Cut 35B microtome blades for soft tissues, low profile 50 pcs Bio/Cut 35B

Bio/Cut 35N microtome blades for thin sections, low profile 50 pcs Bio/Cut 35N

Bio/Cut 35R microtome blades for hard tissues, low profile 50 pcs Bio/Cut 35R

Bio/Cut 35S microtome blades for hard tissues, low profile 50 pcs Bio/Cut 35S

Histoblades Duraedge histology blades - low profile

HistoBlades Duraedge blades with Teflon coating - low profile.

50 pcs 7223

Histoblades Encore histology blades - low profile

Extremely sharp new generation Teflon coated blades - low profile.

50 pcs 7380



Contact and ordering information

BioGnost Ltd.,
Medjugorska 59, Zagreb
10040 Croatia (Europe)

WORKING HOURS
Monday through Friday
8:00 to 16:00 (CET)

e-mail: sales@biognost.hr
website: www.biognost.com
phone: +385 1 2409 997
fax: + 385 1 2404 039



Technical support

Technical questions can be answered by emailing our sales representatives at sandra.ban@biognost.hr, pavle.loncar@biognost.hr, viktoria.liber@biognost.hr and andrej.spehar@biognost.hr.

They will provide you with any MSDS or product information required, also available on our website at www.biognost.com.

Frequently Asked Questions (FAQ) section is available on our website.



Shipping information

Shipping and handling charges are prepaid and added to the invoice. If requested, we can use your account number or carrier of your choice, and also ship on your third-party freight accounts (Incoterms 2020). Charges vary according to destination, weight and content, and are available upon request.

Our employees are trained to classify, package, mark and label products identified as dangerous goods according to ADR, IMDG and IATA to ensure safe transport.



Index

4% Formaldehyde	8
4% Formaldehyde, alcoholic solution	51
4% NB Formaldehyde	8
10% Formaldehyde	9
10% NB Formaldehyde	9

A

Acetic acid, 0.5% solution	53
Acetic acid, 1% solution	53
Acetic acid, Colloidal Iron (A) solution	53
Acetic acid, Colloidal Iron (B) solution	53
Acetic buffer, Grimelius	51
Acid alcohol	24
Acid alcohols	51
Acid solution	51
Acridine Orange, C.I. 46005	58
Acridine Orange solution	51
A.F.O.G. kit	41
A.F.O.G. reagent	51
Alcian Blue 8 GX, C.I. 74240	58
Alcian Blue - P.A.S. kit	41
Alcian Blue pH 1.0 kit	41
Alcian Blue pH 2.5 kit	41
Alcian Blue solution 1%	51
Alcian Blue solution pH 1.0	51
Alcian Blue solution pH 2.5	51
Alcian Green	58
Alcian Yellow, C.I. 12840	58
Alcian Yellow Toluidine Blue kit	41
Alcoholic Eosin Y solutions	25
Alkaline alcohol, solution	51
Amido Black 10B, C.I. 20470	58
Ammonium iron sulfate, solution	52
Ammonium water	52
Aniline, alcoholic solution	52
Aniline Blue, C.I. 42780	58
Aniline Blue reagent	52
Aniline Blue, solution	52
Auramine O, C.I. 41000	58
Azan reagent	52
Azan Trichrome kit	41
Azocarmine G, C.I. 50085	58
Azocarmine, solution	52
Azophloxine, C.I. 18050	58
Azure A, C.I. 52005	58
Azure B, C.I. 52010	58
Azure C, C.I. 52002	58
Azure II	58
Azure II Eosinate	59

B

Basic Activation buffer	52
Benedict's reagent	52
Biebrich Scarlet-Acid Fuchsin reagent	52
Biebrich Scarlet, C.I. 26905	59
Bielschowsky kit	42
BioClear	13
BioClear New	13
Bio/Cut microtome blades - high profile	64
Bio/Cut microtome blades - low profile	64
Bio-Diff 1 reagent	28
Bio-Diff 1 RTU reagent	29
Bio-Diff 2 reagent	28
Bio-Diff 2 RTU reagent	29
Bio-Diff 3 reagent	28
Bio-Diff 3 RTU reagent	29
Bio-Diff kit	28
Bio-Diff RTU kit	29, 42
BioFix	9

BioFix C	9
BioFix GL	9
BioFix HD	9
BioForm S	9
BioGram 4 kit	34
BioGram ECO kit	35
BioGram Histo kit	42
BioMount	18
BioMount Aqua	19, 33
BioMount C	19
BioMount DPX	17
BioMount DPX High	17
BioMount DPX Low	17
BioMount DPX New	17
BioMount High	18
BioMount M	18
BioMount New	18
BioMount New Low	18
Biopsy pads for cassettes, blue	63
BioSchiff reagent	52
BioWax 52/54	14
BioWax 56/58	14
BioWax Blue	14
BioWax Plus 56/58	14
Bismarck Brown Y, C.I. 21000	59
Bluing reagent	24
Borax, solution	52
Bouin's solution	10
Brilliant Cresyl Blue, C.I. 51010	59
Brilliant Cresyl Blue, solution	31
Brilliant Green Certified, C.I. 42040	59
Buffer solution, pH 6.4	30
Buffer solution, pH 6.8	30
Buffer solution, pH 7.2	29
Buffer tablets pH 6.8	30
Buffer tablets pH 7.2	29

C

Carbol Fuchsin Gram 10x concentrated	35
Carmine, C.I. 75470	59
Carminic acid, C.I. 75470	59
Celestine Blue, C.I. 51050	59
Chromotrope 2R / Aniline Blue reagent	53
Chromotrope 2R, C.I. 16570	59
Chromotrope 2R / Fast Green F.C.F. reagent	53
CitoFix	10
CitoSpray	10
CitoSpray solution	10
Clean Lens, solution 1	21
Clean Lens, solution 2	21
Clean Lens, solution 3	21
Colloidal Iron kit	42
Colloidal Iron, stock solution	53
Congo Red, C.I. 22120	59
Congo Red Highman kit	42
Congo Red Puchtler kit	42
Congo Red Puchtler reagent	53
Congo Red reagent	53
Cresyl Violet Acetate Certified	59
Cresyl Violet, solution	53
CryoFix Gel	15
Cryo Spray	15
Crystal Violet, C.I. 42555	59
CytoGnost kit	31

D

DeHematin	9
DeMold solution	15
DeMold Spray	15
Developer stock solution	53
Differentiation reagents based on acetic acid	53
Differentiation reagents based on hydrochloric acid	53
Duoline marker pen	7

LeukoGnost PAS	32	Orange G, C.I. 16230	60
LeukoGnost PERLS	32	Orange II, C.I. 15510	60
LeukoGnost PLUS	33	Orange II reagent, Pap 2B	26
LeukoGnost SPE	32	Orcein, C.I. (ed.1) 1242	60
LeukoGnost SPE kit	45	Orcein kit	47
LeukoGnost SPENSE	32	Orcein reagent	55
Lid for cassettes	62	OsteoFast 1	11
Light Green S.F., C.I. 42095	60	OsteoFast 2	11
Light Green S.F. reagent	54	OsteoSens	11
Liquid Blocker marker pen	7	Oxalic acid solutions	55
Lithium carbonate, Luxol	54		
Lithium carbonate, solution	54	P	
Lugol's C solution	35	Paraldehyde Fuchsin kit	47
Lugol's solution, acid	35	Paraldehyde Fuchsin reagent	55
Lugol's solution, M.I.F.	54	P.A.S. Diastase kit	40, 47
Lugol's solution (Schiller's test)	35	P.A.S. kit	47
Luxol Fast Blue kit	45	Periodic acid, solution	55
Luxol Fast Blue MBSN, C.I. 74180	60	Perls-Van Gieson kit	48
Luxol Fast Blue, solution	54	Phloxine B, C.I. 45410	61
		Phosphomolybdic acid, 1% solution	55
M		Phosphomolybdic acid, 5% solution	55
Malachite Green, C.I. 42000	60	Phosphotungstic acid, 1% solution	55
Mallory Trichrome kit	45	Phosphotungstic acid, 5% solution	55
Martius Scarlet Blue (MSB) kit	46	Picric acid in acetone, solution	55
Martius Yellow C.I. 10315	60	Picro-Sirius Red kit	48
Martius Yellow, solution	54	Picro-Sirius Red reagent	55
Masson Fontana kit	46	Polypropylene containers for transporting samples	65
Masson-Goldner Trichrome kit	46	Polypropylene dissecting boards	63
Masson Trichrome kit	46	Ponceau S, C.I. 27195	61
Mayer's mucicarmine	54	Ponceau S, solution	55
May-Gruenwald	60	Potassium hexacyanoferrate, Colloidal Iron solution	55
May-Gruenwald solution	30	Potassium hexacyanoferrate, solution	56
Metal histology mold	63	Potassium hydroxide, solution	56
Metanil Yellow, C.I. 13065	60	Potassium hydroxide, solution for fungi	39
Metanil Yellow, solution	54	Potassium permanganate	61
Methenamine, solution	55	Potassium permanganate, solutions	56
Methyl Blue Eosin, solution	55	Propylene glycol	56
Methylene Blue 10X concentrated	36	Propylene glycol, 85%	56
Methylene Blue, C.I. 52015	60	P.T.A. - P.M.A. reagent	56
Methylene Blue Loeffler's reagent	36	Pyronine B, C.I. 45010	61
Methylene Green, C.I. 52020	60		
Methyl Green, C.I. 42590	60	R	
Methyl Violet 2B, C.I. 42535	60	Rapid Pap staining kit	27
Methyl Violet 10B, C.I. 42555	60	Reagent for differentiation in Gomori Trichrome kit	53
M.F. reagent	54	Reagent for differentiation in Verhoeff kit	56
Micro Clean	14	Reducing buffer, solution	56
Microtome oil	65	Resorcine Fuchsin reagent	56
M.I.F. kit	45	Resorcinol, C.I. 76505	61
Movat kit	46	Reticulin Contrast kit	48
Mucicarmine kit	46	Reticulin kit	48
Myco-Mount Blue	39	Rhodamine B, C.I. 45170	61
Myco-Mount Red	39	Rhodanine	61
		Rhodanine kit	48
		Rhodanine reagent	56
		Rose Bengal, C.I. 45440	61
N			
NB Formaldehyde concentrated 5x	9	S	
Neutral Red, C.I. 50040	60	Saccommanno's fixative	10
Neutral Red reagent	55	Safranin O, C.I. 50240	61
Nigrosine, isotonic solution	55	Safranin O kit	48
Nigrosine WS, C.I. 50420	60	Safranin O reagent	56
Nuclear Fast Red, C.I. 60760	60	Scott's solution	24
Nuclear Fast Red (Kernechtrot) reagent	55	Shorr's solution	27
		Silver nitrate solutions	56
O		Sirius Red, C.I. 35780	61
OG-6 reagent, Pap 2A	26	Sodium acetate, solution	56
OG-EA31 Pap reagent	27	Sodium hydroxide, solution	56
OG-EA50 Pap reagent	27	Sodium metabisulfite solutions	56
Oil Red O, C.I. 26125	60	Sodium tetraborate, solution	56
Oil Red O kit	47	Sodium thiosulfate solution	56
Oil Red O reagent	55	Sperm-Diff RTU kit	49
Orange G, 1% solution	55		
Orange G / Aniline Blue reagent	55		

Sudan Black B, C.I. 26150	61
Sudan Black B ECO	32
Sudan Black B kit	32
Sudan Black B lipid kit	49
Sudan Black B reagent	57
Sudan Black B, solution	56
Sudan III, C.I. 26100	61
Sudan IV, C.I. 26105	61
Sulfuric acid solutions	57
Sulphosalicylic acid, 20% solution	57

T

TB Armand reagent	37
TB Auramine O reagent	38
TB Auramine-Rhodamine reagent	38
TB Carbol Fuchsin reagent	36, 37
TB Decolorizer	36, 37
TB Decolorizer Fluorescent	38
TB-Fuchsin reagent	36
TB Malachite Green reagent	37
TB Permanganate reagent	38
TB-Stain Auramine O kit	38
TB-Stain Cold kit	37
TB-Stain ECO kit	36
TB-Stain Fluorescent kit	38
TB-Stain Histo kit	49
TB-Stain Hot kit	36
TB-Stain Quick kit	37
ThromboGnost solution	31
Tissue marking dyes	7
Tissue marking dyes, 5 colour set	7
Tissue marking dyes, 7 colour set	7
Tissue Soft	15
Toluidine Blue O, C.I. 52040	61
Toluidine Blue, solution	57
Trichloroacetic acid, solution	57
Tuerk's solution	31

U

UriGnost S kit	49
UriGnost SM kit	49
UriGnost SM reagent	57
UriGnost S reagent	57

V

Van Gieson Trichrome kit	49
Verhoeff kit	50
Victoria Blue B, C.I. 44045	61
Victoria Blue R, C.I. 42563	61
VitroGnost® Plus Ultra adhesive microscope slides	62
Von Kossa kit	50

W

Warthin Starry kit	50
Water Blue, C.I. 42755	61
Weigert-Van Gieson kit	50
Wright-Giemsa solution	31
Wright-Giemsa, solution for reptiles	31
Wright's solution	31
Wright's Stain	61

X

Xylidine Ponceau, C.I. 16150	61
------------------------------	----

Be a part of
our global network
of satisfied customers

BioGnost Ltd.

For development, manufacture and
marketing of pharmaceuticals

Medjugorska 59, 10040 Zagreb, Croatia, EU
phone: +385 1 2409 997, fax: +385 1 2404 039
sales@biognost.hr • www.biognost.com

