

Hygiene

Citrosteril

Heat Disinfectant for Haemodialysis Machines
with Recirculation



 **FRESENIUS
MEDICAL CARE**
Fresenius Medical Care Deutschland GmbH
61352 Bad Homburg • Else-Kröner-Str. 1 • Germany


I.V. Nikolay Kapustin


I.V. Oleksandr Zaika

Citrosteril

Heat Disinfectant for Haemodialysis Machines with Recirculation


Disinfectant for chemo-thermal disinfection of haemodialysis machines¹.

- pH value 1.7 to 2.0
- excellent removal of limescale
- disinfection and decalcification in one process
- active ingredients composed of natural substances
- biodegradable
- odourless
- free from colouring additives

Action

The synergistic effect of its components makes Citrosteril a potent disinfectant solution.

Citrosteril at 84°C has a broad spectrum of micro-biocidal activity and works bactericidal and virucidal² including HBV/HCV/HIV.

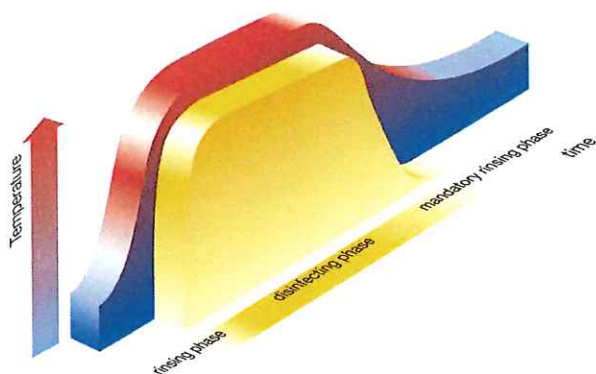

Specification

100 g Citrosteril contains:

21 g citric acid 1-hydrate; lactic acid, malic acid

Ordering information

Unit	Language combination	Art. No.
1 x 5 L (single canister)	multilingual	F00005157
2 x 5 L (carton)	multilingual	F00005158
6 x 2 L (carton)	multilingual	508 536 1

Literature:

1. Solbach W, Universität zu Lübeck; Verification of the sporicidal efficacy of the product Citrosteril at 85°C in the Fresenius 5008 dialysis machine, 20.12.2002
2. Labor Dr. Merk & Kollegen, Ochsenhausen: Antiviral efficacy of Citrosteril against bovine Parvovirus, 9.9.2005

Further information is available on request.



**FRESENIUS
MEDICAL CARE**

Head office: Fresenius Medical Care Deutschland GmbH · 61346 Bad Homburg v.d.H. · Germany
Phone: +49 (0) 6172-609-0 · Fax: +49 (0) 6172-609-2191
www.fmc-ag.com