



[About us](#)



[Products ▼](#)

[Technologies](#)

[Video Library](#)


[News](#)

[Contact](#)

ATTENTION

Products shown may not be available in all countries. Indication claims may vary dependent on geography.

[Your privacy](#)

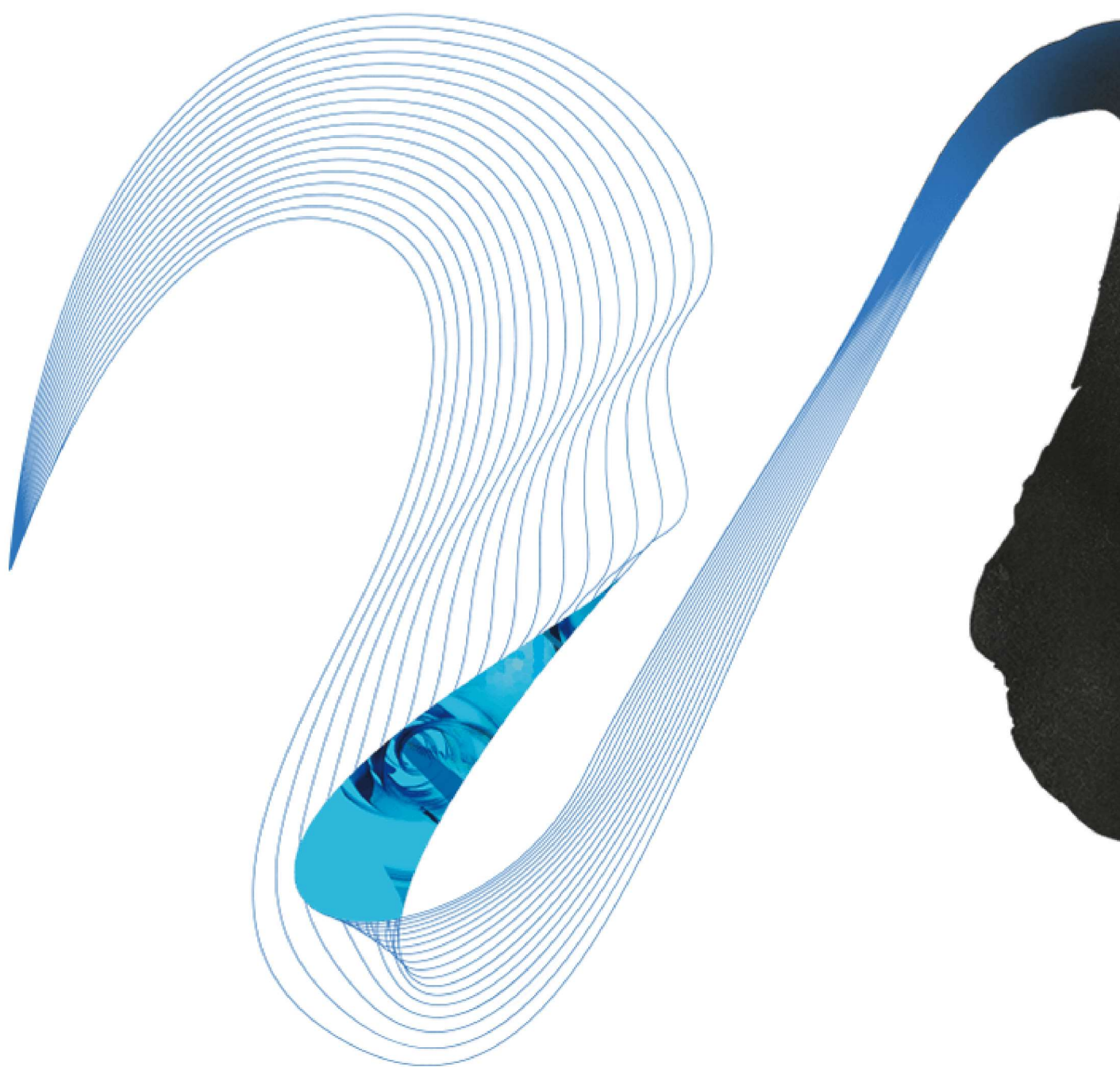


Coronary

ASAHI Corsair Pro

[Your privacy](#)

Delivering excellence through evolution and tradition.

By minimizing the stiffness difference between the tip and the shaft, ASAHI Corsair Pro will be able to track tortuous vessels, which improves its ability to cross challenging anatomy.

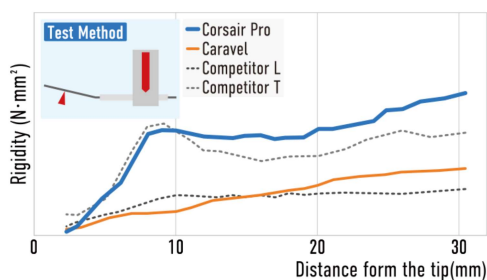


■ Tracking ability

- Optimization of catheter stiffness

Compare to Corsair, the Corsair Pro's stiffness profile has a more linear increase which improves the catheter flexibility and trackability for enhances performance in tortuous anatomy.

■ Balanced joint design

Gradual balance of stiffness from

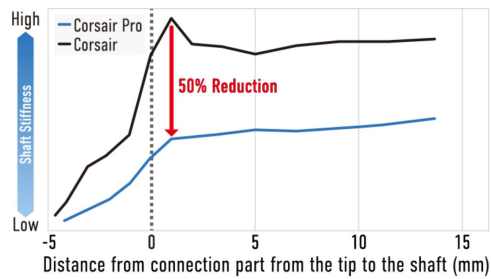

the tip to the joint of the shaft
area enhances
tracking ability.



■ Kink resistance at proximal
shaft

• Spiral Protector

Flexibility of the proximal hub is
improved to protect the
proximal shaft when bent during
catheter advancement or
withdrawal.



■ Torque ability, Pushability and
Support

• SHINKA-Shaft

ASAHI brand proprietary braiding
pattern, which consists of 8
thinner and 2 thicker wires
braided together.

Tungsten braiding enhances
torque, better push and support.



■ Limited rotation

Do not rotate the catheter in the
same direction, clockwise or

counterclockwise, more than 10
times.

* All the data were obtained by company standardized test, which may differ from industry
standardized tests.

* All the data do not prove that all devices have exactly the same performance with the samples
used for these tests.

Technology




Product Spec

Tip Outer 0.42mm (1.3Fr)
Diameter

Hydrophilic60cm
Coating
Length

Usable 135cm / 150cm
Length

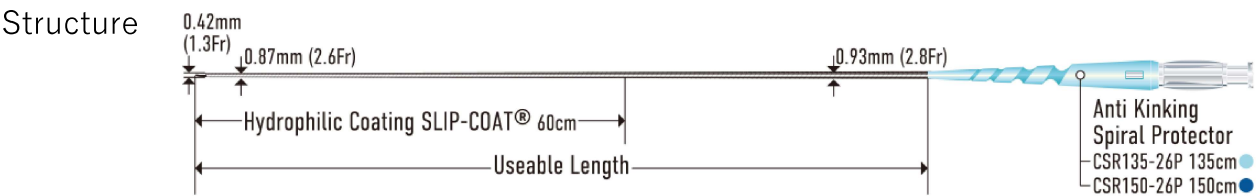

Category [#Microcatheter](#) [#Complex \(/en/products-coronary/?category1=Guidewire&category2=Complex\)](#)

Order Information

.CO

Inquiry about Asahi Intecc medical products

Your privacy

[Contact Form \(/contact/en/?
item=ASAHI Corsair Pro\)](/contact/en/?item=ASAHI Corsair Pro) ➡

SHARE THIS PAGE



[HOME](#) > ASAHI Corsair Pro



Products



[Coronary](/en/products-coronary/)
**[\(/en/products-
coronary/\)](/en/products-coronary/)**

[Peripheral](/en/products-evt/)
[\(/en/products-evt/\)](/en/products-evt/)

[Neuro](/en/products-neuro/)
**[\(/en/products-
neuro/\)](/en/products-neuro/)**

[Interventional](/en/products-ivr/)
[Radiology](/en/products-ivr/)
[\(/en/products-ivr/\)](/en/products-ivr/)

[Technologies](/en/technologies/)
[\(/en/technologies/\)](/en/technologies/)

[Video Library](/en/video-library/)
[\(/en/video-library/\)](/en/video-library/)

[News](/en/news/)
[\(/en/news/\)](/en/news/)


[Contact](/contact/en/)
[\(/contact/en/\)](/contact/en/)

Your privacy

[About us \(/en/about-us/\)](/en/about-us/)

[Privacy&Cookie Policy \(/en/policy/\)](/en/policy/)

 [Corporate site \(http://www.asahi-intecc.co.jp/en/\)](http://www.asahi-intecc.co.jp/en/)
[INTECC EUROPE B.V. \(https://asahi-intecc.eu\)](https://asahi-intecc.eu)

 [ASAHI](#)

© ASAHI
INTECC CO.,
LTD.