
HEMOSTASIS

ClearEndoclip



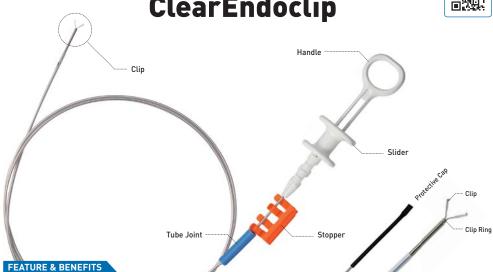

(1) 1:1 Rotation

• Clip rotation with the handle at 1:1 ratio. (Tip: Hold the tube joint and Rotate the handle)

(2) Installed Reopen function

• Reopen and close for correct positioning prior to deployment. [Caution: Open and close up to five times prior to deployment]

(3) Improved safety

• The thickness of the clip is thinner than other products, and it can be tightly deployed to the lesion.





SPECIFICATION

Model	Sheath O.D	Clip Angle	Clip Lengt	Type of Rotation	Working Length	Number of Use	Validity of Sterilzation
FM-CLD11-090(1650) FM-CLD11-090(2300)	¢2.6mm	90°	11mm	Rotational	1,650mm 2,300mm	Single Use Only	3 years

ERCP

Fine-Grab Basket



(1) Standard Four-wire Basket

• These products are made of 4 stainless steel wires in a hexagon shape which is the most popular and suitable for the requirementin the Bile and Pancreatic Ducts.

(2) Precision in Configuration

• 4 wires are arranged in high precision. So this configuration helps operator removes stones fast and safety.

(3) Maintaining the Symmetry

- These products are designed to maintain the symmetry even after the difficult stone manipulation.
- This special and strong structure could make users feel com fortable during ERCP.

SPECIFICATION

Model	Basket Opening Size	Basket Shape & Material	Sheath 0.D	Working Length	Number of Use	Validity of Sterilzation
FM-EG0002 FM-EG0003 FM-EG0004	15mm 20mm 25mm	Hexagon / 4 Wires Stainless Steel	¢2.40mm	2,000mm (+/- 500)	Single Use Only	3 years

20 Imagine What a Woderful difference You have 21