Cooperated with Partners from 101 Countries





BSS MEDICAL SUPPLY CO., LIMITED





⋈ info@bssmedical.com















BSS MEDICAL SUPPLY CO., LIMITED


Better Signal Solution Partner

Why BSS?

- Qifeng Zhou, invented the conductive gel technology for BSS;
- ●FDA, CE, ISO13485, Japanese Registration, RoHS, FSC in EU, SFDA, and dozens of patents;
- The daily capacity is 1.5 million pcs by the advanced intelligent automatic production lines;
- Partners from 101 countries, and products quality approved by thousands of international hospitals;
- •Full products combination meets variety testing purpose: monitoring, stress test, treadmill, holter, MRI, ICG, ABR;

Professor Qifeng Zhou

Academician of Chinese Academy of Sciences (1999) The chairman of IUPAC (2018-2020) President of Peking University (2008-2013)

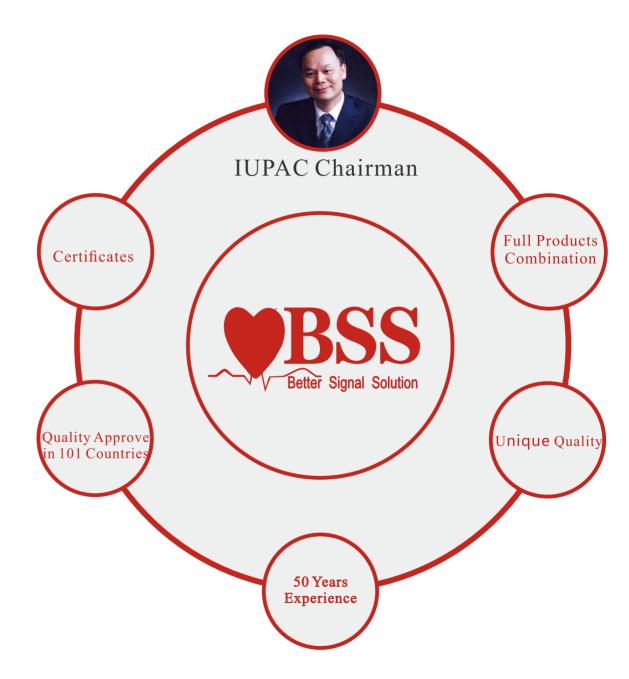





ECG Electrode









Hoter



Resting Test





Neonatal



MRI/X-ray

Better Signal Solution Partner









Wet / Liquid Gel

Wet Gel electrodes could fast catch the ECG signal. They give dependable quick tracing with easy application and removal.

Recommended for: ·EMS ·Emergency



| Item# | Size(mm) | Description | Package |
|--------|----------|---------------|-----------------|
| B4346W | 43*46 | Adult/Oval | 50/pk, 1000/ cs |
| B5050W | 50*50 | Adult / Round | 50/pk, 1000/ cs |

B5050W

ABR Test / EEG

To perform an ABR test, electrodes are placed on the surface of the skin at the patient's head. Good signal quality is achieved by preparation of the skin, improving the conductivity. An electrical potential with very low amplitudes is amplified and recorded, while the ABR system delivers acoustical stimuli to the ears.



Recommended for: ·ABR Test ·EEG Test ·Non-invasive EEG

| Item# | Size(mm) | Description | Package |
|-------|----------|-------------|-----------------|
| B2840 | 28*40 | ABR/EEG | 30/ pk, 1200/cs |

B2840

Pre-wired Neonatal Electrode

The pre-wired neonatal electrode are consist of a conductive adhesive hydro-gel, a Ag/AgCl plated sensing element, a soft cloth or PE Foam substrate with the size of 22mm and a 60 cm flexible tinned copper or carbon lead-wire terminating in a standard DIN connector.





| Item# | Size(mm) | Description | Package |
|-------|----------|---------------|---------------|
| BF22 | 22*23 | Cloth/ Square | 1/pk, 200/cs |
| BY22 | 22*22 | Cloth/ Round | 1/pk, 200/cs |
| FF22 | 22*23 | Foam/ Squre | 1/pk, 200/cs |
| FY22 | 22*22 | Foam/Round | 1/ pk, 200/cs |


Indocyanine green(ICG)Test

Noninvasive ICG clearance test for estimating hepatic blood flow During halothane and isoflurane anaesthesia.

TEB (ICG) electrodes are designed to work with Thoracic Electrical Bioimpedance (TEB) devices..

Recommended for:

ICU/CCU · Emergency · Cardiology · Anesthesiology

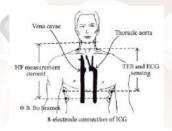

| 2 | 10 | | |
|---|--------------------------------------|-----|-----|
| ê | 00 | e e | ê ê |
| | Name and Address of the Owner, where | No. | |

B4090

| Item# | Size(mm) | Description | Package | |
|--------|----------|-------------|-------------|--|
| B4045B | 40*45 | Foam | 4/pk,100/cs | |
| B4090B | 40*90 | Foam | 4/pk,100/cs | |



B2334

Resting Tab

Tab is used for Resting test to be attached to alligator clip connectors. It has the advantages of fast application and removal from patient's skin.

| Item# | Size(mm) | Description | Package |
|-------|----------|--------------------|-----------------|
| B2334 | 23*34mm | Ag/AgCl coated PET | 100/pk 10000/cs |
| B2034 | 20*34mm | Ag/AgCl coated PET | 100/pk 10000/cs |
| B2634 | 26*34mm | Ag/AgCl coated PET | 100/pk 10000/cs |

OEM/ODM/ Private label solution

BSS also manufacture wide variety of OEM accounts on a private label basis. We offer the options of Gel, Snap, Backing, Private Label, Product shape, Pouch and Carton Design Solutions.

