

DENTAL

Bien-Air® Lubricare 2® Dental Handpiece Care + Maintenance System




LUBRICARE 2

Quickly perform dental handpiece maintenance (cleaning and lubrication) for up to four handpieces at a time with the Lubricare 2 system, an essential addition to your instrument processing workflow. The clean air-blow system ensures thorough cleaning and lubrication inside the handpieces and allows for consistent, repeatable high standards without the intra-operator variability that may occur during manual cleaning and lubrication.

SPRAYNET® + LUBRIFLUID®

Use Spraynet and Lubrifluid for cleaning and lubrication to maintain impeccably hygienic instruments with an unrivaled service life.

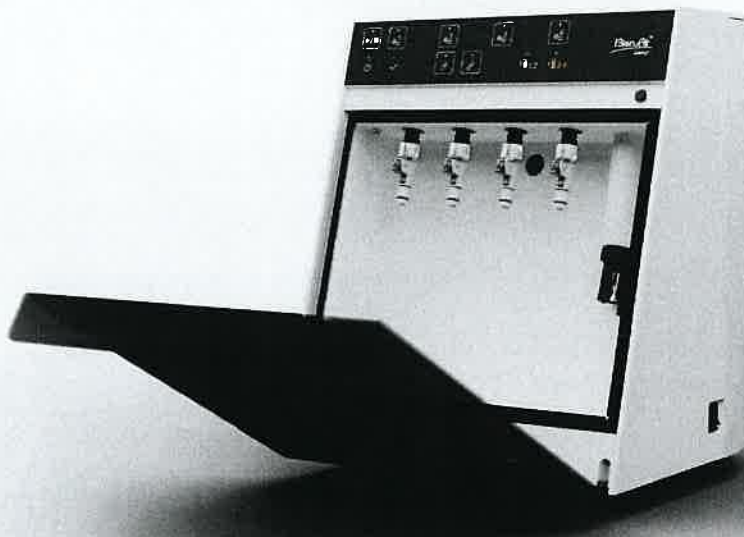

FAST, EFFECTIVE COST-SAVING MAINTENANCE



01 FAST

Clean a turbine in approximately 20 seconds and a contra-angle handpiece in 40 seconds with preprogrammed operating modes that enable automatic maintenance in the Lubricare 2 system.

02 EFFECTIVE

Achieve best cleaning and lubrication outcomes with separate dispensing of cleaner and lubricant for more precise application than systems that dispense the solutions together. Reduce 99.99% of the microbial load with Spraynet cleaning solution.

03 COST-SAVING

Reduce maintenance oil consumption with a more consistent process compared to manual maintenance procedures.

ENHANCE YOUR INSTRUMENT PROCESSING WORKFLOW



LUBRICARE 2® SPECIFICATIONS



Compatibility

E-type attachments
(Additional couplings may be required)

Cycle Duration

20 s turbine, 40 s contra angle

Capacity

4 handpiece lanes (2 cleaning, 2 lubricating) and 1 chuck lane

Rating

AC 100 V–240 V
50–60 Hz

Power Consumption

10 VA–25 VA

Air Flow Rate

60 NL/min–100 NL/min
(0.3 MPa–0.5 MPa)

Air Pressure

43.5 psi–72.5 psi (0.3 MPa–0.5 MPa)
58 psi (0.4 MPa) recommended

Weight

Approx. 17 lb (8 kg)

Outer Dimensions

11.81" L x 11.81" W x 14.37" H
(300 mm x 300 mm x 365 mm)

Bien-Air, Lubricare 2, Spraynet and Lubrifluid are not trademarks of Midmark Corporation and are the sole property of their owner

For more information, contact your Midmark dealer or call 1 800 MIDMARK. Outside the USA, call 1 937 526 3662 or visit our website at midmark.com

© 2022 Midmark Corporation, Miamisburg, Ohio USA
Products subject to improvement changes without notice. Litho in USA
007-10414-00 Rev A1 (2/22)

