

MEDIVATORS™ ISA™

Endoscope Reprocessor
Washer-Disinfecter



EVOLUTIONARY PROCESS, REVOLUTIONARY SYSTEM

The MEDIVATORS ISA Endoscope Reprocessor is a rapid, efficient and flexible reprocessing system designed for endoscopy units demanding the most advanced technology and the highest quality standards.



FAST & SIMPLE OPERATION FOR IMPROVED EFFICIENCY

- Full cycle in 20 minutes with channel monitoring
- Simple, user-friendly handling and operation
- Insertion of endoscopes and operator information through RFID system - helps to avoid errors and saves time



MEDIVATORS ISA Endoscope Reprocessor The Fast And Flexible Reprocessing System

COMPACT & ERGONOMIC DESIGN FOR MAXIMUM USER COMFORT

- Large and easy to load basin accommodates all types of flexible endoscopes
- Intuitive interface assures easy and fast access to all the cycle information
- Compact size contributes to easy installation, operation and maintenance



SAFETY FOR OPERATORS, PATIENTS & ENDOSCOPES

- Hands-free operation and RFID-based technology for improved infection control
- Full traceability through printed and electronic cycle documentation
- Compliant to EN ISO 15883 parts 1, 4 and 5. Tested with multiple test soils and in accordance with German cleaning standards



Foot pedal allows hands-free door lid opening and closing

THE RIGHT FORMULA FOR THE BEST PERFORMANCE & RESULT

- Dedicated room temperature chemicals compatible with endoscope material
- Provides low temperature reprocessing at 20° C to 40° C
- Improved water and chemical consumption
- Designed specifically for use with MEDIVATORS ISA Endoscope Reprocessor



PART OF OUR COMPLETE CIRCLE OF GI REPROCESSING PRODUCTS

Efficient and effective high-level disinfection.

Incorporating STERIS products into every phase of your endoscope reprocessing will help improve department efficiency, reduce the risk of cross-contamination and provide reliable outcomes.

PRODUCT SPECIFICATIONS

Dimensions

700 mm (l) x 1400 mm (h) x 650 mm (d) - monitor included

Weight

75 kg

Electrical Requirements

100 V-240 V, 50-60Hz

Operating Room Temperature

5° - 40° C

Humidity

20 - 80%

Reprocessing Time

20 minutes (full cycle), 12 minutes only (HLD cycle),
25 minutes (double cleaning cycle)

Process Capacity


1 endoscope per cycle

High Level Disinfectant/Sterilant

ISASPOR™ SINGLE SHOT High-Level Disinfectant and Sterilant,
two parts, Peracetic acid and ISAZONE™ Molecule

Detergent

ISACLEAN™ AER Detergent and Decontaminant, Multienzymatic
PROTEAZONE™ PLUS Detergent and Cleaner, Enzymatic

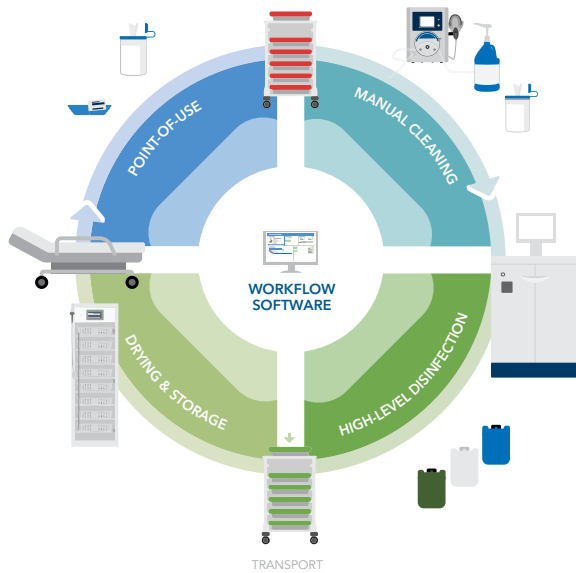

ORDERING INFORMATION

| PART NUMBER | DESCRIPTION | QUANTITY |
|-------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------|-------------------------|
| MEDISA3 | MEDIVATORS ISA Endoscope Reprocessor Washer-Disinfector using ISASPOR SINGLE SHOT High-Level Disinfectant and Sterilant 5L and ISACLEAN AER Detergent 1L | 1 x AER |
| MEDISA4 | MEDIVATORS ISA Endoscope Reprocessor Washer-Disinfector using ISASPOR SINGLE SHOT High-Level Disinfectant and Sterilant 10L and ISACLEAN AER Detergent 1L | 1 x AER |
| MEDISA6 | MEDIVATORS ISA Endoscope Reprocessor Washer-Disinfector using ISASPOR SINGLE SHOT High-Level Disinfectant and Sterilant 5L and ISACLEAN AER Detergent 5L | 1 x AER |
| MEDISA7 | MEDIVATORS ISA Endoscope Reprocessor Washer-Disinfector using ISASPOR SINGLE SHOT High-Level Disinfectant and Sterilant 10L and ISACLEAN AER Detergent 5L | 1 x AER |
| MEDISA8 | MEDIVATORS ISA Endoscope Reprocessor Washer-Disinfector using ISASPOR SINGLE SHOT High-Level Disinfectant and Sterilant 10L and PROTEAZONE PLUS Detergent 5L | 1 x AER |
| MEDISA9 | MEDIVATORS ISA Endoscope Reprocessor Washer-Disinfector using ISASPOR SINGLE SHOT High-Level Disinfectant and Sterilant 5L and PROTEAZONE PLUS Detergent 5L | 1 x AER |
| 500187 | ISASPOR SINGLE SHOT High-Level Disinfectant and Sterilant | 10L Part A + 10L Part B |
| 500188 | ISASPOR SINGLE SHOT High-Level Disinfectant and Sterilant | 5L Part A + 5L Part B |
| 4005635 | ISACLEAN AER Detergent and Decontaminant (EMEA 1) | 12 x 1L |
| 4005636 | ISACLEAN AER Detergent and Decontaminant (EMEA 2) | 12 x 1L |
| 4005637 | ISACLEAN AER Detergent and Decontaminant (EMEA 1) | 2 x 5L |
| 4006082 | ISACLEAN AER Detergent and Decontaminant (EMEA 2) | 2 x 5L |
| 500199 | PROTEAZONE PLUS Detergent and Cleaner | 2 x 5L |
| 11061031 | MEDIVATORS ISA Water Filters Kit (0,1µm + 0,45µm) | 1 |
| ISA4400 | ISA Air Filter | 1 |
| 4008833 | MEDIVATORS ISA Water Filter (0,1µm) | 1 |
| 4008832 | MEDIVATORS ISA Water Filter (0,45µm) | 1 |
| ISAROTS2 | ISA Printer Roll | 1 |

Contact your STERIS representative for more information or for the detergent and disinfectant combination available in your country.