


ERCP FineTome

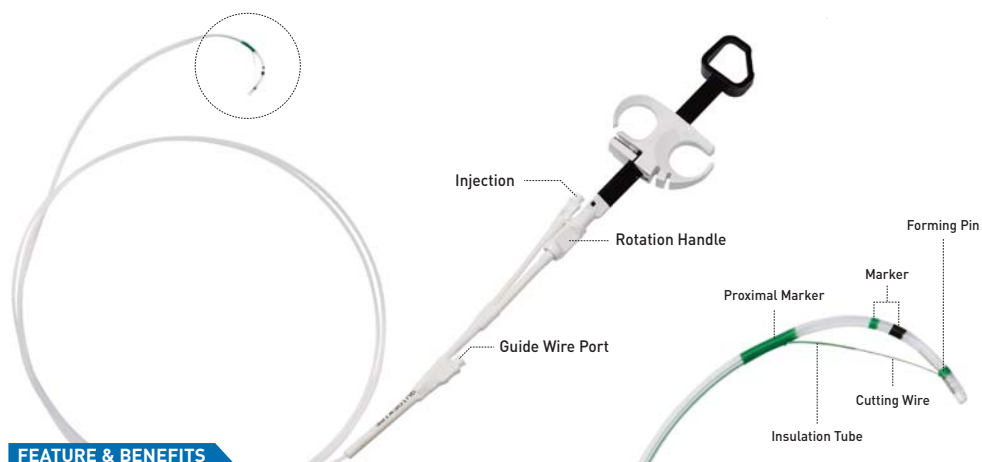

FEATURE & BENEFITS

(1) Enhancing safety by coating

- Using coating on the proximal end of the cutting wire, minimizes the damage to the surrounding tissue of targeted lesion and reduces the possibility of contact between conductive cutting wire and endoscope.

(2) Triple lumen sheath

- It will be used as guidewire, Cutting wire and injection of contrast media and the excellent configuration of 3 lumens facilitates efficient orientation of cutting wire for easy sphincterotomy.

(3) Tapered tip for smooth cannulation

- The distal end of sheath is tapered for easy cannulation into the papilla. It is also used for stricture and minor papilla.

(4) Pre-Curved distal for convenient cutting

- The distal end of this product is pre-curved. The construction of distal end and pre-curved shape facilitates easy positioning and insertion into the papilla.

SPECIFICATION

Model	Lumen	Tip Type	Distal Tip Size	Tip Length	Cutting Wire Length	Insulation Tube	Working Length	Number of Use	Validity of Sterilization
FM-EF0001	Triple Lumen	Tapered	4.4Fr	5mm	25mm	Yes	2,200mm	Single Use Only	3 years
FM-EF0006				5mm	25mm	No			

EUS ClearTip

