


# BD Stent

SX-ELLA Stent Esophageal Degradable BD



# BD Stent

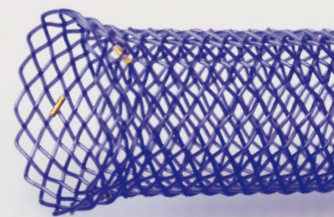

Biodegradation

Unique delivery system

Excellent flexibility

Radiopaque markers



## Basic description

### SX-ELLA Stent Esophageal Degradable BD – BD Stent

BD Stent is designed for treatment of benign esophageal strictures (peptic, anastomotic, caustic) and achalasia resistant to standard therapy.

#### Features / Benefits

- ☑ PULL delivery system / Splittable olive tip / Automatic detachment while stent starts releasing / Prevention of olive entrapment / Easy insertion and safe removal of delivery system
- ☑ Radiopaque markers / Accurate stent positioning
- ☑ Anti-migration design / Dual flared design / Reduced migration rates
- ☑ Integrity and radial force of stent is maintained for 6 – 8 weeks following implantation / Extended period of dilation compared to conventional methods
- ☑ Disintegration of stent occurs after 11 – 12 weeks following implantation / No need for removal
- ☑ Simple compression procedure / Compression tool available

#### Additional description

Stent is supposed to be compressed by compression tool into delivery system just prior to implantation  
5.9 mm / 9.4 mm (18 F / 28 F) delivery system, available in standard length of 75 cm  
Recommended guide wire: 0.035" (0.89 mm) ultra stiff / Minimum length 220 cm

Olive tip / Sheath




Splittable olive tip




Sheath handle / Pusher




Compression tool

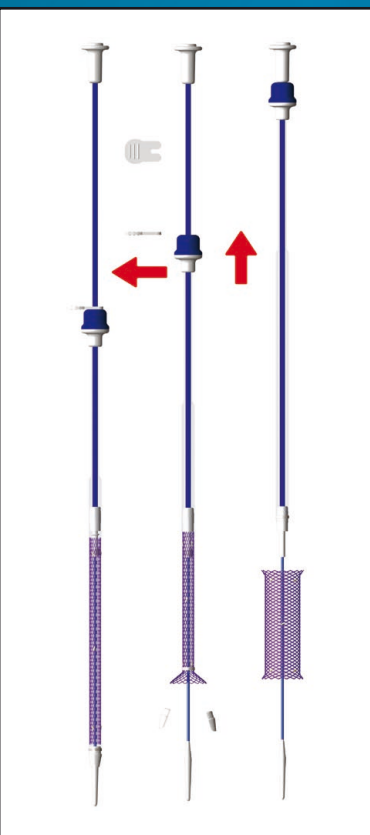

2 weeks post implantation

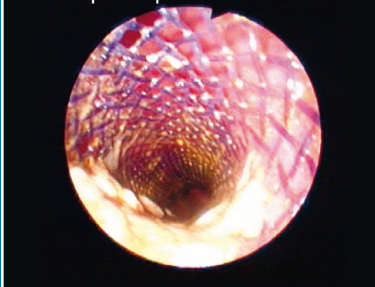

12 weeks post implantation

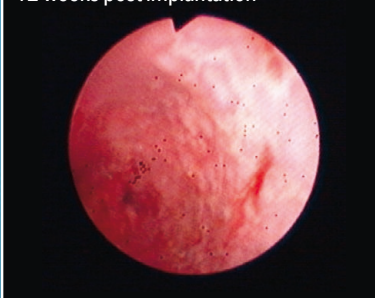

#### Available sizes

| REF No.              | Stent body diameter [mm] | Stent flares diameter [mm] | Stent length [mm] |
|----------------------|--------------------------|----------------------------|-------------------|
| 019-10A-23/18/23-060 | 18                       | 23                         | 60                |
| 019-10A-23/18/23-080 |                          |                            | 80                |
| 019-10A-23/18/23-100 |                          |                            | 100               |
| 019-10A-25/20/25-060 | 20                       | 25                         | 60                |
| 019-10A-25/20/25-080 |                          |                            | 80                |
| 019-10A-25/20/25-100 |                          |                            | 100               |
| 019-10A-28/23/28-060 | 23                       | 28                         | 60                |
| 019-10A-28/23/28-080 |                          |                            | 80                |
| 019-10A-28/23/28-100 |                          |                            | 100               |
| 019-10A-31/25/31-060 | 25                       | 31                         | 60                |
| 019-10A-31/25/31-080 |                          |                            | 80                |
| 019-10A-31/25/31-100 |                          |                            | 100               |
| 019-10A-31/25/31-135 |                          |                            | 135               |

Distributor:

Manufacturer:

**ELLA-CS, s.r.o.**

Milady Horakove 504, 500 06 Hradec Kralove  
Czech Republic

Tel.: +420 495 279 111 Fax: +420 495 265 655  
E-mail: info@ellacs.eu URL: www.ellacs.eu