



TRANSDUCERS

For Ultrasound Systems



CANON MEDICAL SYSTEMS CORPORATION

<https://global.medical.canon>

©Canon Medical Systems Corporation 2019-2025. All rights reserved.
Design and specifications are subject to change without notice.
MCAUS0366EAD 2025-01 CMSC/Produced in Japan

Canon Medical Systems Corporation meets internationally recognized
standards for Quality Management System ISO 9001, ISO 13485.

Canon Medical Systems Corporation meets the
Environmental Management System standard ISO 14001.

Verza, Infiniti Plus, AccuSITE, Ultra-Pro II and BX2 are trademarks of CIVCO Medical Solutions.

Disclaimer: Some transducers and related products presented in this brochure may not be
commercially available on all systems shown or may require the purchase of additional op-
tions. Indication for use depends on the regulatory requirements of each country. Please
contact your local Canon Medical representative for details.

Made For life

TRANSDUCERS



Aplio i700

iBeam Forming



Multi-Sync Pulser

Able to drive simultaneously pure waveforms with high precision creating a fine beam focus in the depth direction.




Multi-Beam Receiver

Enable to receive multiple simultaneous signals at high speed across a wide signal area, resulting in improved lateral and temporal resolution combined with higher frame rates.



Multi-Harmonic Compounding

Merges signals from individual beams with overlapping data from adjacent beams, virtually transforming each of them into thin, straight homogeneous beams, creating higher spatial resolution, contrast resolution and less artifacts.

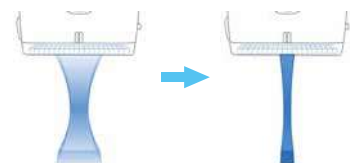

iBeam Slicing



intelligent Dynamic Micro-Slice (iDMS)

Creates a thin slice beam with continuous focus from near field to far field by controlling the main and side elements separately.



intelligent
DMS


Aplio i800

Volume Matrix for Radiology



Real time 3D transducer **Matrix array**

Provides high speed and precise volume imaging, and is compatible with SMI, CEUS, biopsy and full range of sensors and adapters

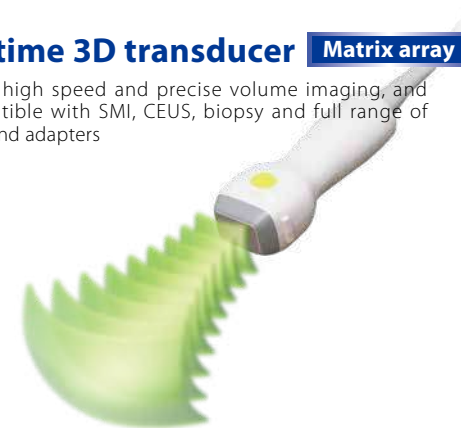

Volume Matrix

Aplio i900 **CV**

Volume Matrix for Cardiology




3D TEE (Transesophageal) transducer **Matrix array**

Acquires high resolution volume data in a short time, with high resolution, high volume rate and multi-plane display.



Trans-thoracic real time 3D transducer

Provides wide area of real time 3D images and 3D color Doppler images.

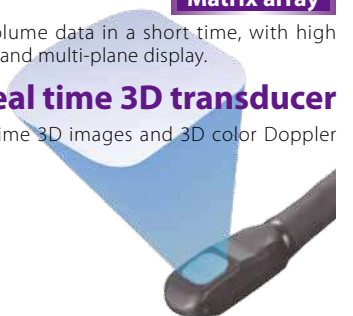

Volume Matrix

Aplio i-series

Prism Edition

Aplio i900 / Prism Edition v8.5 **CV**
Aplio i800 / Prism Edition v8.5
Aplio i700 / Prism Edition v8.5



Sector

PST-25BT

1	2.5	1.7~4.7
2	Approx. 90°	-
3	Cardiac, pediatric, abdominal, adult cephalic*, neonatal cephalic*	
4	680-106	680-086 610-608
5	-	-
6	-	-
7	-	-

5S1

PST-28BT

1	3.0	1.5~6.0
2	Approx. 120°	-
3	Cardiac, pediatric, abdominal, adult cephalic*, neonatal cephalic*	
4	-	-
5	-	-
6	UAFS-009A	-
7	-	-

Single Crystal

6S1

PSI-28VX

1	3.0	1.8~5.0
2	Approx. 120°	-
3	Cardiac, pediatric, abdominal, adult cephalic*, neonatal cephalic*	
4	-	-
5	-	-
6	-	-
7	-	-

Single Crystal

i6SVX1

PST-30BT

1	3.0	1.7~5.2
2	Approx. 90°	-
3	Cardiac, abdominal, adult cephalic*, neonatal cephalic*, pediatric	
4	680-106	680-086 610-608
5	-	-
6	-	-
7	-	-

5S2

PSI-30BX

1	3.0	1.5~6.0
2	Approx. 120°	-
3	Cardiac, pediatric, abdominal, adult cephalic*, neonatal cephalic*, thoracic/pleural	
4	-	-
5	-	-
6	UAFS-009A	-
7	-	-

Single Crystal

i6SX1

PSI-30VX

1	3.0	1.8~5.2
2	Approx. 120°	-
3	Cardiac, pediatric, abdominal, adult cephalic*, neonatal cephalic*	
4	680-156	680-157 610-1500/1502
5	-	-
6	UAFS-009A	-
7	-	-

Single Crystal

i6SVX1

PSI-40VX

1	4.0	1.8~6.0
2	Approx. 120°	-
3	Abdominal, fetal, pediatric, cardiac	
4	680-156	680-157 610-1500/1502
5	-	-
6	UAFS-009A	-
7	-	-

Single Crystal

i6SVX2

PST-50BT

1	5.0	3.0~8.2
2	Approx. 90°	-
3	Cardiac, pediatric, abdominal, adult cephalic*, neonatal cephalic*	
4	-	-
5	-	-
6	-	-
7	-	-

6S3

PSI-50VX

1	5.0	2.2~6.0
2	Approx. 90°	-
3	Cardiac, fetal, pediatric	
4	-	-
5	-	-
6	UAFS-009A	-
7	-	-

Single Crystal

i7SVX2

PST-65BT

1	7.0	3.5~12.0
2	Approx. 120°	-
3	Cardiac, pediatric	
4	-	-
5	-	-
6	-	-
7	-	-

12S4

PSI-70BT

1	7.0	3.5~9.0
2	Approx. 120°	-
3	Cardiac, pediatric, abdominal, adult cephalic*, neonatal cephalic*	
4	-	-
5	-	-
6	-	-
7	-	-

i10S4

PVT-375BT

1	3.5	1.5~6.0
2	Approx. 70°	Approx. 50 mm
3	Abdominal, fetal, pediatric	
4	680-105	680-085 610-608
5	-	-
6	Included with UIFR-AI900A	
7	-	-

6C1

PVT-375SC

1	3.5	1.5~6.0
2	Approx. 70°	Approx. 50 mm
3	Abdominal, fetal, pediatric	
4	680-105	680-085 610-608
5	-	-
6	Included with UIFR-AI900A	
7	-	-

Single Crystal

6Cs1

PVT-382BT

1	3.5	2.0~5.0
2	Approx. 80°	Approx. 20 mm
3	Abdominal, fetal, pediatric	
4	680-111	680-110 680-112
5	UAGV-027A	-
6	UAFS-001A	-
7	-	-

6MC1

PVI-475BT

1	4.0	1.8~6.4
2	Approx. 70°	Approx. 50 mm
3	Abdominal, fetal, pediatric	
4	680-145	680-144 610-608
5	680-155	680-154 610-1500
6	UAFS-007A (included with UIFR-AI900A)	-
7	-	-

Single Crystal

i8C1

PVI-475BX

1	4.0	1.8~6.4
2	Approx. 70°	Approx. 50 mm
3	Abdominal, fetal, pediatric, thoracic/pleural	
4	680-143	680-142 610-608
5	680-153	680-152 610-1500
6	UAFS-007A (included with UIFR-AI900A)	-
7	-	-

Single Crystal

i8CX1

PVT-482BT

1	4.0	1.8~6.0
2	Approx. 80°	Approx. 20 mm
3	Abdominal, fetal, pediatric	
4	-	-
5	-	-
6	UAFS-001A	-
7	-	-

Single Crystal

8MC1

PVI-482BX

1	4.0	1.8~6.0
2	Approx. 80°	Approx. 20 mm
3	Abdominal, fetal, pediatric	
4	-	-
5	-	-
6	UAFS-001A	-
7	-	-

Single Crystal

i8MCX1

PVT-574BT

1	5.0	2.0~9.5
2	Approx. 65°	Approx. 40 mm
3	Abdominal, fetal, pediatric	
4	680-168	680-169 610-1500
5	-	-
6	UAFS-008A	-
7	-	-

Single Crystal

10C1

PVI-574BX

1	5.0	2.0~9.5
2	Approx. 65°	Approx. 40 mm
3	Abdominal, fetal, pediatric	
4	680-166	680-167 610-1500
5	-	-
6	UAFS-008A	-
7	-	-

Single Crystal

i10CX1

PVT-674BT

1	6.0	3.5~9.7
2	Approx. 65°	Approx. 40 mm
3	Abdominal, fetal, pediatric	
4	-	-
5	UAGV-028A	-
6	-	-
7	-	-

10C3

PVT-712BT

1	7.0	4.3~11.0
2	Approx. 100°	Approx. 15 mm
3	Neonatal cephalic*, abdominal, pediatric	
4	680-127	680-126 610-608
5	-	-
6	-	-
7	-	-

11MC4

PLI-605BX

1	6.0	2.8~9.0
2	Approx. 46 mm	-
3	Abdominal, small organs, peripheral vascular, musculoskeletal	
4	680-174	680-175 610-1500
5	680-176	680-177 610-1017 610-1018 610-1019
6	UAFS-008A	-
7	-	-

Single Crystal

i9LX2

PLT-704SBT

1	7.5	4.0~10.5
2	Approx. 38 mm	-
3	Abdominal, small organs, peripheral vascular, musculoskeletal	
4	680-107	680-087 610-608
5	-	-
6	-	-
7	-	-

11L4

PLT-705BT

1	7.0	3.0~11.5
2	Approx. 45 mm	-
3	Abdominal, small organs, peripheral vascular, musculoskeletal	
4	680-129	680-128 610-608
5	-	-
6	-	-
7	-	-

11L3

PLI-705BX

1	7.0	3.5~11.5
2	Approx. 46 mm	-
3	Abdominal, peripheral vascular, small organs, musculoskeletal, thoracic/pleural	
4	680-150	680-151 610-1500
5	680-164	680-165 610-1017 610-1018 610-1019
6	UAFS-008A	-
7	-	-

i11LX3

PLT-1005BT

1	10.0	3.8~14.0
2	Approx. 58 mm	-
3	Abdominal, peripheral vascular, small organs, musculoskeletal, pediatric	
4	680-109	680-090 610-608
5	680-125	680-124 610-1231 610-1235
6	UAFS-004A	-
7	Water bag UAWB-022A	-

14L5

PLT-1202BT

1	12.0	4.5~17.0
2	Approx. 25.6 mm	-
3	Peripheral vascular, small organs, intraoperative, musculoskeletal, pediatric	
4	-	-
5	-	-
6	UAFS-010A	-
7	-	-

17LH7

PLT-1204BT

1	12.0	4.5~18.1
2	Approx. 38 mm	-
3	Abdominal, peripheral vascular, small organs, musculoskeletal	
4	680-102	680-088 610-608
5	680-123	680-122 610-1231 610-1235
6	680-170	680-171 610-1017 610-1018 610-1019
7	UACK-004A	-

18L7

1	Frequency (MHz)	Displayed range of freq. (MHz)
2	Field of view / angle	Curvature
3	Indication for use	
4	Starter Kit*	Bracket
5	Reusable biopsy adapter	Needle guide
6	Fusion adapter	
7	Water bag / Sono coupler	

*1 Includes the bracket and the sterile procedure kit.

*2 Not available in the USA and Canada.

*3 Described as "TCD" on Exam types.

*4 Described as "Neonatal head" on Exam types.

Aplio i-series

Prism Edition

Aplio i900 / Prism Edition v8.5 **CV**

Aplio i800 / Prism Edition v8.5

Aplio i700 / Prism Edition v8.5



Linear

PLI-1205BX

1	12.0	4.0~18.3
2	Approx. 46 mm	-
3	Abdominal, peripheral vascular, small organs, musculoskeletal, pediatric	
4	680-139 680-149 680-132	610-608 610-1500 610-1017 610-1018 610-1019
5	-	-
6	UAFS-008A	-
7	UACK-004A	-

PLI-2002BT

1	20.0	8.8~22.0
2	Approx. 25.6 mm	-
3	Peripheral vascular, small organs, intraoperative, musculoskeletal	
4	-	-
5	-	-
6	UAFS-010A	-
7	-	-

PLI-2004BX

1	20.0	8.8~24.0
2	Approx. 41 mm	-
3	Peripheral vascular, small organs, musculoskeletal, pediatric	
4	680-137 680-147 680-134	610-608 610-1500 610-1017 610-1018 610-1019
5	-	-
6	UAFS-008A	-
7	UACK-004A	-

PLI-3003BX

1	30.0	10.0~33.0
2	Approx. 31 mm	-
3	Musculoskeletal	
4	-	-
5	-	-
6	UAFS-008A	-
7	-	-

PVT-781VT

1	7.0	3.6~11.0
2	Approx. 180°	Approx. 10 mm
3	Abdominal, transrectal, transvaginal	
4	-	680-121
5	UAGV-035A	-
6	UAFS-003A	-
7	-	-

PVT-781VTE

1	7.0	3.6~11.0
2	Approx. 180°	Approx. 10 mm
3	Abdominal, transrectal, transvaginal	
4	-	680-121
5	UAGV-035A	-
6	UAFS-006A	-
7	-	-

PVL-715RST

1	7.5	4.5~9.0
2	Approx. 170°	Approx. 10 mm
3	Abdominal, transrectal	
4	-	-
5	UAGL-001AHA UAGL-002AHA	-
6	UAFS-005A	-
7	-	-

PVT-770RT

1	7.0	4.7~10.0
2	Approx. 170°	Approx. 10 mm
3	Transrectal, transvaginal	
4	-	680-093 680-092
5	UAGV-023A	-
6	-	-
7	-	-

PLT-308BTP*2

1	3.75	2.0~5.5
2	Approx. 85 mm	-
3	Abdominal	
4	-	-
5	UAGL-004AHA (included with transducer)	-
6	-	-
7	-	-

PVT-350BTP*2

1	3.5	2.0~4.5
2	Approx. 70°	Approx. 50 mm
3	Abdominal	
4	-	-
5	UAGV-030A (included with transducer)	-
6	UAFS-002A	-
7	-	-

PVI-450BXP*2

1	4.0	1.8~6.2
2	Approx. 70°	Approx. 50 mm
3	Abdominal, fetal, pediatric	
4	-	680-172
5	UAGV-037A UAGV-038A	-
6	UAFS-011A	-
7	-	-

PVT-745BTV

1	7.0	3.2~10.0
2	Approx. 40°	Approx. 45 mm
3	Abdominal, small organs, intraoperative	
4	-	-
5	UAGV-029A	-
6	-	-
7	-	-

PVT-745BTF

1	7.0	3.0~10.0
2	Approx. 40°	Approx. 45 mm
3	Abdominal, small organs, intraoperative	
4	-	-
5	-	-
6	-	-
7	-	-

PVT-745BTH

1	7.0	3.0~10.0
2	Approx. 40°	Approx. 45 mm
3	Abdominal, small organs, intraoperative	
4	-	-
5	UAGV-034A	-
6	-	-
7	-	-

PLT-705BTF

1	7.0	3.8~8.4
2	Approx. 46 mm	-
3	Abdominal, intraoperative	
4	-	-
5	UAGV020A	-
6	-	-
7	-	-

PLT-705BTH

1	7.0	3.8~8.4
2	Approx. 46 mm	-
3	Abdominal, intraoperative	
4	-	-
5	UAGV020A	-
6	-	-
7	-	-

PET-508MA

1	5.0	2.5~7.0
2	Approx. 90°	-
3	Cardiac (Manual TEE) (transesophageal)	
4	-	-
5	-	-
6	-	-
7	-	-

PET-512MB

1	5.0	3.0~8.5
2	Approx. 90°	-
3	Cardiac (Manual TEE) (transesophageal)	
4	-	-
5	-	-
6	-	-
7	-	-

PET-512MD

1	5.0	3.0~8.5
2	Approx. 90°	-
3	Cardiac (Motor drive TEE) (transesophageal)	
4	-	-
5	-	-
6	-	-
7	-	-

PEI-512VX

1	5.0	1.8~6.6
2	Approx. 90°	-
3	Cardiac (Electric scan 3D TEE) (transesophageal)	
4	-	-
5	-	-
6	-	-
7	-	-

PET-609MA*2

1	6.0	3.6~8.0
2	Approx. 90°	-
3	Cardiac (Manual TEE) (transesophageal)	
4	-	-
5	-	-
6	-	-
7	-	-

PET-835LA

1	8.0	3.5~9.0
2	Approx. 30°	-
3	Abdominal, laparoscopic	
4	-	-
5	-	-
6	-	-
7	-	-

PVT-675MVL

1	6.0	2.5~7.0
2	Approx. 66°	Approx. 50 mm
3	Abdominal, fetal, pediatric	
4	-	-
5	-	-
6	-	-
7	-	-

PVT-675MVS

1	6.0	2.5~7.5
2	Approx. 70°	Approx. 48 mm
3	Abdominal, fetal, pediatric	
4	-	-
5	-	-
6	-	-
7	-	-

PVT-681MVL

1	6.0	3.6~11.0
2	Approx. 180°	Approx. 10 mm
3	Abdominal, transvaginal, transrectal	
4	-	680-131
5	-	-
6	-	-
7	-	-

PC-20M

1	2.0	-
2	-	-
3	Cardiac, peripheral vascular, pediatric	
4	-	-
5	-	-
6	-	-
7	-	-

PC-50M

1	5.0	-
2	-	-
3	Cardiac, peripheral vascular, pediatric	
4	-	-
5	-	-
6	-	-
7	-	-

1	Frequency (MHz)	Displayed range of freq. (MHz)
2	Field of view / angle	Curvature
3	Indication for use	
4	Starter Kit*	Bracket Needle guide
5	Reusable biopsy adapter	
6	Fusion adapter	
7	Water bag / Sono coupler	

*1 Includes the bracket and the sterile procedure kit.
 *2 Not available in the USA and Canada.
 *3 Described as "TCD" on Exam types.
 *4 Described as "Neonatal head" on Exam types.

Aplio beyond^{V2.0}



Sector

PSU-28BT

1	3.0	1.5~6.0
2	Approx. 90°	-
3	Cardiac, abdominal, adult cephalic ^{*5} , neonatal cephalic ^{*6} , thoracic/pleural, pediatric	
4	-	-
5	-	-
6	-	-
7	-	-

Single Crystal
6S1

PSU-30BT

1	3.0	1.7~5.2
2	Approx. 90°	-
3	Cardiac, abdominal, adult cephalic ^{*5} , neonatal cephalic ^{*6} , thoracic/pleural, pediatric	
4	680-106	680-086 610-608
5	-	-
6	-	-
7	-	-

5S2

PSU-50BT

1	5.0	3.0~8.2
2	Approx. 90°	-
3	Cardiac, abdominal, adult cephalic ^{*5} , neonatal cephalic ^{*6} , pediatric	
4	-	-
5	-	-
6	-	-
7	-	-

6S3

PSU-70BT

1	7.0	3.5~9.0
2	Approx. 90°	-
3	Cardiac, abdominal ^{*4} , adult cephalic ^{*5} , neonatal cephalic ^{*6} , pediatric	
4	-	-
5	-	-
6	-	-
7	-	-

10S4

Convex

PVU-375BT

1	3.5	1.5~6.1
2	Approx. 70°	Approx. 50 mm
3	Abdominal, fetal, pediatric, thoracic/pleural, peripheral vascular	
4	680-105	680-085 610-608
5	-	-
6	Included with UHBE-MAGEN	
7	-	-

6C1

PVU-382BT

1	3.5	1.8~5.0
2	Approx. 80°	Approx. 20 mm
3	Abdominal, fetal, pediatric, thoracic/pleural, peripheral vascular	
4	680-111	680-110 680-112
5	UAGV-027A	
6	UAFS-001A	
7	-	-

6MC1

PVU-475BT

1	4.0	1.8~6.4
2	Approx. 70°	Approx. 50 mm
3	Abdominal, fetal, pediatric, thoracic/pleural, peripheral vascular	
4	680-155	680-154 610-1500
5	-	-
6	UHBE-MAGEN or UAFS-007A	
7	-	-

Single Crystal
8C1

PVU-475BTW

1	4.0	1.8~6.4
2	Approx. 70°	Approx. 50 mm
3	Abdominal, fetal, pediatric, thoracic/pleural, peripheral vascular	
4	680-153	680-152 610-1500
5	-	-
6	UHBE-MAGEN or UAFS-007A	
7	-	-

Single Crystal
8CW1

Convex

PVU-574BT

1	5.0	2.0~9.7
2	Approx. 65°	Approx. 40 mm
3	Abdominal, fetal, pediatric, thoracic/pleural	
4	680-168	680-169 610-1500
5	-	-
6	UAFS-008A	
7	-	-

Single Crystal
10C1

PVU-674BT

1	6.0	3.5~9.7
2	Approx. 65°	Approx. 40 mm
3	Abdominal, pediatric, peripheral vascular, fetal	
4	-	-
5	UAGV-028A	
6	-	-
7	-	-

10C3

PVU-712BTE

1	7.0	3.3~11.0
2	Approx. 100°	Approx. 15 mm
3	Abdominal, neonatal cephalic ^{*6} , pediatric	
4	-	-
5	-	-
6	-	-
7	-	-

11MC3

Linear

PLU-704BT

1	7.0	3.0~8.5
2	Approx. 38 mm	-
3	Abdominal ^{*4} , small organs, peripheral vascular, musculoskeletal, thoracic/pleural	
4	680-107	680-087 610-608
5	-	-
6	-	-
7	Water bag UAWB-021A	

11L4

Linear

PLU-805BT

1	8.0	3.0~12.0
2	Approx. 48.8 mm	-
3	Peripheral vascular, small organs, musculoskeletal, thoracic/pleural, abdominal, pediatric ^{*3}	
4	680-184	680-185 610-1500
5	-	-
6	680-186	680-187 610-1018 610-1019
7	-	-

12L4

PLU-1005BT

1	10.0	3.8~14.5
2	Approx. 58 mm	-
3	Small organs, peripheral vascular, musculoskeletal, abdominal, pediatric ^{*3}	
4	680-109	680-090 610-608
5	680-125	680-124 610-1231 610-1235
6	-	-
7	Water bag UAWB-022A	

14L5

PLU-1202BT

1	12.0	4.5~17.0
2	Approx. 25.6 mm	-
3	Peripheral vascular, small organs, intraoperative, musculoskeletal	
4	-	-
5	-	-
6	UAFS-010A	
7	-	-

17LH7

PLU-1204BT

1	12.0	4.5~18.0
2	Approx. 38 mm	-
3	Small organs, peripheral vascular, musculoskeletal	
4	680-102	680-088 610-608
5	680-123	680-122 610-1231 610-1235
6	680-170	680-171 610-1017 610-1018 610-1019
7	-	-

18L7

Linear

PLU-1204BX

1	12.0	4.5~18.0
2	Approx. 38 mm	-
3	Peripheral vascular, small organs, musculoskeletal	
4	680-102	680-088 610-608
5	680-123	680-122 610-1231 610-1235
6	680-170	680-171 610-1017 610-1018 610-1019
7	-	-

18LX7

PLU-2002BT

1	20.0	8.8~22.0
2	Approx. 25.6 mm	-
3	Peripheral vascular, small organs, musculoskeletal, intraoperative	
4	-	-
5	-	-
6	UACK-010A	
7	-	-

22LH8

Endocavity

PVU-781VTE

1	7.0	3.6~11.0
2	Approx. 180°	Approx. 10 mm
3	Transvaginal, transrectal	
4	-	-
5	UAGV-035A	
6	UAFS-006A	
7	-	-

11C3

Endorectal

PVL-715RS

1	7.5	4.5~9.0
2	Approx. 170°	Approx. 10 mm
3	Transrectal	
4	-	-
5	UAGL-001AHA UAGL-002AHA	
6	UAFS-005A	
7	-	-

11CL4

PLU-308BTP^{*2}

1	3.75	2.0~5.5
2	Approx. 85 mm	-
3	Abdominal	
4	-	-
5	UAGL-004AHA	
6	-	-
7	-	-

6LP3

PVU-350BTP^{*2}

1	3.5	2.0~4.5
2	Approx. 70°	Approx. 50 mm
3	Abdominal	
4	-	-
5	UAGV-030A	
6	UAFS-002A	
7	-	-

6CP1

1	Frequency (MHz)	Displayed range of freq. (MHz)
2	Field of view / angle	Curvature
3	Indication for use	
4	Starter Kit ^{*1}	Bracket Needle guide
5	Reusable biopsy adapter	
6	Fusion adapter	
7	Water bag / Sono coupler	

^{*1} Includes the bracket and the sterile procedure kit.
^{*2} Not available in the USA and Canada.
^{*3} Only available in the USA.
^{*4} Only applicable to MDR transducer.
^{*5} Described as "TCD" on Exam types.
^{*6} Described as "Neonatal head" on Exam types.

Aplio beyond^{V2.0}

Aplio a-series

Aplio a550 ^{V6.5}
Aplio a450 ^{V6.5}

Aplio a ^{V6.5}



Intraoperative

	1	7.0	3.0~10.0
	2	Approx. 40°	Approx. 45 mm
	3	Abdominal, small organs, intraoperative	
	4	-	-
	5	-	-
	6	-	-
	7	-	-

	1	7.0	3.8~8.4
	2	Approx. 46 mm	-
	3	Abdominal, intraoperative	
	4	-	-
	5	UAGV020A	
	6	-	-
	7	-	-

	1	6.0	3.6~11.0
	2	Approx. 180°	Approx. 10 mm
	3	Transvaginal, transrectal	
	4	-	680-131
	5	-	-
	6	-	-
	7	-	-

	1	2.5	1.7~4.7
	2	Approx. 90°	-
	3	Cardiac, pediatric, abdominal, adult cephalic ^{*4} , neonatal cephalic ^{*5}	
	4	680-106	680-086 610-608
	5	-	-
	6	-	-
	7	-	-

	1	7.0	3.5~12.0
	2	Approx. 120°	-
	3	Cardiac, pediatric	
	4	-	-
	5	-	-
	6	-	-
	7	-	-

	1	3.5	1.8~5.0
	2	Approx. 80°	Approx. 20 mm
	3	Abdominal, fetal, pediatric	
	4	680-111	680-110 680-112
	5	UAGV-027A	
	6	UAFS-001A	
	7	-	-

	1	7.0	3.0~10.0
	2	Approx. 40°	Approx. 45 mm
	3	Abdominal, small organs, intraoperative	
	4	-	-
	5	UAGV-034A	
	6	-	-
	7	-	-

	1	5.0	3.0~8.5
	2	Approx. 90°	-
	3	Cardiac (Motor drive TEE) (transesophageal)	
	4	-	-
	5	-	-
	6	-	-
	7	-	-

	1	2.0	-
	2	-	-
	3	Cardiac, peripheral vascular, pediatric	
	4	-	-
	5	-	-
	6	-	-
	7	-	-

	1	3.0	1.7~6.0
	2	Approx. 120°	-
	3	Cardiac, pediatric, abdominal, adult cephalic ^{*4} , neonatal cephalic ^{*5} , thoracic/pleural	
	4	-	-
	5	-	-
	6	UAFS-009A	
	7	-	-

	1	7.0	3.5~9.0
	2	Approx. 120°	-
	3	Cardiac, pediatric, abdominal, adult cephalic ^{*4} , neonatal cephalic ^{*5}	
	4	-	-
	5	-	-
	6	-	-
	7	-	-

	1	4.0	1.8~6.2
	2	Approx. 70°	Approx. 50 mm
	3	Abdominal, fetal, pediatric, thoracic/pleural	
	4	680-155	680-154 610-1500
	5	-	-
	6	UAFS-007A (Included with UIFR-AA550A)	
	7	-	-

	1	7.0	3.2~10.0
	2	Approx. 40°	Approx. 45 mm
	3	Abdominal, small organs, intraoperative	
	4	-	-
	5	UAGV-029A	
	6	-	-
	7	-	-

	1	8.0	3.5~9.0
	2	Approx. 30°	-
	3	Abdominal, laparoscopic	
	4	-	-
	5	-	-
	6	-	-
	7	-	-

	1	5.0	-
	2	-	-
	3	Cardiac, peripheral vascular, pediatric	
	4	-	-
	5	-	-
	6	-	-
	7	-	-

	1	3.0	1.7~5.2
	2	Approx. 90°	-
	3	Cardiac, abdominal, adult cephalic ^{*4} , neonatal cephalic ^{*5} , pediatric	
	4	680-106	680-086 610-608
	5	-	-
	6	-	-
	7	-	-

	1	3.5	1.5~6.1
	2	Approx. 70°	Approx. 50 mm
	3	Abdominal, fetal, pediatric	
	4	680-105	680-085 610-608
	5	-	-
	6	Included with UIFR-AA550A	
	7	-	-

	1	4.0	1.8~6.0
	2	Approx. 80°	Approx. 20 mm
	3	Abdominal, fetal, pediatric	
	4	-	610-1500/1502
	5	-	-
	6	UAFS-001A	
	7	-	-

	1	7.0	3.8~8.4
	2	Approx. 46 mm	-
	3	Abdominal, intraoperative	
	4	-	-
	5	UAGV020A	
	6	-	-
	7	-	-

	1	6.0	2.5~7.5
	2	Approx. 70°	Approx. 40 mm
	3	Abdominal, fetal, pediatric	
	4	-	-
	5	-	-
	6	-	-
	7	-	-

1	Frequency (MHz)		Displayed range of freq. (MHz)	
2	Field of view / angle		Curvature	
3	Indication for use			
4	Starter Kit ^{a1}	Bracket	Needle guide	
5	Reusable biopsy adapter			
6	Fusion adapter			
7	Water bag / Sono coupler			

^{*1} Includes the bracket and the sterile procedure kit.
^{*2} Not available in the USA and Canada.
^{*3} Only available in the USA.
^{*4} Only applicable to MDR transducer.
^{*5} Described as "TCD" on Exam types.
^{*6} Described as "Neonatal head" on Exam types.

	1	5.0	3.0~8.2
	2	Approx. 90°	-
	3	Cardiac, pediatric, abdominal, adult cephalic ^{*4} , neonatal cephalic ^{*5}	
	4	-	-
	5	-	-
	6	-	-
	7	-	-

	1	3.5	1.5~6.0
	2	Approx. 70°	Approx. 50 mm
	3	Abdominal, fetal, pediatric	
	4	680-105	680-085 610-608
	5	-	-
	6	included with UIFR-AA550A	
	7	-	-

1	Frequency (MHz)		Displayed range of freq. (MHz)	
2	Field of view / angle		Curvature	
3	Indication for use			
4	Starter Kit ^{a1}	Bracket	Needle guide	
5	Reusable biopsy adapter			
6	Fusion adapter			
7	Water bag / Sono coupler			

^{*1} Includes the bracket and the sterile procedure kit.
^{*2} Not available in the USA and Canada.
^{*3} Not compatible with Aplio a450.
^{*4} Described as "TCD" on Exam types.
^{*5} Described as "Neonatal head" on Exam types.

Aplio a-series

Aplio a550 V6.5
Aplio a450 V6.5

Aplio a V6.5



Convex

PVT-574BT

1	5.0	2.0~9.7
2	Approx. 65°	Approx. 40 mm
3	Abdominal, fetal, pediatric	
4	680-168	680-169 610-1500
5	-	-
6	UAFS-008A	
7	-	-

Single Crystal

10C1

PVT-674BT

1	6.0	3.5~9.7
2	Approx. 65°	Approx. 40 mm
3	Abdominal, fetal, pediatric	
4	-	-
5	UAGV-028A	
6	-	-
7	-	-

10C3

PVT-712BT

1	7.0	4.3~11.0
2	Approx. 100°	Approx. 15 mm
3	Neonatal cephalic*, abdominal, pediatric	
4	680-127	680-126 610-608
5	-	-
6	-	-
7	-	-

11MC4

Linear

PLT-704SBT

1	7.5	4.0~10.0
2	Approx. 38 mm	-
3	Abdominal, small organs, peripheral vascular, musculoskeletal, thoracic/pleural	
4	680-107	680-087 610-608
5	-	-
6	-	-
7	-	-

11L4

PLT-705BT

1	7.0	3.0~8.5
2	Approx. 45 mm	-
3	Cardiac, pediatric, abdominal, adult cephalic*, neonatal cephalic*, thoracic/pleural	
4	680-129	680-128 610-608
5	-	-
6	-	-
7	-	-

11L3

PLT-1005BT

1	10.0	3.8~14.0
2	Approx. 58 mm	-
3	Abdominal, peripheral vascular, small organs, musculoskeletal, pediatric	
4	680-109	680-090 610-608
4	680-125	680-124 610-1231 610-1235
5	-	-
6	UAFS-004A	
7	Water bag UAWB-022A	

14L5

PLT-1202BT

1	12.0	4.5~17.0
2	Approx. 25.6 mm	-
3	Peripheral vascular, small organs, intraoperative, musculoskeletal, pediatric	
4	-	-
5	-	-
6	UAFS-010A	
7	-	-

17LH7

PLT-1204BT

1	12.0	4.5~18.0
2	Approx. 38 mm	-
3	Abdominal, peripheral vascular, small organs, musculoskeletal	
4	680-102	680-088 610-608
4	680-123	680-122 610-1231 610-1235
4	680-170	680-171 610-1017 610-1018 610-1019
5	-	-
6	-	-
7	UACK-004A	

18L7

Linear

PLT-1204BX*3

1	12.0	4.5~18.0
2	Approx. 38 mm	-
3	Peripheral vascular, small organs, musculoskeletal	
4	680-102	680-088 610-608
4	680-123	680-122 610-1231 610-1017 610-1018 610-1019
5	-	-
6	-	-
7	-	-

18LX7

PLI-2002BT

1	20.0	8.8~22.0
2	Approx. 25.6 mm	-
3	Small organs, peripheral vascular, musculoskeletal, abdominal	
4	-	-
5	-	-
6	UAFS-010A	
7	-	-

i22LH8

Endocavity

PVT-781VT

1	7.0	3.6~11.0
2	Approx. 180°	Approx. 10 mm
3	Abdominal, transrectal, transvaginal	
4	-	-
5	UAGV-035A	
6	UAFS-003A	
7	-	-

11C3

PVT-781VTE

1	7.0	3.6~11.0
2	Approx. 180°	Approx. 10 mm
3	Abdominal, transrectal, transvaginal	
4	-	-
5	UAGV-035A	
6	UAFS-006A	
7	-	-

11C3

Endorectal

PVL-715RT

1	7.5	4.5~9.0
2	Approx. 170°	Approx. 10 mm
3	Abdominal, transrectal	
4	-	-
5	UAGL-001AHA UAGL-002AHA	
6	UAFS-005A	
7	-	-

11CL4

PVT-770RT

1	7.0	4.7~10.0
2	Approx. 170°	Approx. 10 mm
3	Transrectal, transvaginal	
4	-	-
5	UAGV-023A	
6	-	-
7	-	-

10C5

Biopsy

PLT-308BTP*2

1	3.75	2.0~5.5
2	Approx. 85 mm	-
3	Abdominal	
4	-	-
5	UAGL-004AHA (included with transducer)	
6	-	-
7	-	-

6LP3

PVT-350BTP*2

1	3.5	2.0~4.5
2	Approx. 70°	Approx. 50 mm
3	Abdominal	
4	-	-
5	UAGV-030A (included with transducer)	
6	UAFS-002A	
7	-	-

6CP1

Intraoperative

PVT-745BTV

1	7.0	3.2~10.0
2	Approx. 40°	Approx. 45 mm
3	Abdominal, small organs, intraoperative	
4	-	-
5	UAGV-029A	
6	-	-
7	-	-

11CI4

PVT-745BTF

1	7.0	3.0~10.0
2	Approx. 40°	Approx. 45 mm
3	Abdominal, small organs, intraoperative	
4	-	-
5	-	-
6	-	-
7	-	-

11CI4

PVT-745BTH

1	7.0	3.0~10.0
2	Approx. 40°	Approx. 45 mm
3	Abdominal, small organs, intraoperative	
4	-	-
5	UAGV-034A	
6	-	-
7	-	-

11CI4

PLT-705BTF

1	7.0	3.8~8.4
2	Approx. 46 mm	-
3	Abdominal, intraoperative	
4	-	-
5	UAGV020A	
6	-	-
7	-	-

11LI4

Intraoperative

PLT-705BTH

1	7.0	3.8~8.4
2	Approx. 46 mm	-
3	Abdominal, Intraoperative	
4	-	-
5	UAGV020A	
6	-	-
7	-	-

11LI4

Laparoscopic

PET-835LA

1	8.0	3.5~9.0
2	Approx. 30°	-
3	Abdominal, laparoscopic	
4	-	-
5	-	-
6	-	-
7	-	-

12CI4

Transesophagus

PET-508MA

1	5.0	2.5~7.0
2	Approx. 90°	-
3	Cardiac (Manual TEE) (transesophageal)	
4	-	-
5	-	-
6	-	-
7	-	-

7S3

1	Frequency (MHz)		Displayed range of freq. (MHz)
2	Field of view / angle		Curvature
3	Indication for use		
4	Starter Kit*	Bracket	Needle guide
5	Reusable biopsy adapter		
6	Fusion adapter		
7	Water bag / Sono coupler		

*1 Includes the bracket and the sterile procedure kit.
*2 Not available in the USA and Canada.
*3 Not compatible with Aplio a450.
*4 Described as "TCD" on Exam types.
*5 Described as "Neonatal head" on Exam types.

Aplio a-series

Aplio a550 V6.5
Aplio a450 V6.5

Aplio a V6.5

Aplio me V2.0



Transesophagus	PET-512MB *2		
		1 5.0 3.0~8.5	2 Approx. 90° -
		3 Cardiac (Manual TEE) (transesophageal)	4 -
	8S2	5 -	6 -
PET-512MD			7 -
		1 5.0 3.0~8.5	2 Approx. 90° -
		3 Cardiac (Motor drive TEE) (transesophageal)	4 - - -
	8SM2	5 -	6 -
		7 -	
PET-609MA *2			
		1 6.0 3.6~8.0	2 Approx. 90° -
		3 Cardiac (Manual TEE) (transesophageal)	4 -
	8S3	5 -	6 -
		7 -	
Mechanical 4D	PVT-675MVS		
		1 6.0 2.5~7.5	2 Approx. 70° Approx. 48 mm
		3 Abdominal, fetal, pediatric	4 - - -
	9CV2	5 -	6 -
PVT-681MVL			7 -
		1 6.0 3.6~11.0	2 Approx. 180° Approx. 10 mm
		3 Abdominal, transvaginal, transrectal	4 - - 680-131
	11CV3	5 -	6 -
		7 -	
Pencil	PC-20M		
		1 2.0 -	2 - -
		3 Cardiac, peripheral vascular, pediatric	4 - - -
	P2	5 -	6 -
PC-50M			7 -
Mechanical 4D		1 6.0 2.5~7.0	2 Approx. 66° Approx. 50 mm
		3 Abdominal, fetal, pediatric	4 - - -
	9CV2	5 -	6 -
		7 -	

1	Frequency (MHz)	Displayed range of freq. (MHz)
2	Field of view / angle	Curvature
3	Indication for use	
4	Starter Kit*1	Bracket Needle guide
5	Reusable biopsy adapter	
6	Fusion adapter	
7	Water bag / Sono coupler	

*1 Includes the bracket and the sterile procedure kit.
*2 Not available in the USA and Canada.
*3 Not compatible with Aplio a450.
*4 Described as "TCD" on Exam types.
*5 Described as "Neonatal head" on Exam types.

Sector	PSU-28BT		
		1 3.0 1.5~6.0	2 Approx. 90° -
		3 Cardiac, abdominal, adult cephalic*5, neonatal cephalic*6, thoracic/pleural, pediatric	4 - - -
	6S1	5 -	6 -
PSU-30BT			7 -
		1 3.0 1.7~5.2	2 Approx. 90° -
		3 Cardiac, abdominal, adult cephalic*5, neonatal cephalic*6, thoracic/pleural, pediatric	4 680-106 680-086 610-608
	5S2	5 -	6 -
		7 -	
	PSU-50BT		
		1 5.0 3.0~8.2	2 Approx. 90° -
		3 Cardiac, abdominal, adult cephalic*5, neonatal cephalic*6, pediatric	4 - - -
	6S3	5 -	6 -
PSU-70BT			7 -
		1 7.0 3.5~9.0	2 Approx. 90° -
		3 Cardiac, abdominal*4, adult cephalic*5, neonatal cephalic*6, pediatric	4 - - -
	10S4	5 -	6 -
		7 -	
Convex	PVU-375BT		
		1 3.5 1.5~6.1	2 Approx. 70° Approx. 50 mm
		3 Abdominal, fetal, pediatric, thoracic/pleural, peripheral vascular	4 680-105 680-085 610-608
	6C1	5 -	6 -
PVU-382BT			7 -
		1 3.5 1.8~5.0	2 Approx. 80° Approx. 20 mm
		3 Abdominal, fetal, pediatric	4 680-111 680-110 680-112
	6MC1	5 UAGV-027A	6 -
		7 -	
Convex	PVU-674BT		
		1 6.0 3.5~9.7	2 Approx. 65° Approx. 40 mm
		3 Abdominal, pediatric, peripheral vascular, fetal	4 - - -
	10C3	5 UAGV-028A	6 -
PVU-712BTE			7 -
		1 7.0 3.3~11.0	2 Approx. 100° Approx. 15 mm
		3 Abdominal, neonatal cephalic*6, pediatric	4 -
	11MC3	5 -	6 -
		7 -	
Linear	PLU-704BT		
		1 7.0 3.0~8.5	2 Approx. 38 mm -
		3 Abdominal*4, small organs, peripheral vascular, musculoskeletal, thoracic/pleural	4 680-107 680-087 610-608
	11L4	5 -	6 -
PVU-475BT			7 Water bag UAWB-021A
		1 4.0 1.8~6.4	2 Approx. 70° Approx. 50 mm
		3 Abdominal, fetal, pediatric, thoracic/pleural, peripheral vascular	4 680-155 680-154 610-1500
	8C1	5 -	6 -
		7 -	
	PVU-574BT		
		1 5.0 2.0~9.7	2 Approx. 65° Approx. 40 mm
		3 Abdominal, fetal, pediatric, thoracic/pleural	4 680-168 680-169 610-1500
	10C1	5 -	6 -
PVU-574BT			7 -
		1 5.0 2.0~9.7	2 Approx. 65° Approx. 40 mm
		3 Abdominal, fetal, pediatric, thoracic/pleural	4 680-168 680-169 610-1500
	10C1	5 -	6 -
		7 -	

1	Frequency (MHz)	Displayed range of freq. (MHz)
2	Field of view / angle	Curvature
3	Indication for use	
4	Starter Kit*1	Bracket Needle guide
5	Reusable biopsy adapter	
6	Fusion adapter	
7	Water bag / Sono coupler	

*1 Includes the bracket and the sterile procedure kit.
*2 Not available in the USA and Canada.
*3 Only available in the USA.
*4 Only applicable to MDR transducer.
*5 Described as "TCD" on Exam types.
*6 Described as "Neonatal head" on Exam types.

Aplio me

V2.0



Linear

PLU-805BT



1	8.0	3.0~12.0
2	Approx. 48.8 mm	-
3	Peripheral vascular, small organs, musculoskeletal, thoracic/pleural, abdominal, pediatric*3	
4	680-184	680-185 610-1500
	680-186	680-187 610-1017 610-1018 610-1019
5	-	-
6	-	-
7	UACK-004A	

12L4

PLU-1005BT



1	10.0	3.8~14.5
2	Approx. 58 mm	-
3	Small organs, peripheral vascular, musculoskeletal, abdominal, pediatric*3	
4	680-109	680-090 610-608
	680-125	680-124 610-1231 610-1235
5	-	-
6	-	-
7	Water bag UAWB-022A	

14L5

PLU-1202BT



1	12.0	4.5~17.0
2	Approx. 25.6 mm	-
3	Peripheral vascular, small organs, intraoperative, musculoskeletal	
4	-	-
5	-	-
6	-	-
7	-	-

17LH7

PLU-1204BT



1	12.0	4.5~18.0
2	Approx. 38 mm	-
3	Small organs, peripheral vascular, musculoskeletal	
4	680-102	680-088 610-608
	680-123	680-122 610-1231 610-1235
	680-170	680-171 610-1017 610-1018 610-1019
5	-	-
6	-	-
7	UACK-004A	

18L7

Linear

PLU-1204BX



1	12.0	4.5~18.0
2	Approx. 38 mm	-
3	Peripheral vascular, small organs, musculoskeletal	
4	680-102	680-088 610-608
	680-123	680-122 610-1231 610-1235
	680-170	680-171 610-1017 610-1018 610-1019
5	-	-
6	-	-
7	-	-

18LX7

PLU-2002BT



1	20.0	8.8~22.0
2	Approx. 25.6 mm	-
3	Peripheral vascular, small organs, musculoskeletal, intraoperative	
4	-	-
5	-	-
6	-	-
7	-	-

22LH8

Endocavity

PVU-781VTE



1	7.0	3.6~11.0
2	Approx. 180°	Approx. 10 mm
3	Transvaginal, transrectal	
4	-	-
5	UAGV-035A	
6	-	-
7	-	-

11C3

Endorectal

PVL-715RS



1	7.5	4.5~9.0
2	Approx. 170°	Approx. 10 mm
3	Transrectal	
4	-	-
5	UAGL-001AHA UAGL-002AHA	
6	-	-
7	-	-

11CL4

Biopsy

PLU-308BTP^{*2}



1	3.75	2.0~5.5
2	Approx. 85 mm	-
3	Abdominal	
4	-	-
5	UAGL-004AHA	
6	-	-
7	-	-

6LP3

PVU-350BTP^{*2}



1	3.5	2.0~4.5
2	Approx. 70°	Approx. 50 mm
3	Abdominal	
4	-	-
5	UAGV-030A	
6	-	-
7	-	-

6CP1

Intraoperative

PVU-745BTF



1	7.0	3.0~10.0
2	Approx. 40°	Approx. 45 mm
3	Abdominal, small organs, intraoperative	
4	-	-
5	-	-
6	-	-
7	-	-

11CI4

PVU-745BTH



1	7.0	3.0~10.0
2	Approx. 40°	Approx. 45 mm
3	Abdominal, small organs, intraoperative	
4	-	-
5	UAGV-034A	
6	-	-
7	-	-

11CI4

Intraoperative

PVU-745BTV



1	7.0	3.2~10.0
2	Approx. 40°	Approx. 45 mm
3	Abdominal, small organs, intraoperative	
4	-	-
5	UAGV-029A	
6	-	-
7	-	-

11CI4

PLU-705BTF



1	7.0	3.8~8.4
2	Approx. 46 mm	-
3	Abdominal, intraoperative	
4	-	-
5	UAGV020A	
6	-	-
7	-	-

11LI4

PLU-705BTH



1	7.0	3.8~8.4
2	Approx. 46 mm	-
3	Abdominal, intraoperative	
4	-	-
5	UAGV020A	
6	-	-
7	-	-

11LI4

Transesophagus

PET-512MD



1	5.0	3.0~8.5
2	Approx. 90°	-
3	Cardiac (Motor drive TEE) (transesophageal)	
4	-	-
5	-	-
6	-	-
7	-	-

8SM2

Laparoscopic

PET-835LA



1	8.0	3.5~9.0
2	Approx. 30°	-
3	Abdominal, laparoscopic	
4	-	-
5	-	-
6	-	-
7	-	-

12CI4

Mechanical 4D

PVU-674MVS



1	6.0	2.5~7.5
2	Approx. 70°	Approx. 40 mm
3	Abdominal, fetal, pediatric	
4	-	-
5	-	-
6	-	-
7	-	-

9CV2

PVU-681MVL



1	6.0	3.6~11.0
2	Approx. 180°	Approx. 10 mm
3	Transvaginal, transrectal	
4	-	-
5	-	-
6	-	-
7	-	-

11CV3

Pencil

PC-20M



1	2.0	-
2	-	-
3	Cardiac, peripheral vascular, pediatric	
4	-	-
5	-	-
6	-	-
7	-	-

P2

PC-50M



1	5.0	-
2	-	-
3	Cardiac, peripheral vascular, pediatric	
4	-	-
5	-	-
6	-	-
7	-	-

P5

1	Frequency (MHz)	Displayed range of freq. (MHz)	
2	Field of view / angle	Curvature	
3	Indication for use		
4	Starter Kit* ¹	Bracket	Needle guide
5	Reusable biopsy adapter		
6	Fusion adapter		
7	Water bag / Sono coupler		

*¹ Includes the bracket and the sterile procedure kit.
*² Not available in the USA and Canada.
*³ Only available in the USA.
*⁴ Only applicable to MDR transducer.
*⁵ Described as "TCD" on Exam types.
*⁶ Described as "Neonatal head" on Exam types.

Aplio flex V2.0

Aplio go V2.0



Sector

PSC-25LT



1	2.5	1.8~4.2
2	Approx. 90°	-
3	Cardiac, pediatric, abdominal, adult cephalic ² , neonatal cephalic ³ , thoracic/pleural	
4	-	-
5	-	-
6	-	-
7	-	-

c5S1

PSU-50BT



1	5.0	2.8~6.2
2	Approx. 90°	-
3	Cardiac, pediatric, neonatal cephalic ³	
4	-	-
5	-	-
6	-	-
7	-	-

6S3

Convex

PVC-375LT



1	3.5	1.8~6.0
2	Approx. 70°	Approx. 50 mm
3	Abdominal, fetal, pediatric, thoracic/pleural	
4	680-143	680-142 610-608
5	680-153	680-152 610-1500
6	-	-
7	-	-

c6C1

PVU-382BT



1	3.5	1.8~6.0
2	Approx. 80°	Approx. 20 mm
3	Abdominal, fetal, pediatric	
4	680-111	680-110 610-112
5	UAGV-027A	
6	-	-
7	-	-

6MC1

Convex

PVU-712BTE



1	7.0	3.3~11.0
2	Approx. 100°	Approx. 15 mm
3	Neonatal cephalic ³ , abdominal, pediatric	
4	-	-
5	-	-
6	-	-
7	-	-

11MC3

Linear

PLC-704LT



1	7.0	3.0~11.0
2	Approx. 38 mm	-
3	Peripheral vascular, small organs, musculoskeletal, thoracic/pleural	
4	680-137 680-136 610-608	680-147 680-146 610-1500
5	680-134 680-135 610-1017 610-1018 610-1019	
6	-	-
7	UACK-004A	

c11L4

PLC-1004LT



1	10.0	4.4~14.0
2	Approx. 38 mm	-
3	Peripheral vascular, small organs, musculoskeletal	
4	680-137 680-136 610-608	680-147 680-146 610-1500
5	680-134 680-135 610-1017 610-1018 610-1019	
6	-	-
7	UACK-004A	

c14L5

PLU-1204BT



1	12.0	4.4~18.0
2	Approx. 38 mm	-
3	Small organs, peripheral vascular, musculoskeletal	
4	680-102 680-088 610-608	680-123 680-122 610-1231 610-1235
5	680-170 680-171 610-1017 610-1018 610-1019	
6	-	-
7	UACK-004A	

18L7

Endocavity

PVC-781VLT



1	7.0	3.6~11.0
2	Approx. 180°	Approx. 10 mm
3	Transvaginal, transrectal, abdominal	
4	-	680-121
5	UAGV-035A	
6	-	-
7	-	-

c11C3

Mechanical 4D

PVU-674MVS



1	6.0	1.8~7.5
2	Approx. 70°	Approx. 40 mm
3	Abdominal, fetal, pediatric	
4	-	-
5	-	-
6	-	-
7	-	-

9CV2

PVU-681MVL



1	6.0	3.6~11.0
2	Approx. 180°	Approx. 10 mm
3	Transvaginal, transrectal	
4	-	680-131
5	-	-
6	-	-
7	-	-

11CV3

1	Frequency (MHz)	Displayed range of freq. (MHz)
2	Field of view / angle	Curvature
3	Indication for use	
4	Starter Kit ¹	Bracket Needle guide
5	Reusable biopsy adapter	
6	Fusion adapter	
7	Water bag / Sono coupler	

¹ Includes the bracket and the sterile procedure kit.
² Described as "TCD" on Exam types.
³ Described as "Neonatal head" on Exam types.



Better diagnostics start here

A complete set of lightweight transducers with proven ergonomic design, excellent clinical versatility and outstanding ease of use allows you to work efficiently and comfortably.



Disposable biopsy adapters

Disposable type needle guides			
680-092	680-121	680-131	680-172
			
1 Replacement Kit	1 Disposable Endocavity Needle Guide	1 Disposable Endocavity Needle Guide	1 Integrated Needle Guide
2 PVT-770RT	2 PVT-781VT, PVT-781VTE PVU-781VTE, PVC-781VLT	2 PVT-681MVL PVU-681MVL	2 PVI-450BXP
3 16~18	3 16~18	3 16~18	3 11,12,13,14,15,16,17,18,20,21,22
1 Product name			
2 Compatible transducer			
3 Applicable gauge size			

Needle guides (disposable parts)			
A 610-608	B 610-1500/1502*1	C1 610-1017	C2 610-1018
			
1 Ultra-Pro II™ Needle Guide Replacement Kit	1 Verza™ Guidance System Replacement Kit	1 AccuSITE™ Needle Guide Replacement Kit	1 AccuSITE™ Needle Guide Replacement Kit
2 8.5F, 14, 15, 16, 17, 18, 20, 21, 22, 23	2 7F, 8F, 9F, 10F, 12F, 14, 15, 16, 17, 18, 20, 21, 22, 25	2 18	2 20
3 -	3 -	3 0.5, 1.0, 1.5, 2.0, 2.5, 3.0, 3.5	3 0.5, 1.0, 1.5, 2.0, 2.5, 3.0, 3.5
D1 610-1231	D2 610-1235	E 680-112	
			
1 Infiniti Plus™ Needle Guide	1 Infiniti Plus™ Needle Guide	1 BX2™ Needle Guide	
2 18, 20, 21, 22, 25	2 12, 14, 16, 18	2 8.5F, 13, 14, 15, 16, 17, 18, 20, 21, 22, 23	
3 -	3 -	3 -	
Puncture angle using 610-1500/1502			
PSI-30VX			
PSI-40VX			
PVI-475BT			
PVI-475BX			
PVT-475BT			
PVU-475BT			
PVC-375LT			
PLI-705BX			
PLI-1205BX			
PLI-2004BX			
PLC-704LT			
PLC-1004LT			
PVT-482BT			
PVI-482BX			
PVI-574BX			
PVT-574BT			
PVU-574BT			
PLI-605BX			
PLU-805BT			
43°, 55°, 65°, 74°, 80°			
35°, 47°, 57°, 66°, 72°			
33°, 45°, 55°, 64°, 70°			
27°, 39°, 49°, 58°, 64°			
44°, 56°, 66°, 75°, 81°			
36°, 48°, 58°, 67°, 73°			
10°, 22°, 32°, 41°, 47°			
*1 610-1500 : Transducer Cover Size 14cmx91.5cm 610-1502 : Transducer Cover Size 8.9cmx91.5cm			
Biopsy brackets and needle guides manufactured by CIVCO Medical Solutions can be used together with Canon transducers. Use disposable biopsy adapters and needle guides when using these biopsy brackets. For details of the biopsy adapters and needle guides that can be used with each biopsy bracket, contact CIVCO Medical Solutions.			
1 Product name			
2 Applicable gauge size (F-French Size)			
3 Puncture depth(cm)			

Needle guides descriptions



Maximum versatility with a five-angle approach for improved anatomical access while also featuring an expanded gauge range. Its reduced blind zone provides high stability for improved accuracy. The needle angle is adjustable even when the transducer cover is on.



More flexibility and advanced features for increased productivity during central line placement procedures. This intuitive guide can be set for 7 different depths and is ready for use after just one click. The quick-release function ensures easy needle detachment.




A wide variety of gauge sizes offer multiple angle capabilities for shallow and deep access. Instrument positioning in the scan plane provides a safe and efficient mean for ultrasound-guided procedures.

Reusable biopsy adapters
Sono couplers

Reusable biopsy adapters					
UAGL-001AHA	UAGL-002AHA	UAGL-004AHA*	UAGL-015A	UAGV020A	UAGV-023A
					
1 PVL-715RS PVL-715RST	1 PVL-715RS PVL-715RST	1 PLT-308BTP PLU-308BTP	1 PLT-1204MV	1 PLT-705BTF, PLT-705BTH PLU-705BTF, PLU-705BTH	1 PVT-770RT
2 16	2 18	2 11~22 (Replaceable)	2 14~23 (Adjustable)	2 14~22	2 18, 19
UAGV-027A	UAGV-028A	UAGV-029A	UAGV-030A*	UAGV-032A	UAGV-033A
					
1 PVT-382BT PVU-382BT	1 PVT-674BT PVU-674BT	1 PVT-745BTV PVU-745BTV	1 PVT-350BTP PVU-350BTP	1 PVT-382MV	1 PVT-375MV
2 13~22 (Adjustable)	2 13~22 (Adjustable)	2 14~23 (Adjustable)	2 13~23 (Replaceable)	2 12~22 (Adjustable)	2 12~22 (Adjustable)
UAGV-034A	UAGV-035A	UAGV-037A	UAGV-038A		
					
1 PVT-745BTH PVU-745BTH	1 PVT-781VT, PVT-781VTE PVU-781VTE, PVC-781VLT	1 PVI-450BXP	1 PVI-450BXP		
2 14~23 (Adjustable)	2 15~22 (Replaceable)	2 11~22 (Replaceable)	2 11~22 (Replaceable)		
1 Compatible transducer					
2 Acceptable needle gauge size					
* Included with transducer					

Sono couplers

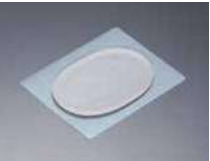Water bags (reusable)		Sono coupler (disposable)	
UAWB-021A	UAWB-022A	UACK-004A	
			
1 PLU-704BT	1 PLT-1005BT PLU-1005BT	1 PLI-1205BX, PLI-2004BX PLT-1204BT, PLU-1204BT PLU-805BT	
1 Compatible transducer			